



DDMT Manual

Version: 30/06/2021

This is an “in progress document” - incomplete



Contents

i.	System overview	2
ii.	DDMT overview	3
iii.	Some keywords	4
iv.	Setting up the Daily Diary logger for use	6
v.	Powering up the Daily Diary	11
vi.	Quick DDMT commands	12
0.	First run of DDMT - setting the graphics window dimensions	13
1.	Importing and Exporting Data with DDMT	15
	Extracting Daily Diary data from an SD card	15
	Loading data files into the software	15
	Setting the Working Directory	16
	The <i>COL</i> (column) data format for loading other manufacturer's acc data	17
	The data file, broken down by split number	18
	Exporting data / calibration data	18
	i. Save data in current split	18
	ii. Save data from one event to another in file	19
	iii. Save data by <i>Split</i> number	19
	What to save	21
	Raw (binary) data files	22
	Under-sampling to "thin" the data	22
	Exporting by marked events	22
	Merging DD raw files	23
	Merging csv / tab delimited files	24
2.	Main Graphs (2D) Windows	25
	Time indicator bars at the top and bottom of the main drawing windows	25
	Displaying data in "under-sampled mode" for improved application performance	26



Displaying data in “Visual averaging mode”	27
Two <i>2D graphing windows</i>	27
Navigating through a large dataset	28
Customising colour schemes of 2D graphed data	29
Horizontal guidelines	30
Values overlay	31
Highlighting variables	32
Adding items to the values overlay	32
Other display options	32
Adjustments to the visual vertical amplitude / offsets / controls of 2D graphs	35
i. Acceleration controls	35
ii. Magnetometry controls	36
iii. Temperature controls	37
iv. Pressure controls	37
v. Metrics – VeDBA/VeSBA/ODBA/VeSBA+ controls	38
vi. Metrics – Pitch / Roll / Heading controls	39
Display simple 2D compass on graphing window	39
vii. Metric – Variance controls	40
viii. GPS controls	41
ix. Metrics - Angular /Absolute angular velocity controls	42
x. Light / Battery	43
xi. Metric - Cumulative heading	43
xii. Analogue	43
xiii. Activity / SD card resets	44
xiv. Gridlines (for acceleration, geomagnetism, and metrics)	44
xv. Bitwise – multiple channels embedded within the one)	45
Zoom panel	46
Animal image to visually show animal orientation with respect to surge/sway/heave	47



Slide function for animation	48
3. Time / date, and acceleration / geomagnetism data corrections	49
Daily Diary on-board Real Time Clock – Setting the start time for your dataset	49
Magnetometry hard/soft iron corrections	50
Acceleration offsets and scaling controls	52
TDO (time, date, acceleration/magnetometer offset) files	55
Further optional data adjustments	56
Rotation correction	56
4. Differential Channels	58
Differential channel creation	58
Displaying the differential events values on the values overlay	59
Turning point finder	60
5. Marked Events and Bitwise Layers	63
Marked Events / controls	63
Bitwise Layers formation	66
Bitwise Layers in the Behaviour Builder	67
Bitwise Layers Export and Import of Global Bitwise Layer data	68
6. Bookmarks	70
Manually creating bookmarks on the 2D graphing windows	70
Categories	70
Renaming categories	71
Adjusting the vertical height of the bookmark bands	71
Save/load .bmk files	72
Export .raw files	72
Convert Bookmarks to Marked Events	72
Delete bookmarks	73
Additional bookmark options	73
Behaviour ethogram files	74



Multisession (data export)	77
Auto-create bookmarks	78
Auto-create bookmarks in split for Marked Events	79
7. Tilt-Compensated Compass	80
Setting the DD orientation and exporting magnetic heading	81
Pitch and Roll	83
8. Fast Fourier Transform (FFT)	84
9. Importing GPS Data (Pre-Loading)	87
Syncing GPS data to Daily Diary data	88
GPS straight-line interpolation	90
Bracketing	90
10. Behaviour Builder	92
Using Behaviour Builder to create Marked Events	93
Behaviour builder pre-smoothing	96
Behaviour Builder function – Variance	96
Gradient function	97
Blindspot function	97
Differential channels in the Behaviour Builder	98
Categories	98
Summary	98
Time series – Advanced Behaviour Builder	99
First example of Time series Behaviour builder	100
Second example of Time series Behaviour builder	102
Extend VALID to n events before next element	104
11. 3D Visualisation of Data	105
Enable 3D, and syncing with changes in data view on 2D graph window	105
Show 3D help text	106



Create a visualisation	107
Visualisation controls	108

Visualisation controls by type:

XY plot	113
XYZ plot	115
XYZ + Height plot	117
1D/2D histogram	119
Spherical histogram	121
Orientation Sphere	123
Rose plot	125
3D histogram	126
Point compass	128
Dead reckoning plot	129
Using the dead reckoning algorithm	137
Colourmap generator	139



System overview



Standard / original Daily Diary with connectors for battery/external sensors. Additional light sensor




Elongated Daily Diary with no connectors, power switch on/off. No light sensor. Battery / external sensors soldered to board.



Thumb Daily Diary with no connectors, no power switch on/off. No light sensor. Battery soldered to the board.

The Daily Diary is typically powered by a lithium cell (rechargeable or non-rechargeable), 3.6 – 4.2V; the cell should be capable of supplying peak current up to ~ 100 mA for the MicroSD card. This supply is regulated to 3.3 V, powering the CPU, MicroSD data storage, and the sensors Temp/Pressure/Acc/Mag (+Light on some variants). The Daily Diary monitors, and records, battery voltage, by using a potential divider to divide the incoming battery supply voltage by two, which is then read on one of the analogue inputs by the CPU and compared to the 3.3 V regulated supply.



The 32 kHz crystal is present to tick-over the Real Time Clock (RTC), and to regulate timing of sensor data collection.

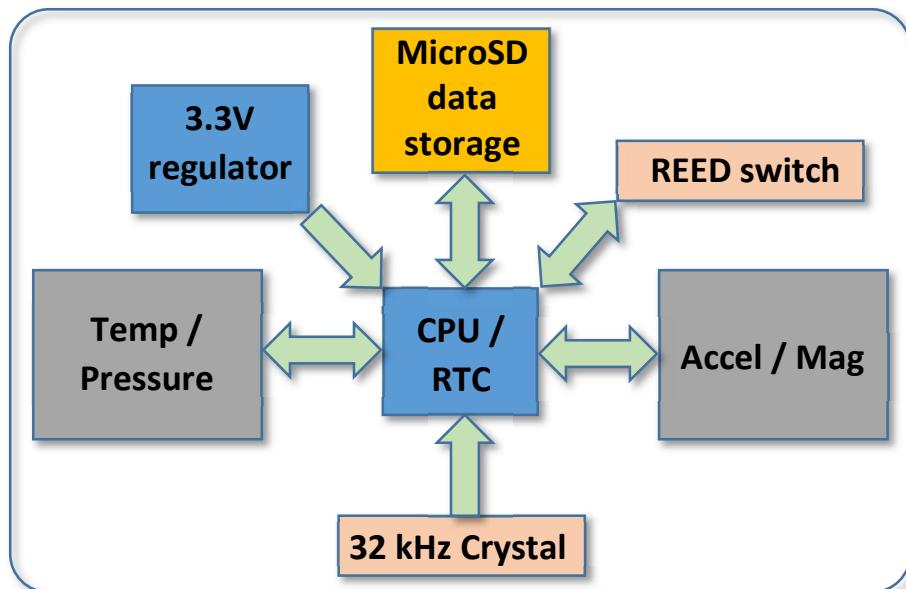

Typically, the accelerometer is the primary sensor, and so, if the user selects 40 Hz as the sensor collection frequency, acceleration X, Y, and Z axes data will be collected at this frequency. Magnetometry, used more for direction information, doesn't need to be at such a high frequency, and so is normally collected every ~4th acceleration data point. Temperature / Pressure are less important still, and due to their inability to sense quick changes in their environment (due to housings dampening their response time), their data is collected every 8th acceleration data point. The frequency ratio of acceleration to magnetometry to temperature/pressure is therefore approximately 40:13:4:4.



DDMT overview

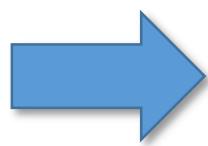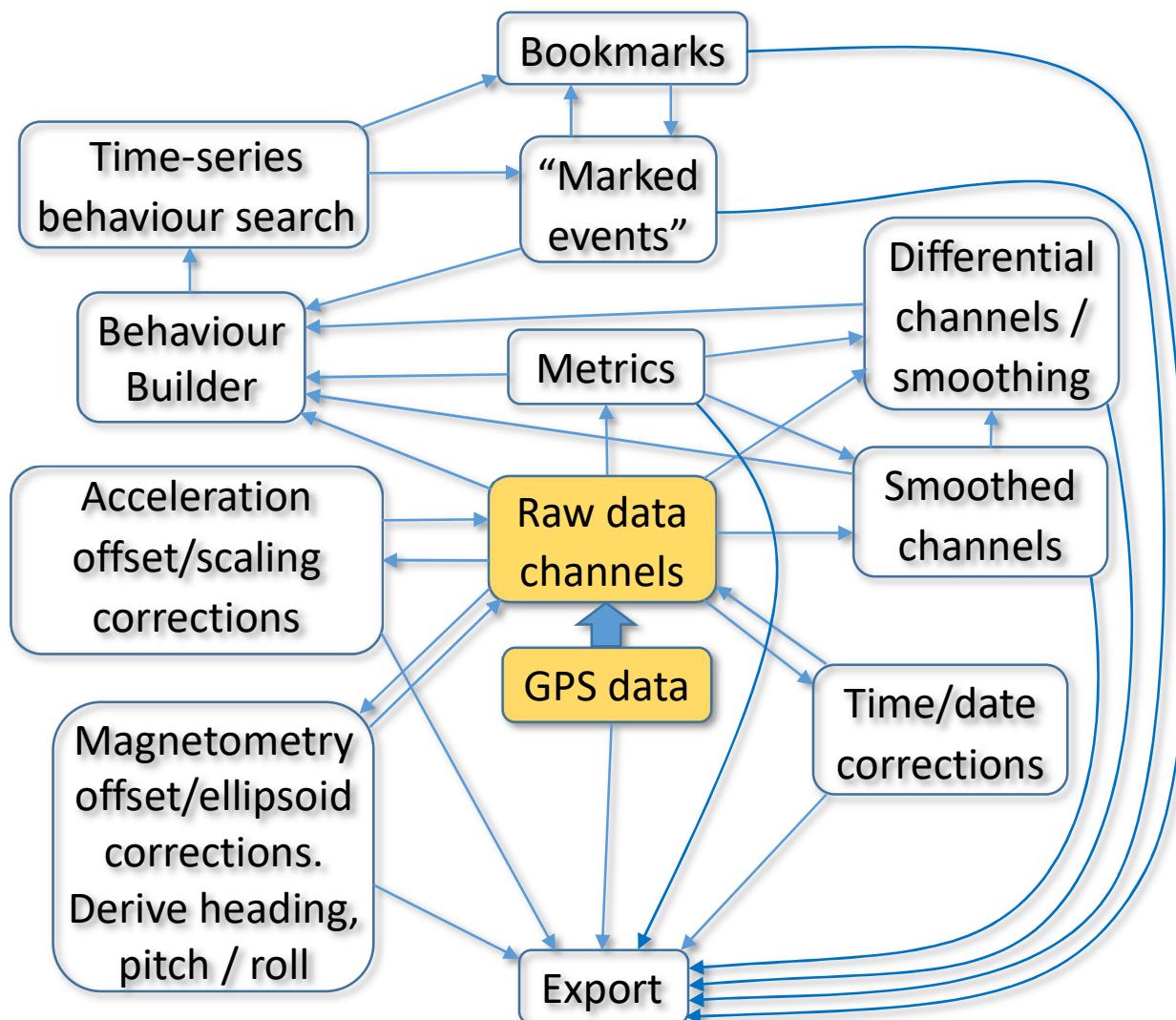

2D time graphs
3D plots

- Offsets corrects (acceleration / magnetometry) semi-auto/guided
- Time corrections (manual)
- Smoothed channels derived from primary raw data channels (user can select level of smoothing)
- Metrics created from both raw and smoothed channel data
- Smoothed channels of metrics can also be created (user can select level of smoothing)
- Differentials can be created from any channels / all (user can selectively control level of smoothing if required)
- All resulting channel data made available to the visualisation, search functions, and available for immediate export
- Any changes to level of smoothing or differential stepping has a ripple effect through the system to update visualisations



Some keywords

Event

In this software / manual, an **Event** is a data point in time where one or more sensor values exist within a large data array. For a logger capturing data at perhaps 40 Hz, this is the equivalent to 40 **Events** per second.

Split

A subset of the data file. The data that is currently in memory. The split size is chosen by the user prior to loading any data. The split size is then fixed for that session. Example below shows a data file containing 10,000,000 events, while the user has chosen to load a split size of 1,000,000 events



Marked Event


is where an **Event** has been marked with an integer value between 1 and 9 (inclusive) either as the result of a search function, or manually marked by the user. **Marked Events** are persistent for a session throughout the whole file i.e. moving to another split and returning, any **Marked Events** will still be present/in-view. The purpose of **Marked Events** is to:

- i. highlight data to be saved
- ii. to be used as part of a Boolean search function
- iii. or simply as a visual indicator of the meaning of a part of the data
- iv. to be used as a filter for displayed data on the 3D visualisation system
- v. to be used as a filter for the Dead Reckoning algorithm, built into the 3D visualisation system

Marked Events from a given session can be exported as a file and reloaded at a later session; see [Exporting Global Marked Events](#)

An **Event** can only be marked with a single **Marked Event** value (1-9). If it is marked again, the old **Marked Event** value will be replaced with the new value


Marked Events are discussed in later chapters where relevant

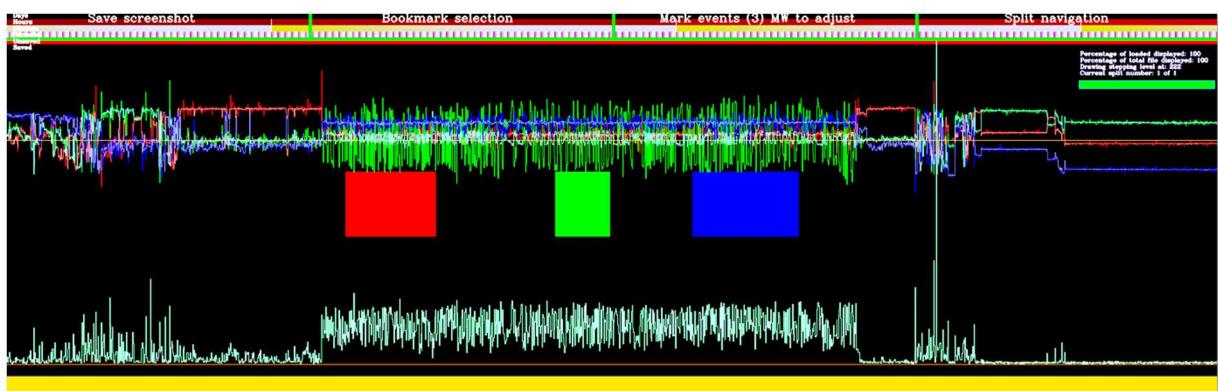

Marked events shown as highlighted regions in red, green, and blue (Marked Events values 1,2, and 3)

**Bookmark**


is a highlighted section of multiple events that has been noted either as the result of a search function, or manually noted by the user. **Bookmarks** can be assigned any 1 of 14 different “categories”; categories initially named Category A to Category N, the names of which can be updated by the user to give better meaning. **Bookmarks** from a given session can be exported as a file and reloaded at a later session. Unlike **Marked Events**, **Bookmarks** can overlap, although, on the display this will not be immediately apparent. **Bookmarks** can be used to export multiple sections of data as a single file, filtered by category if required, with numerous statistics associated with every **Bookmark**. **Bookmarks** are discussed in-depth in Chapter 6.

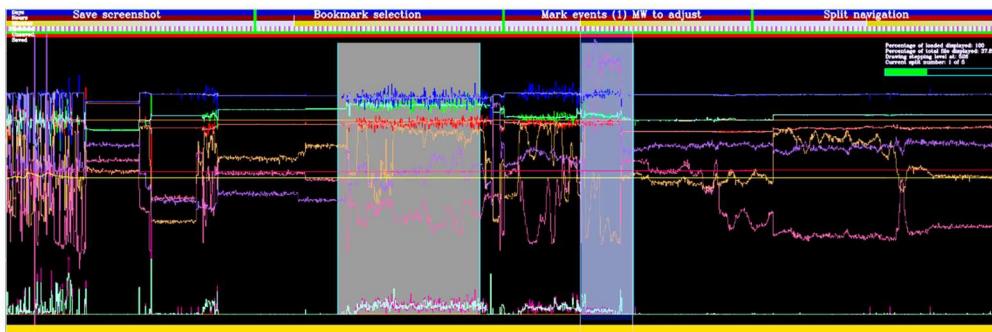

Example of a couple of bookmarks

Differential


is simply the rate of change of a channel of data. In the figure below, the bottom trace is **Acceleration Y** (smoothed), while the top trace is the differential of this, calculated, per point n , as $(x_{n+y/2} - x_{n-y/2})$ where x is the value at point n and y is *stepping difference* or the number of points over which to measure the rate of change i.e. for point n , the differential is calculated using half the stepping distance before/after

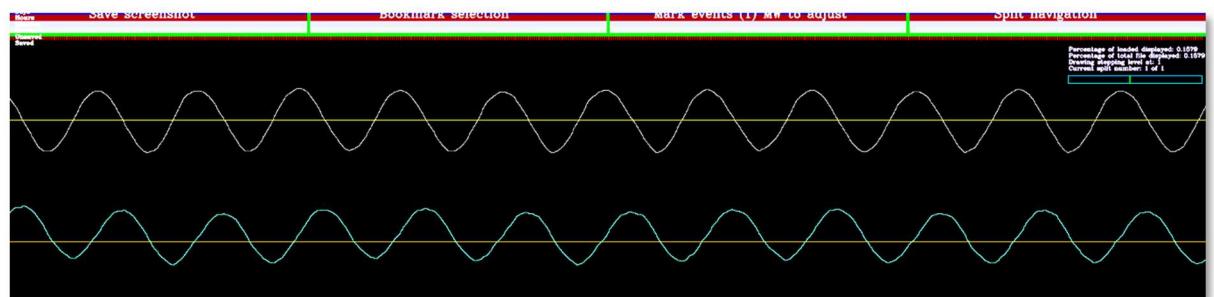

Example of a differential

Layers

are multiple bits of information ('1' or '0') packed into a single channel named **Bitwise**. Information may be embedded within a file read in, or layers within the **Bitwise** channel can be set '1' by copying (or merging) **Marked Events** across, or cleared '0'. Single or multiple **Layers** may be searched for using the **Behaviour Builder** search engine to generate **Marked Events / Bookmarks** to export data, or to act as a filter for the various 3D visualisations. The **2D graphing window** displays the layers vertically separated

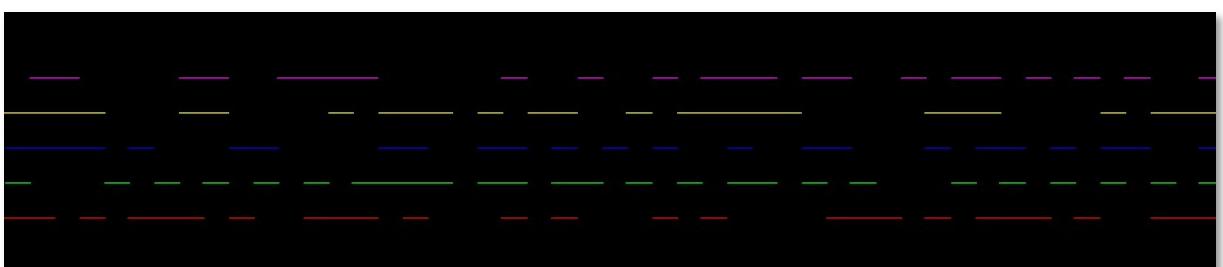




Setting up the Daily Diary logger for use

There are two types of storage with the Daily Diaries from Wildbyte Technologies. One type with on-board flash memory where an app is used to configure the device for deployment and to extract the data post-retrieval, and the other that uses MicroSD flash cards for large storage capacity, up to 32 GB.

It is recommended that a new MicroSD card at least be quick formatted before first use

2 GB MicroSD card logger

Start the 2GB application:

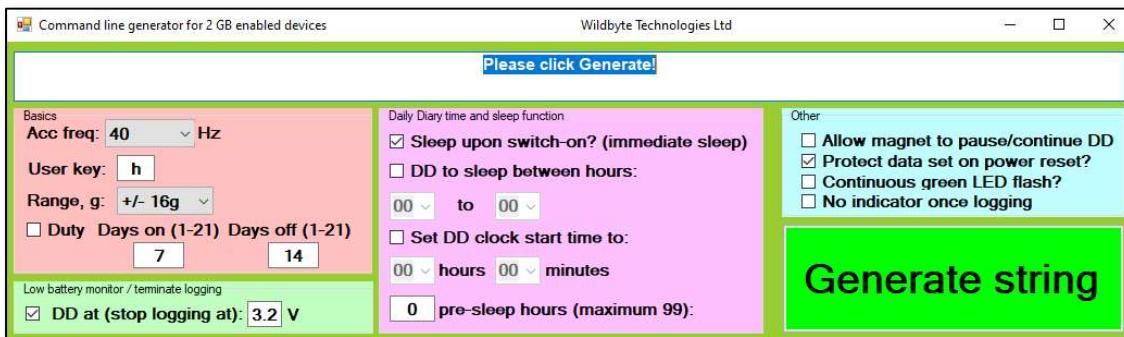

Select options (or just take the default – 40 Hz acceleration (13 Hz magnetometry, 4 Hz temperature and pressure), initial sleep (wake with a magnet)) and click **Generate string** to obtain the command string e.g.

AMTPXng,F40,R0144,Kh,s,ul,G060,s0000,Q030,I030,S300,b,t0000,P,S,g0000,gsd0000,grc0000,Dsd1985,T,I,S5,A30,t000,B20,T,p,p 00,R16,p00,d0000,t0000000000000000

On the MicroSD card (formatted, no files present), in Windows, right click in the empty root directory, and select **New text document**, and name it text. Open this text document and copy/paste the above command string into the file and save/close it. Use the **Safely remove hardware** option in Windows to eject the MicroSD from the laptop/computer to ensure the text document was saved correctly to the card. Insert into the Daily Diary and apply power. The options in the app interface are discussed below:

(Note that some older firmware versions may not have all these options present)

Basics

Acc freq This selects the logging frequency of the accelerometer (5, 10, 20, 30, 40, or 80 Hz). The accelerometer, magnetometer, temperature, and pressure sensors log with the ratio 40:13:4:4 (approximately), so if the accelerometer is set here to 10 Hz, then the magnetometer will sample at ~6 Hz, and the temperature/pressure sensors will log at ~ 2 Hz

User key All data on the MicroSD card is stored in sectors, with each sector holding information on timing, data, and also this key. When formatting in Windows, the default is to perform a Quick Format, that simply only erases the directory structure and not the actual data on the card. If one were to log for 2 hours, and then restart the logger and log for only one hour, then the data on the card would be read as a 2 hour log, with the first hour of data being the second logging session, and the second hour of data would be the partially overwritten first logging session. To get around this, one can either perform a Full Format, or simply change the **User key** to a character (any ASCII character on the keyboard i.e., a-z, A-Z, 0-9 etc. If the first logging session were with a key = 'h', and the second logging session with a key = 'g', then when the data is read into DDMT, it would see that there was a key change, and terminate reading the data at the end of the first hour with the key = 'g'

Range, g The accelerometer can be set to a range of +/- 2, 4, 7, or 16 g, with 2 g having the highest resolution, typically sub-millig, depending on the sensor on the board itself

Duty (cycle) With **Duty** the logger can be instructed to log for **X** days on and **Y** days off, switching at midnight



Low battery monitor / terminate logging

DD at X V

Daily Diaries are equipped with a simple voltage monitor circuit. When the voltage drops below the set point defined here, it will terminate logging, and stop writing data to the Micro SD card. This helps to prevent corruption of the SD card due to insufficient current as flash cards can require the order of “up to” 100 mA burst current to store data, but typically a lot less

Daily Diary time and sleep function

Sleep upon switch-on? (immediate sleep)

When the logger is first powered up, it will read the command string that has been placed by the user on the SD card (discussed later). If the device is set to sleep, it will do nothing else, and immediately enter sleep. The user must either swipe a magnet by the device to wake it (yellow LED flash), or with the square logger (with connectors and a button), press the button to begin logging (if a “pre-sleep” is not also configured)

DD to sleep between hours

The Daily Diary can be configured to sleep between specific hours (24 hr clock). If set to 04 and 08, the logger would sleep at 04:00, and continue logging from 08:00

Set DD clock start time to

The Daily Diary can be configured to set its time to the data defined here (hrs:mins). This time will not be stored in the clock until woken from sleep (if configured). The idea is that the clock can be configured to 10:00, and the device powered up at 09:45. When the user see their watch switch to 10:00, they would swipe a magnet by the device to set the clock and begin logging

Pre-sleep hours

The Daily Diary can be configured to sleep for **X** hours before logging begins. This might be for a reason such as not logging until the animal has recovered from any stress related to its capture/release

Allow magnet to pause/continue DD

Allows a magnet to swipe past the device and toggle the logging process on/off. The on-board clock (RTC time/date) continues to “tick” even when logging is suspended

Protect data set on power reset

If this option is selected and previously recorded data is present on the card, it will alert the user with a short and long yellow (/amber) LED flash indicating old data present. The Daily Diary will not initiate the logging sequence until a magnet is brought to its REED switch (magnetic switch) such that the yellow (/amber) LED illuminates continuously for 3 seconds and is then removed. For the square logger with a button present, this can instead be depressed for 3 seconds and released to acknowledge overwrite. Bear in mind that overwriting data on an SD card may result in “data overhang” from previous recording(s). It is ultimately better, for a true deployment, to either fully format the card, or to change the **User key** to a character that has not been used before with the SD card in question

Continuous green LED flash

If this is selected then, instead of the green LED intermittent flashing (indicates logging is ongoing) switching to dim flashing after a few minutes, it will continue to flash brightly the whole time. This will of course impact on battery life, but might be more useful when testing and not likely to influence animal behaviour in the wild (at night)


No indicator once logging

This will switch off the green LED from flashing at all when logging



8-32 GB MicroSD card logger

Start the 32GB application:

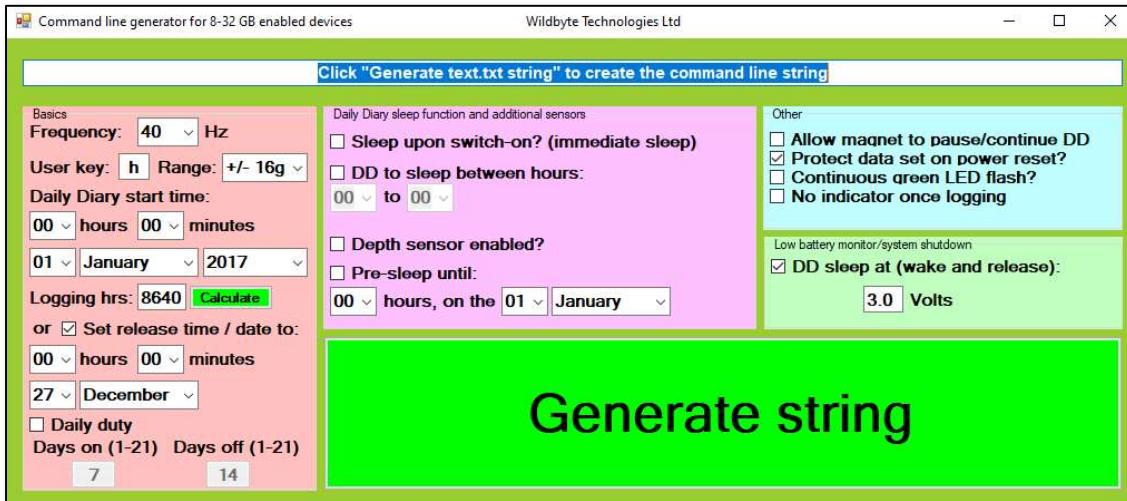

Select options (or just take the default – 40 Hz acceleration (13 Hz magnetometry, 4 Hz temperature and pressure), initial sleep (wake with a magnet)) and click **Generate string** to obtain the command string e.g.

AMTPXng,F40,R8640,Kh,s,ul,s0000,b,T000001012017,R000027120000,P,S,Dsd1861,T,I,B20,t,p,p00,R16,p000000,d0000

On the MicroSD card (formatted, no files present), in Windows, right click in the empty root directory, and select **New text document**, and name it “text”. Open this “text.txt” document and copy/paste the above command string into the file and save/close it. Use the **Safely remove hardware** option in Windows to eject the MicroSD from the laptop/computer to ensure the text document was actually saved correctly to the card. Insert into the Daily Diary and apply power. The options in the app interface are discussed below:

(Note that some older firmware versions may not have all of these options present)

Basics

Frequency This selects the logging frequency of the accelerometer (5, 10, 20, 30, 40, or 80 Hz). The accelerometer, magnetometer, temperature, and pressure sensors log with the ratio 40:13:4:4 (approximately), so if the accelerometer is set here to 10 Hz, then the magnetometer will sample at ~6 Hz, and the temperature/pressure sensors will log at ~ 2 Hz

User key All data on the MicroSD card is stored in sectors, with each sector holding information on timing, data, and also this key. When formatting in Windows, the default is to perform a Quick Format, that simply only erases the directory structure and not the actual data on the card. If one were to log for 2 hours, and then restart the logger and log for only one hour, then the data on the card would be read as a 2 hour log, with the first hour of data being the second logging session, and the second hour of data would be the partially overwritten first logging session. To get around this, one can either perform a Full Format, or simply change the **User key** to a character (any ASCII character on the keyboard i.e. a-z, A-Z, 0-9 etc. If the first logging session were with a key = ‘h’, and the second logging session with a key = ‘g’, then when the data is read into DDMT, it would see that there was a key change, and terminate reading the data at the end of the first hour with the key = ‘g’

Range, g The accelerometer can be set to a range of +/- 2, 4, 7, or 16 g, with 2 g having the highest resolution, typically sub- milli’ g, depending on the sensor on the board itself

Daily Diary start time

The Daily Diary can be configured to set its time to the data defined here (hrs:mins). This time will not be stored in the clock until woken from sleep (if configured). The idea is that the clock can be configured to 10:00, and the device powered up at 09:45. When the user see their watch switch to 10:00, they would swipe a magnet by the device to set the clock and begin logging. As this is the larger capacity device, the user can also set the date here.



This is not absolutely necessary as this can be set later in DDMT, but using this option fixes the date within the datafile itself

Logging hrs If the logger is only required to log for **X** hours, then enter the hours here, and click calculate, and the app will determine the stop time. This might be for testing where battery capacity is limited and required for multiple tests

Set release date/time to

Or enter manually the date/time here for it to terminate the logging

Daily duty With **Duty** the logger can be instructed to log for **X** days on and **Y** days off, switching at midnight

Daily Diary sleep function and additional sensor

Sleep upon switch-on?

When the logger is first powered up, it will read the command string that has been placed by the user on the SD card (discussed later). The device will first set the start time and date and if the device is set to sleep, it will do nothing else, and immediately enter sleep. The user must either swipe a magnet by the device to wake it (yellow LED flash), or with the square logger (with connectors and a button), press the button to begin logging

DD to sleep between hour

The Daily Diary can be configured to sleep between specific hours (24 hr clock). If set to 04 and 08, the logger would sleep at 04:00, and continue logging from 08:00

Depth sensor enabled

If a depth sensor (Keller) is connected to the logger, this needs to be enabled for the logger to include it in its sensor data log

Pre-sleep until The Daily Diary can be configured to wait until a specific hour, on a specific date, before it starts the logging process



Other

Allow magnet to pause/continue DD

Allows a magnet to swipe past the device and toggle the logging process on/off. The on-board clock (RTC time/date) continues to “tick” even when logging is suspended

Protect data set on power reset

If this option is selected and previously recorded data is present on the card, it will alert the user with a short and long yellow (/amber) LED flash indicating old data present. The Daily Diary will not initiate the logging sequence until a magnet is brought to its REED switch (magnetic switch) such that the yellow (/amber) LED illuminates continuously for 3 seconds and is then removed. For the square logger with a button present, this can instead be depressed for 3 seconds and released to acknowledge overwrite. Bear in mind that overwriting data on an SD card may result in “data overhang” from previous recording(s). It is ultimately better, for a true deployment, to either fully format the card, or to change the **User key** to a character that has not been used before with the SD card in question

Continuous green LED flash

If this is selected then, instead of the green LED intermittent flashing (indicates logging is ongoing) switching to dim flashing after a few minutes, it will continue to flash brightly the whole time. This will of course impact on battery life, but might be more useful when testing and not likely to influence animal behaviour in the wild (at night)

No indicator once logging

This will switch off the green LED from flashing at all when logging

Low battery monitor / terminate logging

DD at X V(volts) Daily Diaries are equipped with a simple voltage monitor circuit. When the voltage drops below the set point defined here, it will terminate logging, and stop writing data to the Micro SD card. This helps to prevent corruption of the SD card due to insufficient current as flash cards can require the order of “up to” 100 mA burst current to store data, but typically a lot less




Powering up the device

When power is applied to the Daily Diary, it will quickly read the command string within the text.txt file on the SD card and configures itself. If no initial sleep is enabled, then it will flash a few times, and if all is well, it will initiate logging. When it's logging data to the SD card, it will flash the green LED once per ~40 data points, so if logging acceleration at 40 Hz, this will flash at approximately 1 Hz. Some firmware will cause a bright green flash for the first 8-10 minutes and then switch to a very dim flash. A darkened room may be required to see this dim flash. If an error occurs on bootup such as no SD card, or the incorrect type (2 GB instead of 8-32 GB or vice versa), or it finds an error in the command string to stop its proper interpretation, it will flash the red LED. The number of flashes (repeated with pause) is indicative of the actual error

Sleep mode	A brief yellow(/amber) flash every 2-3 seconds
Logging mode	Brief green flash every 1-3 seconds, depending on the sampling frequency
Error	2 red flashes Can't find the SD card
	3 red flashes No text.txt file found on the SD card. If present on the card, then perform a quick format and recreate the text file on the card. A full format of any SD card prior to a proper deployment is recommended, rather than a quick format
	5 red flashes Command string sequence error
	6 red flashes Low battery – this is shown only once (non-repeated), so connecting a low battery, it begins logging, and then flashing the red LED 6 times and going to sleep (intermittent yellow (/amber) LED tells you low battery i.e. below the DD sleep at X Volts setting

To stop the device logging, either pop out the SD card (Thumb DD, just pull the card, other devices click in to eject), switch off power (Elongated 'long' DD), or just pull the battery connector (square logger).

It is **HIGHLY** recommended that when taking the SD card out of the logger and putting it into a PC/laptop to extract the data, that the SD card be inserted into an SD card adapter with a "Lock" set to lock. This prevents Microsoft Windows from overwriting any data on the SD card as the Daily Diary writes the data into "free space" on the card.

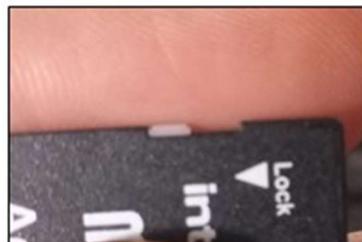

Please see chapter 1 on extracting data from the SD card using DDMT.



Quick DDMT commands

A selection of controls for use on the **2D graphing windows**

These controls are always a combination of keyboard buttons Shift, Ctrl, Alt, and mouse Left, Right, and Centre buttons, including the mouse wheel

What you want to do	How to do it
Bring up data point values overlay and probe values at specific time points	Click middle mouse wheel button Left and Right mouse click for “left / right white lines” to get data point values/ move mouse around to update overlay
Move “left / right white lines” to probe data points	Move mouse while holding down L or R click on 2D graphing window
Zoom into specific area within the data	Place left / right white lines around data of interest and then with Ctrl, Shift & Alt held down together, L click
Zoom out to view all data	Use L / R mouse click to drag the orange bar along the bottom of the drawing out
Zoom out a little	Use mouse wheel (stroke down) in the centre of the 2D graphing window to zoom in or out a little
Scroll through data	When zoomed in to some data within the split, hold down Alt + L click, and drag mouse left / right across the vertical centre of the 2D graphing window to drag the bottom orange bar
Jump quickly through data	When ‘Jump by half screen?’ is ticked on the display options tab, by holding down Ctrl and clicking L / R on the 2D graphing window , the data position will jump by half a screen to the side
Move position of the values overlay	Hold down Shift & L mouse button to move
Highlight (yellow) parameter(s) on the values overlay	Hold down Ctrl & Alt and move mouse over first (parameters description) column in the values overlay. Items will change colour to green. L mouse click to toggle highlight to yellow
Toggle display of “zoomed in” window of drawing	Ctrl, Shift & Alt held down together, then R mouse click
Move zoomed window	Hold down shift and use R mouse click to move
Adjust zoom level of zoom panel	Left Alt and mouse wheel
Moving to the next split	Left and right click simultaneously in the fourth panel across the top of the drawing; labelled “Split navigation”
Bookmarking data	Left and Right click to set parameters of the data you want on the drawing, then Left and right click in the second panel across the top



0. First run of DDMT - setting the graphics window dimensions

When first running DDMT, the user is presented with the control panel shown in Figure 0.1.

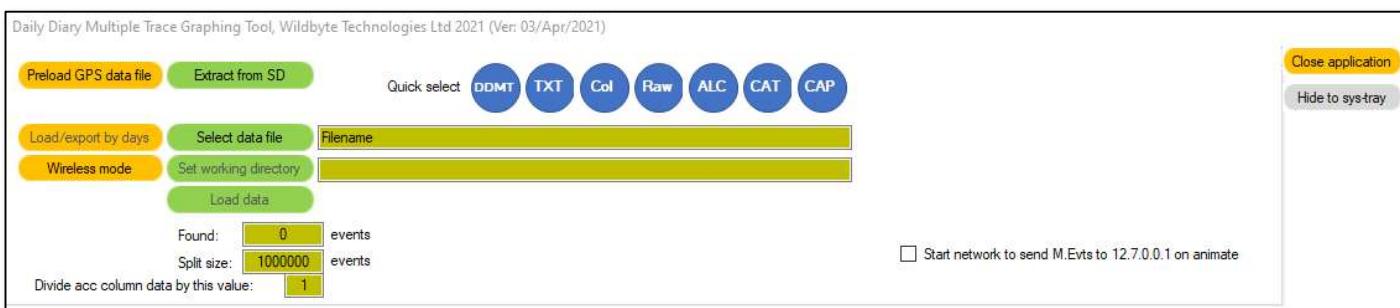


Figure 0.1 Initial control panel – setting the graphics window dimensions

By default, the **Main graphs** window dimensions will default to less than the screen resolution of the user's desktop PC / laptop resolution.

Note that a great many laptops will default to a non-standard scaling of the screen. For instance, a laptop having a full HD resolution of 1920x1080, might then be rescaled by Windows to 125%, or even 150%, and so a window having dimensions 1900x600 would be "blown up" to nearly 3000x1500. A great deal of the window content will be off the edge of the display as a result. To resolve this, right-click on Windows' Desktop, select **Display settings**, and change the **Scale and layout** to 100%. It might not be necessary to alter the **Display resolution**; Figure 0.2.

The **2D graphing windows** can be resized manually by the user by grabbing the window corners as usual in Windows.

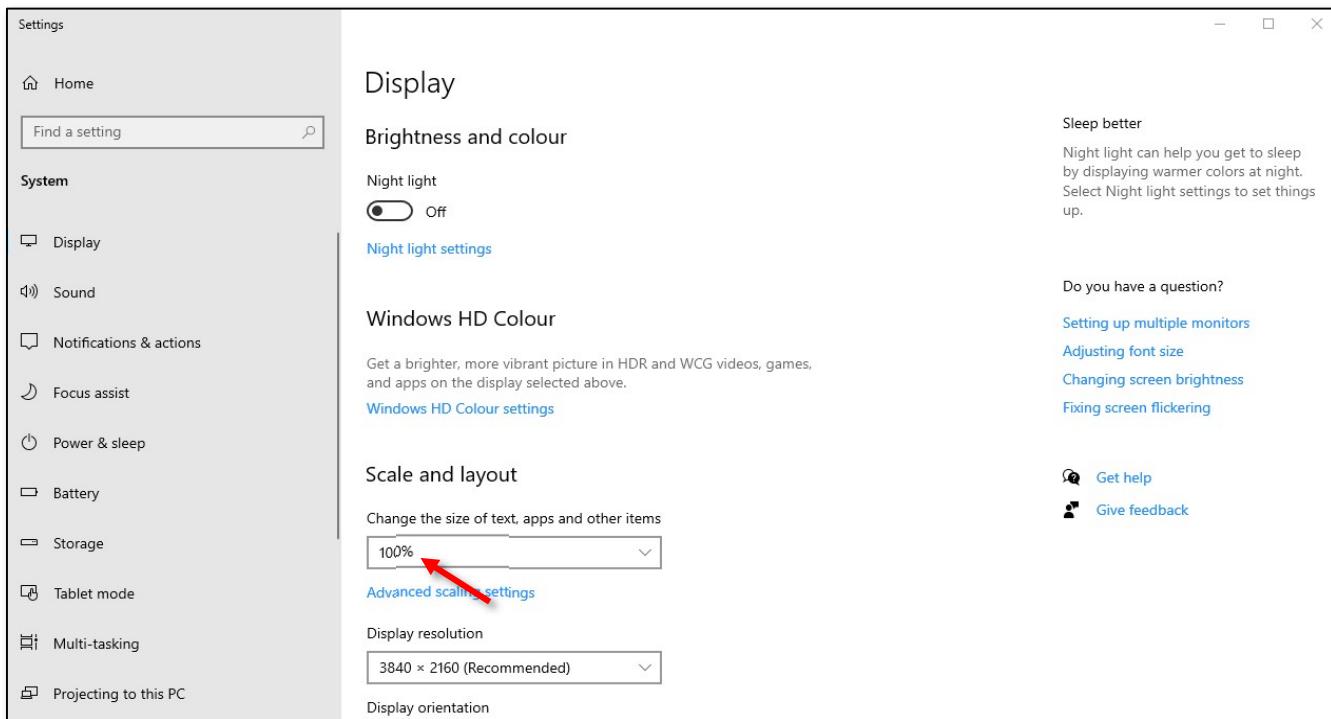


Figure 0.2 Adjusting Windows Desktop Display to 100% scaling



If, either on a PC or laptop, after loading a data file into DDMT, you don't see any of the windows on your screen, it may be because at some earlier time you had a second external monitor connected, and Windows has created the graphics windows for DDMT on that screen space (that no longer exists!). Hovering the mouse over DDMT on the taskbar, one might see greyed windows such as:




To bring the windows onto the current display, hover the mouse over each of the grey windows one at a time, right-click the mouse, left click **Move**, and tap one of the keyboard cursor keys. This will lock the window to the mouse pointer. You can now move the mouse and bring the window into view. Repeat for all.

See "Method 2" on this Microsoft link: <https://support.microsoft.com/en-us/topic/some-programs-are-displayed-off-the-screen-b2a915b7-4705-c455-febf-8ddb6e7f0ae4>

Common error:

Should you, for the first time you run DDMT, experience the following error, or similar:

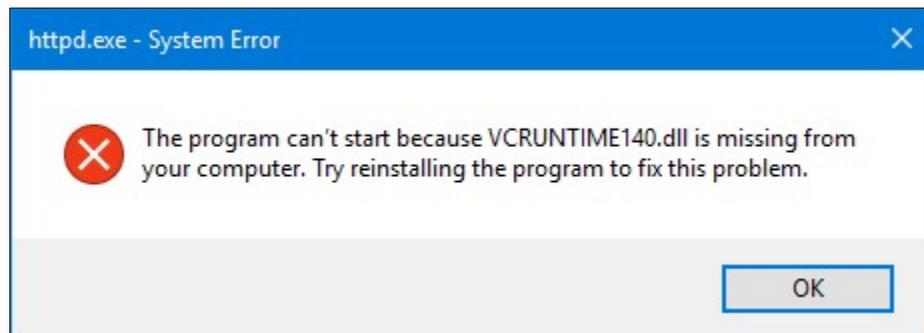

Please install the file named **vcredistx64.exe** contained within the DDMT zip file, to install the runtime libraries required for DDMT to function. Note that a restart of your computer may be required to complete this installation.



1. Importing and Exporting Data with DDMT

Extracting Daily Diary data from an SD card

With the current Daily Diary firmware, it is necessary to use Daily Diary Multiple Trace software (hereafter DDMT) to first extract the collected data from the MicroSD card.

When DDMT is first started, the user is presented with the control panel, Figure 1.1. To begin, click the button labelled **Extract from SD card**. This will initially present a small box, into which the name of the file that will be created with the Daily Diary data to be extracted should be entered. After **Set** is clicked, the user will then be asked to provide first the location of the SD card using a browser box, and second the destination directory for the new file containing the Daily Diary data. DDMT will then proceed to pull the data off of the SD card and store in 2 GB files ion the user's destination folder.

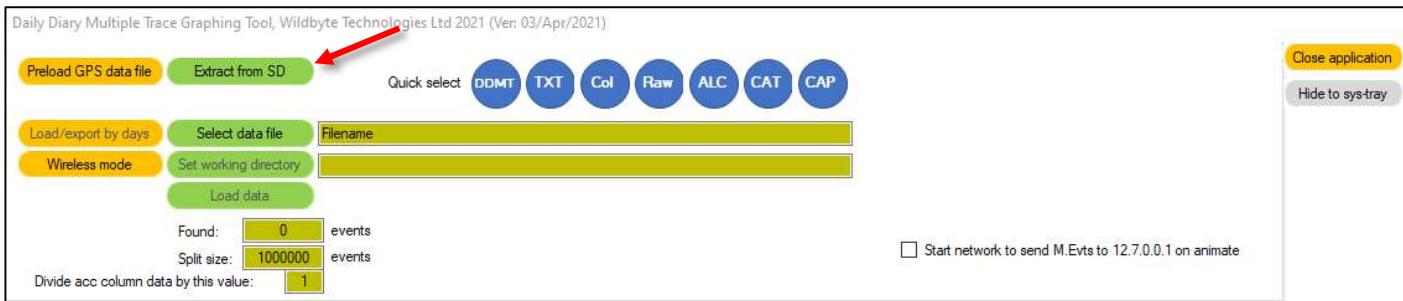


Figure 1.1 Initial control panel for data extraction from MicroSD cards (from a Daily Diary logging session)

Once DDMT informs the user that the data have been extracted, the user is then free to select and load this data file by clicking **Select data file**.

Loading data files into the software

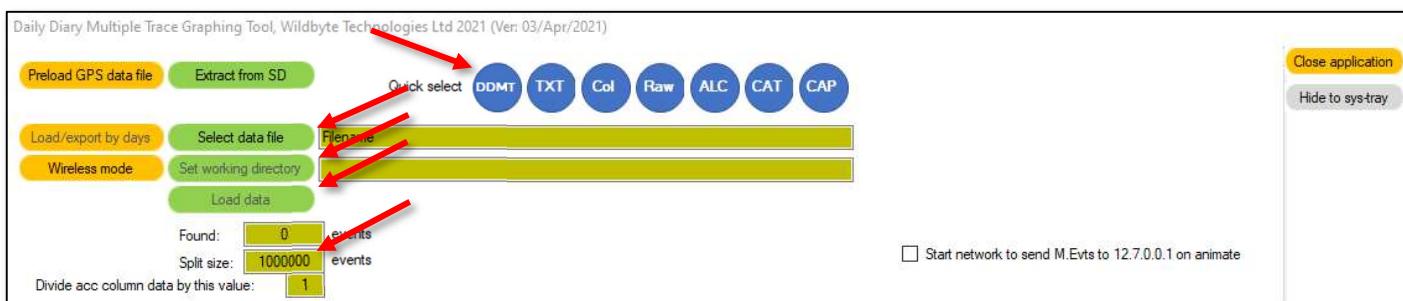


Figure 1.2 Click **Select data file** to locate the data to load into the application, followed by **Set working directory** (if required; Desktop is the default), and finally **Load data** to load some, or all, of the data into memory

There are several different file types that can be loaded into the software.

TXT	(extension *.txt) data file created by the Daily Diary, binary format (not human readable)
RAW	(extension *.raw) created within the software, binary format
COL	(extension *.txt) tri-axial acceleration data from other devices, ascii format
DDMT	(extension *.ddmt) various data channels (acc, mag, etc.) from other devices, binary format

Additionally, the user can pre-load other data structures before loading sensor data.

Preload GPS data file preload a GPS data file into memory; ascii/CSV format, discussed in Chapter 9



Select data file

a browser window will appear giving access to your files; user can manually select the data type from the above types (*.txt, *.raw, *.txt, *.ddmt, and other legacy types). Alternatively, select the file type using the blue buttons as indicated in Figure 1.2. Select the text file you would like to load either directly from your PC/laptop/storage device.

The application will read through the file and determine how many ‘events’ (data points) are recorded in the file (40 Hz = 40 ‘events’ (data points) per second; each event may contain one or more sensor values). The total number of events and the number of events to load will be displayed in the **Found** and **Load** boxes respectively.

Set working directory

to indicate where to store the outputs, e.g. data files or screenshots, from the software. The user’s **Desktop** is selected by default.

Note, some versions of Windows do not allow DDMT to write to the user’s Desktop, the default Working Directory. If you find no file found on export, please select an alternative Working Directory i.e. an alternative drive or an external drive

Load data

to display your first ‘split’ of data, Figure 1.3. A split is defined as either all of the data if the user is able to load all in one go, or a subset of the whole. A series of windows will now appear. The window named **drawing** will contain a series of graphs. **Drawing 2** is simply a second window to display the same, or other, data – think of it as more “floor space”.



To load a data file into DDMT:

1. **Select data file**
2. Change the file type to that required (*.txt, *.raw, *.ddmt etc.)
3. Set **Working directory** if the default (User’s Desktop) is not preferred
4. Adjust the **Split size** from 1,000,000 events to something larger (keep < 3,000,000 else the software will lag)
5. Click **Load data**

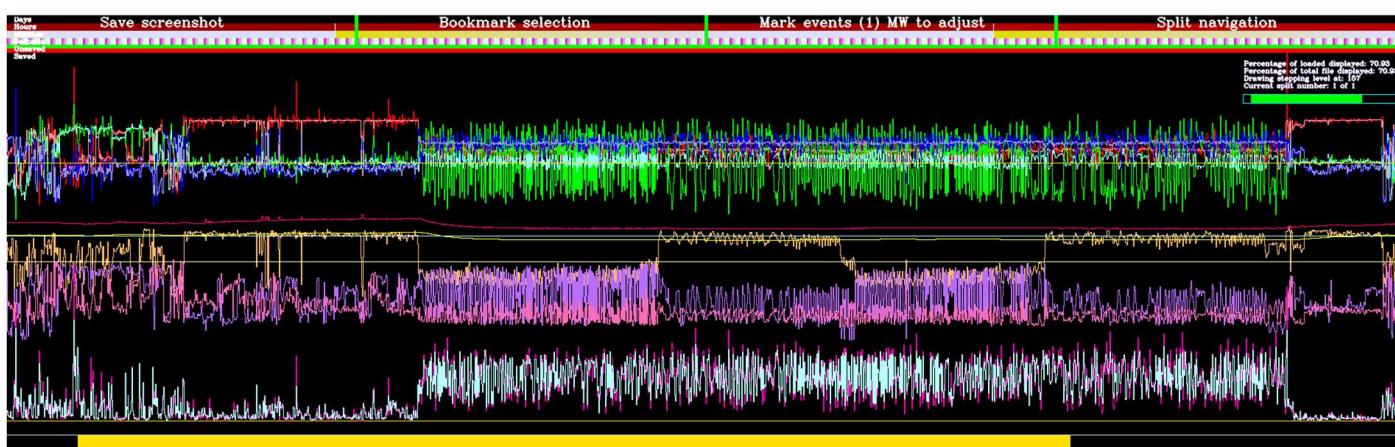


Figure 1.3 Main graphs window (2D)

By default, the top red/green/blue graphs represent acceleration X/Y/Z, the middle red and yellow represent pressure and temperature respectively, the lower brown/purple/pink represent magnetometer X/Y/Z, and finally at the bottom is VeDBA (a measure of Vectorial Dynamic Body Acceleration) in pink and VeDBA smoothed over the top in pale green. The temperature and



pressure vertical scales and positions will automatically adjust depending on the data to be displayed in order to maintain a relatively central position.

Close application to close the software at any point; located on the right-hand side of the control panel. This performs an orderly shutdown of the application.

Hide to Sys-tray will minimise the main control panel to the system tray. This is useful when it is set to “stay on top” but is required to be temporarily out of the way.

The **col** (Column) file format for importing triaxial accelerometry into DDMT

Triaxial accelerometry data can be loaded into DDMT from any source, not just Wildbyte's Daily Diaries. Simply create a file containing your 3 acceleration channels in the order X, Y, and Z. The delimiter must be “tab”, and there should be no header row. The file should have a file extension .txt.

As an example, here is some data in Excel; Figure 1.4.

	A	B	C
1	0.292606	0.206278	0.47943
2	0.079364	0.657229	0.454833
3	0.813417	0.255477	0.695476
4	0.949998	0.48214	0.002423
5	0.709955	0.763539	0.554591
6	0.94285	0.582487	0.788854
7	0.963425	0.489661	0.151246
8	0.652498	0.059643	0.71832
9	0.686657	0.667858	0.41335
10	0.637207	0.961429	0.513011
11	0.479828	0.84392	0.202608

Figure 1.4 Triaxial accelerometry in Excel. No header row present

Select **File -> Save As**, browse to the destination folder, and change the file type to “Text (Tab delimited) (*.Txt)”.

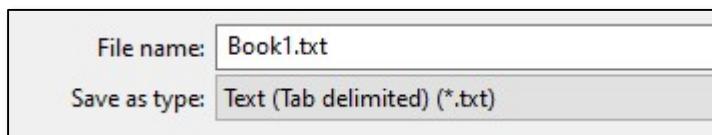


Figure 1.5 Export the data as file-type “Text (Tab delimited) (*.Txt)”

Previously, it was required that the extension .txt be adjusted to .col. This is no longer the case; it must now remain as .txt. Selecting **Col** as the file type on the initial control panel instructs DDMT of the incoming file format (*tab delimited, tri-axial only*).



The data file, broken down by split number

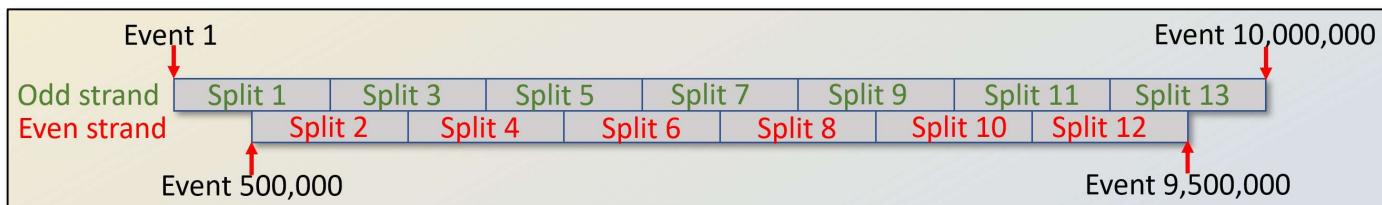


Figure 1.6 Data file “split” structure

Some definitions:

- “Event” A single instance in time containing all sensor values i.e. data point
“Split” A slice of the whole data file

Number of splits for a data file = $\sim ((\text{data points}) / (\text{split size})) \times 2 - 1$

In the above example, there are a total of 10,000,000 events (data points) contained within the file. In this case, the user loaded the default number of events for a “split” as 1,000,000 events. The *odd* splits span the total file, while the *even* splits are offset at the midway points of the *odd* splits, allowing the user some continuity as they work through the data.

Exporting data / calibration data

Data can be saved in two main formats, either as a raw format which can be loaded back into the software for further visualisation/algorithmic searches, or as tab delimited columns, which cannot be reloaded into DDMT, but can be transferred for analysis in other software packages. Below is the **File output/Save** tab and an explanation of each export option.

The export data function DDMT has now grown to allow a high level of control over what channels are output.

The user can change the default filename from “Filename”, to which are several characters that define the type of export, followed by a numerical suffix that is incremented after the execution of each save to help generate unique filenames.

By ticking **Omit headers**, when exporting ASCII data i.e. human readable format, headers will be removed from the first row.

The data to be exported is determined by how the user selects the range. There is a control in the centre that allows the user to switch between 3 types of data selection:

1. Save data in current split

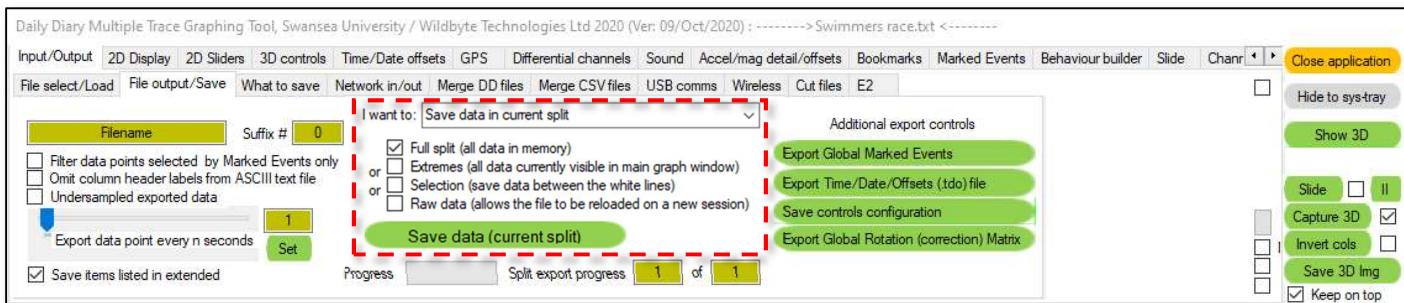


Figure 1.7a Save data in current split

This allows the user to export data from the current split, quickly and simply, as either ASCII for reading into another application such as Excel, MatLab, or R, or as a binary format, .raw, that can be reread into DDMT during another session.



Full split	will export all the data from the current split you have open
Extremes	will export all the data-points currently on display in the 2D graphing window .
Selection	will export all the data-points between the left and right white lines (left/right mouse buttons)
Raw data	creates a file that can be reloaded into DDMT (discussed a few pages down)

2. Save data from one event to another (total event number, not the local split event number) i.e. from any event within one split to any event number within another split

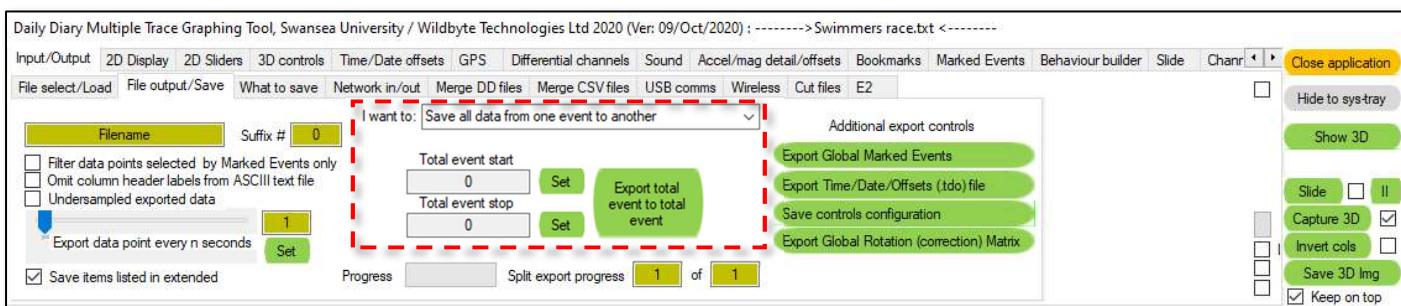


Figure 1.7b Save data from one event to another.

DDMT will move through the splits and export all data between the start and end points. The user must place the left white line at the start and end points and click **Set** to define these, and then click the button labelled "Export total event to total event"

3. Save data as per the split list provided

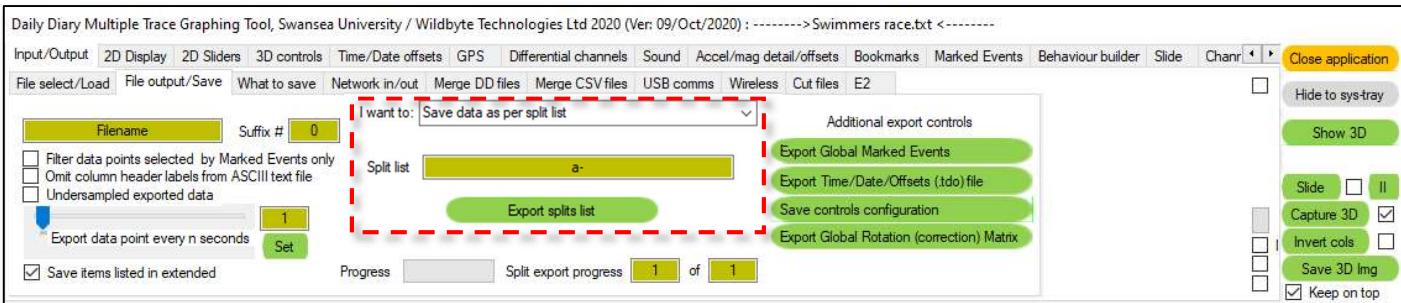


Figure 1.7c Save data by split number as dictated by the list default "a-"

Export splits list

Specific splits can be selected for export using the syntax illustrated in the examples below. Note that while letters below are given in uppercase, these commands are not case sensitive, i.e. E4-10 and e4-10 generate the same outputs.

- E4-10** even-numbered splits only, from 4 to 10 inclusive, i.e. splits 4, 6, 8 and 10 are exported
- O3-10** odd-numbered splits only, from 3 to 10 inclusive, i.e. splits 3, 5, 7 and 9 are exported
- E-** all even-numbered splits in the data file
- O-** all odd-numbered splits in the data file
- all odd-numbered splits in the data file
- 3** split 3 only
- 3-** all splits from split 3 to the end of the data file, both odd- and even-numbered
- 7** splits 1 to 7 inclusive, both odd- and even-numbered



- E-7** even-numbered splits from 1 to 7, i.e. splits 2, 4 and 6 are exported
- O-8** odd-numbered splits from 1 to 8, i.e. splits 1, 3, 5 and 7 are exported

Additional syntax

- C or c** signifies the current split, i.e. the data split currently loaded in DDMT
- +** used with the current split command above and followed by a number exports the current split and the specified number of following splits (inclusive)
- c-** exports all *odd* or all *even* splits (depending on whether the current split is *odd* or *even*) from the current split to the end. Note if exporting *even* splits, it will miss the last section of data which is contained in the *odd* splits only which are staggered by half a split
- ec-** exports all *even* splits from the current split to the end of the data file. If the current split is *odd* it will not be included
- oc-** exports all *odd* splits from the current split to the end of the data file. If the current split is *even* it will not be included
- ec-9** exports all *even* splits from the current split to split 9. If the current split is 3 then splits 4, 6 and 8 are exported
- c** exports the current split
- c+2** the current split and the next 2 splits; the system will look to see if the current split is on the *odd* or *even* strand and export the current split and the number of splits listed
- ec+3** forces DDMT to *even* strand. If the current split is *odd*, it will be skipped and only the next 3 *even* splits are exported
- e-c** all *even* splits from the start of the data file (i.e. 2 onwards) to the current split (or the one before the current if the current split is *odd*)
- c** either all *even* or all *odd* splits until the current split. Whether *odd* or *even* splits are exported is determined by whether the current split is *odd* or *even*
- 3-c** here, the current split has the priority of the strand. If the current split is *even*, split 3 is excluded, but 4, 6, 8, etc. to the current split will export
- E3-c** here, the *even* strand is forced so that 4, 6, 8, etc. to the current split are exported (unless the current split is *odd*, in which case, it will be the split before the current)

The user can also choose to export alternating splits, where DDMT saves data from all splits (therefore alternating between odd and even splits) but excludes the first and last 25% of data in each individual split. This exports all data without overlap and avoids zero values in the smoothed channels at split boundaries. ***This function cannot be combined with other search strings.***

- a-** all splits from start to finish, dropping the first and last 25% of each split. This command is the default option
- a-7** from the start of the data file to the end of split 7, dropping the first and last 25% of each split
- a3-7** from split 3 to split, dropping the first and last 25% of each split
- a7-** from split 7 to the final split, dropping the first and last 25% of each split



The first split defined here will include the first 25% of data and the final defined split will include the last 25% of data, so no data is lost from the beginning or end of the total data file.

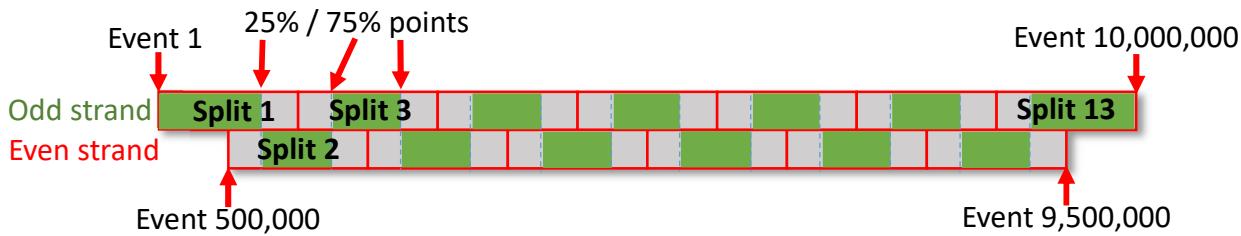


Figure 1.8 Using the “a” prefix will make the export function switch (alternate) between the odd and even strand at the 25% / 75% points to ensure that it doesn’t export the left and right extremes of each split where smoothed data is zero (due to the centred window averaging function)

What to save (tab)

When exporting as tabulated data (**Raw data** unticked), this determines which sensors and metrics are exported. The **General / Primary channels** tab contains the timing and raw sensor data, while **Additional / metrics channels** displays all Daily Diary channels; dead reckoning coordinate data can also be exported (see chapter 10).

If **Save items listed in extended** is unticked, none of the selections within **What to save** will be exported, regardless of their being selected or not; a quick way of omitting multiple selections.

Daily Diary Multiple Trace Graphing Tool, Swansea University / Wildbyte Technologies Ltd 2020 (Ver: 09/Oct/2020) : ----->Swimmers race.txt <-----

Input/Output 2D Display 2D Sliders 3D controls Time/Date offsets GPS Differential channels Sound Accel/mag detail/offsets Bookmarks Marked Events Behaviour builder Slide Chanr ▾ Close application

File select/Load File output/Save What to save Network in/out Merge DD files Merge CSV files USB comms Wireless Cut files E2

General / Primary channels Additional / metrics channels

Counter Time (HH:MM:SS)
 Local event number Time (HH:MM:SS.ddd)
 Total event number Total time
 Behaviour Decimal seconds
 Accel X/Y/Z Analogue
 Mag X/Y/Z Light level (issue 3C+)
 Temperature Battery V (issue 3C+)
 Pressure Time (HH MM SS)
 Date (DD:MM:YYYY)

Accel X/Y/Z Sm Altitude Changed Mag smoothed Dead Reck Long/Lat/Alt Marked Events
 Accel vectorial sum and + VeDBA Mag Vec-Sum AAV Differentials
 Accel sm vectorial sum and + VeDBA smoothed Mag Vec-Sum smooth GPS fix Long /Lat/Alt FFT (x4)
 Dyn Accel X/Y/Z VeSBA Median Mag X/Y/Z Fly_Wire Smoothed FFT sweep
 Accel X/Y/Z Submean VeSBA+ Angular velocity Speed mag X/Y/Z FFT zero crossings
 Smoothed Pressure ODBA Normalised Mag Speed Vel Smoothed FFT zero c (manual)
 Median Pressure Analogue smoothed Magnetic heading Vel m/s Separate DD MM and YY
 Pressure Changed Activity Pitch/Roll angles Vel/VeDBA
 Altitude Smoothed Cumul Hdg

Select all Deselect all

Close application

Hide to sys-tray

Show 3D

Slide Capture 3D

Invert cols Save 3D Img

Keep on top

Figure 1.9 What to save General / Primary channels / Save items listed in Additional / metrics channels sub-tabs on the control panel



Raw (binary) data files

Sections of data can be exported in *.raw* format (binary files, as opposed to ASCII – human readable) and reloaded into the software later. This might be useful where the user is cutting out useful sections of the data, and thus discarding useless / unwanted data. Note that *.raw* data files can be merged into a single large *.raw* file.

Note that the **Export splits list** function outlined in the previous section can also be used to export *.raw* files, simply select the **Raw data** checkbox on the **File output/Save** sub-tab.

Please ensure that any time corrections and magnetometer / accelerometer corrections have been applied to the data prior to exporting as *.raw* files, as these correction functions are disabled when loading in *.raw* files.

Exporting with Under-sampling to “thin” the data

Exported data may be ‘thinned’ by selecting the **Undersample exported data** checkbox on the **Input / Output File output / Save** sub-tab. The degree of under-sampling is determined using the **Export data point every n seconds** slider, or by entering a value in the box and clicking **Set**. Allowed values are “every 0.1 s”, “every 0.2 s”, “every 0.5 s”, or every second from 1 to 297. To sample at sub-second levels, it is easier to enter 0.1, 0.2, or 0.5 into the box and clicking **Set**, but one can use the slider, moving it to the left, and then left-clicking in the space just to the right of the slider to job it up to 0.2 or 0.5.

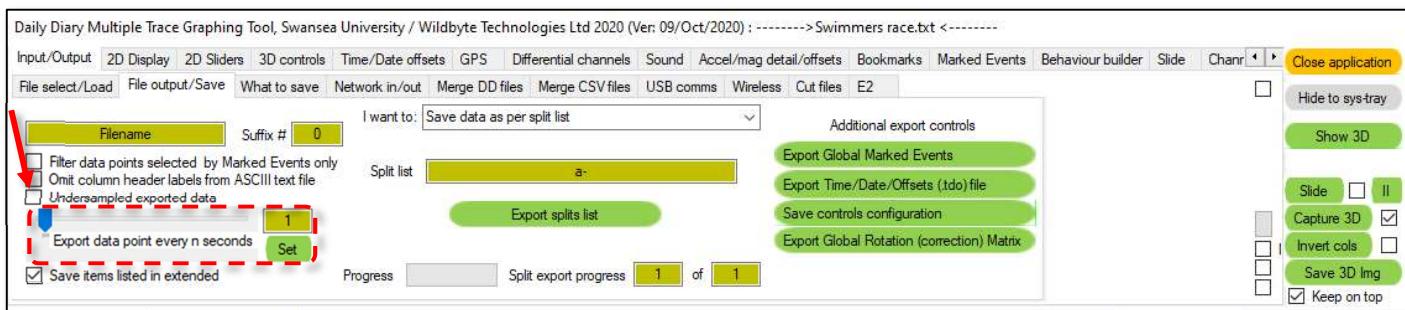


Figure 1.10 File output/Save sub-tab with under-sampling controls indicated

Exporting data by Marked Events

Marked Events are literally just that, they are events that have been given a marker (value 1-9) either manually by the user, by converting **Bookmarks** to **Marked Events** or through the result of the Behaviour Builder/Time series Search functions. Manually marking data is discussed in Chapter 10.

It is also possible to export **Marked Events** as defined by behaviour builder (see Chapter 10) by selecting the **Export behaviour builder marked events only** checkbox in the **Input / Output File output/Save** sub-tab. In doing so, behaviours, or simply sections of identified data, will export along with the **Marked Event** value. Possibly useful for further analysis in other software packages such as R, or MatLab.

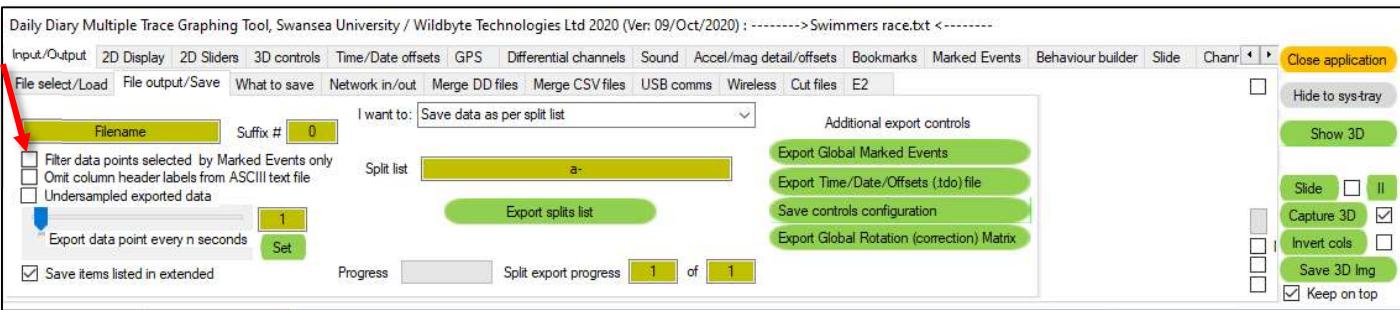


Figure 1.11 File output/Save sub-tab with export marked events option indicated



Merging raw files

Exported .raw files (e.g. perhaps with each containing data of a specific behaviour or time period) can be merged into a single contiguous file, and then reloaded into DDMT, as long as they are from the same original .txt file.

Note that in order to merge files, one must already have “any” file loaded in memory to access the **Merge file** tab.

On the **Merge files** sub-tab, select files with the **Merge multiple .raw files** button. This will open a window allowing you to choose the files you wish to merge. The selected files will then be listed in the box below (*All files to be merged must reside within the same directory*).

Next use the **Set working directory** button to determine where the new merged file should be saved.

Finally, select the **Merge the files** button.

A dialogue box stating the name of the new file will appear when the software has finished merging the files. The software will indicate if the selected files are standard .raw files, or if they are bookmarked .raw files. Standard .raw files are simply selections of the data saved in .raw format, while bookmarked .raw files are exported through a separate process using the functions within the **Bookmarks** tab. Bookmarked .raw files and standard .raw files cannot be merged. A bookmarked .raw file can be converted to a standard .raw file simply by loading it and then using the standard **Save** function with the **Raw data** checkbox selected from the list on the **File output / Save** tab.

The file merging function will not affect any data currently in memory. Merged standard .raw files will result in a standard .raw file, while merged bookmarked .raw files will result in a larger bookmarked .raw file, with all the bookmarks from the two files in their expected positions.

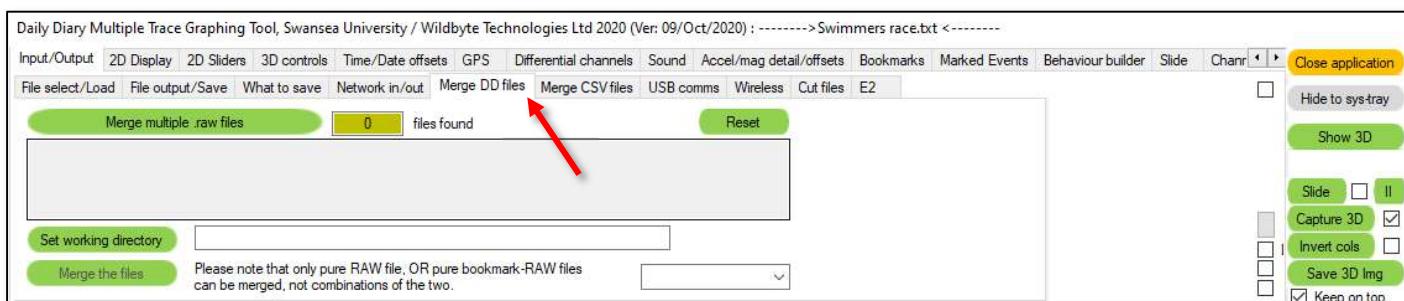


Figure 1.12 Merge DD files tab on the control panel

Note that merged .raw files maintain their timing data (both date and time) so there will be “time jumps” (on the **Values overlay**) as the user slides through the data.



Merging csv / tab delimited files

When exporting data from DDMT, often is the case that this results in many long files of either tab delimited, or csv files. The majority of the export functions work on a per-split basis with the export/save data function and so, if the user requires that the final result is a single file of human-readable data, these multiple files will need to be merged.

Files may be merged using the tab **Merge CSV/Tab files** (includes the tab-delimited format).

Clicking **Select files to merge** opens a browser dialog. Hold down *ctrl* to select multiple files and click *Open* to apply these file names into the box below the button. Note that the files are added into the box in the order that they will be merged, starting with the first at the top. To change the order of the files added, use Windows Explorer column header controls to order them such that the first file is at the top of the file selection dialog prior to clicking *Open*. Within the Windows file dialog, there is a dropdown box at the bottom to switch between csv and tab (txt) formats.

Set the working directory, where the merged files will come together in a new file, appropriately labelled, and click **Merge csv/tab files**. While it is merging the files, a small counter at the bottom left will increment as it steps through the list of files.

A tickbox labelled **Omit columns headers from all but first file** should ordinarily be ticked.

The **Reset** button clears the list of files in the box.

Note that it is assumed that all files will contain the same number of columns; the function simply reads each line of data, and immediately exports these into the final file.

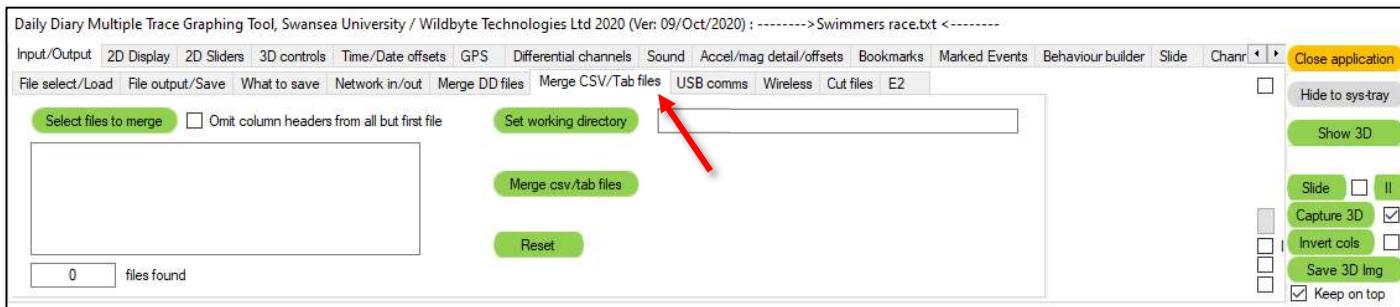


Figure 1.13 Merge csv/tab files tab on the control panel



2. Main graphs (2D) windows

Time indicator bars at the top and bottom of the main drawing windows

Graphically, data from these tags can be difficult to understand at times, more so when one cannot easily discern the time scale of changing events. The coloured bars at the top of the screen give an indication of evolving time from the Daily Diary data (Figure 2.1).

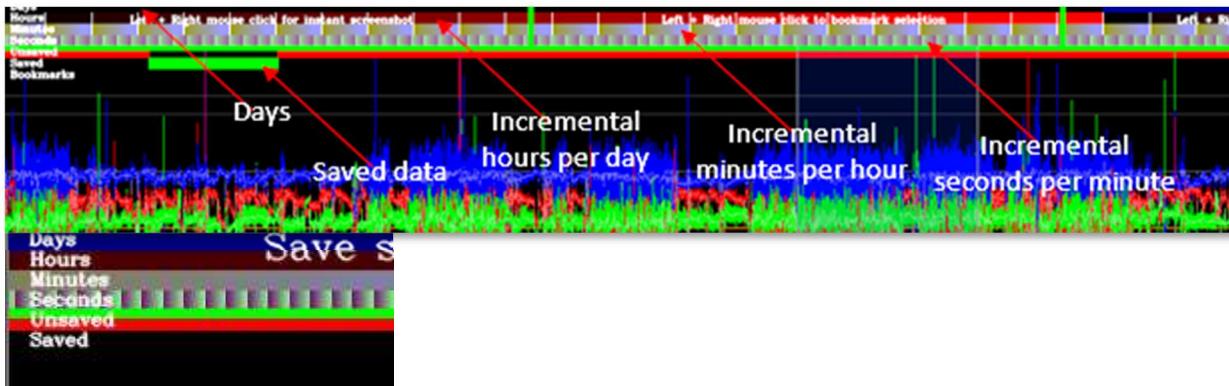


Figure 2.1 Explanation of the top indicator bars on the **drawing** window

The coloured bars, as indicated above, represent Days, Hours, Minutes, Seconds, and whether sections of data have been saved or not. The red bar level with the white **Unsaved** label to the left of the **drawing** window indicates nothing so far has been saved; the green line below (which intersects the red line above it) shows this section of data has been exported/saved. The **Incremental hours per day** bar changes colour from black to bright red with increasing hour value, returning to black at midnight, while the **Incremental minutes per hour** bar changes from green to mauve through the hour, and underneath likewise for **Incremental seconds per minute** (Figure 2.1).

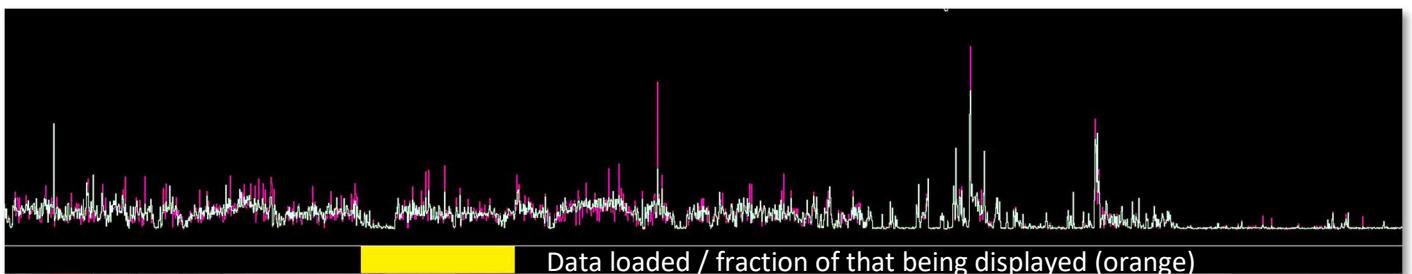


Figure 2.2 Bottom indicator bar on **drawing** window

The orange bar at the bottom of the screen represents the fraction of the current split on display (a full width bar represents the full split). To zoom in on a section of the data, simply left-click to define the left **Extreme**, and right-click to define the right **Extreme** on this orange bar, which will rescale accordingly in size. The position of the bar corresponds to the position of the data currently in-view as part of the currently loaded split.

Percentages of displayed data as a function of that loaded and the whole dataset are displayed in the top right corner. The percentage of data loaded from the current split and total file are given, as well as the drawing stepping level and current split number. These values can be removed from view by deselecting **Display percentage visuals on 2D graph window** under **Display / Display options**. The above bar indicator will make more sense once the user has explored the data a little.



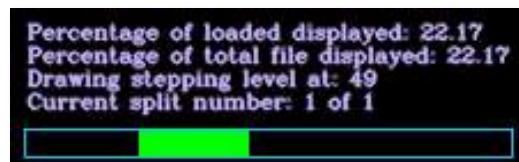


Figure 2.3 Percentage loaded/displayed indicators

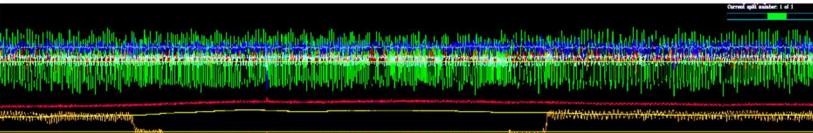


Displaying data in “under-sampled mode” for improved application performance

The data drawn in the two **Drawing windows** is, by default, under-sampled automatically to improve application performance. Most display monitors these days are typically 2k pixels across (full HD) or 4k pixels across (4k UHD), and so, if there are 1,000,000 events / data points in the current split and the user has zoomed out enough to display all 1,000,000 of these for several different channels of data, this would cause a severe performance problem. What the application does is to look at how many events it has to display, and divides this down by the horizontal resolution of the display, and uses this value as the “stepping” value through the data. This under-sampling can be disabled on the **Display** tab by deselecting “Dynamic undersampling”:

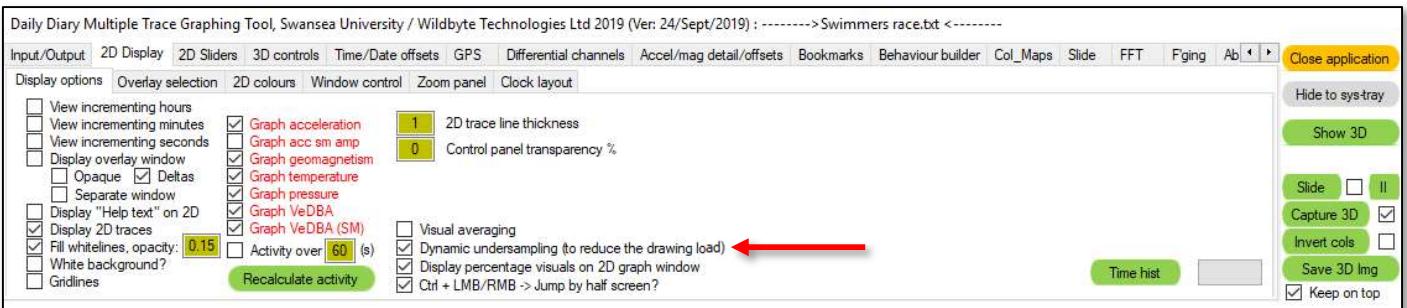

Left: Under-sampled data view; more responsive, but key peaks may not be displayed

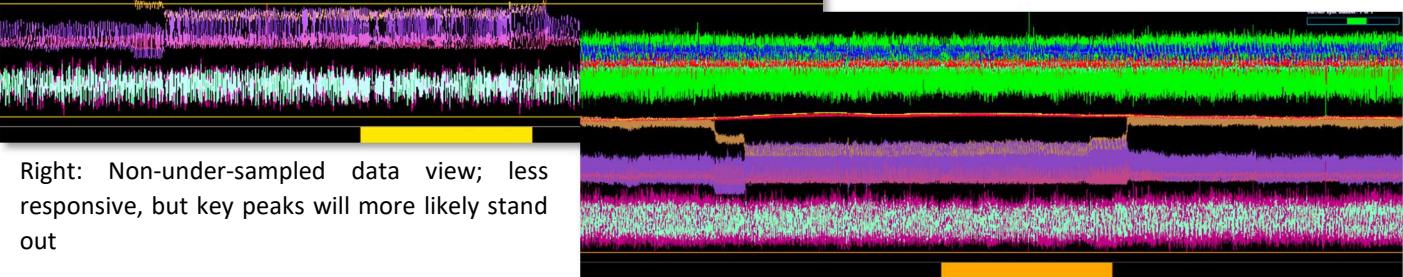

Right: Non-under-sampled data view; less responsive, but key peaks will more likely stand out

Figure 2.4 Under-sampled vs non under-sampled drawing modes on the two **drawing windows**



Displaying data in “Visual averaging mode”

The alternate drawing mode is **Visual Averaging**. To use this drawing mode, select **Visual averaging enable**. When displaying averaged data, a slider, **Degree of averaging**, will be available to set the degree of smoothing of the currently visible channel data. This defines a centred sliding average window width that visually smooths the data set; Note that this does not in any way affect any exported data; it is a visual effect only.

Two 2D graphing windows

When considering multiple data channels, it can be helpful to use two graphing windows. By default, selected data channels are displayed in the main graphing window, titled **drawing**. An additional graphing window, **drawing2**, is blank by default but can be useful for showing more channels at larger amplitudes across two screens. To adjust these settings, select the tab labelled **Window control**, a sub-tab of **Display**. This allows the user to set which data channels are displayed in each of the two **drawing windows**. Note that the two windows represent the same time period and adjusting the data on view in one graphing window automatically updates the other to match. Figure 2.5 shows an example where the second graphing window has been used to display the differential channels (see chapter 4) corresponding to Acceleration X, Y and Z.

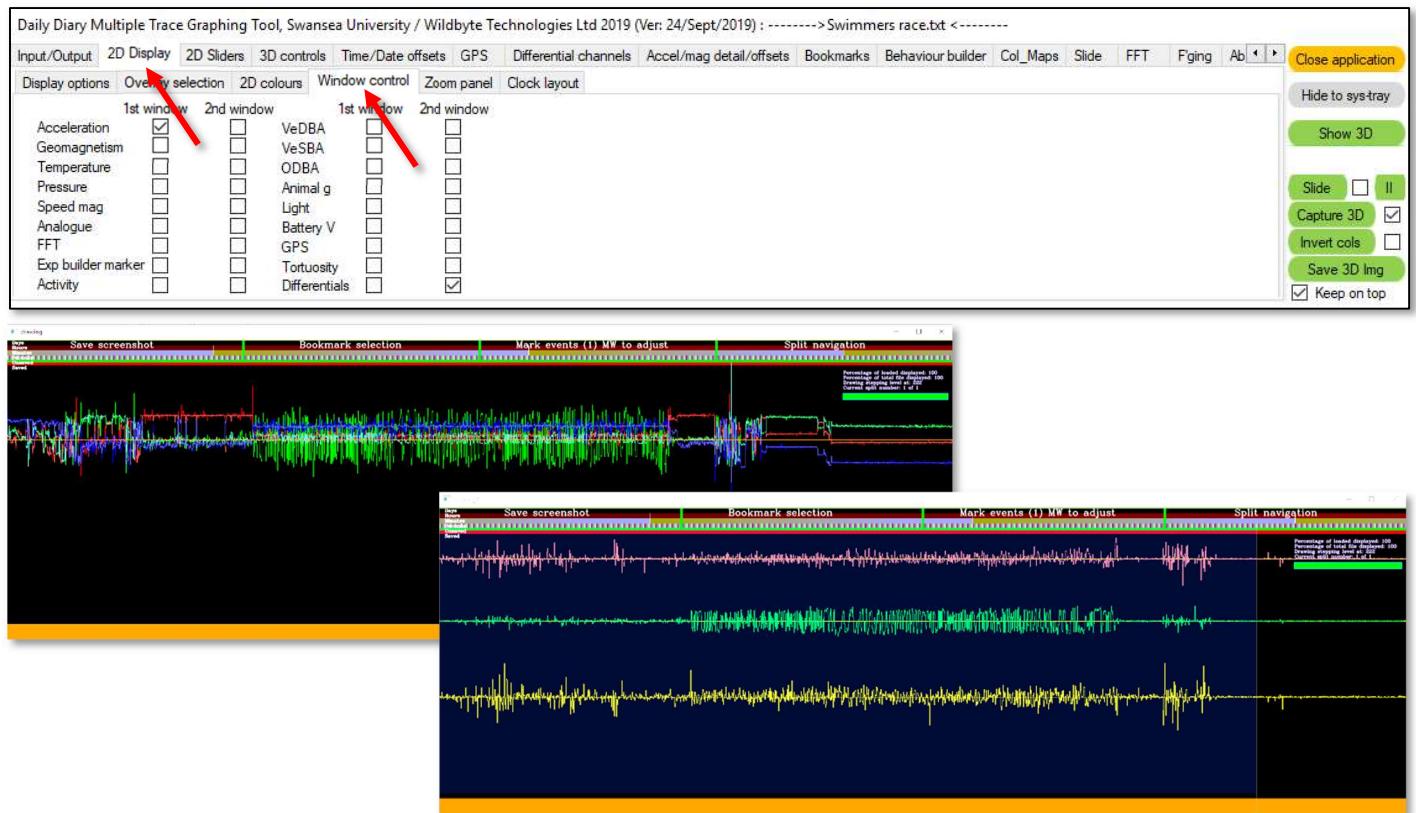


Figure 2.5 The settings on the are set to show acceleration (X, Y, and Z) on **drawing** and the differentials of these on **drawing2**



Navigating through a large dataset

To move to another split, simultaneously left and right mouse-click on the fourth box along the top of the drawing (**Left + Right mouse / Left CTRL + Left mouse Navigate through data**).

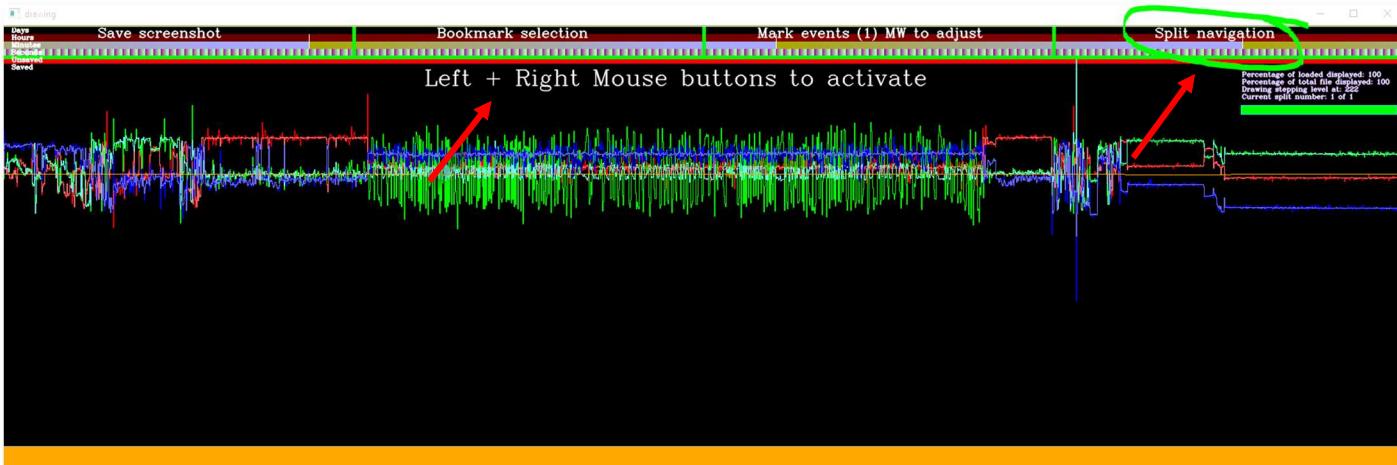


Figure 2.6 Bringing up the split navigation menu

A sub-box will then appear in the centre of the drawing; if this doesn't appear, hold down the left and right mouse-click on the fourth box again for a couple of seconds, and then release. The top row of the sub-box represents the full set of data, while in the bottom row, each unit represents half of each of the two top units it touches (Figure 3.6). These units are “splits” of data. The splits are alternately numbered, with odd numbers on top and even numbers on the bottom. The number of splits present in the file is determined by the number of events loaded when the file is first opened (refer back to Chapter 1 – Exporting Data).

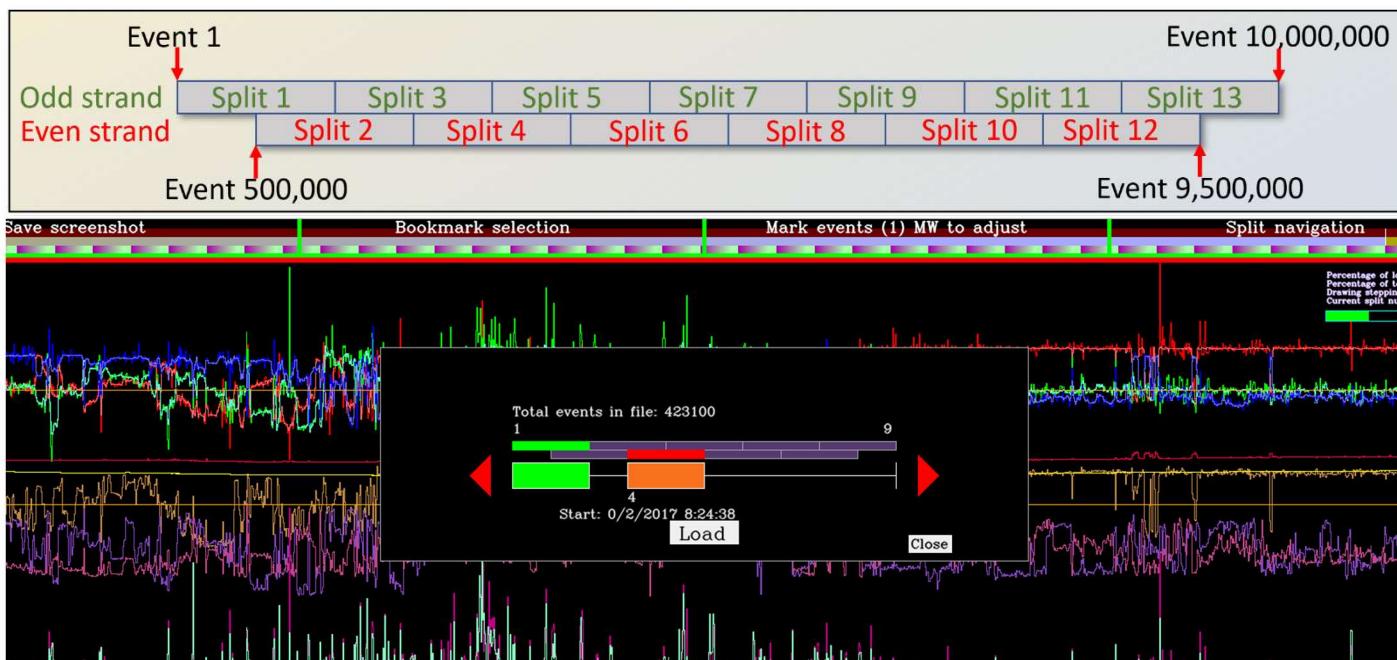


Figure 2.7 The split navigation menu

The currently loaded split is coloured green. To select a new split to load, left mouse click and drag the orange block along the panel until you reach the split you would like to load. The start time and date of the split at the position of the orange box selected is given below the bar for reference. Alternatively, especially for files containing a large number of splits, the two red arrows either side of the navigation menu can be used to increment/decrement the split number. Click **Load** to load and process this split for display.



Generating custom colour schemes for 2D graphed data

The default colour scheme of the graphing windows can be modified for ease of access for users with colour blindness. All data channel colours can be customised fully. The third sub-tab of the **Display** tab is **2D colours**. The controls shown on this tab enable the user to select a variable from the drop-down menu (Acc X, for example), and adjust the red/green/blue colour sliders to adjust its final colour. The rectangle beneath represents the resulting colour mix. New colour schemes can be saved as a .csv file and loaded again later or reset to the default using the green button commands under the drop-down menu.

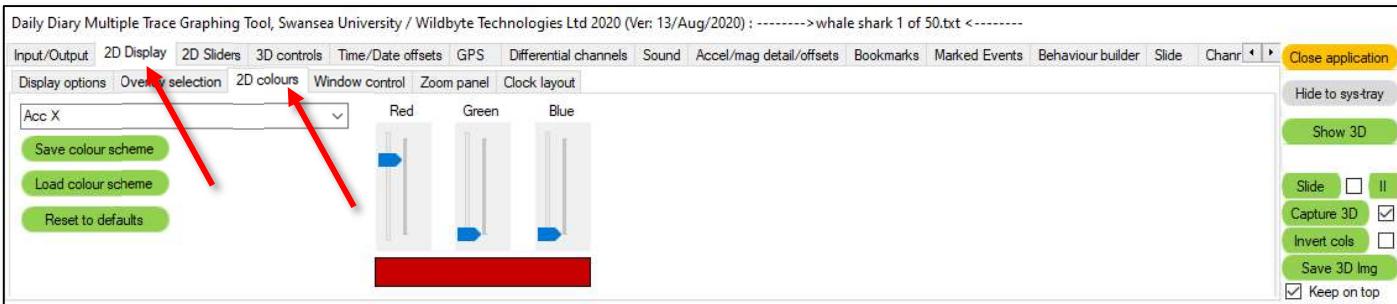


Figure 2.8 Adjusting the colour of traces and indicators within the main drawing

Figure 2.9 shows the effect of adjusting the Acceleration X/Y/Z colour balance from the standard red/green/blue. This can also be useful for screenshots, especially when captured using a white background.

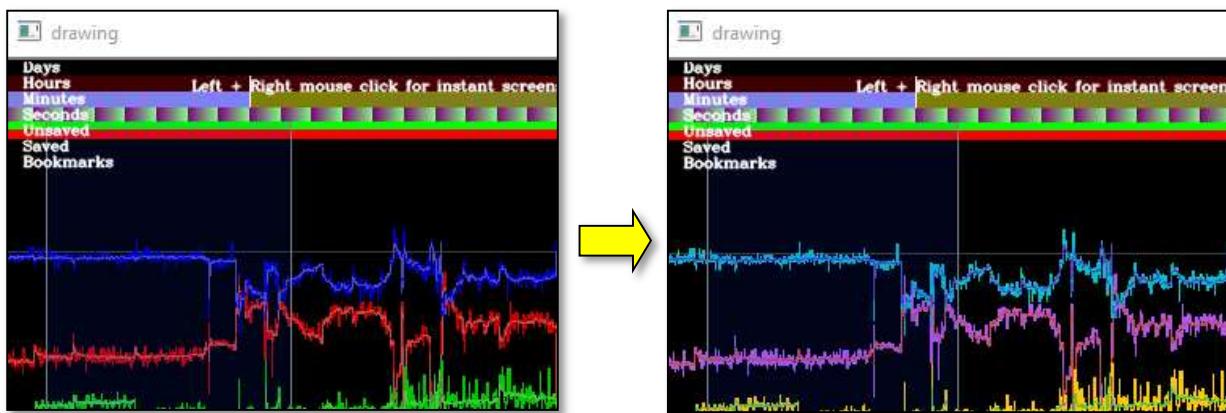


Figure 2.9 Colours for Acceleration X/Y/Z adjusted from standard red/green/blue



Horizontal guidelines

On the **2D main windows**, by holding down **Ctrl** and **Alt**, and then left or right mouse clicking, one can produce 2 horizontal lines (different colours). These are merely for comparing the heights of data points on the display. They are also used for defining the vertical height of the **Bookmarks** and **Marked Events** bands.

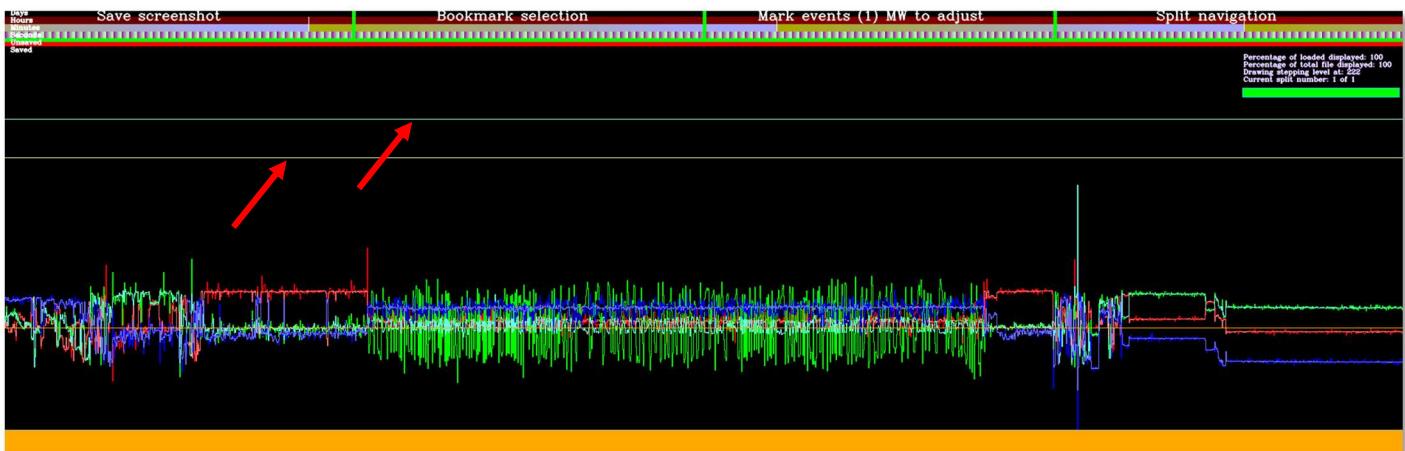


Figure 2.10 Horizontal lines for data comparison / levelling

It is also possible to add more horizontal lines that have the only purpose of data level comparison. To add a new horizontal line, hold down the **Alt** key; a coloured band will appear on the left and right sides of the **2D graphing windows**. Right-clicking in the left side band will generate a new horizontal line that can be adjusted up/down at any time (while holding down the **Alt** key and right-clicking the mouse). Holding down **Alt** and right clicking in the left band again will create another, and another etc. horizontal band. Holding down the **Alt** key and right-clicking in the band on the right side of the window will clear all these supplemental horizontal lines.

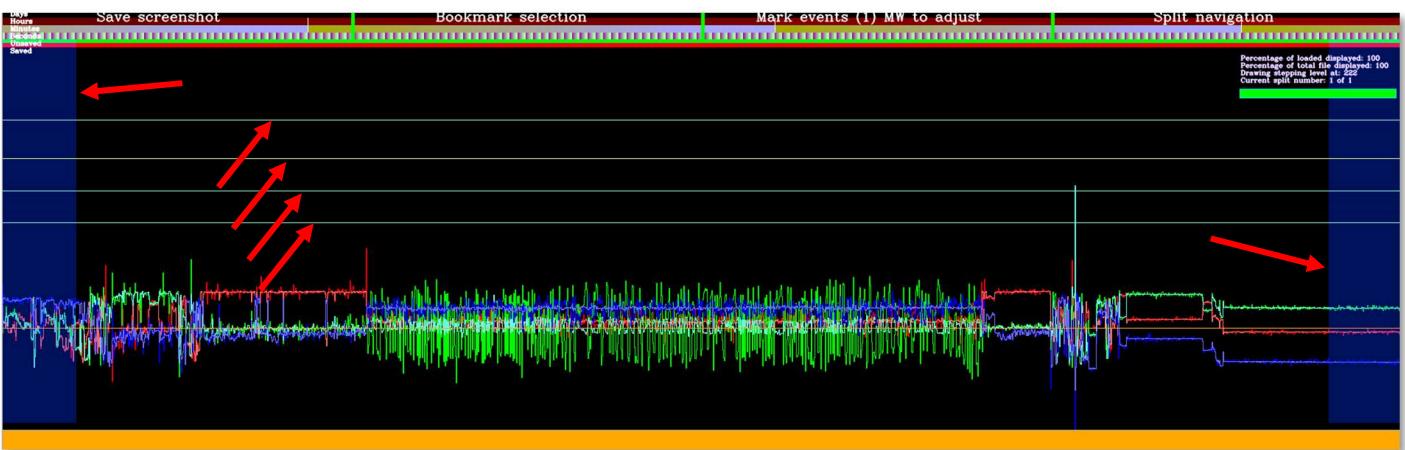


Figure 2.11 Additional horizontal lines for data comparison



Values overlay

To view, in real-time, numerical values from the data displayed in DDMT, select the **Display** tab, and **Display options** sub-tab on the control panel, the default tab after first loading a data file. Using this tab, you can control a lot of what is displayed on the **2D graphing window**.

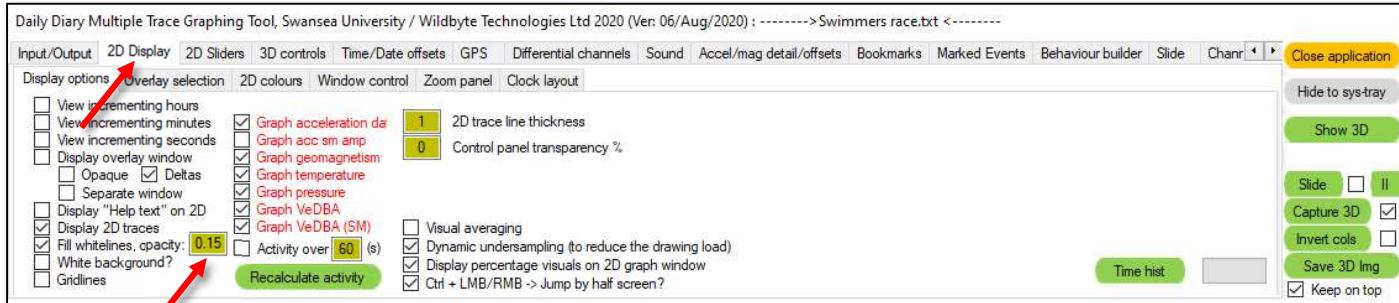


Figure 2.12 Display/Display options tab

On the left side of the window you will see a tick box named **Display overlay window**; check this box and a list of names and values should appear on the drawing. An alternative, quicker, method to bring the overlay window up is to depress the mouse-wheel button. To move the overlay window, hold down shift on your keyboard and drag the overlay to where you would like it using a left mouse click. Additionally, by selecting **Opaque** under **Display overlay window** on the control panel, the overlay is put on a dark background for ease of viewing. If using the mouse wheel to bring up the overlay window, it appears on a dark background by default. The overlay can also be displayed in an external window by selecting the **Separate window** tick-box in the **Display options** tab, under **Display overlay window**.

Left and right click respectively on the drawing to select two vertical parameters to compare data from, example in Figure 2.13. The event data along the left and right vertical white lines will be displayed in the left- and right-hand column on the overlay. Additionally, a third column will be displayed; these delta values are the difference of the left and right column. This column can be removed from the overlay by deselecting **Deltas**.

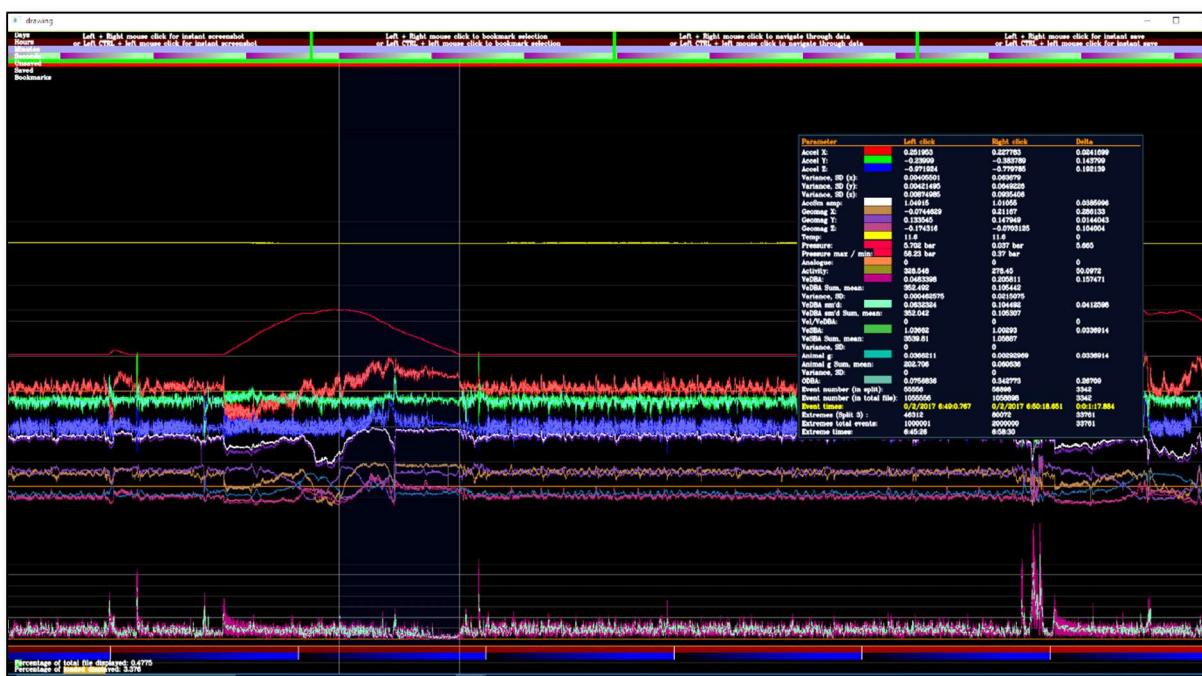


Figure 2.13 Main drawing with the opaque overlay and left and right mouse line position parameters

Variance and Standard Deviation for the acceleration (raw acceleration, not smoothed), VeDBA, and VeDBA sum/ VeDBA mean etc. relate to the data between the left and right white lines. These values therefore will change when the left or right white lines are adjusted in position.



Highlighting variables

Variables of interest on the overlay window can be highlighted for ease of viewing. Items in the list can have their display colour toggled from white to yellow for viewing contrast. This can be done by holding down **Left ctrl + Left alt** and running the mouse over the left side of the overlay on the parameter labels and left-clicking on one of the parameter names; the parameter name will increase in size and turn green, and toggle yellow or white with a left mouse click. The **Event times** parameter is yellow and cannot be unhighlighted (Figure 2.14).

Parameter	Left click	Right click	Delta
Acceleration X (raw):	0.171875	0.355469	0.183594
Acceleration Y (raw):	0.746094	-0.164062	0.910156
Acceleration Z (raw):	0.738281	0.453125	0.285156
Acc Variance, SD (x):	0.0371921	0.192852	
Acc Variance, SD (y):	0.733651	0.856534	
Acc Variance, SD (z):	0.0341892	0.184903	
Accel smoothed vectorial sum:	0.61792	0.602051	0.0158691
Accel smoothed vectorial sum+:	1.38208	1.39795	0.0158691
Accel sm vectorial sum+ (sum, mean):	24994.5	1.38489	
Accel vectorial sum:	1.06348	0.598633	0.337891
Accel vectorial sum+:	1.06348	1.40137	0.464844
Accel vectorial sum+ (sum, mean):	21396.1	1.18551	

Figure 2.14 Highlight items on the overlay

Adding items to the values overlay

Variables which are not of interest can be removed from the overlay window list entirely. The parameters displayed on the overlay are customisable, check or uncheck your choices using the list on the **Display/Overlay selection** sub-tab (Figure 2.15).

Daily Diary Multiple Trace Graphing Tool, Swansea University / Wildbyte Technologies Ltd 2020 (Ver: 06/Aug/2020) : ----->Swimmers race.txt <-----

Input/Output 2D Display 2D Sliders 3D controls Time/Date offsets GPS Differential channels Sound Accel/mag detail/offsets Bookmarks Marked Events Behaviour builder Slide Chanr. ▾ ▾

Display options Overlay selection 2D colours Window control Zoom panel Clock layout

Overlay display to include:

<input type="checkbox"/> Acceleration X/Y/Z	<input type="checkbox"/> Geomag	<input type="checkbox"/> Magnetic heading	<input type="checkbox"/> VeDBA	<input type="checkbox"/> Battery voltage	<input type="checkbox"/> Display DR data and markers
<input type="checkbox"/> Acceleration XYZ smoothed	<input type="checkbox"/> Geomag normalised	<input type="checkbox"/> Cumulative Heading %	<input type="checkbox"/> VeSBA	<input type="checkbox"/> Light level (au)	<input type="checkbox"/> GPS Long/Lat/Alt
<input type="checkbox"/> Acceleration XYZ mean dynamic/smooth	<input type="checkbox"/> Geomag smoothed	<input type="checkbox"/> Absolute angular velocity	<input type="checkbox"/> VeSBA+	<input type="checkbox"/> GPS mean velocity m/s	
<input type="checkbox"/> Acceleration smoothed for compass use	<input type="checkbox"/> Vec_sum	<input type="checkbox"/> Angular velocity	<input type="checkbox"/> ODBA	<input type="checkbox"/> FFT freq/Period	<input type="checkbox"/> GPS ground speed m/s
<input type="checkbox"/> Acceleration raw / smoothed vec sum (and vec sum plus)		<input type="checkbox"/> Pitch/Roll angles	<input type="checkbox"/> ODBA Smooth	<input type="checkbox"/> FFT Search, Hz/Amplitude	
<input type="checkbox"/> Temperature	<input type="checkbox"/> Variance of selected channel			<input type="checkbox"/> Analogue	
<input type="checkbox"/> Pressure <input type="checkbox"/> Metres?		<input type="checkbox"/> Marked Events		<input type="checkbox"/> Analogue smoothed	
<input type="checkbox"/> Pressure sm'd	<input type="checkbox"/> Altitude (pressure)?			<input type="checkbox"/> Activity	
				<input type="checkbox"/> Side time bars	

Close application Hide to system tray Show 3D

Slide II Capture 3D Invert cols Save 3D Img Keep on top

Figure 2.15 Overlay display options

Other display options

The **2D Display / Display options** tab (Figure 2.14) has some additional functions that allow the user to customise their view in the main graphing windows.

The **View incrementing hours / minutes / seconds** tick-boxes can be selected to create a shaded area behind the data to illustrate the changing time. These synchronise with the top indicator bars. In the example below, **Seconds** has been ticked and each incremental unit of shading from black to green represents one minute (Figure 2.17).

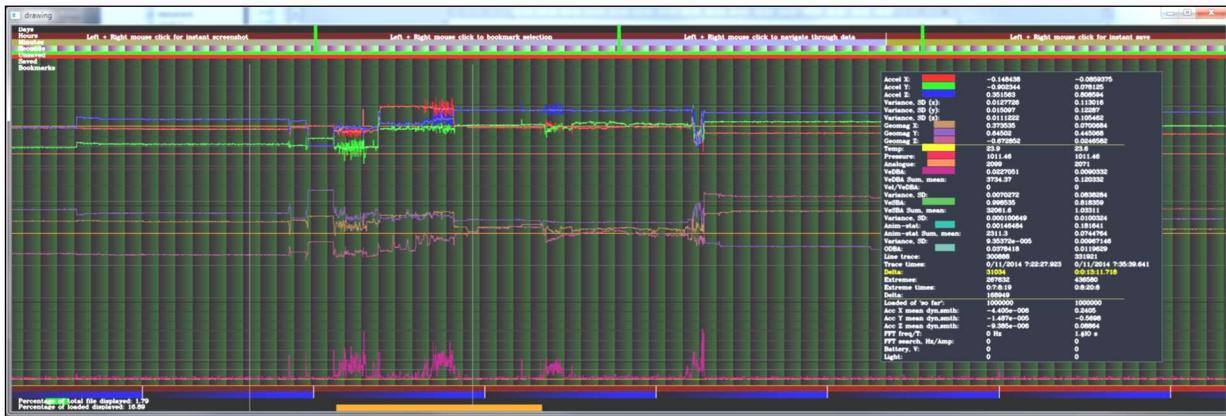


Figure 2.16 Background time progression indications of seconds, minutes, or hours

Display “Help text” on 2D brings up text in the top left corner of the drawing that describes various mouse/key combinations to navigate through the data (Figure 2.17).

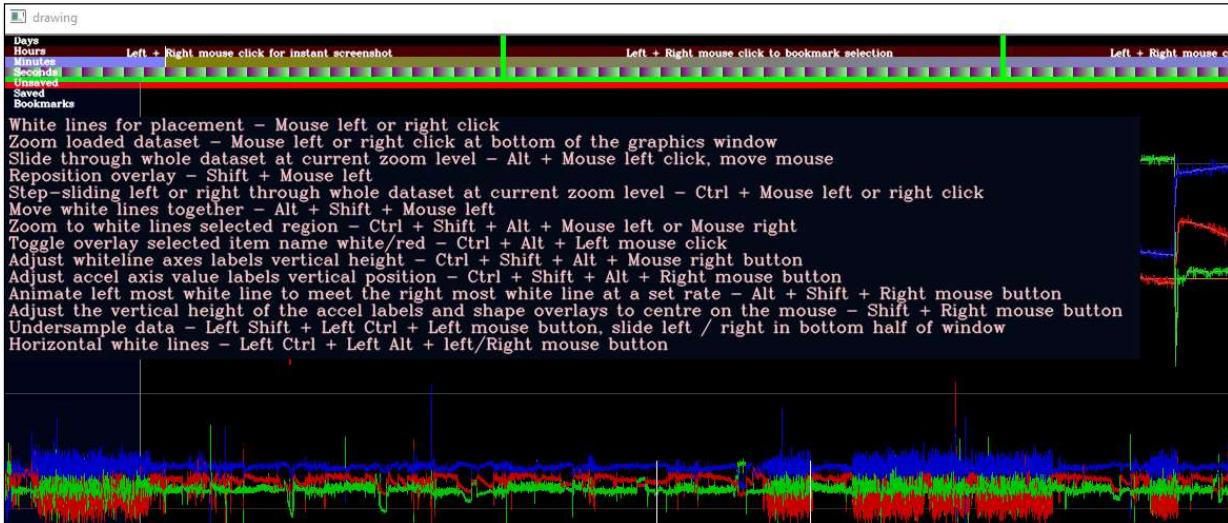


Figure 2.17 Displaying the controls help text

Display 2D traces is ticked by default and tells the application to draw the graphs. If unticked, the graphs are not drawn, and the drawing will appear blank. This is useful when in 3D mode as it lessens the load on the CPU, allowing it to devote more power to the 3D environment. Note that the neighbouring column allows the user to remove specific data channels (acceleration, geomagnetism, etc) from the drawing (see red text) or add in a measure of activity over a customisable time period.

**Fill white-lines, opacity**

(range 0 to 1.0) changes the opacity of the shaded area between the two white lines. This value is 0.15 by default but can be adjusted in the adjacent text box (Figure 2.12).

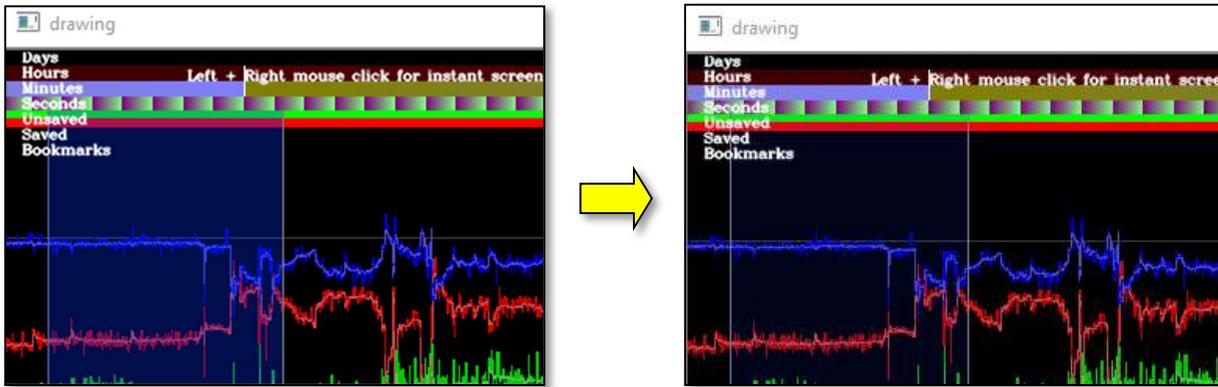


Figure 2.18 Adjusting the opacity of the rectangular fill between the left/right white lines

White background?

changes the background colour of the drawing window from black to white.

Gridlines redirects to **2D Sliders/Gridlines** (discussed earlier in this chapter).

2D trace line thickness

increases the line thickness of the graphed data on the 2D graphing window. This is often useful when using screenshots from DDMT in presentations.

Control panel transparency, %

allows the user to change the transparency of the above control panel; particularly useful when a small-screened laptop is in use. The maximum is limited to 60% so there will always be 'some' visibility.

If the limits of accelerometry or magnetometry are required, this can be enabled by heading to the **Accel/mag detail/offsets** tab, then **Accel/Mag trace limits** sub-tab. Once the **View max/min limits** box is ticked any data between the white lines (left and right click) will be scanned to determine the maximum and minimum values for the acceleration and magnetism channels (Figure 2.19).



Figure 2.19 Determining maximum and minimum acceleration and magnetism values between the white lines



Adjustments to the visual vertical amplitude / offsets / controls of 2D graphs

By selecting the **2D Sliders** tab, you are presented with a tab **Slider controls grouped by purpose** which represents the various sensor channels from the Daily Diary. A control on this tab allows the user to switch between the different sensors and alter the visual amplitude and vertical offset of graphed data within the drawing window, along with controlling the level of smoothing for these channels. The area highlighted in Figure 2.20 will switch by action of the users selection of the combo on the left side to show the relevant controls.

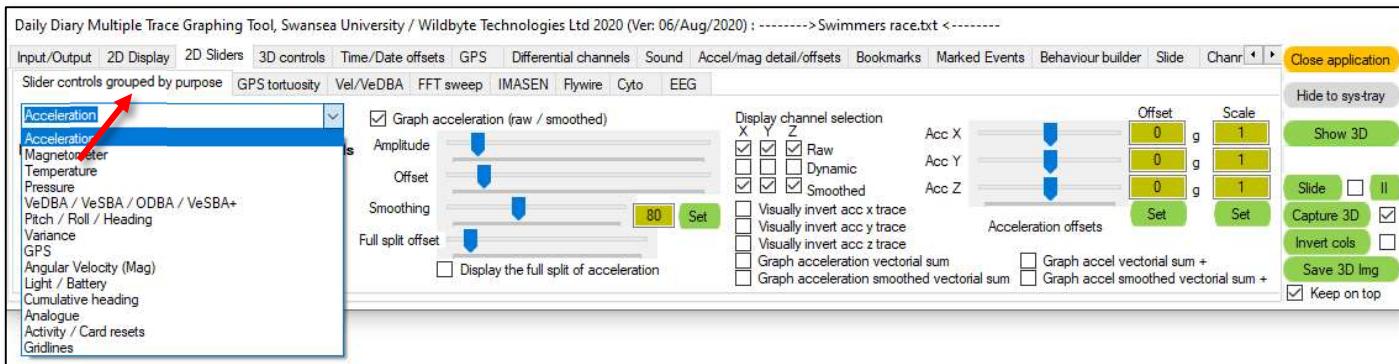


Figure 2.20 Controls list for the various sensors and other metrics

The 14 sets of controls represent **Acceleration**, **Magnetometry**, **Temperature**, **Pressure**, **VeDBA / VeSBA**, **ODBA**, **VeSBA+**, **Pitch / Roll / Heading**, **Vairance**, **GPS**, **Angular Velocity (Mag)**, **Light / Battery**, **Cumulative Heading**, **Analogue**, **Activity / Card resets**, and **Gridlines**.

A discussion of each of these 14 sets of controls is discussed below.

→ Acceleration controls

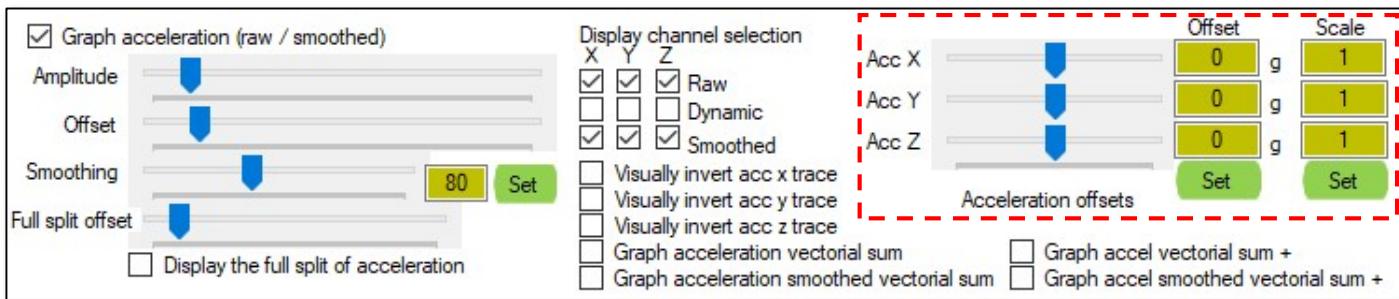


Figure 2.21 Acceleration controls

**(Checkbox) Display the full split of acceleration**

Switches on a display of the full split's acceleration, along with markers to show the section of data currently being viewed at zoom.

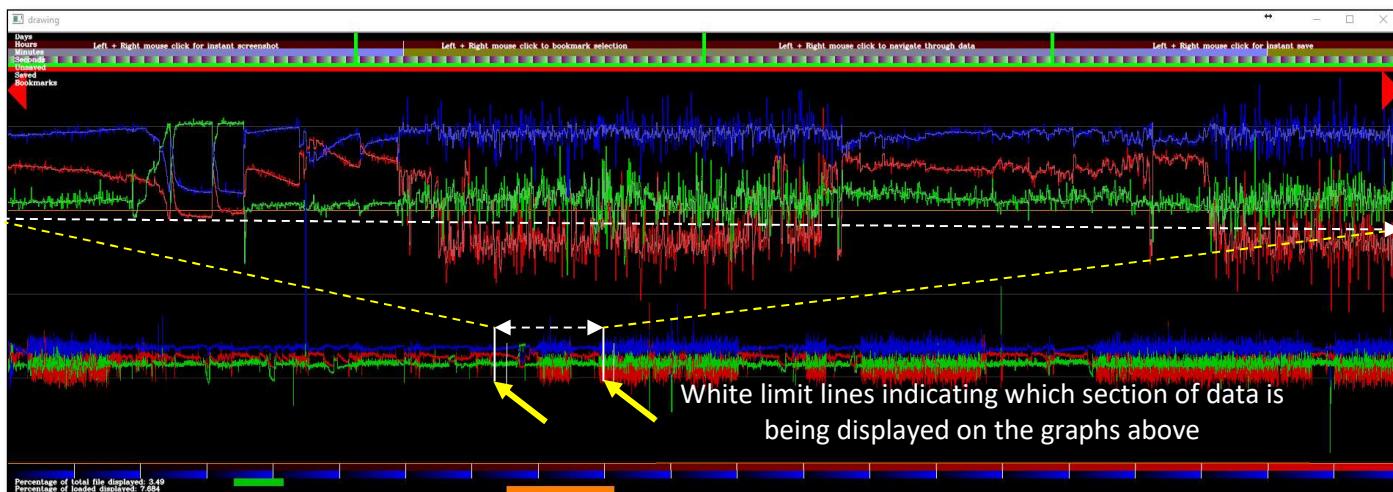


Figure 2.22 Displaying the full acceleration split with the zoomed portion

- Amplitude, and Offset** adjusts the amplitude and offset of the 2D graphed acceleration data
Smoothing adjusts the centre averaged smoothing window of the acceleration data
Full Split Offset adjusts the vertical offset of the optional Full split (selectable below)
Graph acceleration (raw / smoothed) enables/disables the display of the acceleration traces. Fine control of which axes are displayed is achieved through the checkboxes in the centre
Graph acc sm amp enables/disables the display of the smoothed acceleration amplitude
(This is the square root of the sum of the squares of the smoothed accel X/Y/Z traces. At rest, this would be 1.0 g)
Visually invert acc x/y/z traces visually invert the vertical trace of the acceleration traces

Graph acceleration (smoothed) vectorial sum
display the graph of $\sqrt{Acc_x^2 + Acc_y^2 + Acc_z^2}$

Graph accel (smoothed) vectorial sum +
display the graph of $1 - \sqrt{Acc_x^2 + Acc_y^2 + Acc_z^2} + 1$

The controls in the area highlighted in a dashed red box – offsets and scale are discussed in Chapter 3.

→ Magnetometry controls

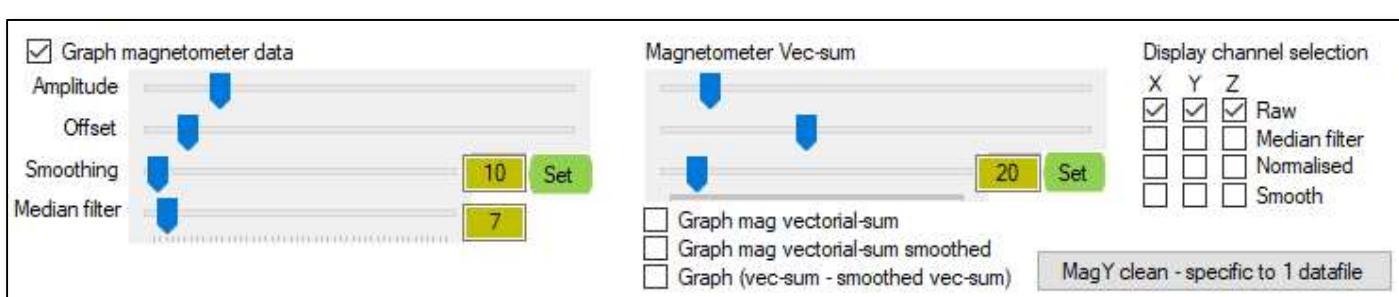


Figure 2.23 Magnetometry controls



Amplitude, and Offset	adjusts the amplitude and offset of the 2D graphed magnetometry (also mag vec-sum) data
Smoothing	adjusts the centre averaged smoothing window of the magnetometry data
Graph magnetometry	enables/disables the display of the magnetometry traces. Fine control of which axes are displayed is achieved through the checkboxes on the right side
Graph mag (smoothed) vectorial sum	display the graph of $\sqrt{Mag_x^2 + Mag_y^2 + Mag_z^2}$
Graph (vec-sum - smoothed vec-sum)	display the graph of $\sqrt{(Mag_x^2 - MagSm_x^2) + (Mag_y^2 - MagSm_y^2) + (Mag_z^2 - MagSm_z^2)}$

Normalised channels are simply the individual X, Y, and Z mag channels divided by the magnitude of the XYZ vector

→ Temperature controls

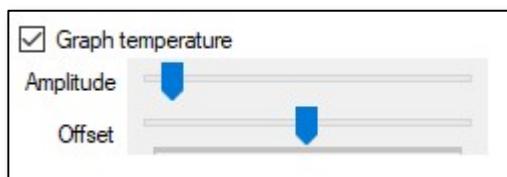


Figure 2.24 Temperature controls

Graph temperature	enables/disables the display of the temperature traces
Amplitude	slider adjusts the amplitude of the graphed temperature data
Offset	slider adjusts the vertical offset of the graphed temperature data

→ Pressure controls



Figure 2.25 Pressure controls

Amplitude	adjusts the amplitude of the graphed pressure data
Offset	adjusts the vertical offset of the graphed pressure data
Smoothing	adjusts the centre averaged smoothing window of the pressure data
Median pressure	enables/disables display of the median filtered pressure data
Smooth pressure	enables/disables display of the smoothed pressure data (the smoothing value is individually set on the Channel Smoothing tab)
Pressure	enables/disables display of the pressure data
Invert visually	enables/disables inversion of any of the above pressure data, to see dives more clearly.
Altitude, metres	adds an altitude channel to the drawing window (default colour is light blue)
Ground pressure	used to set the air pressure at sea level to enable the calculation of altitude directly from barometric pressure sensor readings
Altitude/Press delta step	is used to adjust the altitude and pressure stepping rate



→ Metrics – VeDBA, VeSBA, ODBA, and VeSBA+ controls

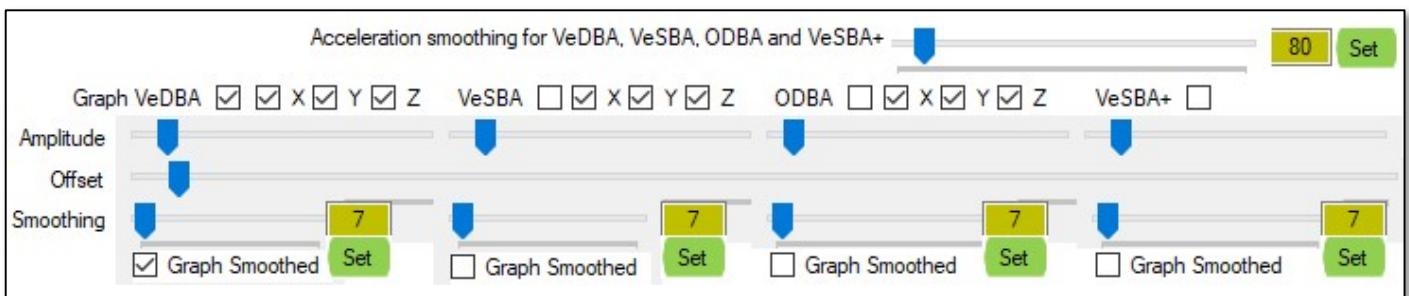


Figure 2.26 Metrics controls

Amplitude

adjusts the amplitude of the graphed metrics data

Offset

adjusts the vertical offset of the graphed metrics data

Smoothing

adjusts the centre averaged smoothing window of the metrics data after they are calculated i.e. VeDBA smoothed, VeSBA smoothed etc.

4 different graphs are controlled here; VeDBA, VeSBA, ODBA, and VeSBA+. These are 4 metrics derived from the 3 raw acceleration channels, and smoothed acceleration – note this is its own level of smoothed acceleration (top slider in Figure 3.27), not the level of acceleration defined on the acceleration controls panel. There is only the one level of smoothing of the raw acceleration channels for all 4 of these metrics.

There are also optional X, Y, and Z tickboxes to allow the first 3 metrics constitute of either 2 or 3 channels. Note that deselecting more than 2 channels will invoke the auto selection of the other 7 unticked channel.

VeDBA (vector of the dynamic body acceleration):

$$VeDBA = \sqrt{((AccX - AccX_{sm})^2 + (AccY - AccY_{sm})^2 + (AccZ - AccZ_{sm})^2)}$$

VeSBA (vector of the static body acceleration):

$$VeSBA = \sqrt{((AccX_{sm})^2 + (AccY_{sm})^2 + (AccZ_{sm})^2)}$$

ODBA (Overall dynamic body acceleration):

$$VeDBA = \sqrt{(AccX - AccX_{sm})^2} + \sqrt{(AccY - AccY_{sm})^2} + \sqrt{(AccZ - AccZ_{sm})^2}$$

VeSBA+ (vector of the static body acceleration):

$$VeSBA+ = (VeSBA \geq 1) ? VeSBA : 2 - VeSBA$$

i.e. if VeSBA is ≥ 1 then $VeSBA+ = VeSBA$, else $VeSBA+ = 2 - VeSBA$

All of these 8 channels (the 4 basic metrics, and their smoothed counterparts) are then available in the **Behaviour Builder** function, and also as axes on the 3D visualisation part of the software.

Each of the 4 primary graphs can be toggled on/off on the 2D drawing display by ticking the **Graph VeDBA**, **Graph VeSBA** etc. tickboxes, while their smoothed channels can also be displayed by ticking the **Graph Smoothed** tickboxes at the bottom of each of the 4 sets of controls.



→ Metrics – Pitch / Roll / Heading controls

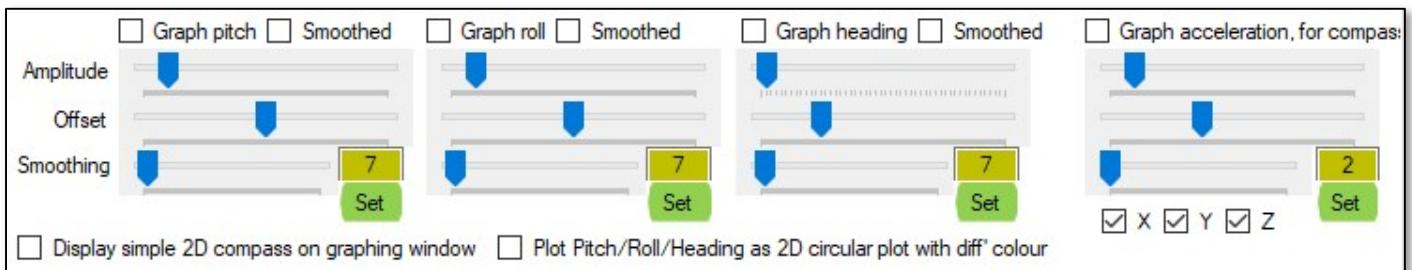


Figure 2.27 Metrics controls

With both triaxial acceleration and triaxial magnetometry channels present, it is possible to calculate heading, pitch, and roll (HPR) of the logger. When the logging device is tilted with respect to the ground, HPR will be miscalculated if this tilt were not considered. Therefore, the acceleration channels are utilised in a well known algorithm known as the **Tilt compensated compass** whereby the pitch and roll of the logger are determined, and the magnetometry channels are rotated back to the horizontal plane, providing a non-skewed measure of HPR. This is discussed later in Chapter 8.

Amplitude

adjusts the amplitude of the graphed pitch, roll, or heading data

Offset

adjusts the vertical offset of the graphed pitch, roll, or heading data

Smoothing

adjusts the centre averaged smoothing window of the pitch, roll, or heading data after they are calculated i.e. pitch, roll, or heading smoothed, pitch, roll, or heading smoothed

The first 3 sets of controls relate to pitch, roll, and heading respectively. Each can be post-smoothed, and either the raw channel, or its smoothed counterpart can be graphed individually.

The 4th set of controls allows the user to view the acceleration used by the tilt-compensated compass algorithm. This is to allow the user to ensure the acceleration channels are not being under/over smoothed, and to view the degree of “roughness” of the trace. The third slider on this 4th set of controls, **Smoothing**, defines the level of smoothing of the raw acceleration channels. It is the smoothed acceleration channel that is then passed onto the **Tilt compensated compass** algorithm where pitch, roll, and heading are derived. The 3 tickboxes for X, Y, and Z are simply to allow the user to switch on/off whichever channels are of interest on the **2D graphing window**.

Display simple 2D compass on graphing window (checkbox)

This displays a simple compass that displays the heading (both raw and mean) at the event selected with the left white line – note that heading is already derived from both smoothed acceleration and smoothed magnetometry channels:

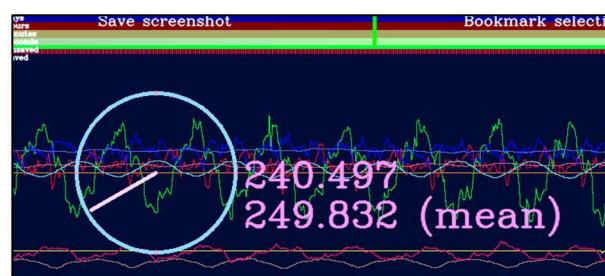
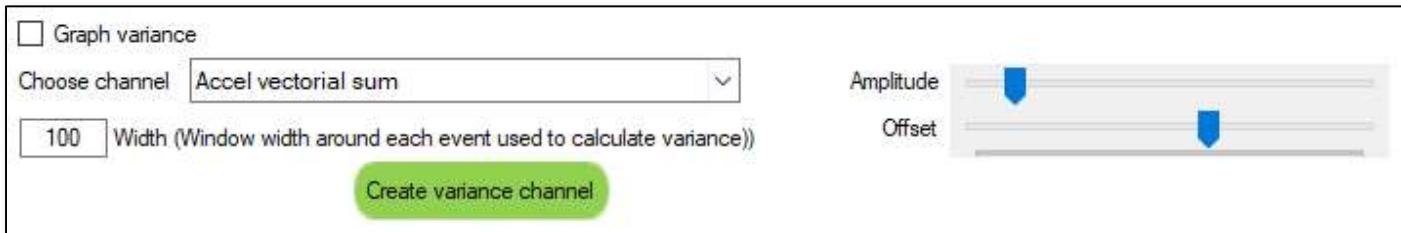


Figure 2.28 Simple compass overlay showing heading as defined at the left white line

Note - By first turning off the values overlay with the middle mouse button, one can hold down left shift and left mouse click to drag this small compass around the **2D graphing window**.

**Plot Pitch/Roll/Heading as 2D circular plot with diff' colour (checkbox)**


This switches the plots on the **2D graphing window** from linear plots of these 3 metrics, to the cosine of. Also, the colour is switched when either positive or negative for pitch and roll, and when either 0-180, >180 with heading.

→ Metric – Variance controls**Figure 2.29 Variance controls**

In order for variance to be used as a channel in the Behaviour Builder, one needs to determine the limits of a particular channel of interest. What this does here is to create a channel named Variance of any other data channel in the list. This Variance channel is calculated based on a centred window, with each point being the variance of a window of data around it. The user can choose the data channel and the window width, finally clicking **Create variance channel** to create that channel.

- Amplitude** adjusts the amplitude of the graphed variance data
Offset adjusts the vertical offset of the graphed variance data

This does not create a channel that can be used anywhere else. It is merely to allow the user to investigate likely limits by observing values present on the display over (enable Variance on the values overlay by opting it on the **2D Display / Overlay selection** tab:

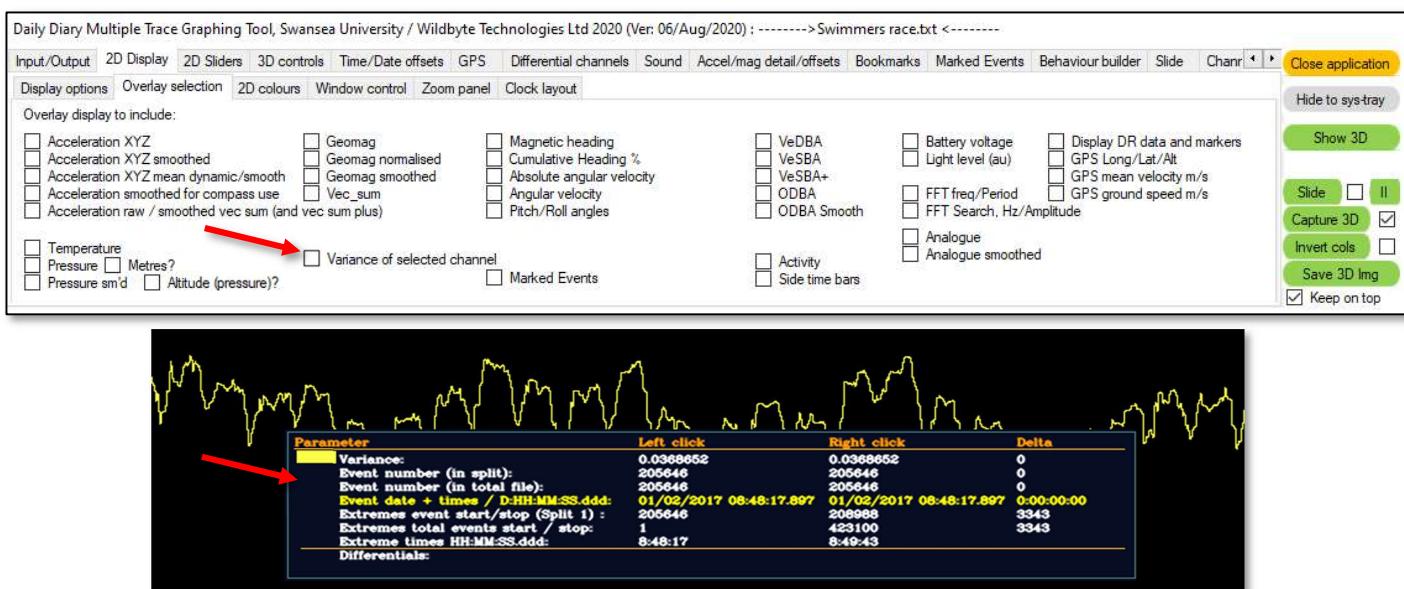**Figure 2.30 Variance on the Values Overlay**



→ GPS controls

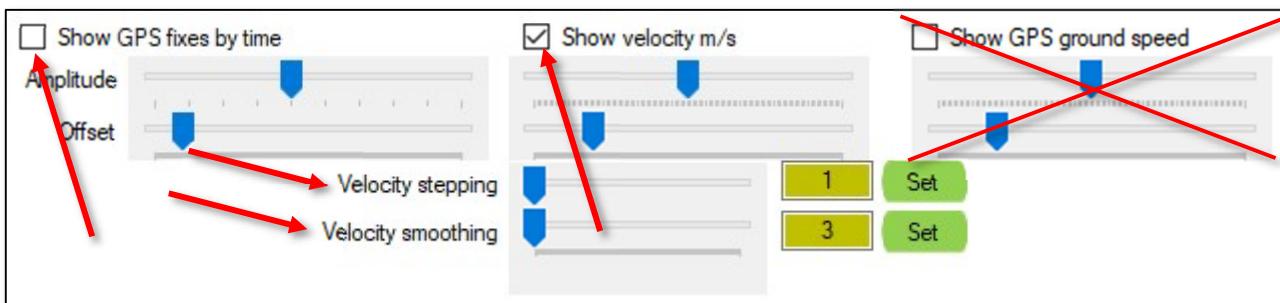


Figure 2.31 GPS fixes and velocity controls

For GPS, markers can be shown on the **2D graphing window** showing where a GPS fix synchronises with a time point in the logger's data.

Amplitude adjusts the amplitude of the graphed variance data

Offset adjusts the vertical offset of the graphed variance data

Velocity stepping adjusts the GPS-GPS counter for determining velocity

Velocity smoothing adjusts the centred average smoothing of the resulting velocity graph



Figure 2.32 Visualising where GPS synchronises with the logger's data

Where a vertical blue line exists, a GPS fix exists in memory with that date and time (to the integer second) value. There are multiple blue lines per GPS fix as acceleration loggers typically log at > 1 Hz. The red dot that precedes each stretch of vertical blue lines represents the first point for a given GPS fix i.e. once per second at 1 Hz GPS.

Secondly, **Show velocity m/s** can be ticked to overlay a graph, Figure 2.33.



Figure 2.33 Visualising GPS velocity (m/s) and mean velocity



→ Metrics - Angular /Absolute angular velocity controls

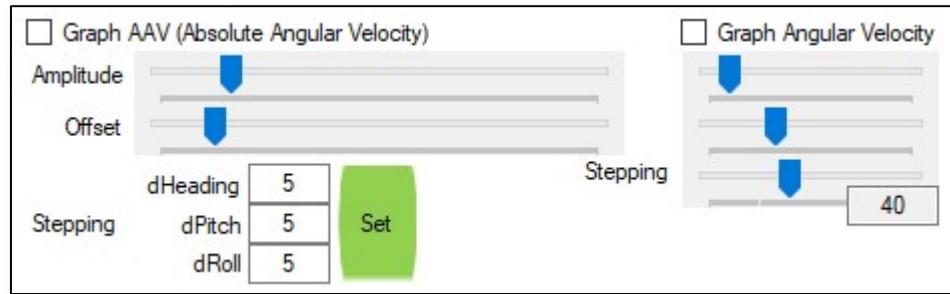


Figure 2.34 Absolute angular velocity controls

Absolute angular velocity is a metric that combines rate of change of 3 channels, heading, pitch, and roll. The rate of change is separately adjustable for each of the 3, each calculated centred around each event.

For example, if *dHeading*, *dPitch*, and *dRoll*, were 3, 5, and 7 respectively, then when the channel is created, event 1000 would be calculated based on the difference of data at events 999 and 1001 for *dHeading*, 998 and 1002 for *dPitch*, and 997 and 1003 for *dRoll*.

The metric is calculated as:

$$AAV = \sqrt{(dHeading_1 - dHeading_2)^2 + (dPitch_1 - dPitch_2)^2 + (dRoll_1 - dRoll_2)^2}$$

The value is stored in degrees.

Angular velocity is a metric that determines the rate of change across the sphere of heading vs pitch. This sphere is the same sphere also described as the **Orientation sphere**, discussed in the 3D visualisation section later in this manual.

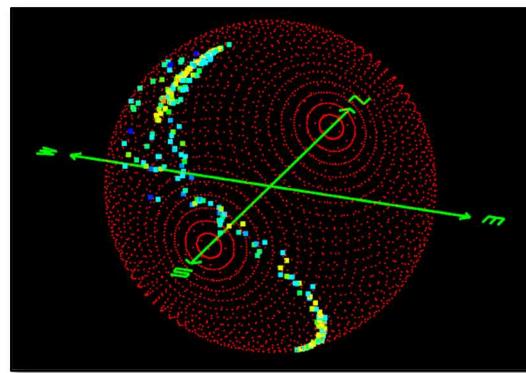


Figure 2.35 Angular velocity controls

A point is located on this surface based first on its heading value, to determine the position around the “equator”, and then rotated vertically up (down) for pitch positive (negative).

The value is stored in degrees. For this alternative angular velocity measure, the user can also adjust the stepping using its 3rd slider. Note that this calculation is based on the data point at the current event and the events *n* steps behind i.e. not a “centred measure”.



For both controls,

- Amplitude** adjusts the amplitude of the graphed angular velocity data
Offset adjusts the vertical offset of the graphed angular velocity data

→ Light / Battery

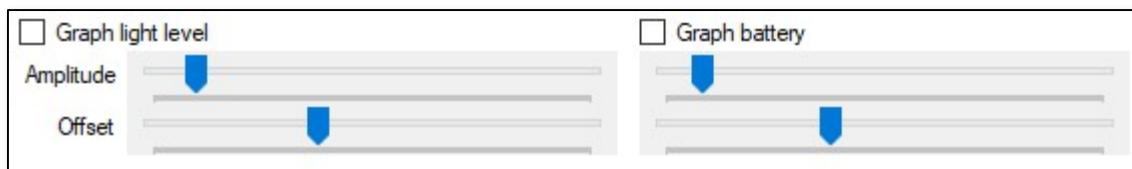


Figure 2.36 Battery and light level controls

These controls simply allow the viewing of battery voltage and light level (the square Daily Diary with connectors).

- Amplitude** adjusts the amplitude of the graphed battery voltage / light level data
Offset adjusts the vertical offset of the graphed battery voltage / light level data

→ Metric - Cumulative heading

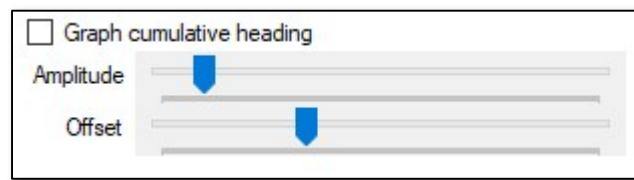


Figure 2.37 Cumulative heading controls

Cumulative heading is still under development.

- Amplitude** adjusts the amplitude of the graphed cumulative heading data
Offset adjusts the vertical offset of the graphed cumulative heading data

→ Analogue

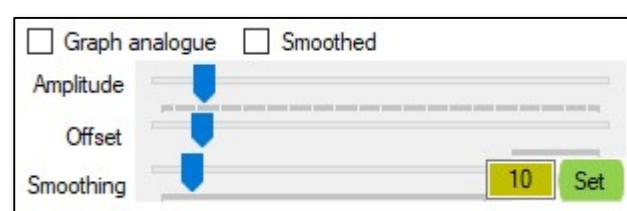


Figure 2.38 Analogue controls

Analogue data, available with some types of loggers.

- Amplitude** adjusts the amplitude of the graphed analogue data
Offset adjusts the vertical offset of the graphed analogue data
Smoothing adjusts the centre averaged smoothing window of the analogue channel



→ Activity / SD card resets



Figure 2.39 Activity / SD card resets controls

Activity is simply the sum of VeDBA per unit time. The user can define the number of seconds over which this is calculated.

The SD card reset line colour is by default green, and when a reset has occurred where power was momentarily lost to the SD card and the logger had paused collecting data (but the real-time clock had continued to increment), there will be a vertical red line. At this point, the user will notice a time jump/step where the logger resumed collecting/storing data to the SD card.

- Amplitude** adjusts the amplitude of the graphed activity data
Offset adjusts the vertical offset of the graphed SD card reset data

→ Gridlines (for acceleration, geomagnetism, and metrics)

The main graphing windows can show graduated lines, customisable by the user on the **2D Sliders / Gridlines** sub-tab.

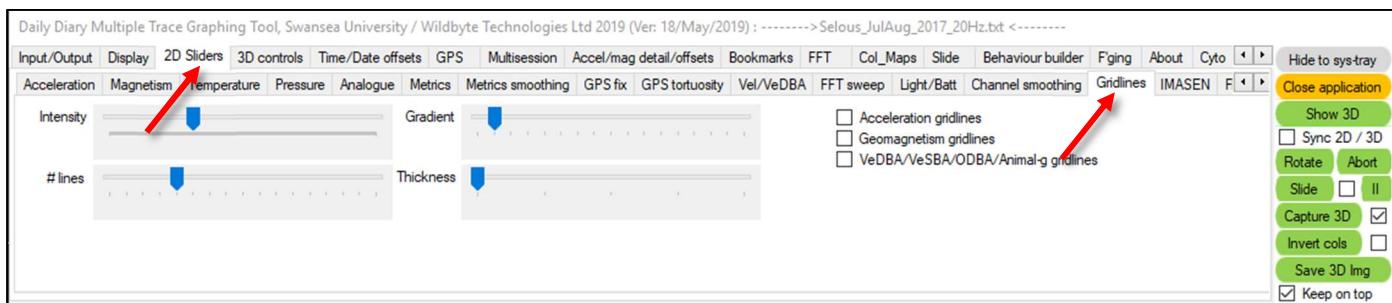


Figure 2.40 The Gridlines menu

- Intensity** adjusts the intensity of the gridlines displayed around the Acceleration, Geomagnetism, and VeDBA traces
lines adjusts the number of positive and negative gridlines visible on the graph. For Acceleration and Geomagnetism traces, the gridlines are +/- 1, +/- 2, +/- 3 g or gauss etc. For the VeDBA trace, the gridlines are spaced positively at 0.2, 0.4, 0.6 etc
Gradient determines how quickly the intensity of the gridlines falls within increasing absolute range value
Thickness determines how thick the drawn lines are. This is useful in presentations etc. to help illustrate the amplitude of your dataset at various time points
Acceleration gridlines shows acceleration gridlines when selected
Geomagnetism gridlines shows geomagnetism when selected
VeDBA/VeSBA/ODBA/Animal-g gridlines shows metric gridlines when selected

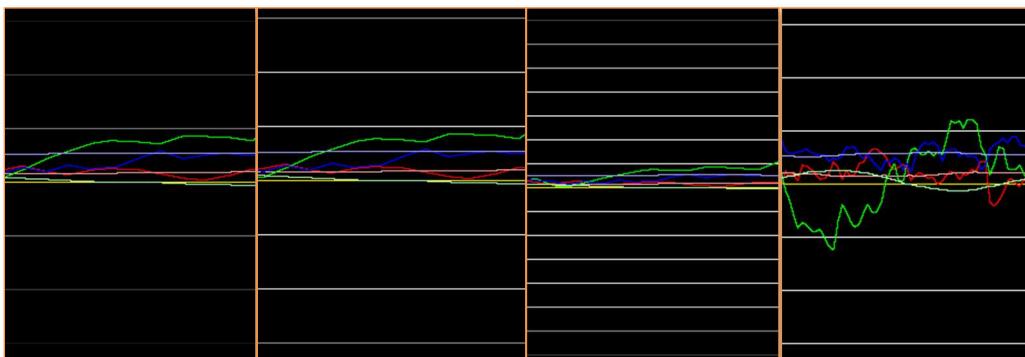


Figure 2.41 The effects of dim, bright, gradient and flat gridlines for the accelerometry trace

If the gridlines do not display even when the checkboxes outlined above on the **Gridlines** sub-tab are selected, ensure that the separate **Gridlines** checkbox on the **Display/Display options** sub-tab has also been selected.

→ Bitwise

The main graphing windows can show the **Bitwise** channel split up into multiple layers. The channel is 64 bits wide, meaning it could potentially hold 64 layers of '1' or '0' per data point, but currently only 16 bits are used. **Marked Events** (the presence of) can be copied into any one of the first 16 **Layers** of the **Bitwise** channel. This channel is global, and so moving to another split and returning will preserve any data within the 16 **Layers**. See chapter 5 **Marked Events and Bitwise Layers** for details.

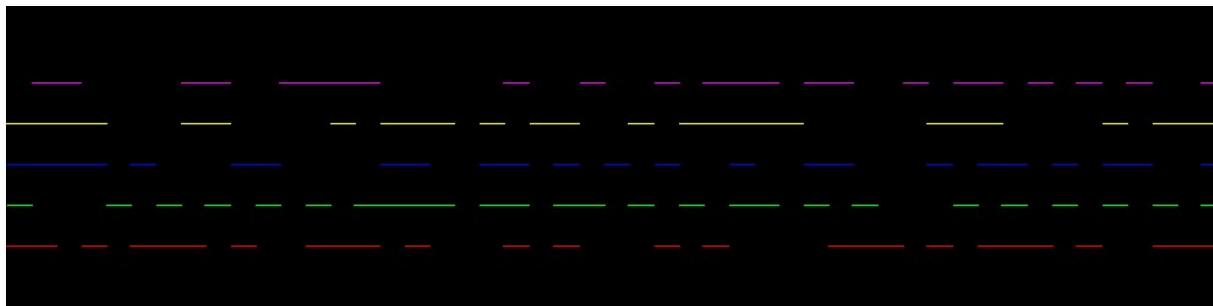
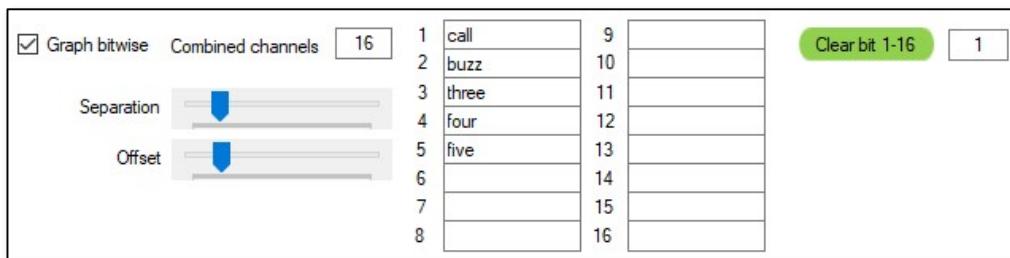


Figure 2.42 Bitwise channel controls

Offset

Separation

Combined channels

Clear bit 1-16

Labels

adjusts the vertical offset of the different bits of the **Bitwise** channel
adjusts the relative separation of the different bits of the **Bitwise** channel
limits how many channels (from #1) are displayed on the **2D graphing window**
(button) will clear (globally / all splits) the chosen bit number, unless the selected bit number is part of a data file (some data file types have **Bitwise** channels embedded within them and therefore cannot be overwritten by clearing, or by copying **Marked Events** to them)

On the **2D graphing windows**, if **Bitwise** is selected on the **Overlay selection** tab, then any bits where the left or right mouse clicks, relevant labels will be displayed on the overlay



Figure 2.43 Bitwise channels on the overlay. Notice that the left and right click labels are different as they highlight different bits on the graph (left = 1 (red)+2 (green)+3(blue)+5(pink), right = 2(green)+5(pink))

The Zoom panel – a quick localised zoom of the data currently in-view

There is a zoom panel available within the graphing window for closer inspection of data. To toggle this panel, hold **Ctrl**, **Shift**, and **Alt** together and right click the mouse. Move the zoom panel across the graphing window by holding down the **left shift** key and dragging with a **right** mouse click. Use the **left Alt** key and mouse wheel to adjust the data width (the zoom level) displayed in the zoom panel.

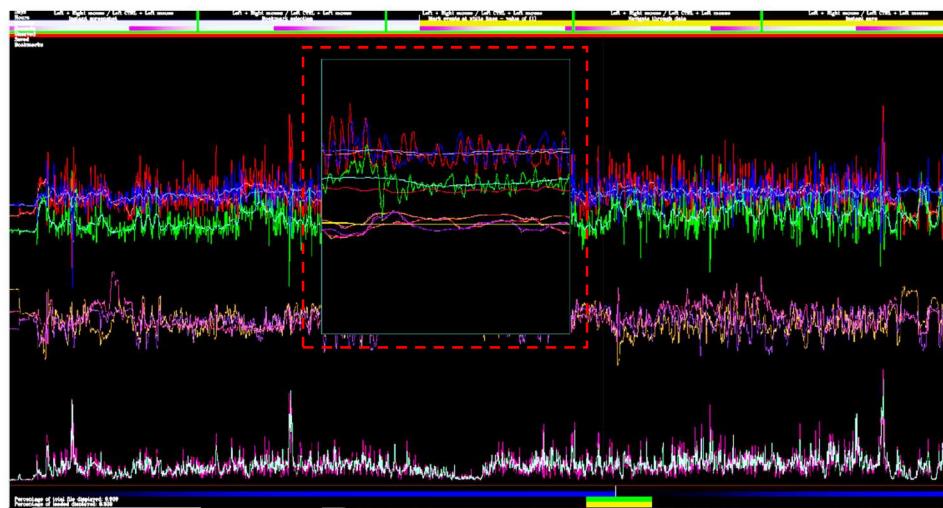


Figure 2.44 Main graphing window with zoom panel enabled

Currently only acceleration (+smoothed), magnetometry (+smoothed), temperature, and pressure can be displayed on the zoom panel. The controls can be found in Figure 2.45 below. These are primarily offsets for the individual data channels, along with a data width slider (also adjustable with key and mouse combination as described above).

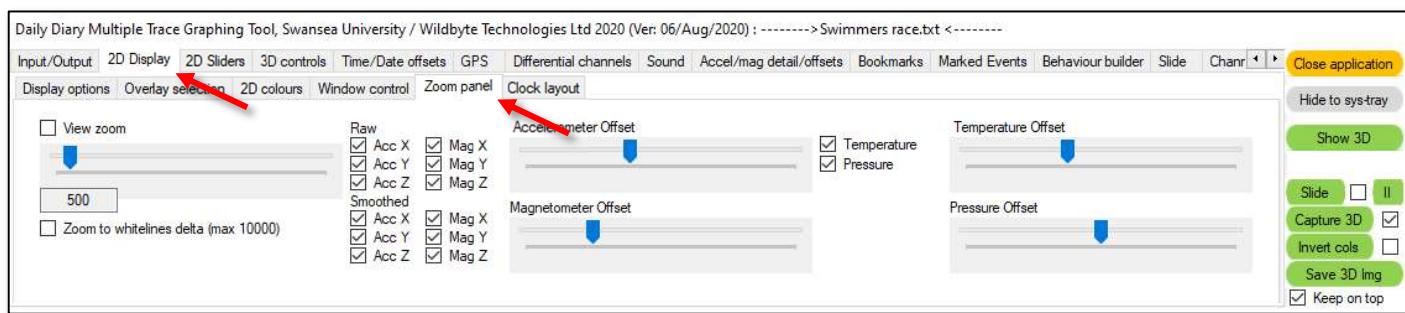


Figure 2.45 The zoom panel controls



Animal image overlay to show animal orientation with respect to surge/sway/heave

To gain a better understanding of animal orientation from the data, the programme allows the user to place three animal drawings down the left and right white lines symbolising the surge, sway and heave acceleration axes. These are automatically oriented as a function of the smoothed acceleration data channels. By using the smoothed channels, we are effectively obtaining the static acceleration as most dynamic acceleration spikes are smoothed out.

Under the **Accel/mag detail/offsets/Animals 2D** tab you will find the controls for the animal drawings. To make them appear the **X**, **Y**, **Z** check boxes on the left of the coloured squares need to be ticked.

The axis labels can be renamed from their default “Surge”, “Sway”, and “Heave”. These labels depend on the orientation of the tag on the animal, so is left to the user to determine the appropriate labelling. Similarly, the check boxes to the right of the coloured squares invert the images, which can also be useful depending on the position of the tag on the animal. The **Axes offsets** (coloured boxes, see Figure 2.46) allow the images to be rotated into the correct resting orientation.

By default, both animal images and angle values are displayed. The angle is a simple *arcsine* function of the value from each of the accelerometer traces (if smoothed, normally within the range -1 to +1; values greater than this are clipped at these limits). This feature is particularly useful for talks/demonstrations.

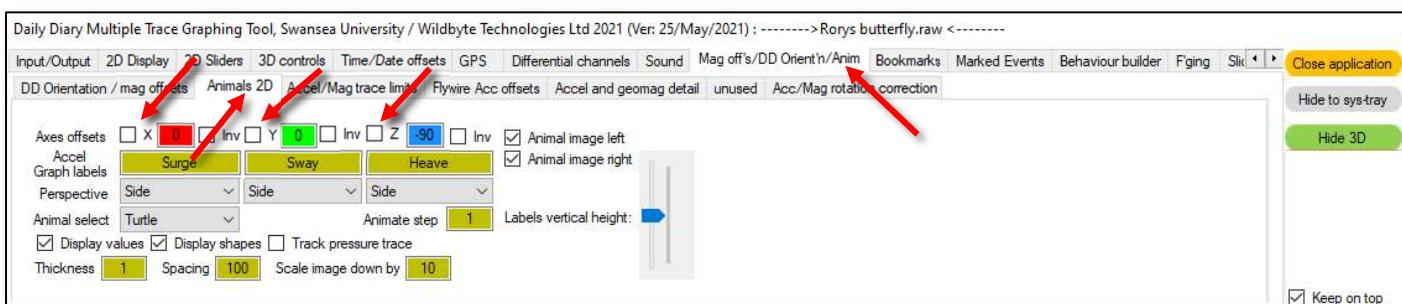


Figure 2.46 Accel/mag detail/offsets/Animals 2D tab

The animal images can be adjusted in size by altering the **Scale down factor** (default value of 10). Larger values shrink the size of the animal images, while smaller values increase their size. When the **Scale down factor** is changed, the image sizes will be updated as soon as the mouse moves over the graphic window. The vertical spacing between the images can also be adjusted; necessary if the images are made significantly larger (default spacing of 100 pixels).

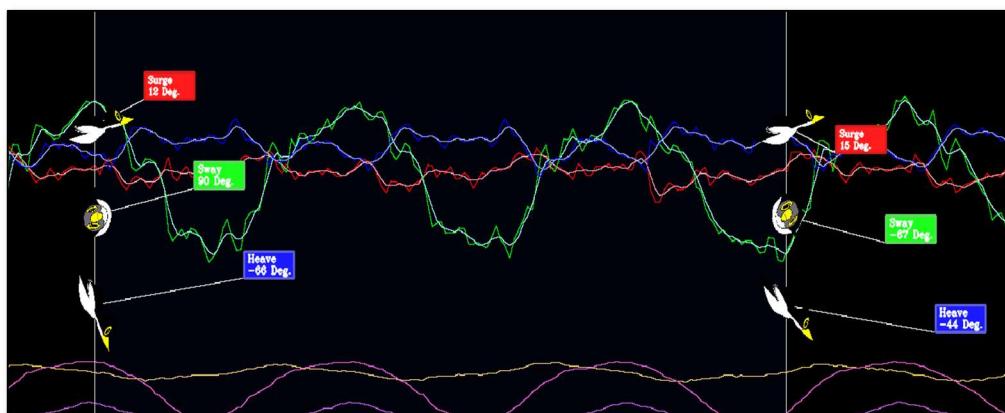


Figure 2.47 Accel/mag detail/offsets/Animals 2D tab, set

Depending on the size of the graphic window, the user might find that one of the animal traces falls below the bottom of the screen. To raise the animals' centres up the screen, hold down **Left shift** and hold down the right mouse button. Mouse movement then determines the vertical position of these images.



Slide function for animation

For discussions/presentations, it is sometimes useful to have the animal and the white line move across the screen with changes to animal posture automatically animated as it progresses. To do this, select start and stop positions with left and right mouse clicks to set the extremes. Then hold down **Left Shift** and **Left Alt** and click the right mouse button in the graphic area. The left white line will now move to the right towards the right white line until they meet and stop. The animal image(s) will automatically assume the orientation determined by the smoothed acceleration channels as the white line travels from left to right. By default, the white line will step one “event” at a time. For a large amount of data points between the two white lines means this might take some time and so an **Animate step** variable is provided within the window, allowing it to jump X events with every frame update (default value of 1).

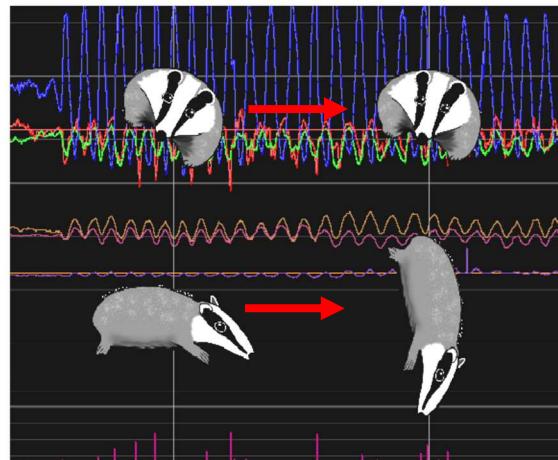


Figure 2.48 Animate function automatically drags the left white line over to meet the right white line, causing the images to rotate according to the smoothed acceleration data at each step

Additionally, if pressure data is a good variable for the animal such as with condors (changes in barometric pressure during circling), or turtles (depth sensor pressure), the animation can also be made to track vertically up and down along the pressure trace by ticking the **Track pressure trace** check box.

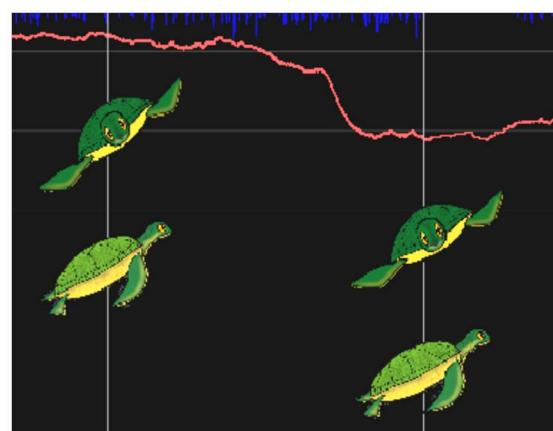


Figure 2.49 Track pressure trace function with the animation



3. Time / date, and acceleration / geomagnetism data corrections

Daily Diary on-board Real Time Clock – Setting the start time for your dataset

As discussed in the **System Overview**, the Daily Diary unit possesses an on-board Real-Time Clock (RTC). When data is stored to the logger's memory, a time stamp, or partial stamp, is stored with that single event, to aid in providing the best timing accuracy knowledge. There are several factors to consider when managing time data in DDMT.

The **Input/output** tab will show the start and finish times of the entire dataset in **Days : Hours : Minutes : Seconds** (Figure 2.1). Note that as the Daily Diary doesn't have a link to actual time i.e. GMT, this is a relative value. The user could have perhaps set the clock within the command string on the memory card, such that upon power up, it assumes a specific time value. Alternatively, DDMT allows the user to input an offset to the time and date attributes to correct the timing aspect of the data. Note that for some data files, this information will either not be displayed, or could be incorrect.

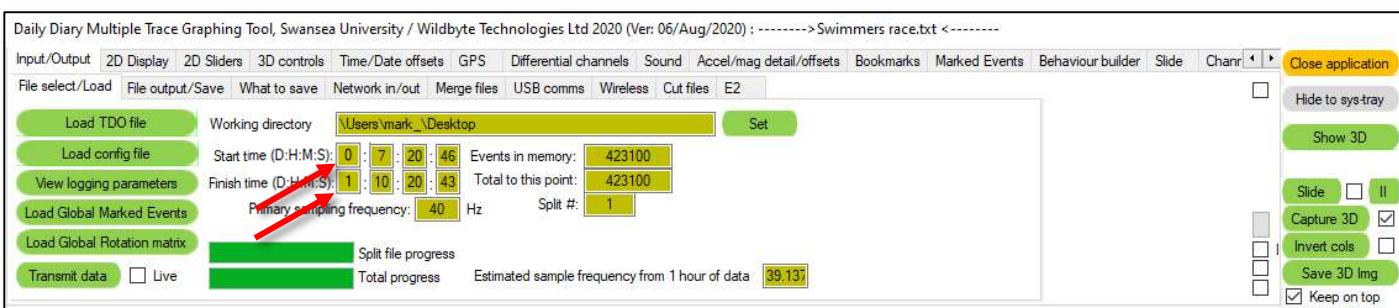


Figure 3.1 The start and finish times of the source text file (DD:HH:MM:SS)

To manually set the correct time for your data click on the **Time / Date offsets** tab along the top of the control panel. The **Daily Diary settings** tab will automatically open and display the screen in Figure 2.2. Knowing your deployment details is crucial to set the correct time and date for your data. Select an event that began at a known time to the second (usually this is the time the Daily Diary was turned on or the start time of device calibration), ensuring that the left white line (displayed when left cursor is clicked on the main graph window) matches this event exactly. Set the date by entering the day / month / year in the top three boxes on the right-hand side of the window (circled). When the correct date is entered, press the bright green **Set start date on Daily Diary data** button on the left-hand side to apply date corrections to the whole file. Next, enter the time of the selected event in the three green boxes under the date using hour / minute / second format. When the correct time has been entered, press the bright green **Set time for DD data set from white line** button to the left-hand side. Make sure the time is correct, as once set the time cannot be changed. If set incorrectly, the user will need to close and reopen the source text file.

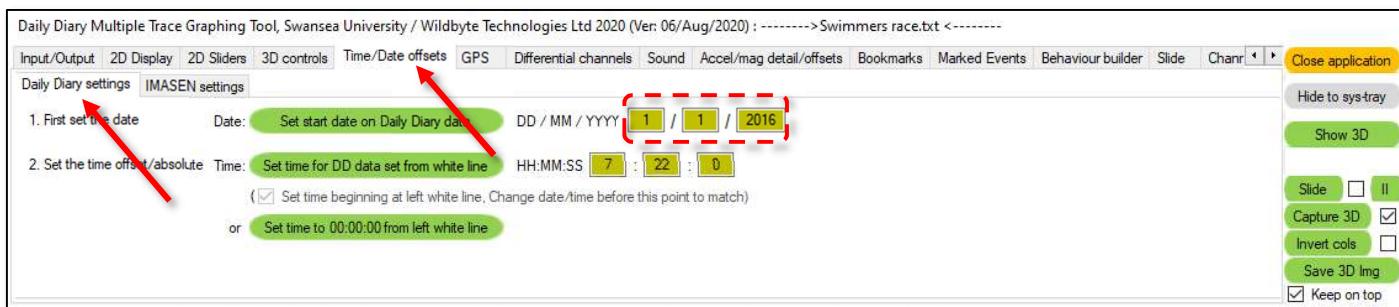


Figure 3.2 The screen displayed when opening the Time/Date offsets tab

Decimal seconds are created by the software when the data file is first analysed / loaded into memory. For 40 Hz data, there will be approximately 40 events per second with the same "seconds" value and so a decimal value is assigned to every data point with the same time / date (to the second). Note that due to the data being collected asynchronously i.e. not "by the RTC", for 40 Hz sampling, there might be 39 events in one second, and 41 events in the next, so bear this in mind when looking at the decimal seconds value; it's just to provide time-ordering of sequential events.



Magnetometry hard/soft iron corrections

In order to perform the magnetometer correction, the user must first enable the 3D side of the software by clicking **Show 3D** on the right side of the control panel. Next find the button labelled **Initiate magnetometer correction algorithm**



Figure 3.3 Initiating magnetometer correction algorithm on the 3D window

Upon clicking this, the following will be shown.

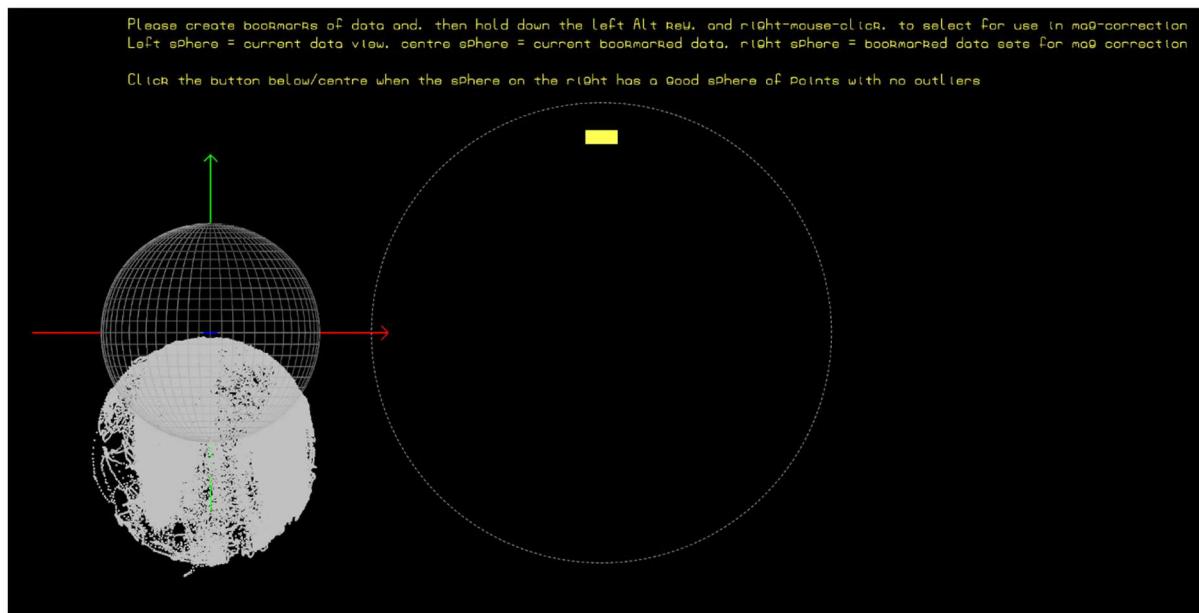


Figure 3.4 First stage of magnetometer offset/ellipsoid correction; selecting data by bookmarks

The data shown in the grey sphere above represents all the data currently on view in the **2D graphing window**. Zooming in on the **2D graphing window** will change the data contained in the above sphere to the left. The aim here is to find magnetometry data that roughly forms either a sphere or an ellipsoid, with no outliers, as is not the case in the data shown above. By zooming in on some data in the **2D graphing window** that shows a sphere with no outliers, the user can then bookmark this data on the **2D graphing window**, as in Figure 3.5 below.

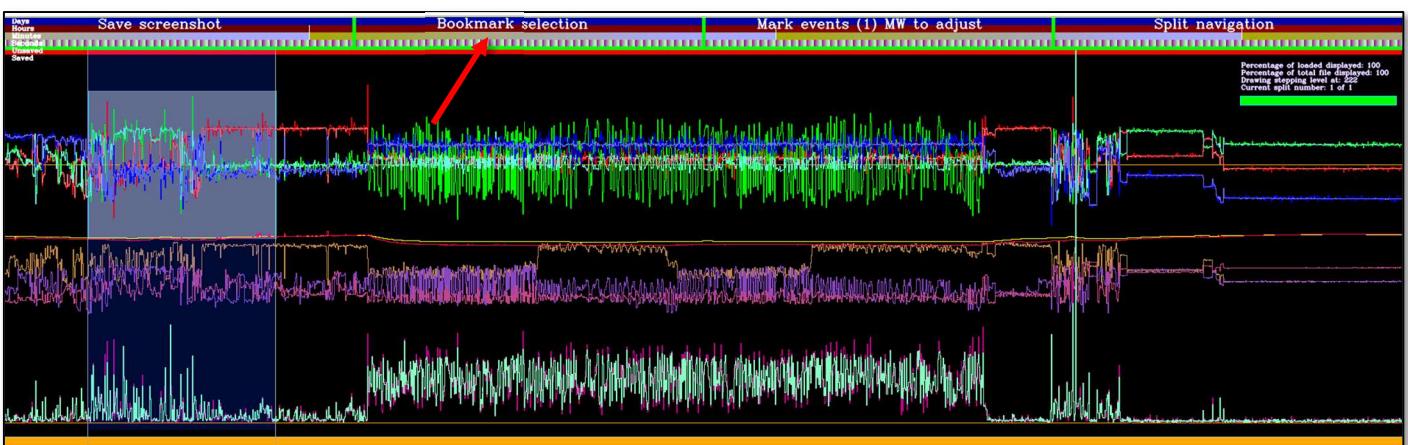


Figure 3.5 Bookmark data by left and right clicking on "Bookmark selection" at the top of the **2D graphing window**, with a section of bookmarked data shown in grey



Note that on the 3D window, there is a large, dashed circle in the centre region of the screen that can be used to drag/roll the sphere of data around (rotation) with the left mouse button.

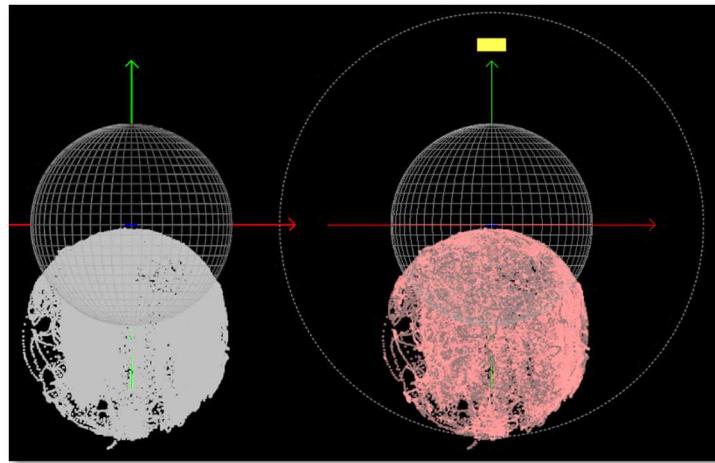


Figure 3.6 Second stage of magnetometer offset/ellipsoid correction; in the left image, a new sphere is now shown, the result of the bookmark that has been created on the left side of the **2D graphing window** in the image on the right

The next step is to create one or more bookmarks that encompass data that define near-spherical data. The bookmarks must exist in the current split. Any data contained in bookmarks (in the current split) will add to the second (centre) sphere. When enough data is present in the centre sphere, and no outliers are visible (checked by left clicking in the dashed circle and rotating the spheres of data), then begin “selecting” the bookmarks that will be used for the magnetometer correction, by holding down the **left Alt** key and **right-clicking** on the bookmarks required to create the sphere of data. As each bookmark is “selected” for inclusion in the magnetometer correction, they will slowly add to the third sphere on the right side; see Figure 3.7.

When in this data selection mode part of the magnetometer correction, the yellow label that appears when the mouse moves over a bookmark on the **2D graphing window** will have a new line stating whether or not that bookmark is “selected” or “not selected” into the 3rd (right side) sphere as part of the magnetometer correction routine:

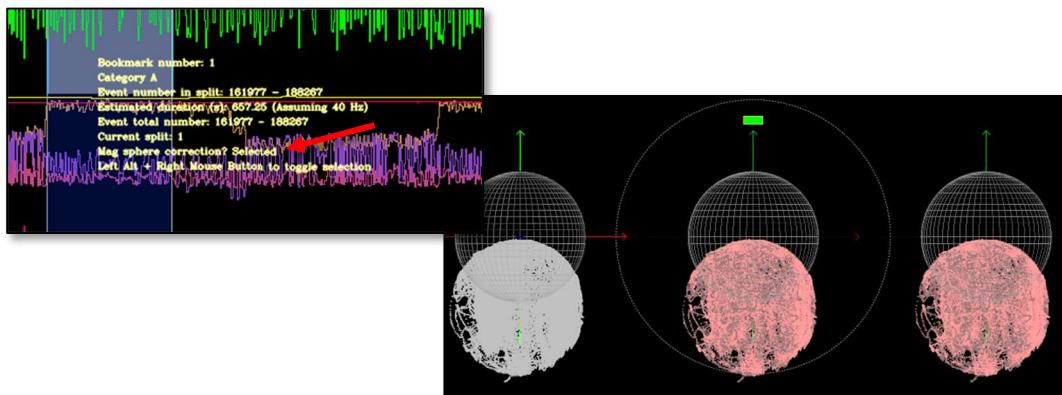


Figure 3.7 Third stage of magnetometer offset/ellipsoid correction; in the left image, the bookmark has been “selected” for inclusion in the magnetometer correction algorithm, resulting in the data being used to form the 3rd sphere on the right-side image

When the 3rd sphere on the right side has enough data forming a sphere or ellipsoid, click the rectangular button centre-top of the 3D window to progress to the next stage.

The 3D environment will now show two spheres on either side of a central square with white and purple cross hairs. The left side sphere represents the data that was in the 3rd sphere on the right side in the previous stage of the correction algorithm. X/Y/Z offsets in the data have now been subtracted from all data, centring the data on the origin. By right-clicking within the centre square, DDMT will run through all possible ellipsoid corrections. As it does this, the square will slowly fill with colour as shown in



Figure 3.8 below. The more data selected for the correction process, the longer this will take. The darkest blue areas mark the best fit (most spherical) correction, while the darkest red areas highlight the worst selections. The purple cross hairs will automatically move to the best fit area found by the software once it has finished looking at the data presented to it. Alternatively, the user can move the white cross hairs to a point of their choosing using **Left Shift** and dragging with a left mouse click. Click the orange button above the square when finished to apply the correction. A text box will appear on screen at the end of this process to verify the magnetic correction is complete.

Note that, with good data selection, the ellipsoidal correction now incorporates the ability to also solve for scalene effects – where an ellipsoid is squashed along the long axis. In this video (<https://www.youtube.com/watch?v=Hu-dWmhqvLk>), pseudo random data, at random orientations, is generated, and immediately corrected back to a sphere (the right side shows the XY, and the XZ perspectives of the corrected ellipsoid).

Note that orange help text appears in the 3D environment throughout the magnetometer correction process to advise/guide the user through the process.

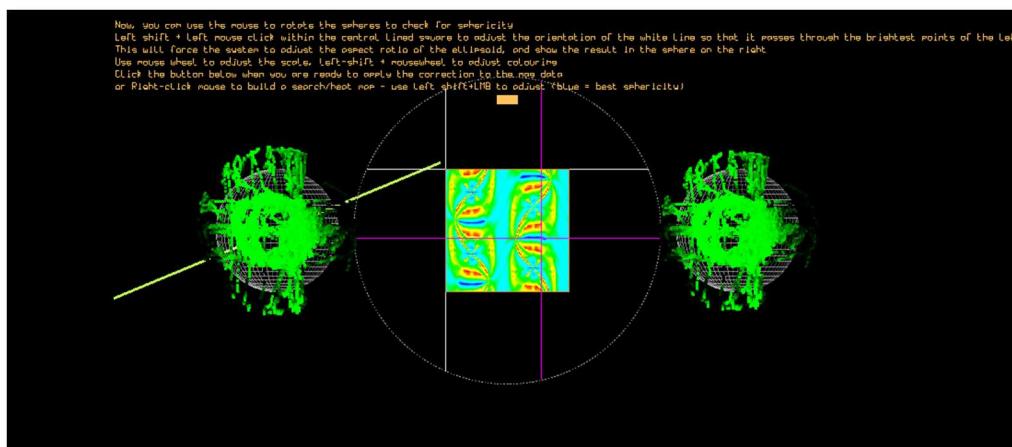


Figure 3.8 The software generates a heat map to aid the user in finding the most spherical magnetometer correction

Acceleration offsets and scaling controls

The acceleration offsets are corrected manually by the user by navigating to the **2D Sliders / Slider controls grouped by purpose** tab and then the scrolling through to **Acceleration** as in Figure 3.9. Click **Show 3D** to open the 3D visualisation window and either click the button near the top labelled “**Create smoothed accel sphere**” or right-click in the centre of the 3D area and create a new XYZ plot (Acceleration by default). If the data you see does not present a suitable sphere of data, move around the split, or progress to a latter/earlier split, ideally one that has ~1g acceleration covering roughly the whole surface of the sphere. Or, if the possibility exists, collect some data similar to that shown in Figure 3.11 below.

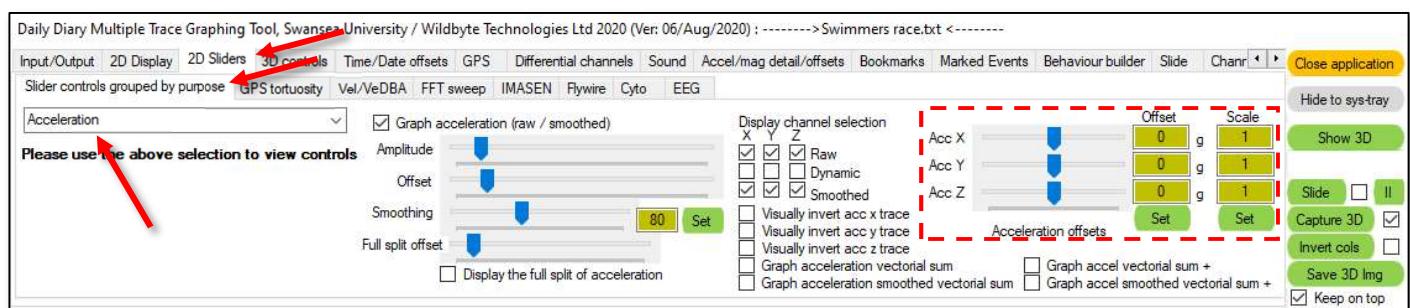


Figure 3.9 The acceleration offsets and scaling control panel



The aim here is to adjust the sliders so that the sphere is centred in all 3 X, Y, and Z axes, unlike in the image below in Figure 3.10. To make the offset more visible you can increase the line intensity of your sphere / data, and manually rotate it on the X axis to 90 degrees – see the 3D visualisation chapter(s) later in this manual.

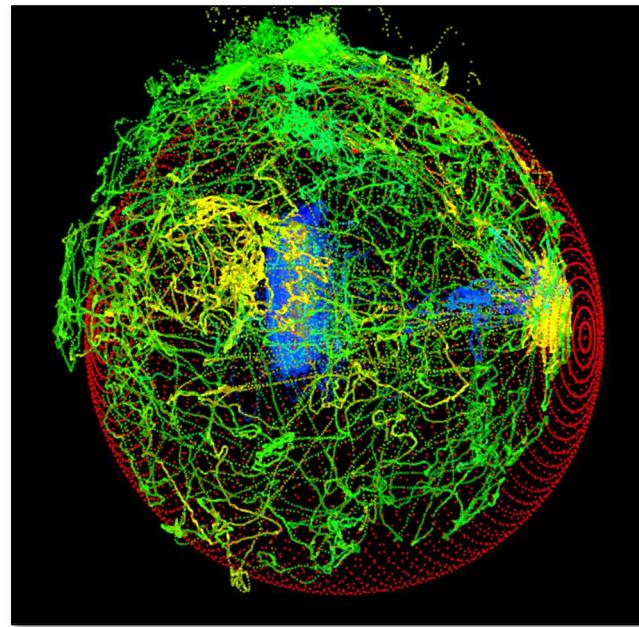


Figure 3.10 The sliders allow the data points to be repositioned so that they sit on the sphere

Once the sphere has been centred in X and Y, pull the sphere down so that the X axis is into/out of the “screen” and the Y axis is vertical on-screen. Then adjust the Z axis to centre the data on the origin.

Once the sphere is centred, It is then time to consider acceleration axes scaling, where each axis is expected to be a minimum/maximum of 1 g at some point on the sphere when the logger is stationary. To do this, it is advised to collect some data using the logger in question (acceleration offsets do not change for a particular logger, they are the result of changes that occur during the PCB manufacture. By holding the logger while it is collecting data with each “face” of the logger pointing vertically away from the floor, and then “very gently” wobble it slightly around in small circles to ensure that a near zero degree offset from the vertical point is definitely captured (for each of the 6x faces), we capture something similar to that in Figure 3.11.

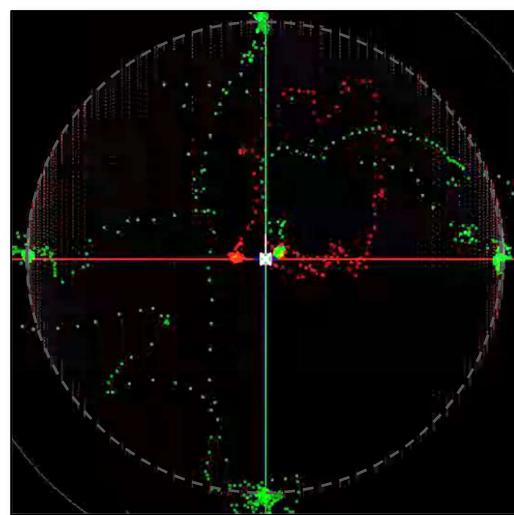


Figure 3.11 6x clusters of acceleration data for use in correcting the acceleration scaling (applicable to correcting offsets also)

Now enter decimal values such as 0.95, or 1.05 etc. to each of the 3 axes scaling boxes shown in Figure 3.9 so that these small clusters of data sit on the 1g sphere above – the sphere has a dashed ring drawn over the top to highlight the location of the digital 1g sphere.



Note that the 3D viewer now includes a small ortho-viewer model at the top right of the screen. Right-clicking on this, will switch/rotate the currently selected visual through different rotations i.e. X-Y, X-Z etc. The ortho-viewer can also be left-clicked and dragged anywhere on the screen, and resized with the mouse-wheel:

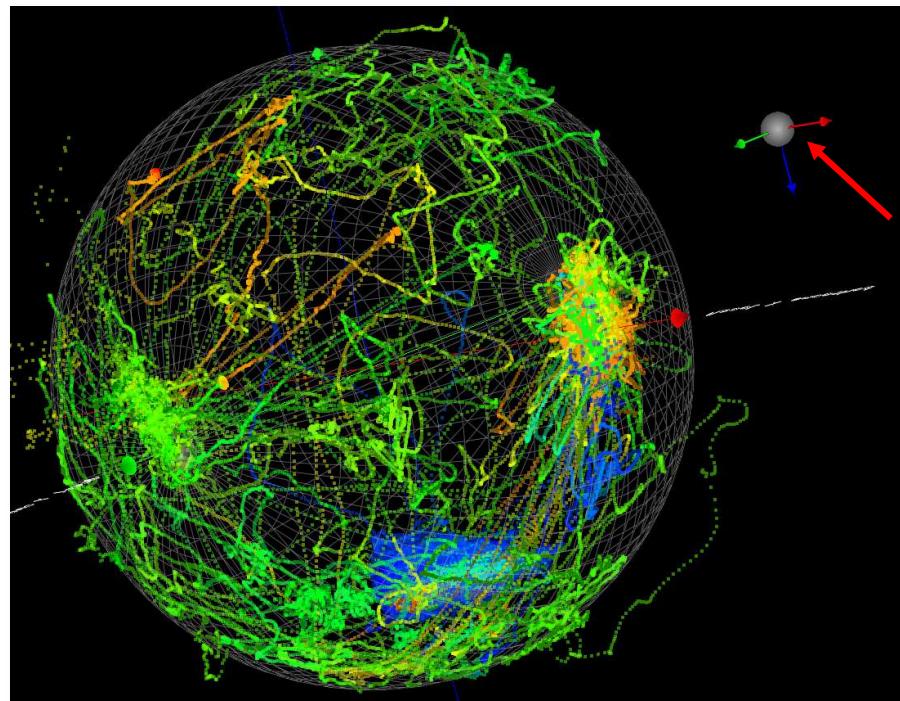


Figure 3.12 The ortho-viewer that helps illustrate the orientation of the current visualisation



TDO (time, date, acceleration/magnetometer offset) files

These time adjustments and acceleration/magnetometry correction settings can then be exported as a **.tdo** configuration file that can be reloaded at a later date. As well as saving the time correction, this can facilitate data sharing and aid users working with previously collected data, and synchronisation with GPS data. See Figure 3.2 of this chapter for correcting the time and date for your data in DDMT.

To export the **.tdo** corrections file, simply click the green button labelled **Export Time/Date/Offsets (.tdo) file** on the **Input/Output / File output/Save** tab. Note that if no magnetometry ellipsoid correction has been performed an error message will alert the user that this must be done before saving a **.tdo** file.

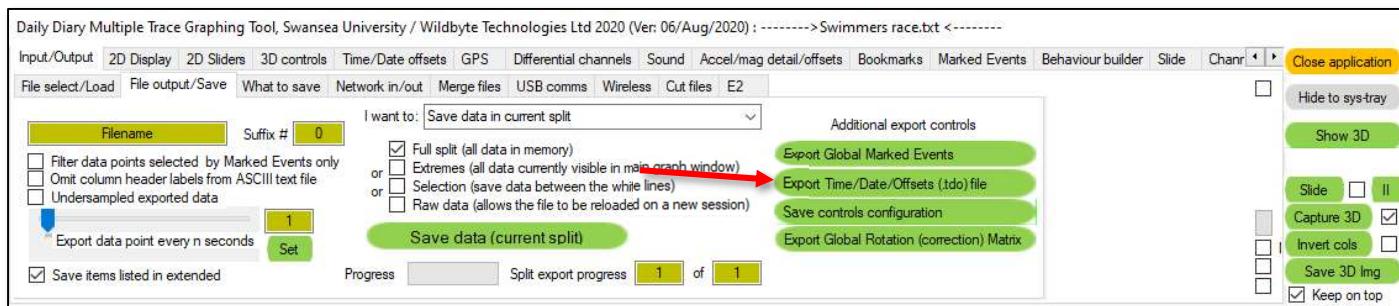


Figure 3.13 The File output/Save sub-tab with the Export Time/Date/Offsets (.tdo) file button indicated.

The next time the file is used, the previously exported **.tdo** file can be reloaded at the start of the session using the **Load TDO file** function on the **Input/Output / File select/Load** sub-tab.

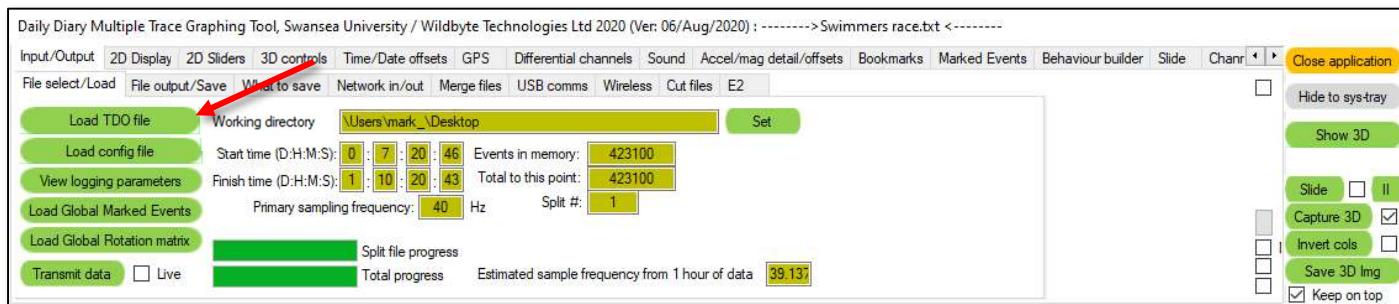


Figure 3.14 File select/Load tab on the control panel

Overall, the **.tdo** file contains:

Acceleration offset correction

Acceleration scale correction

Magnetometry offset correction

Magnetometry spherical/ellipsoid correction

Time / date correction

Device orientation settings (to enable correct calculation of heading/pitch/roll – discussed later)



Further optional data adjustments

Rotation correction

When attaching tags to animals or objects, there is always a chance that the attachment may shift its central position to some different angle. By looking at the data in either 3D or possibly in the **2D graphing window**, this can be determined when the tag is in a resting state. Data can be sectioned and manually rotated on any or all the 3 rotational axes.

First tick **Enable rotation correction controls**, highlight a section of data that will be the **Reference set**, and click the left button labelled **Set**. Next highlight a second section of data, **Adjustment range**, that is considered to require rotation relative to the first highlighted section and click the second button labelled **Set**. Two data sets should appear on the 3D window to separate from the centre. The **Data** scale slider can be adjusted using the mouse wheel to scale the data up.

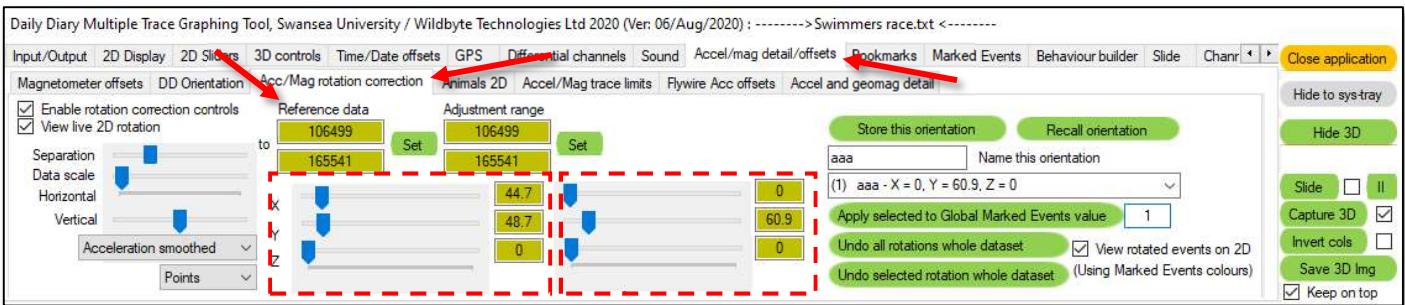


Figure 3.15 Rot' correction control panel

Tick the “View live 2D rotation” and “View rotated events on 2D” tickboxes

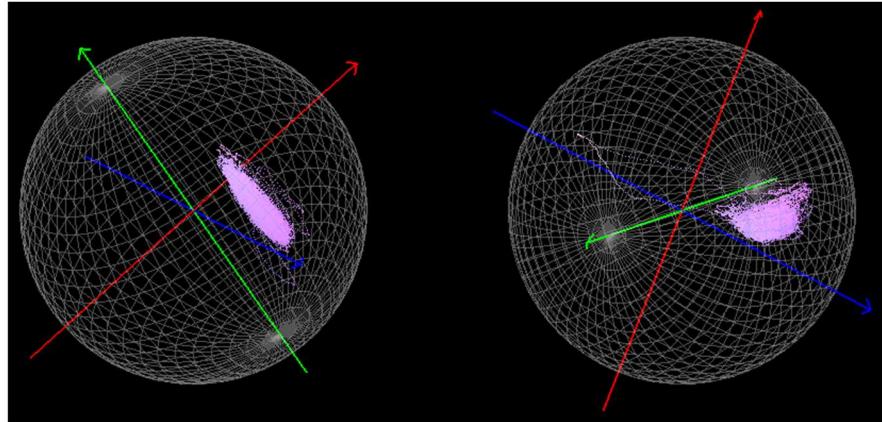


Figure 3.16 The reference sphere on the left and the section of data being adjusted relative to it on the right

The idea here is that a rotation is defined using the second set of XYZ rotation sliders shown in Figure 3.15. Once the user is confident the rotations are correct, these are stored by clicking **Store this rotation** (a name can be assigned first). Currently, no rotation has actually been applied to any data, only the rotation parameters have been defined. To apply it to data, the user must first use **Marked Events** to select the data that requires rotation. So, either sections of data can be selected (throughout the entire file) by either manually selecting data, or using the Behaviour Builder to select data, and then click **Apply selected to Global Marked Events value** and the selected rotation parameters (up to 9 rotations can be stored) will be applied to all data with the **Marked Events** value in the box adjacent to this button. Once the rotation parameters have been applied, the 2D graphing window should show the change, as should all 3D visualisations. Marked Events can now be deleted / reset to something else; this will not affect the rotation parameters that have been put in place. Note that an “event” can only store one set of rotation values (for XYZ) i.e. it is not possible to apply multiple rotations to an event or sections of events. Any attempt to apply a second rotation to an event will result in a message popup pointing out that this is not possible. Rotations can be undone using either of the two buttons **Undo all rotations whole dataset** or **Undo selected rotation whole dataset** i.e. one can completely reset all rotations, or, undo the currently selected rotation across the whole dataset.



To know which data has already been rotated, and with which Marked Event value, DDMT colours the rotated sections with the same colour scheme as is allocated to the **Marked Events** system.

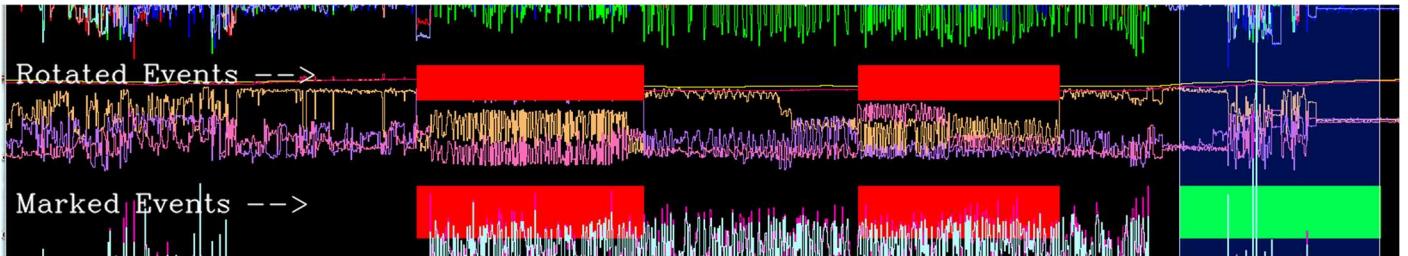


Figure 3.17 Rotation data has been applied to the two sections marked red, while some other data has Marked Events value 2 (green) that has not yet been rotated. It is not possible to rotate any events (data points) that have already been marked in this way without clearing them first.

Note that the **Rotated Events** colour is the colour of the number of the stored parameter set. In the above, the two sets of data are red because they were rotated with the values from the first stored set. With two more sets stored – now 3 in total, and the third selected, and a value of 2 in the box adjacent to the button labelled **Apply selected to Global Marked Events value**, to rotate the green band of data, it now appears blue, as blue is the third default colour for **Marked Events**. See Figure 3.15.

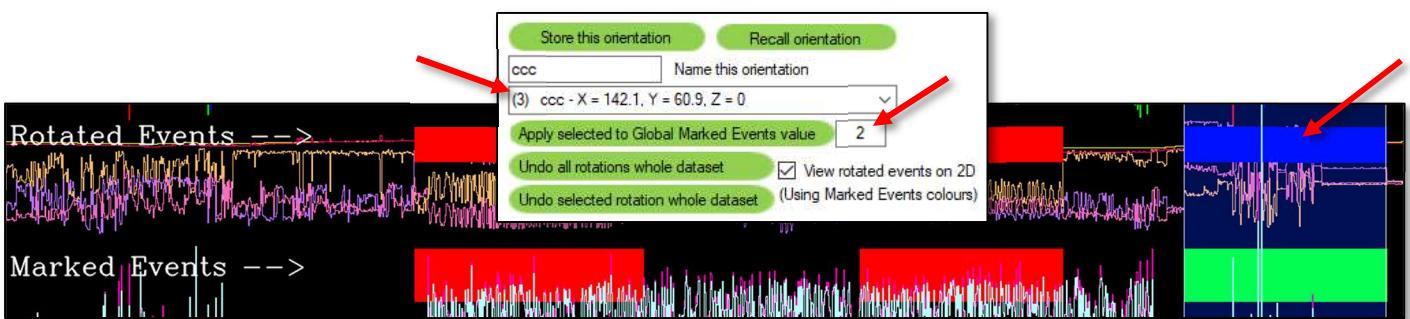


Figure 3.18 The **Marked Events** data in green is rotated as it is valued 2, the same as in the text box above, and the rotation colour is blue, because it is the 3rd stored equation – blue is the **Marked Events** third colour

This sets up a rotation matrix for the entire data file in memory. This rotation matrix can be exported and reimported later i.e. overlaid onto the primary data file (DDMT never makes any changes to the original data file, so it is necessary to load the *.tdo* file for time/date corrections, and acceleration/magnetometry corrections, and now there is the ability to load the rotation correction matrix).

To export the “global rotation correction matrix”, click on the button **Export Global Rotation (correction) Matrix** button on the **Input/Output / File Output/Save** tab.

To import this correction matrix, click on the button **Load Global Rotation Matrix** on the **Input/Output / File select/Load** tab.




4. Differential Channels

Differential channels offer a method of extracting rate of change of various channels. Differential channels stepping ranges, per channel, are user-defined. This is particularly useful when using Behaviour Builder (Chapter 12). Differential channel settings are found on the **Differential channels** tab.

Differential channel creation

The full list of variables for which differential channels can be created is listed under **Available channels**. Select the channel(s) of interest and specify the range using **Differential calculated over range** value.

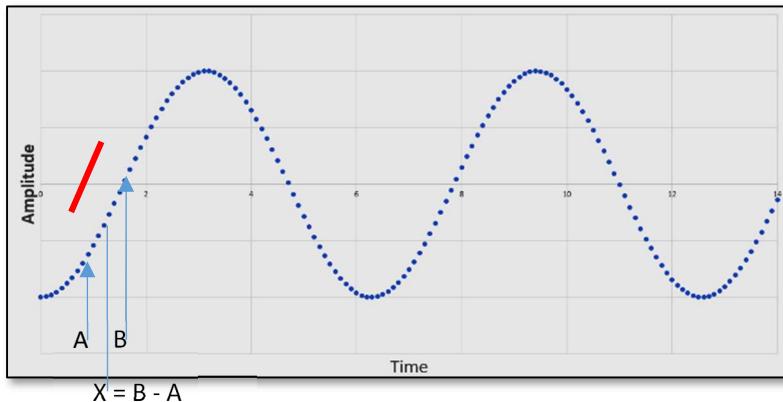

The differential algorithm runs through each point calculating $B-A$, centred around the point with a separation of *Range n*, default 5. If $n = 5$ then B will be 2 events after the current point, and A will be 2 events before. The same channel can be added to the differentials list many times, each with a different *Range* value. In this way, we obtain the gradient (red line) for the point X .

Figure 4.1 Differential channel calculation

Once you have selected the channel(s) of interest and set the range, click the adjacent green **Add** button. The selected channels will now be displayed in the next box **Channels selected for processing**. Click the green **Process channels** button under this box. The channels will now be displayed in the third and final box on the **Differential channels** window and on the 2D graphing window.

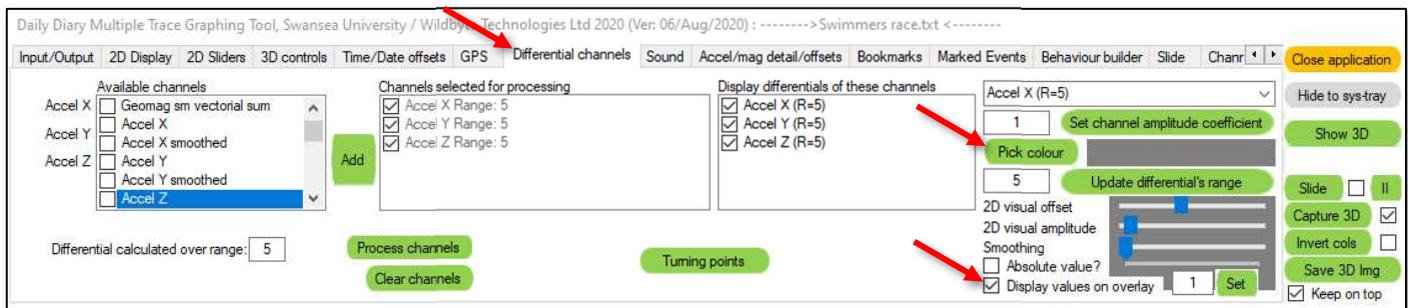


Figure 4.2 Differential channels with Accel X,Y,Z selected and processed with Range values of 5 events (+/- 2 from each event)

Note that **Set channel amplitude coefficient** allows some small or large differential waveforms to be displayed on the **2D graphing windows** with appropriate scaling; if the waveform is too large, perhaps enter a value of 0.1, or even 0.01 here and click the **Set amplitude** button. This is a visual / cosmetic change and does not affect the data itself in any way. The **2D visual offset** and **amplitude** sliders are also just visual / cosmetic adjustments to allow the user to individually adjust the vertical and scale of each of the differentials on the display so that multiple waveforms may be separated vertically.

Smoothing does, however, affect the channels of differential data. When the value is 1 then no smoothing is applied. Each differential channel has its own individual level of smoothing.

A differential channel's **Range** can be altered by manually entering a value in the box on the right side and clicking **Update differential's range** button.

Note that the default colour for a differential channel is grey and if processing multiple differential channels at once all will appear in the same colour on the same location of the 2D graphing window. Select individual channels using the drop-down menu on the right-hand side of the **Differential channels** sub-tab. Use **Pick colour** to update the channel colour and adjust its offset, amplitude and smoothing using the sliders underneath, shown in the same colour as that selected for the differential channel currently



selected. When selected, the **Absolute value?** checkbox will display the magnitude of the differential channel values without regard to its sign (no negative values given).

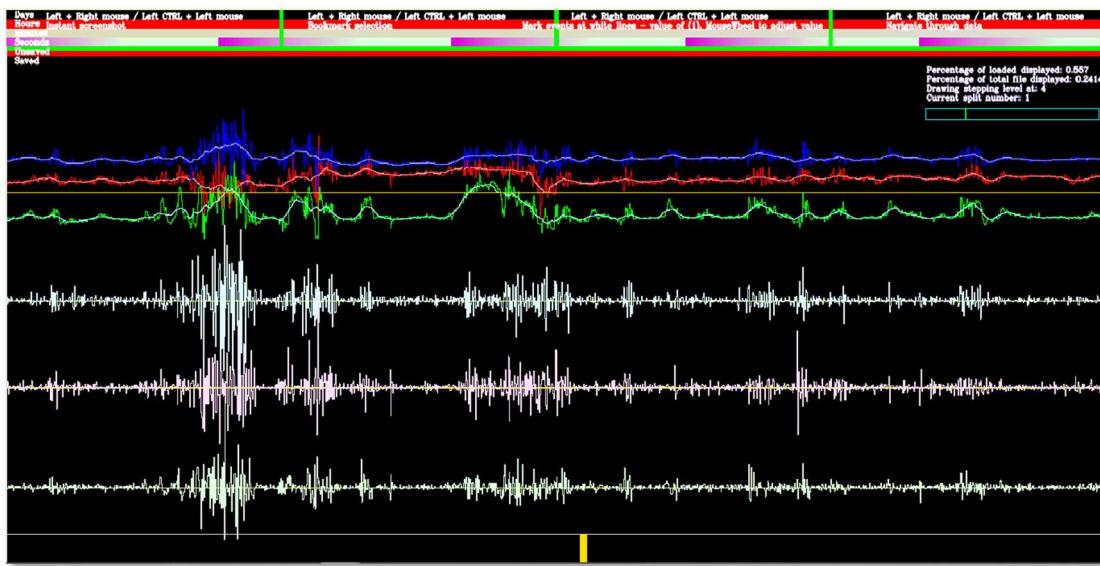


Figure 4.3 Acceleration channels with corresponding differential channels after colour and offset adjustments

Displaying the differential events values on the values overlay

The **Display values on overlay** checkbox on the **Differential channels** sub-tab is selected by default, so unless unchecked the newly created differential channels will be added to the bottom of the display overlay. The overlay will update when colours are customised and include the range (stepping) over which the channels were calculated. Note that the same channel can be processed multiple times at different stepping ranges to explore which is most useful for detecting patterns of interest.

Note that the overlay now includes the *mean* and *absolute(mean)* of the values between the left/right white lines on the 2D graphing window.

Parameter	Left click	Right click	Delta
Acceleration X (raw):	0.373291	-0.161621	0.534912
Acceleration Y (raw):	-0.671631	0.359375	1.03101
Acceleration Z (raw):	0.634521	0.822021	0.1875
Acc Variance, SD (x):	0.0433757	0.208268	
Acc Variance, SD (y):	0.765566	0.874966	
Acc Variance, SD (z):	0.0422289	0.205497	
Event number (in split):	222566	222598	32
Event number (in total file):	222566	222598	32
Event date + times / D:HH:MM:SS.ddd:	01/02/2017 08:55:29.512	01/02/2017 08:55:30.333	0:00:00:00.821
Extremes event start/stop (Split 1) :	221688	223248	1561
Extremes total events start / stop:	1	423100	1561
Extreme times HH:MM:SS.ddd:	8:55:7	8:55:46	
<hr/>			
Differentials:			
Accel X (R=5)	0.078125	0.171875	0.09375
mean, absolute(mean)	-0.05904	0.1821	
Accel Y (R=5)	0.1445	-0.2734	0.418
mean, absolute(mean)	0.1468	0.3343	
Accel Z (R=5)	-0.05859	0.3672	0.4258
mean, absolute(mean)	0.009233	0.1882	
Accel Z (R=10)	-0.4219	0.3984	0.8203
mean, absolute(mean)	-0.03078	0.3321	

Figure 4.4 Differential channels automatically added to the display overlay



Turning point finder

All channels in DDMT can be passed through the **Differential channels** generator, therefore, it is relatively straightforward to isolate turning points in data by honing in on the differential values around zero.

Built into the **Differential channels** algorithm, one can click the **Turning points** button; Figure 4.5.

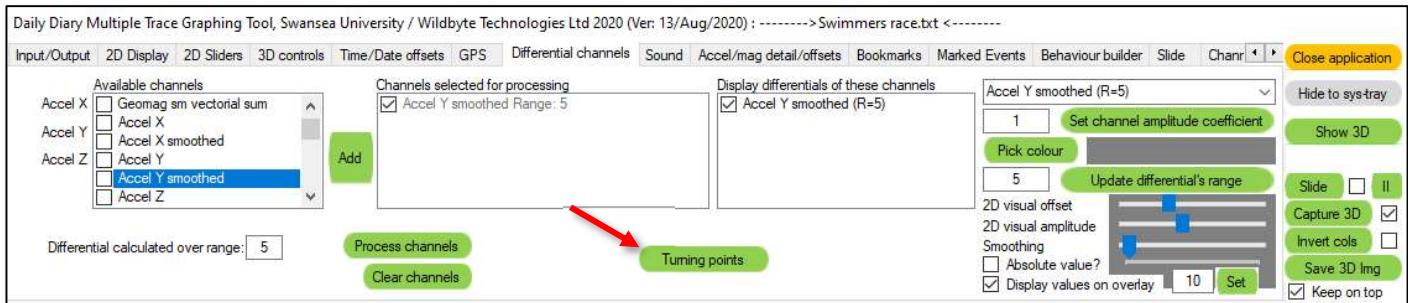


Figure 4.5 Differential channels **Turning points**

This reveals the **Turning points** controls panel; Figure 4.6

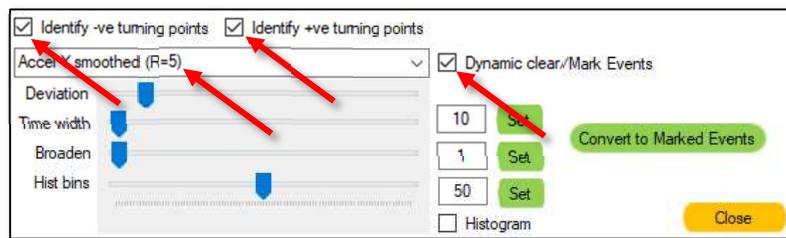


Figure 4.6 Turning points controls

Options are present to select either peaks (+ve turning points) or dips (-ve turning points). Tick **Dynamic clear/Mark Events** so that as a region of data is highlighted on the **2D graphing window** with the left/right white lines, this region is dynamically assessed for peaks/dips.

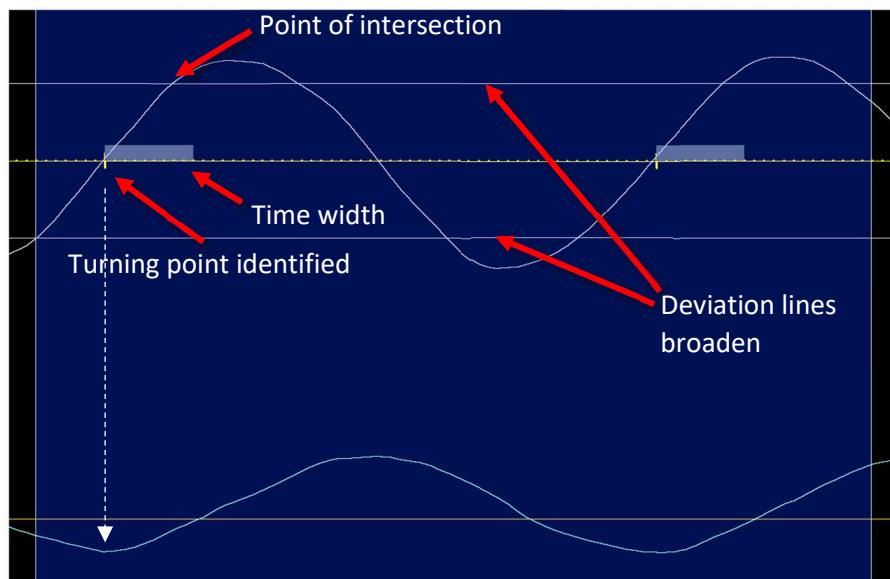


Figure 4.7 Turning points deviation lines



As the **Deviation** slider is adjusted up, the **Deviation** lines widen. The grey patch after the turning point is the **Time width** which can be expanded with the **Time width** slider. As long as the **Differential** line intersects with the **Deviation** line before reaching beyond the **Time width** ends, the turning point will exist (small vertical line; yellow for a dip, purple for a peak).

The **Broaden** slider adjusts the single line indicator at the turning point to multiple points; Figure 4.8.

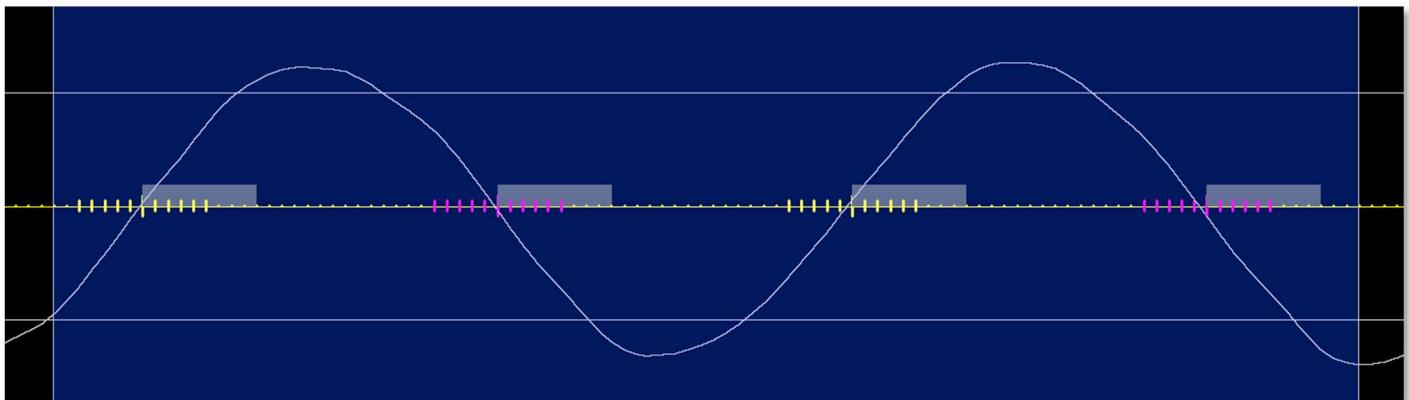


Figure 4.8 Turning points and broadening

Here, in Figure 4.8, we've also selected the positive turning points (purple). The turning points only exist temporarily. To convert them to **Marked Events**, click **Convert to Marked Events**.

Here, Figure 4.9, these markers have been converted to **Marked Events**; value 1 (red) for dips, and value 2 (green) for peaks.

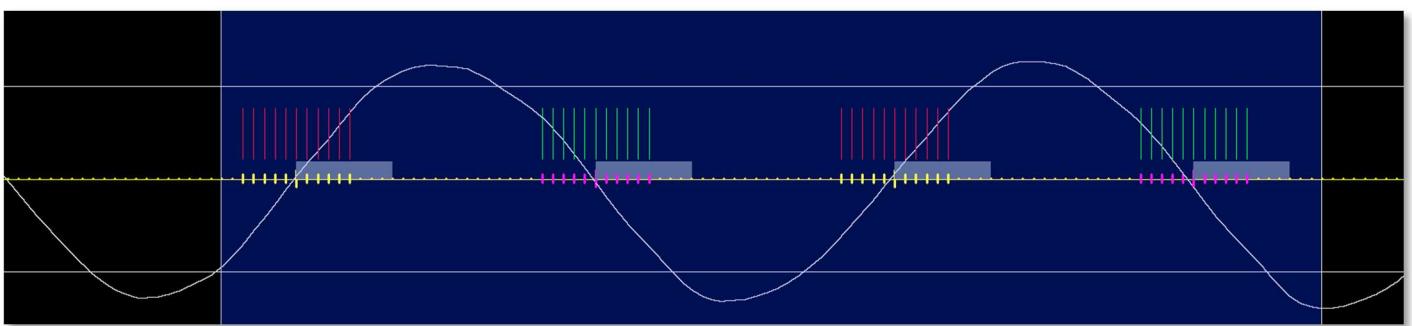


Figure 4.9 Converting Turning points to Marked Events

When a lot of turning points are available, a histogram is available by ticking **Histogram**.

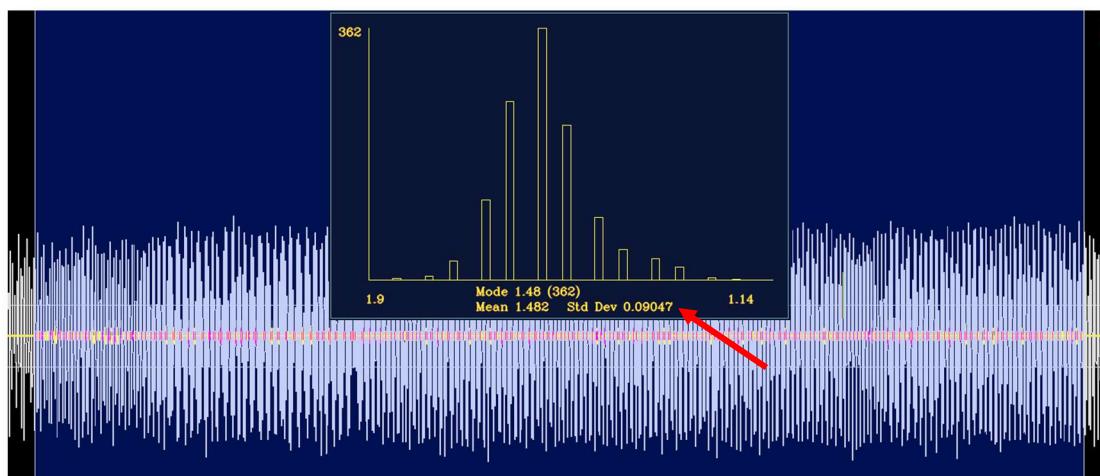


Figure 4.10 Histogram of Turning points showing the peak-peak period Mode, Mean, and Standard Deviation



The histogram is one of timing and illustrates the peak-peak and dip-dip period data for all identified turning points for the selected differential. Figure 4.10 identifies the Mode and Mean around 1.4-1.5. Zooming in on one of the oscillations gives a little over 1.4 s.



Figure 4.11 A quick check on the period of an oscillation to compare against that of the histogram in Figure 4.10

With these turning points converted to **Marked Events**, these data can either be saved directly, or used in conjunction with new expressions in the **Behaviour Builder** to generate **Bookmarks** etc. to export data of interest.



5. Marked Events and Bitwise Layers

A **Marked Event** is where an **Event** has been marked with a value between 1 and 9 (inclusive) either as the result of search function, or manually marked by the user. The purpose of **Marked Events** is to:

- i. highlight data to be saved (Marked Events can be used as a filter in the export process)
- ii. to be used as part of a Boolean search function
- iii. or simply as a visual indicator of the meaning of a part of the data

Marked Events from a given session can be exported as a file and reloaded at a later session.

An **Event** can only be marked with a single **Marked Event** value (1-9). If it is marked again, the old **Marked Event** value will be replaced with the new value.

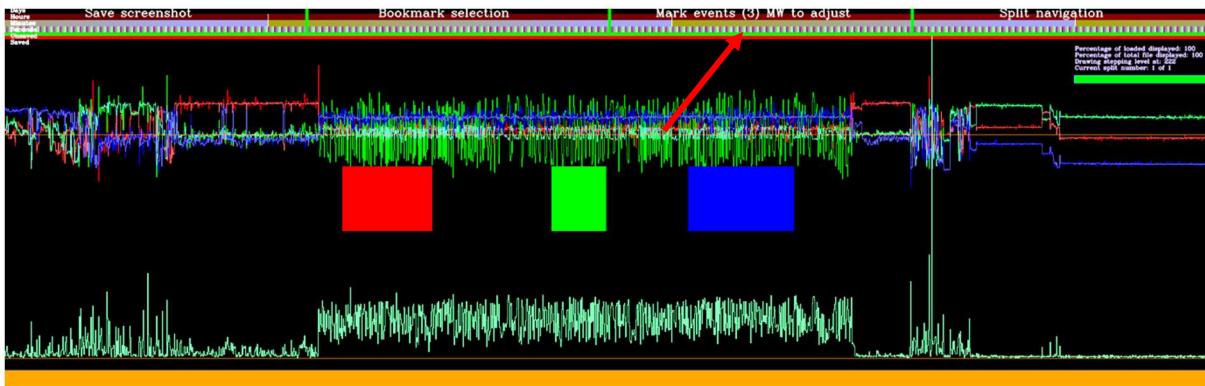


Figure 5.1 Marked Events shown as highlighted regions in red, green, and blue (Marked Events values 1,2, and 3)

Marked Events may be created either manually by the user by placing the left and right white lines on the **2D graphing window** as above and then left and right clicking simultaneously on the area marked in Figure 5.1 labelled “Mark events”. Also, hovering the mouse over this area and rolling the mouse wheel will change the **Marked Events** value between 0 and 9. If the value is set to 0 then the region selected will have any current **Marked Events** cleared, while a value between 1 and 9 will set that area to that value, overwriting any previous values present.

A new method of marking events (**dynamic marking**) has recently been added. By hovering the mouse over the 3rd portion of the top banner on the **2D graphing window**, a new message appears informing the user that holding down the **CTRL** and clicking the middle mouse button toggles on/off dynamic marking events on the window. The first **MMB** click enables this feature and then *left clicked* (+optional dragging the mouse left/right across other data points – slowly) events will be marked with the current **Marked Events** value if < 9, else a value of 1 (red) will be selected. Right clicking in the same way will mark events with the current **Marked Events** value +1 i.e. if 8 is showing, it'll set the **Marked Events** value to 9.

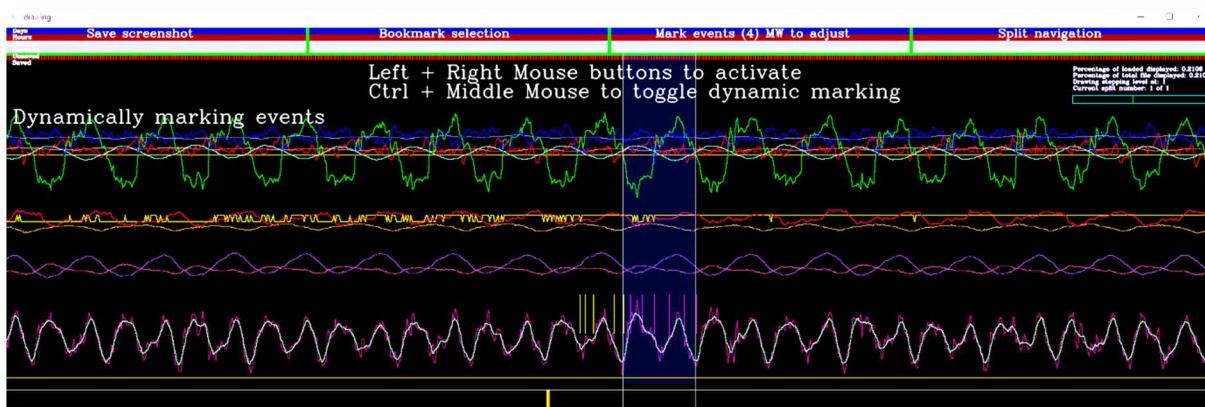


Figure 5.2 Marked Events via ‘Dynamic marking’ (with the mouse left/right buttons)



Marked Events are persistent for the whole session, while moving through various splits, as is the case with **Bookmarks** (discussed in the next chapter).

For more **Marked Events** controls, click on the **Marked events controls** button on the **Behavioural Builder** tab. Here the colours for the **Marked Events** 1-9 can be set, **Marked Events** between the white lines, or throughout the entire data file can be cleared etc. The vertical height of the **Marked Events** bands can also be set by using the two horizontal lines created by holding down **Alt** and **Ctrl** and left or right clicking on the **2D graphing window** defining the vertical spacing of the bands, and clicking **Set visual height of Marked Events** to set it. Finally, **Marked Events** that are close together, and of the same value, can be merged by clicking **Merge Marked Events with proximity**.

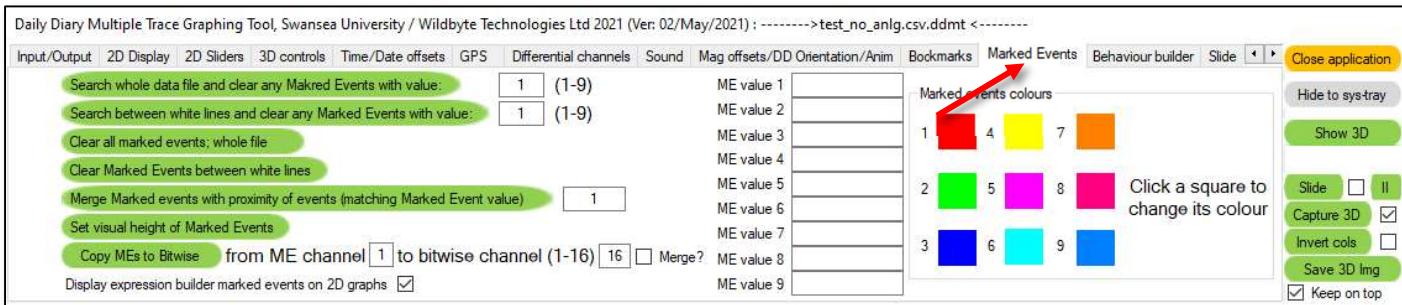


Figure 5.2 Marked Events control panel

There are 9 text boxes for the 9 Marked Events values. In Figure 5.2, these have been marked as “Walking”, “Running” etc. These serve no purpose other than general reference for the user. Some files that have annotations for the data will import labels with the data into these 9 boxes.

The vertical height/position of the coloured banding on the **2D graphing window** for **Marked Events** can be adjusted by holding **left ctrl** and **left alt**, and then left / right clicking to define the vertical position of two reference lines. Clicking the button labelled **Set visual height of Marked Events** then sets the vertical limits for this feature.

Marked Events were initially the output/result from the **Behaviour Builder** function, discussed later in Chapter 10 (**Behaviour Builder**). They can now also be used as part of the **Behaviour Builder** search functions themselves i.e. the results of one or more searches can be used as a variable of subsequent searches.

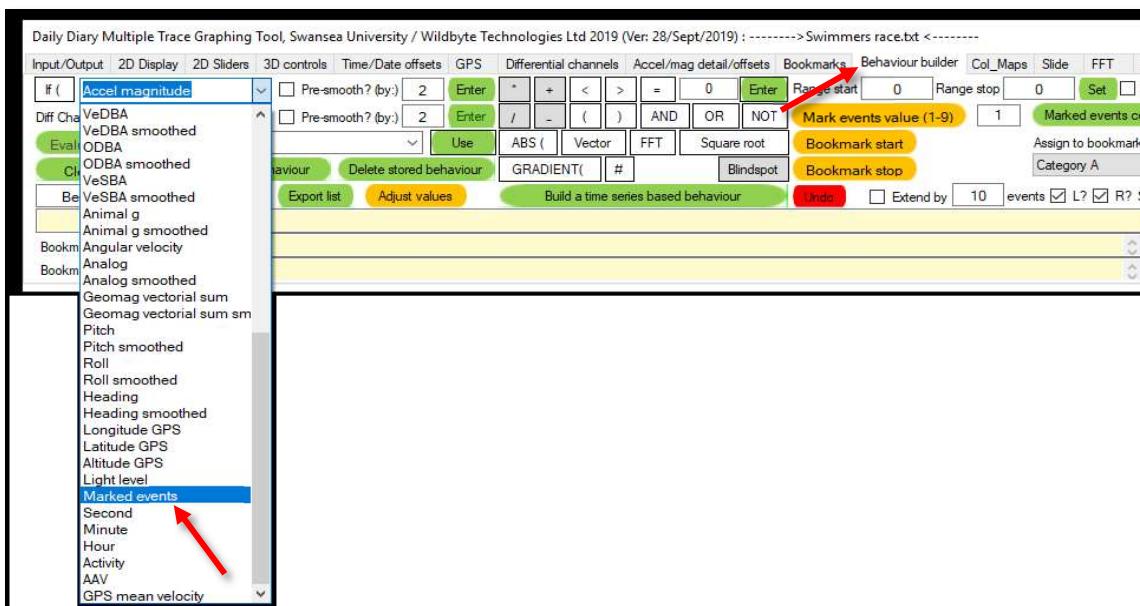


Figure 5.3 Marked Events is one of the variables that can be used for Boolean search expressions i.e. $\text{Marked Events} > 0$, or $\text{Marked Events} = 3$



Marked Events can also be used as a filter for the visualisations. Ordinarily, a visualisation is constructed from the data viewed in the **2D graphing window**. Figure 6.3 shows some data in a **2D graphing window** and the result in a 3D dot plot. The data on the right side has **Marked Events** value 1, red, while the left side has a value of 0.

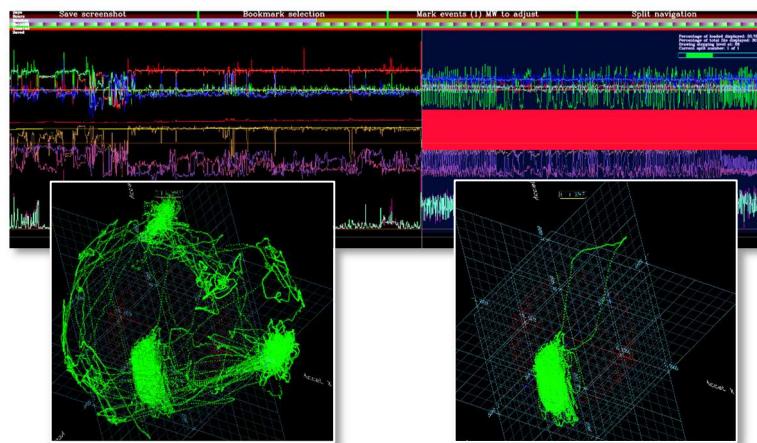


Figure 5.4 **Marked Events** used as a filter to only include specific events in the 3D visualisations

Continuing with the 3D visualisation side, **Marked Events** play an important role in the Dead Reckoning visual, allowing up to 9 different coefficients of movement or stationary behaviours.

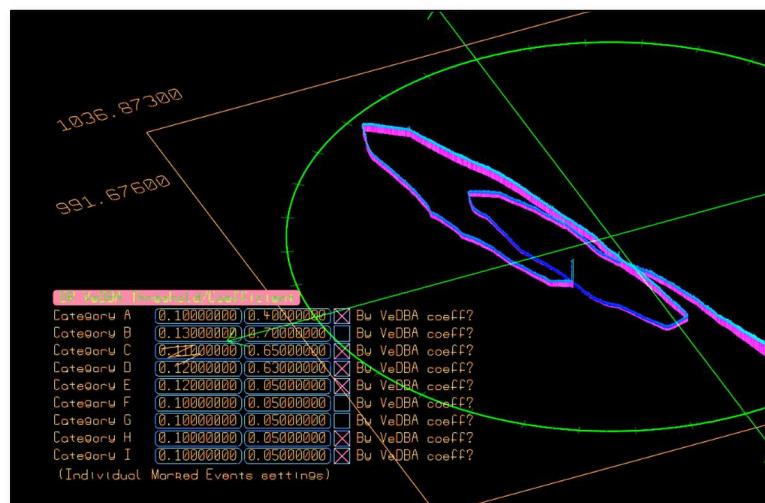


Figure 5.5 **Marked Events** used as different threshold/coefficients for various behaviours in the Dead Reckoning the 3D visualisations. This allows for fine tuning the path

Marked events can also be used to as a filter on the export/save function.

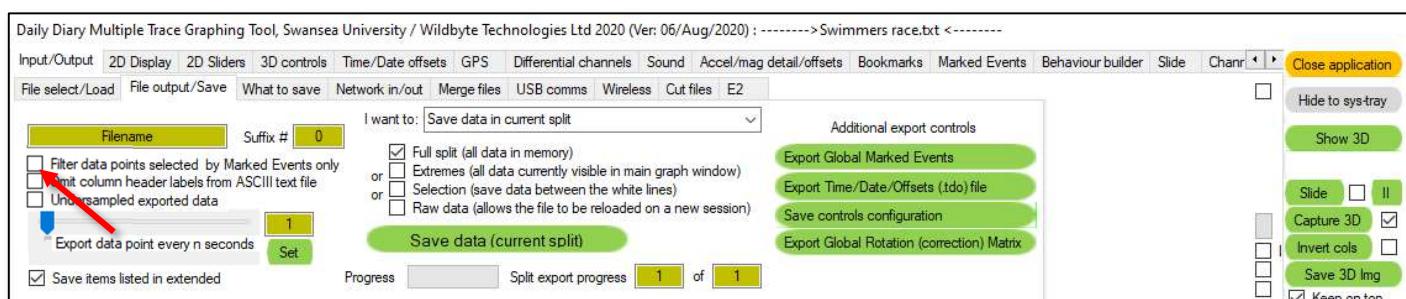


Figure 5.6 **Marked Events** used as a filter to only include specific events in the 3D visualisations

By selecting the tickbox above, only events that have **Marked Events** associated will be exported.



Bitwise Layers is a feature that allows for the storage of multiple **Marked Events** by copying or merging them across to 1 of 16 different bits within the **Bitwise** channel.

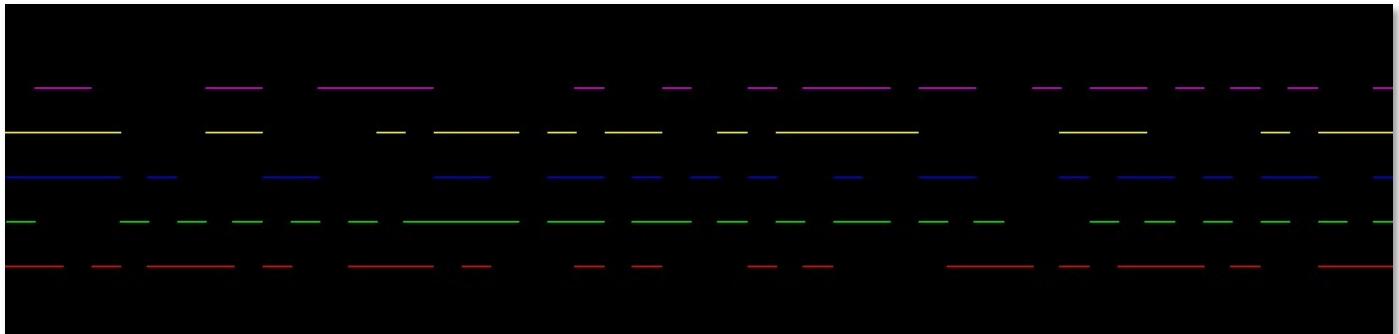


Figure 5.7 Bitwise Layers on the 2D graphing window

If the primary DDMT source file comes from a conversion app, it may contain a **Bitwise** channel with embedded info ('1's and '0's)

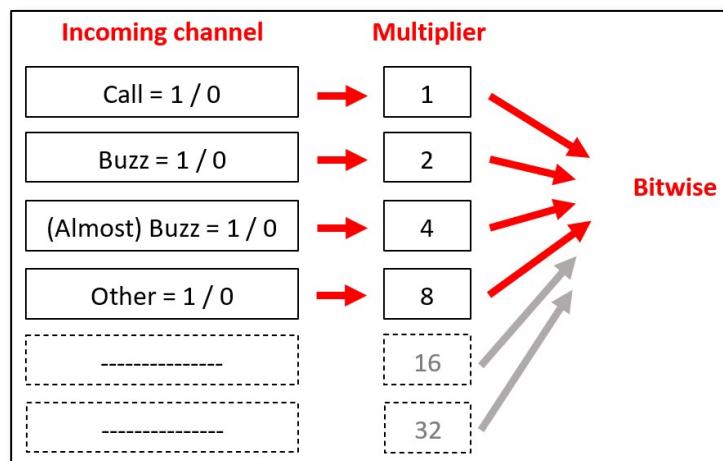


Figure 5.8 Creation of the **Bitwise Layers** using an external app to merge multiple channels into one. Each channel is multiplied by a power of 2 (from 0-15) so that when they combine, no information is lost/merged

Slider controls are available on the **2D sliders** tab that allow the vertical offset on the **2D graphing window** and their relative vertical separation from each other.

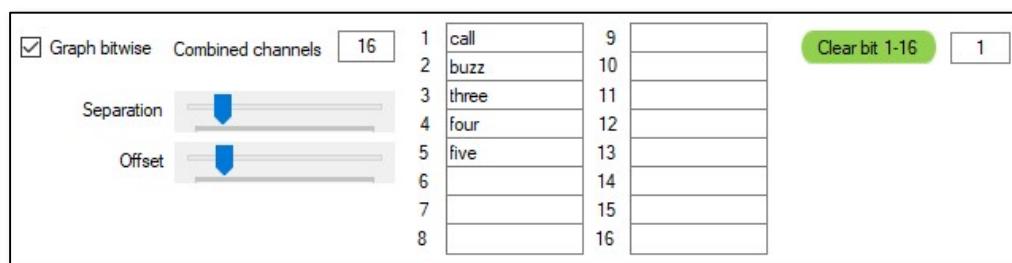


Figure 5.9 Controls for the **Bitwise** channel on the **2D graphing window**

The **Combined channels** box limits how many channels are graphed on the 2D window.

The 16x text boxes are labels for the **Display overlay** on the 2d window. They are also filled (to the number of channels) if the primary DDMT source file has embedded **Bitwise** channel data, from the column headers.

The **Clear bit 1-16** button will clear down and **Layers** in the **Bitwise** channel, although any channels embedded within the source file will not be cleared; they're locked out from any change.



The **Bitwise** channel's **Layers** can be used as part of an expression in the **Behaviour Builder**. For instance,

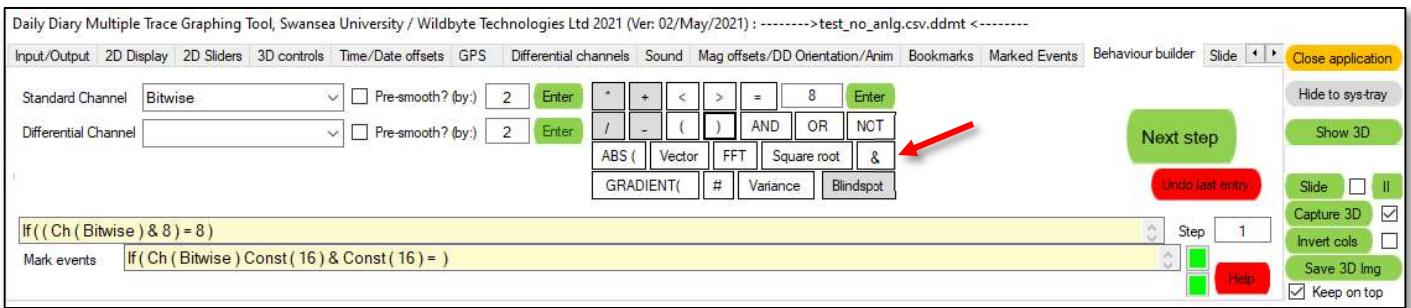


Figure 5.10 The **Bitwise** channel requires the use of the ampersand ‘&’ in order to mask out specific bits within the channel

An example expression using the **Bitwise** channel might be:

If((Ch(Bitwise) & 8) = 8) then Mark Events

would **Mark Events** for any data point that had a ‘1’ in the 4th **Layer** of the **Bitwise** channel (as $2^3 = 8$, while counting in powers from 0), while

If((Ch(Bitwise) & 10) = 10) then Mark Events

would **Mark Events** for any data point that had a ‘1’ in the 2nd and 4th **Layers**. Using this expression would also mark any data points that had other layers. The “& 10” acts as a mask, ignoring the presence or absence of other **Layers** within the **Bitwise** channel.

Note that the ampersand is equated to the value of the power of $2^{(\text{Layer}-1)}$, not the Layer value i.e. to filter out the 4th **Layer**, one would equate to $2^3 = 8$. Layer values 1-16 would therefore be equated to 1, 2, 4, 8, 16, 32, 64, 128 (for layers 1-8), and 256, 512, 1024, 2048, 4096, 8192, 16384, 32768 (for layers 9-16) respectively.

Such an expression, using the **Bitwise** channel, could be combined with other channels such as

If(((Ch(Bitwise) & 24) = 24) AND (Ch(Pressure sm) < 984.23)) then Mark Events

Note that the two parts of the expression are both enclosed in their own bracket structure, and a logical **AND** separating them. This level of bracketing is required for the breakdown/translation (by DDMT) of the expression.

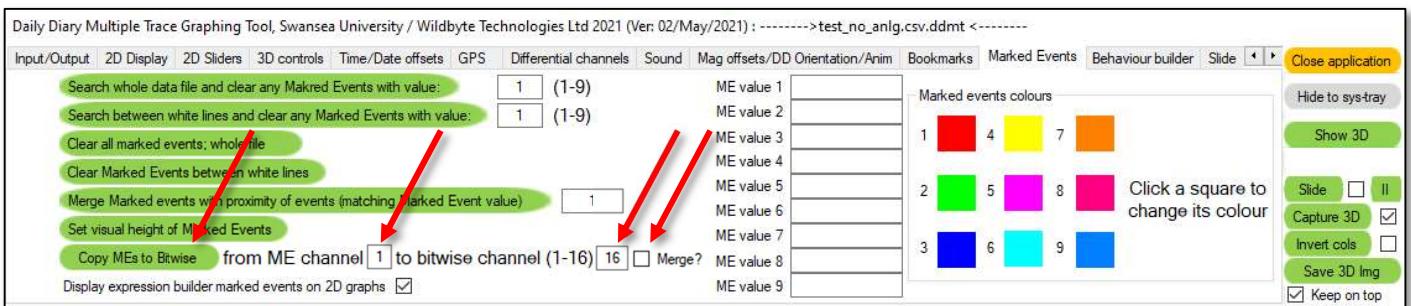


Figure 5.11 Controls for copying/merging **Marked Events** into a **Layer** in the **Bitwise** channel

The presence of **Marked Events** of a particular value can be copied across into any **Layer** within the **Bitwise** channel. In the above (Figure 5.10), clicking **Copy MEs to Bitwise** would copy (overwrite the entire **Layer**) from the **Marked Events** channel values of 1, into **Layer** 16 ($2^{15} = 32768$ / e.g. **If((Ch(Bitwise) & 32768) = 32768) then Mark Events**). “Globally” (across the entire data file) **Marked Events** of value 1 will be entered into **Bitwise Layer** 16. If a Marked Event is ‘0’ or any other value, then the **Bitwise Layer** bit will be set to ‘0’.

If **Merge** was ticked prior to clicking the **Copy MEs to Bitwise** button, then only data points within the **Bitwise Layer** where a **Marked Event** was =1 would be set to ‘1’ in the **Layer**.



The global (across all **Splits**) **Bitwise** channel (containing all 16 **Layers**) can be exported to a file (.gbw) by clicking **Export Global Bitwise channel data** button on the **Input/Output / File output/Save** tab:

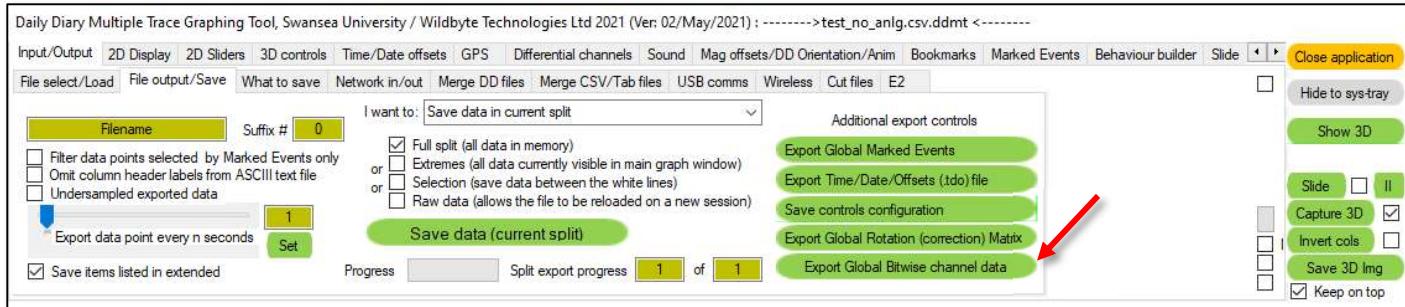


Figure 5.12 Export Global Bitwise channel data button on the File output/Save tab

While a file containing the global **Bitwise** channel data (from a previous session with this primary DDMT source file) can be reimported by clicking the **Load Global Bitwise matrix** button on the **File select/Load** tab.

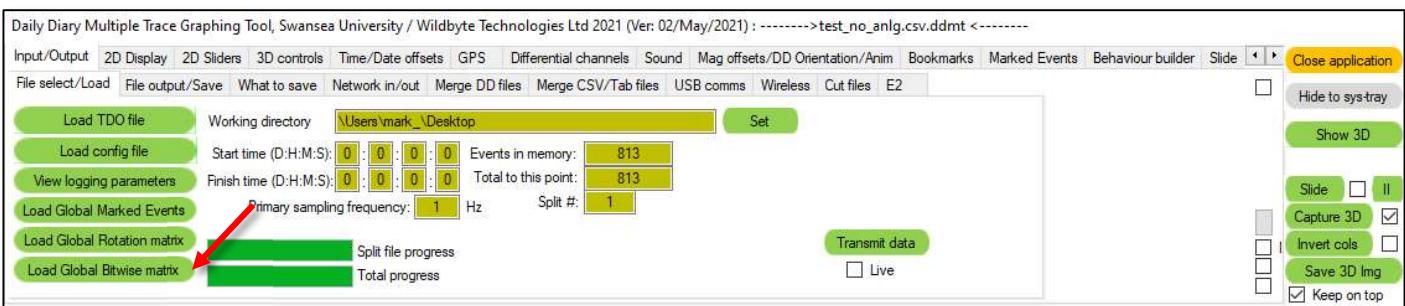


Figure 5.13 Load Global Bitwise matrix on the File select/Load

The presence of '1's within a particular **Layer** can be viewed, along with the **Layer**'s label, by ticking the **Bitwise** checkbox on the **Overlay selection** tab. The **Layer**'s labels will be displayed for both the left and right white mouse lines:



Figure 5.14 Labels are displayed on the Overlay for any layers containing a '1' at the position of the left/right white lines

The **Bitwise** channel information can be exported by selecting the **Bitwise** check box on the **What to save / Additional** tab. Every data point will then show all 16 **Bitwise** channels, with the channel labels as defined on the 2D **Bitwise** controls tab, with the data either being a '1' (present), or a '0' (not present).



1	Event no.	Bitwise 1 Call	Bitwise 2 Buzz	Bitwise 3 Three	Bitwise 4 Five	Bitwise 5	Bitwise 6	Bitwise 7	Bitwise 8
454	165972	1	0	1	1	0	0	0	0
455	165973	1	0	1	1	0	0	0	0
456	165974	1	0	1	1	0	0	0	0
457	165975	1	0	1	1	0	0	0	0
458	165976	1	0	1	1	0	0	0	0
459	165977	1	0	1	1	0	0	0	0
460	165978	1	0	1	1	0	0	0	0
461	165979	1	0	1	1	0	0	0	0
462	165980	1	0	1	1	0	0	0	0
463	165981	1	0	1	1	0	0	0	0
464	165982	1	0	1	1	0	0	0	0
465	165983	1	0	1	1	0	0	0	0
466	165984	0	1	1	0	0	0	0	0
467	165985	0	1	1	0	0	0	0	0
468	165986	0	1	1	0	0	0	0	0
469	165987	0	1	1	0	0	0	0	0
470	165988	0	1	1	0	0	0	0	0
471	165989	0	1	1	0	0	0	0	0
472	165990	0	1	1	0	0	0	0	0
473	165991	0	1	1	0	0	0	0	0

Figure 5.15 Exported Bitwise channel data showing the presence '1' or absence '0' of data in each bit layer



6. Bookmarks

Manually creating bookmarks on the 2D graphing windows

Bookmarks highlight user-defined sections of data which can be saved and reloaded later. Bookmarks can be created manually by placing the two white lines either side of the data, and then simultaneously left and right mouse-clicking in the second box along at the top of the main graphing window where it is labelled **Bookmark Selection**. The selected data are then highlighted with two blue lines denoting the start and end points and the bookmarked area showing a grey background. Collections of bookmarks can be saved to file (.bmk) and reloaded during a later session. As the cursor is moved over a bookmark in the **drawing** window, event limits, category, and bookmark number are displayed beside the cursor in yellow (Figure 7.1). This can be switched off by deselecting the tick box labelled **Show bookmark tags** on the **Bookmarks** tab. Multiple bookmarks may be created throughout the entire data file, within any / all of the splits. Moving from one split to another will refresh the data view and reveal any bookmarks previously created within that split.

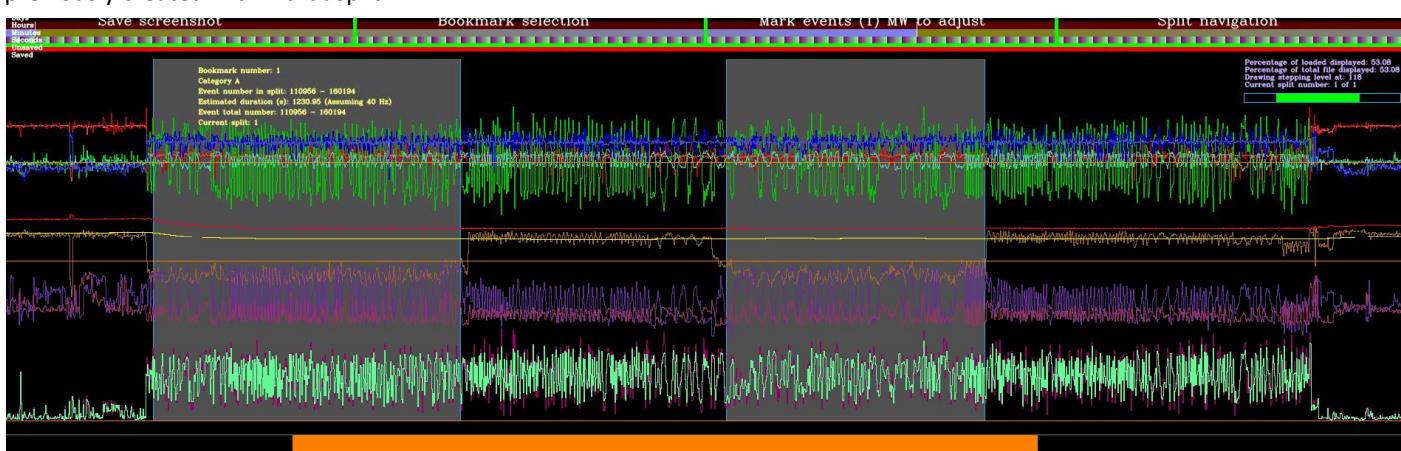


Figure 6.1 Selecting data and bookmarking; bookmark information is displayed on each bookmark as the cursor moves over it showing start/stop event numbers, category name, bookmark time width (s), and split#.

Categories

After selecting the **Bookmarks** tab, a drop-down box will show the bookmarks that have been created this session numbered/ordered chronologically (not necessarily the order in which the bookmarks are created). Note that **Bookmarks** "inherit" the currently selected category. Once a bookmark is selected, its category can be changed to any one of 14 different labelled categories, allowing sections of data to be reclassified. To change the currently selected category, use the second drop-down menu in the centre of the **Bookmarks** tab and then click the adjacent green **Set** button.

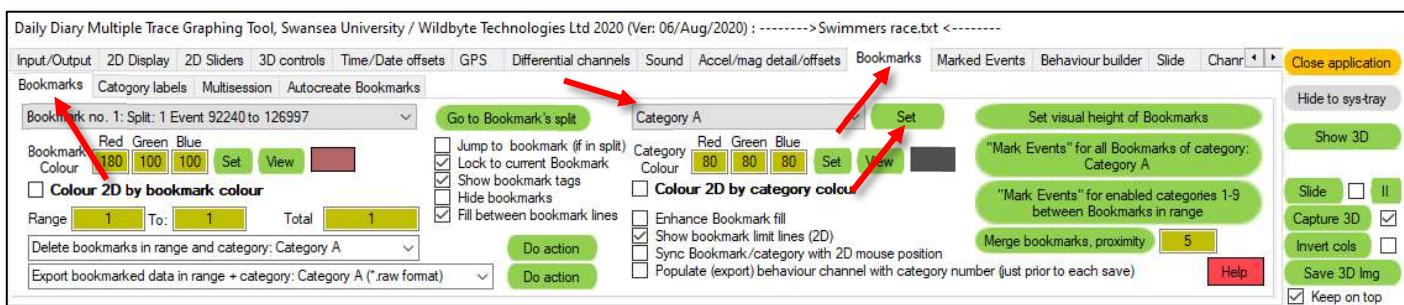


Figure 6.2 The drop-down box for selecting categories is indicated by the red arrow.



Renaming categories

Category names can be changed to anything from the default **Category X**. Category names will also be stored within the exported **.bmk** file. To rename categories, navigate to the **Bookmarks/Category labels** sub-tab. Click on the text box corresponding to the category you would like to rename and type the new name in its place. Click the **Update categories using list** button to save the new labels.

Note that the checkboxes on this tab allow each category to be added or removed from both the 2D and 3D graphing environments. By default, each category label is enabled in both 2D and 3D.

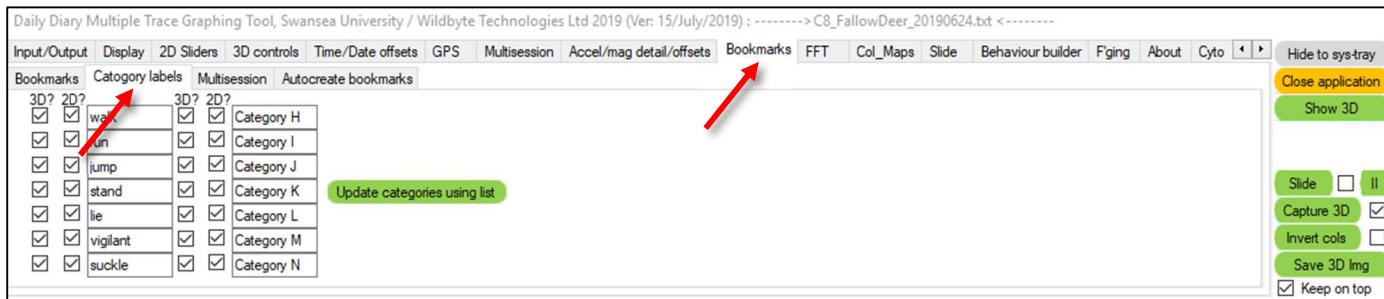


Figure 6.3 Setting new category labels

Adjusting the vertical height of the bookmark bands

By default, the vertical height of the bookmark bands is from the top to the bottom of the **2D graphing windows**, but this can be changed to suit the data on display by holding down **Ctrl** and **Alt** and using the left and right mouse buttons to place the two horizontal lines and clicking **Set visual height of bookmarks**. The placing of these defines the top and bottom limits of the bookmarks.

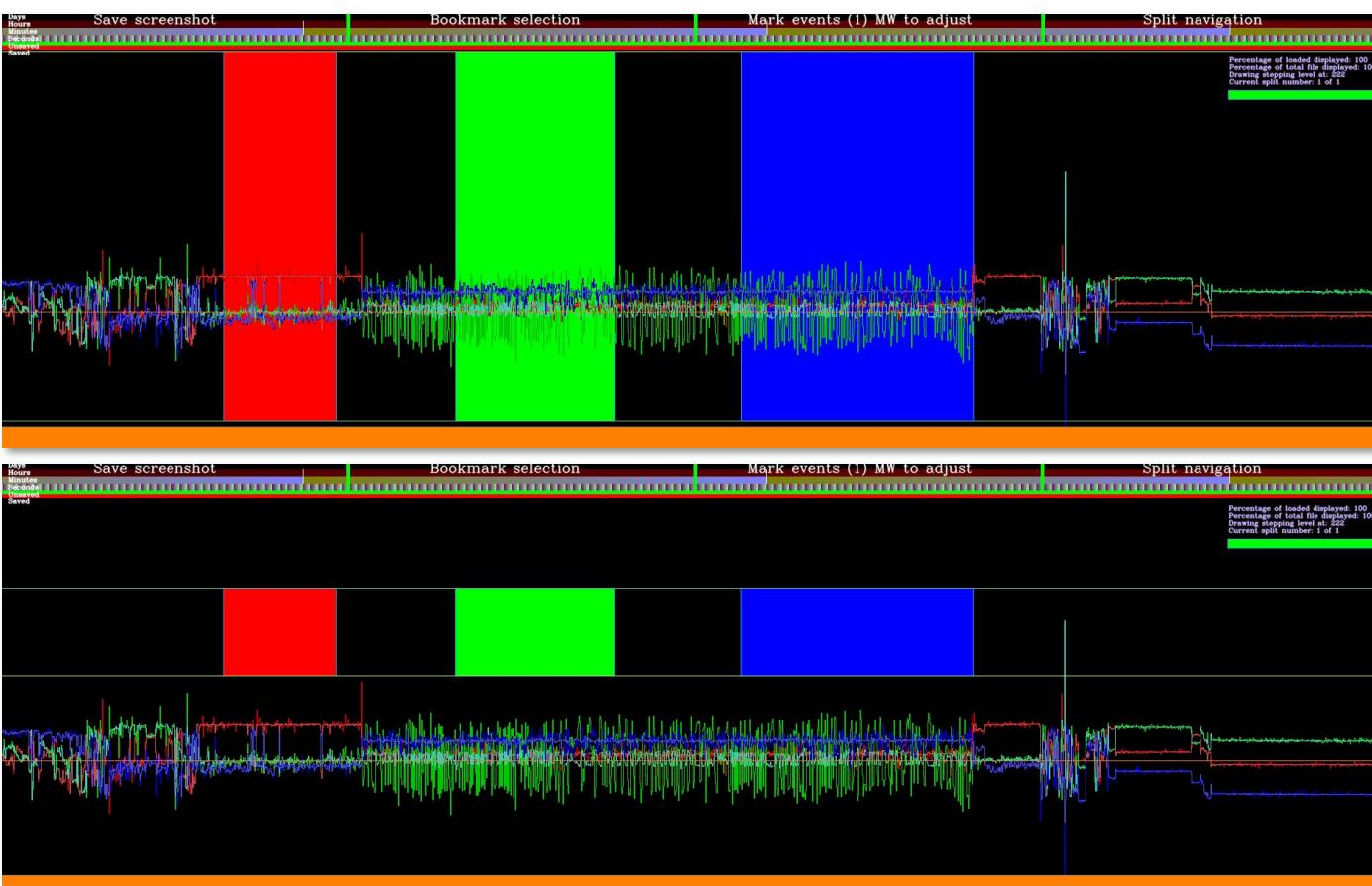


Figure 6.4 Adjusting the top/bottom height of the bookmarks



Save/load .bmk files

Load bookmarks loads previously saved bookmarks into the original data file from which they were exported. **Save bookmarks (.bmk + .txt)** saves the current bookmarks in your DDMT file as both text and .bmk files. The .bmk files enable you to reload previously defined bookmarks into DDMT again during future sessions, while the text files can be used for inspection or analysis in other programmes such as Excel and R.

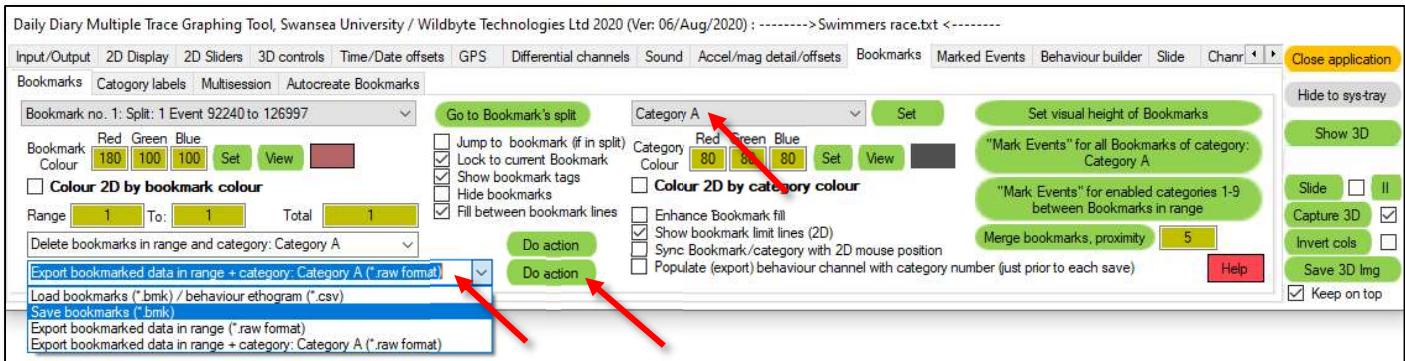


Figure 6.5 Saving and loading .bmk files

Export .raw files

Export bookmarks in range (raw) exports all the data from bookmark X to bookmark Y, where X and Y are the values in the two textboxes above this button, as a single .raw file that can be reloaded into DDMT. Note that the bookmark markers will be maintained, showing the separate sections of data. **Export in range and category (raw)** exports all bookmarked data that matches the category selected in the drop down menu (as indicated in Figure 6.5).

Convert Bookmarks to Marked Events

Use the **Mark events for selected bookmark category** button to quickly mark all events in all bookmarks of the currently selected category. When working with multiple categories, **Mark events for enabled categories 1-9 between bookmarks** can be used to mark all events in all bookmarks, with each category represented by a separate colour.

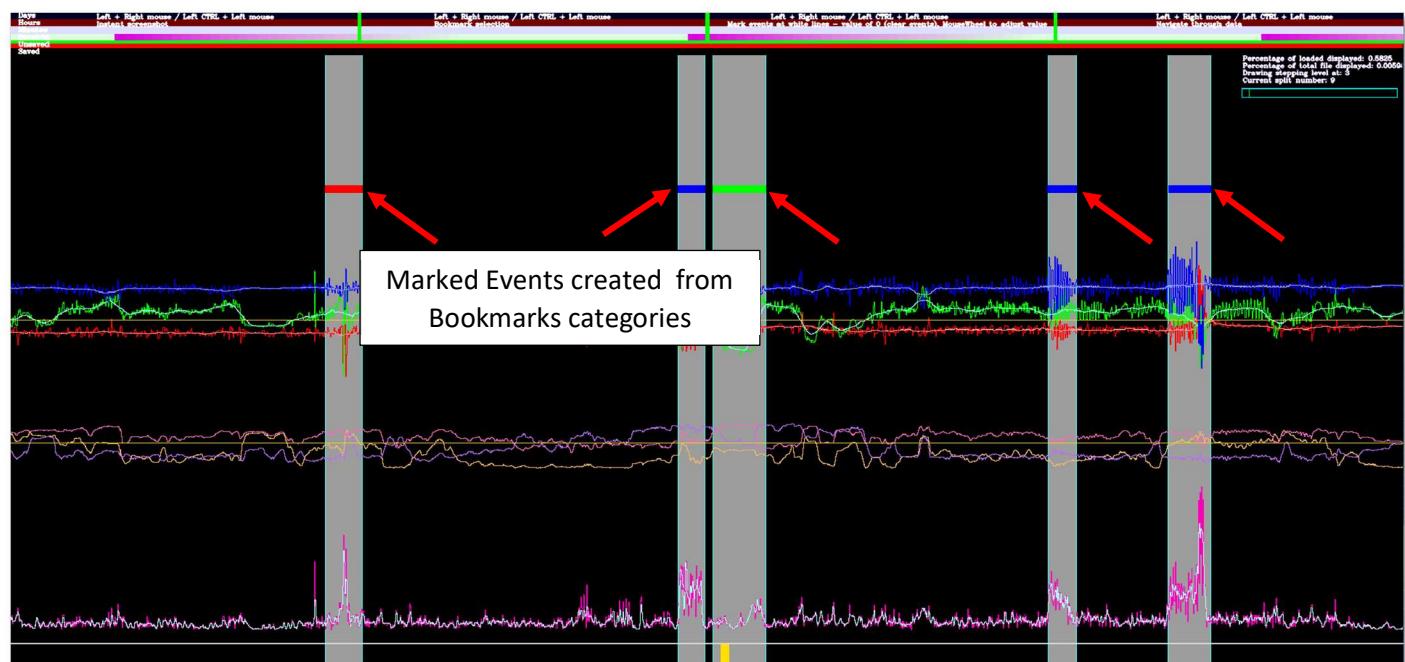


Figure 6.6 Marking events where each bookmark category has a distinct colour



Delete bookmarks

The delete functions are listed in the second column of the **Bookmarks / Bookmarks** tab. **Delete all bookmarks** deletes all bookmarks currently in file while **Delete current** only deletes the currently selected bookmark. Alternatively, the user can delete specific bookmarks by their size, range or category. **Delete bookmarks < size** deletes bookmarks below a certain number of events, which can be set in the adjacent box. **Delete bookmarks in range** deletes bookmarks X to Y in the yellow **Range** boxes. **Delete bookmarks in range / category** deletes bookmarks which are included in the currently set **Range** and match the set category in the adjacent drop-down list.

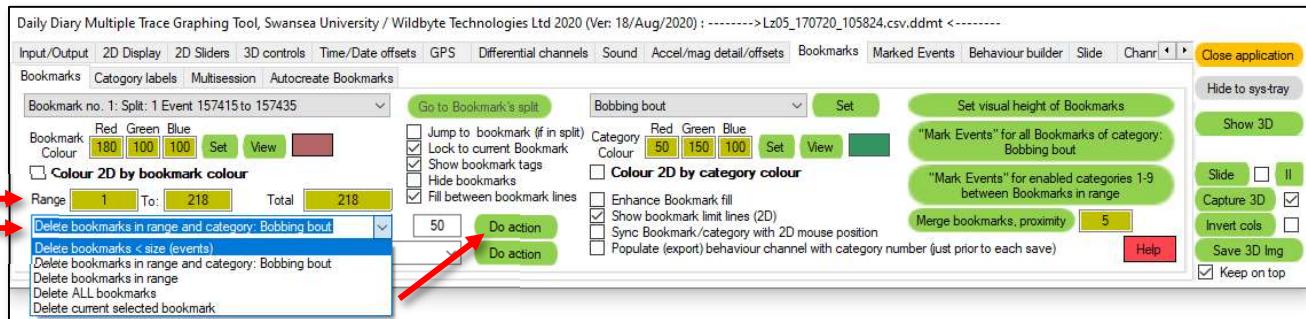


Figure 6.7 Deleting bookmarks in range and category

Additional bookmark options

Merge bookmarks, proximity

merges two bookmarks in close proximity to create one larger bookmark. Proximity is defined as being the number of events between the end of one bookmark, and the start of the next

Colour 2D by bookmark colour

each bookmark can have its own individually assigned colour with RGB value 0-255. Enter these colour values in the boxes and click Set to assign these values

Colour 2D by category colour

colours the bookmark to colours set for that category, which can be customised by adjusting the values in the **Category Colour** boxes, then using **Set** and **View** (separately for each category using the categories drop-down menu)

Lock

Jump to bookmark

enables navigation to whichever bookmark is selected in the bookmark drop-down menu (if within the currently loaded split). This zooms in on the bookmark

Sync bookmark and category with mouse

when selected will update bookmark and category drop-down lists to match any bookmark that the mouse moves over in the 2D graphing windows

Enhance banding

improves the contrast of the bookmarks on the black background

Fill between bookmark lines

controls whether the area between the bookmark limits is shaded or not

Go to bookmark's split

navigates to the split containing the bookmark currently selected in the bookmark drop-down menu

Populate behaviour channel

adds bookmark behaviour categories as a column in exported .txt file data (ASCII format) per event

Show bookmark tags

turns the information about the bookmark visible in yellow writing on and off

Show bookmark limit lines (2D)

toggles the visibility of the bookmarks start/stop boundary lines on / off

Hide bookmarks

allows all bookmarks to be hidden from the 2D "Main graphs" window

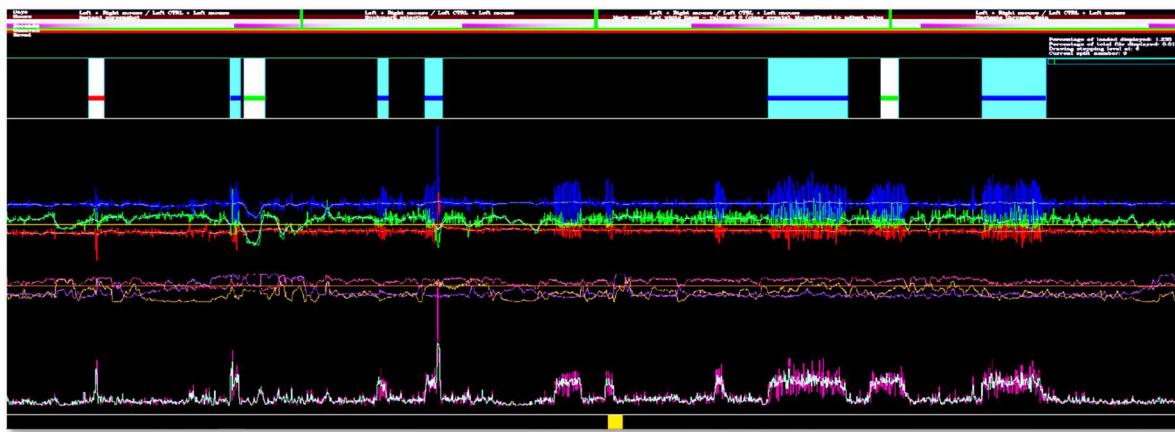


Figure 6.8 Bookmarks with marked events, adjusted visual heights and custom colours by category

Behaviour ethogram files

Behaviour ethogram files can be created and loaded into DDMT to create bookmarks based on lists of behaviours where behavioural observations have been conducted. Lists of observations should include the date, time and duration of the behaviour in seconds and behaviour category as a value. DDMT will then search for a match based on the start time and date, note the duration as width (seconds), and create bookmarks (by category per line) where appropriate.

Behaviour ethogram files only work when the format matches that shown in Figure 6.8 i.e. *DD/MM/YYYY* and *HH:MM:SS*, as DDMT looks to specific cell positions for the information it requires to create bookmarks. The cells shown in grey in Figure 6.8 should be copied exactly in the new file, in terms of both cell content and position. The **Red / Green / Blue** value columns (C-E) can be any value from 0-255 to allow custom colours for each category. The **type** column (A) allows category labels to be named, which are applied adopted by the category name lists in DDMT.

If using **R** to create this list, when you use *write.table*, ensure you add command *command = l.e.. write.table(df, "xxx", quotes=FALSE)* else the date and time column data will have quotes (Excel won't show these) and DDMT will refuse to load the data.

Two different formats exist to the user for creating ethograms.

Either **Date, Start, Duration, and Behaviour** (with **Events width** set to zero)

date	start	duration	behavior	decimal	events width
01/02/2017	08:42:11	55	1	0.15	0
01/02/2017	08:42:21	177	2	0.85	0

Figure 6.9 Bookmarks with marked events, adjusted visual heights and custom colours by category

This will create **Bookmarks** to the second with duration in seconds (events width must be zero for the software to recognise this is the format required). Decimal seconds will be used as best it can. In Figure 6.9, the start times would be taken as being 08:42:11.15, and 08:42:21.85. The duration values can also be decimalised, with carry-over from the start time's decimal value i.e. 08:42:21.55 with duration of 2.5 seconds, would result in the bookmark stop being set at 08:42:24.05 seconds.

Or, **Date, Start, Behaviour, Decimal, and Events width** (with **Events width** being >0, the software recognises that the **Duration** channel is to be ignored and that the width of the **Bookmark** will be in **events** not seconds).

date	start	duration	behavior	decimal	events width
01/02/2017	08:42:11	55	1	0.15	55
01/02/2017	08:42:21	177	2	0.85	177

Figure 6.10 Bookmarks with marked events, adjusted visual heights and custom colours by category

Here, decimal seconds is taken into account. The **Bookmarks** in the above data will start at approximately 08:42:11.15 and 08:42:21.85 s, and will be 55 and 177 events wide respectively.



	A	B	C	D	E	F	G	H	I	J	K
1	type	code	Red (0-255)	Green (0-255)	Blue (0-255)						
2	Playing	1	255	0	0						
3	Walking	2	0	255	0						
4	Standing	3	0	0	255						
5	Trotting	4	255	255	0						
6	Sniffing, walking	5	0	255	255						
7	Vigilant	6	255	0	255						
8	Running	7	128	0	255						
9	Shook	8	0	128	255						
10	Sniffing	9	255	0	128						
11	Rolling on back	10	255	128	0						
12	Social interaction	11	128	80	128						
13	Out of sight	12	128	50	80						
14	Blank	13	110	70	20						
15	Feeding	14	255	64	64						
16	Note that if a value in the events column below > 0 then -duration (seconds) is ignored and the width is calculated from the events value										
17	date	start	duration	behavior	decimal	events	width				
18	01/02/2017	08:42:11	0	1	0.15	177					
19	01/02/2017	08:42:21	0	2	0.85	50					
20											

Figure 6.11 Behaviour ethogram file format, cells in yellow must be included exactly as shown. Do not delete/alter them

As an example, after loading the main daily diary text file and setting the time and date (either manually or using a .tdo file), use the **Load bookmarks** function on the **Bookmarks** tab (Figure 6.12).

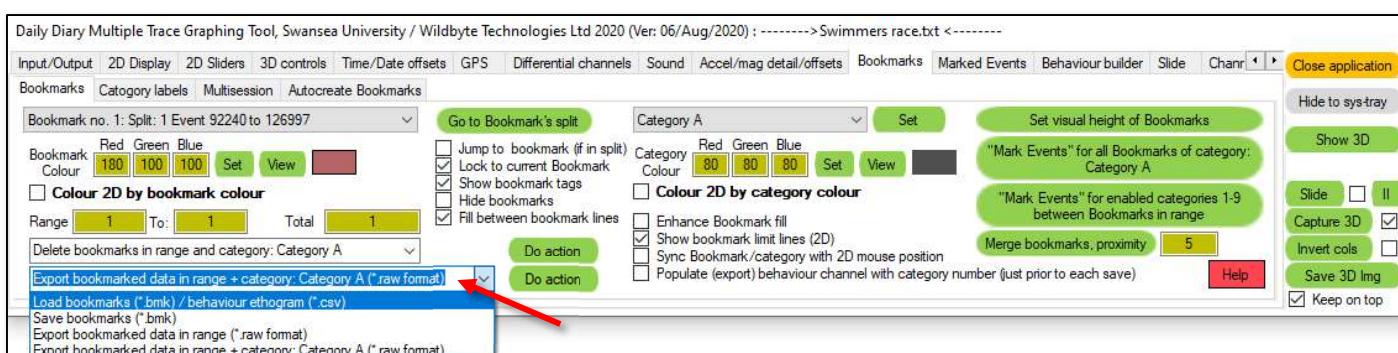


Figure 6.12 Select the first in the list and click **Do action** to load the ethogram

A window will open allowing the user to select the relevant file. Change **Files of Type:** to **Behaviour ethogram files (*.csv)** and select your file (Figure 6.13). After loading the behaviour ethogram file, a pop-up window will confirm the number of successful bookmarks loaded and the number which have failed to load, if applicable. This is when it may be helpful to use the **Console**.

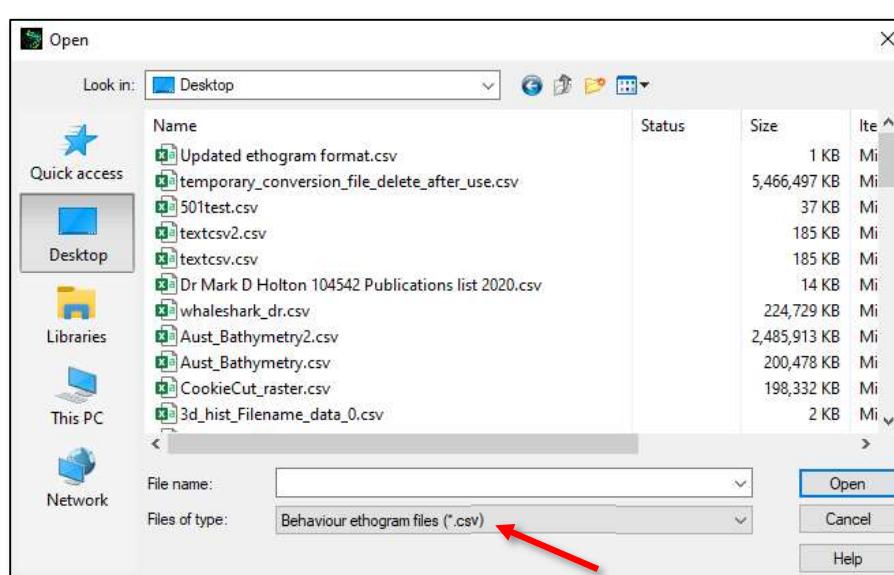


Figure 6.13 Select your file, ensuring you have selected the correct file type



DDMT will search for a time match for each ethogram line, add the duration in seconds and create a bookmark at that point for the noted behaviour category. The blue bookmark limit lines will not appear and the **Colour 2D by category colour** will be enabled by default. The bookmarks will fill the entire screen height initially, but the visual height can be manually adjusted as outlined above, **Adjusting the vertical height of the bookmark bands**. All bookmarks will be automatically added to the bookmark drop-down menu on the **Bookmarks** tab, including split and event numbers.

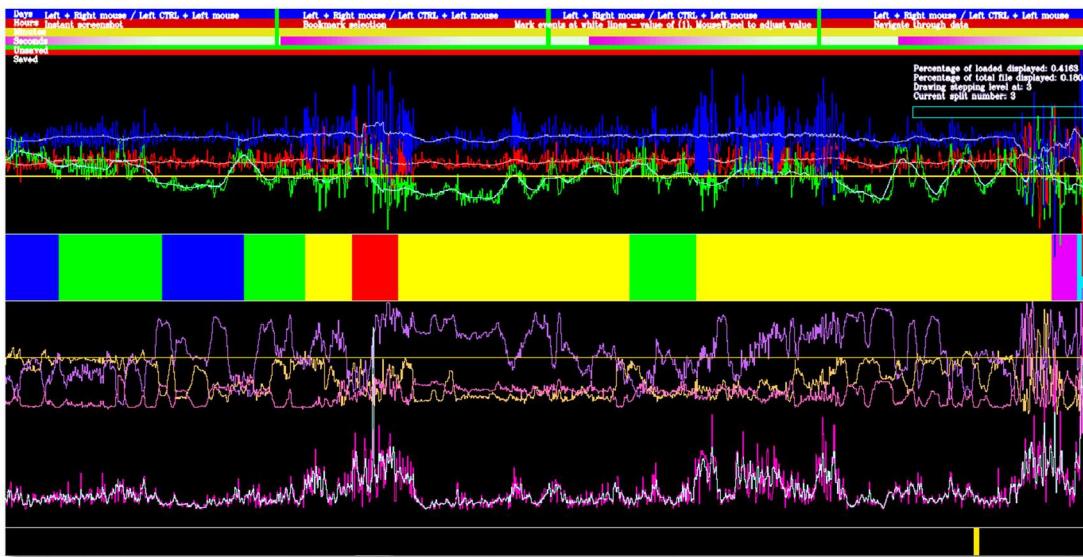


Figure 6.14 Bookmarks created by a behaviour ethogram file, with visual height of the bookmarks adjusted to a central band

Note that Behavioural ethogram files do not need to cover the whole span of time; gaps in the timeline are allowed.



Multisession (*data export*) of bookmarked data

A **Bookmark Multisession** is a file containing ASCII data from one or more **Bookmarks**. The columns contained within this file (partially named Master data file) are primarily dictated by the selection criteria in the **Input / Output / What to save** tab. Accompanying this Master data file is a Header file with a single entry per **Bookmark** of data in the Master data file that denotes where each **Bookmark**'s data starts and how many lines of (events) it has, along with other statistical information per **Bookmark**. This **Multisession** function is a quick way of exporting data that is of relevance/interest, perhaps containing only one type of behaviour, for analysis in another software package.

Multisession is controlled from the **Bookmarks / Multisession** tab where the user can select the parameters of interest (of each Bookmark) to be exported in the Header file in the blue box, Figure 6.15. These parameters include maximum and minimum values, means, medians, standard deviations, variances and ranges, as well as event times and dates. The **Header summaries only (no master data file)** checkbox controls whether an optional Master file is exported along with summaries i.e. sometimes it is not required to generate the Master data file, only the summary Header file.

To begin the multisession, click **Initiate BM multisession (text/column)**. A pop-up window will appear with the message "Bookmark Multi-session initialised." Clicking **OK** will bring up a second pop-up window stating, "Stop BM Multi-session to re-show Close application button." Notice that on the DDMT **Graphing Tool**, the orange **Close application** button is no longer visible and the **Initialise BM multisession (text/column)** button now reads **Stop BM multisession**.

Create bookmarks as normal. Use **Save bookmarks in range to multisession** or **Save bookmarks in range + categ to multisession** to export only bookmarks that have the currently selected Category and in the numbered range shown in the two boxes on the primary **Bookmarks** tab. Remember to click **Stop BM multisession** when finished, after file export is complete.


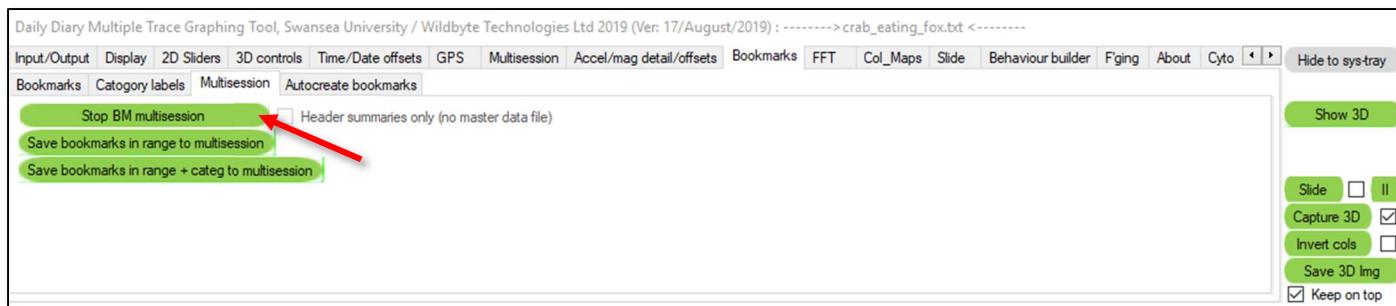


Figure 6.15 Bookmarks/Multisession tab before and after initialising the session

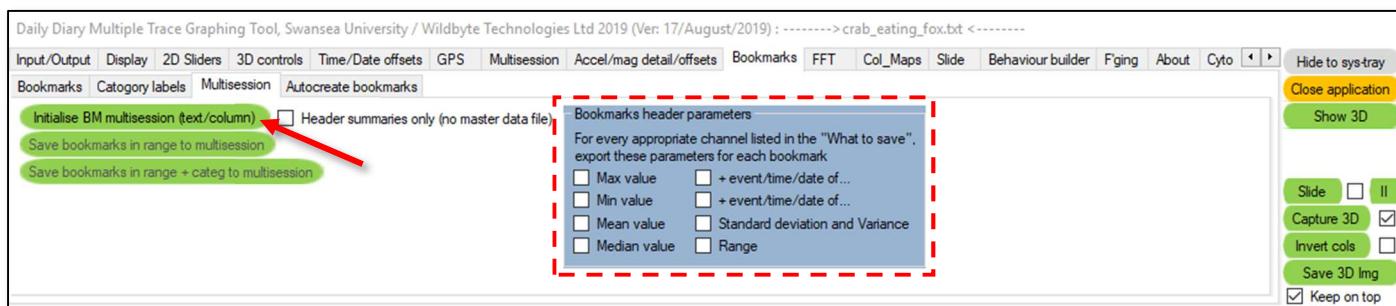

Auto-create bookmarks

Bookmarks can be created automatically for regular time periods throughout the whole file on the **Bookmarks/Autocreate bookmarks** tab. The **Minutes per bookmark** drop-down menu allows the user to select a 1-60 minute time period. Place the left white limit line (left mouse click) at the desired start point and click **Create bookmarks from left white line of this length** to create bookmarks for the specified time period from this start point to the end of the file (bookmarks created across splits in larger files). Alternatively, use **Autocreate bookmarks in split for Marked Events 1-9**.

Note that the time the first bookmark begins at the next boundary of “seconds = 0” from where the left white line has been placed (this was just algorithmically more simple to program!).



Figure 6.16 Autocreate bookmarks tab



Auto-create bookmarks in split for Marked Events

The user may have bands of **Marked events** within the split, or even perhaps in other splits, and might prefer to have these as **Bookmarks** in order to take advantage of the **Multisession** function that provide summary statistics of the data contained within each **Bookmark**. By clicking **Autocreate bookmarks in split for Marked Events 1-9** (for only **Marked events** within the currently loaded **Split**), or **Autocreate bookmarks, whole file, for Marked Events 1-9**, contiguous banks of **Marked Events** (contiguous being every event and of the same **Marked event** value) a **Bookmark** will be generated around them.

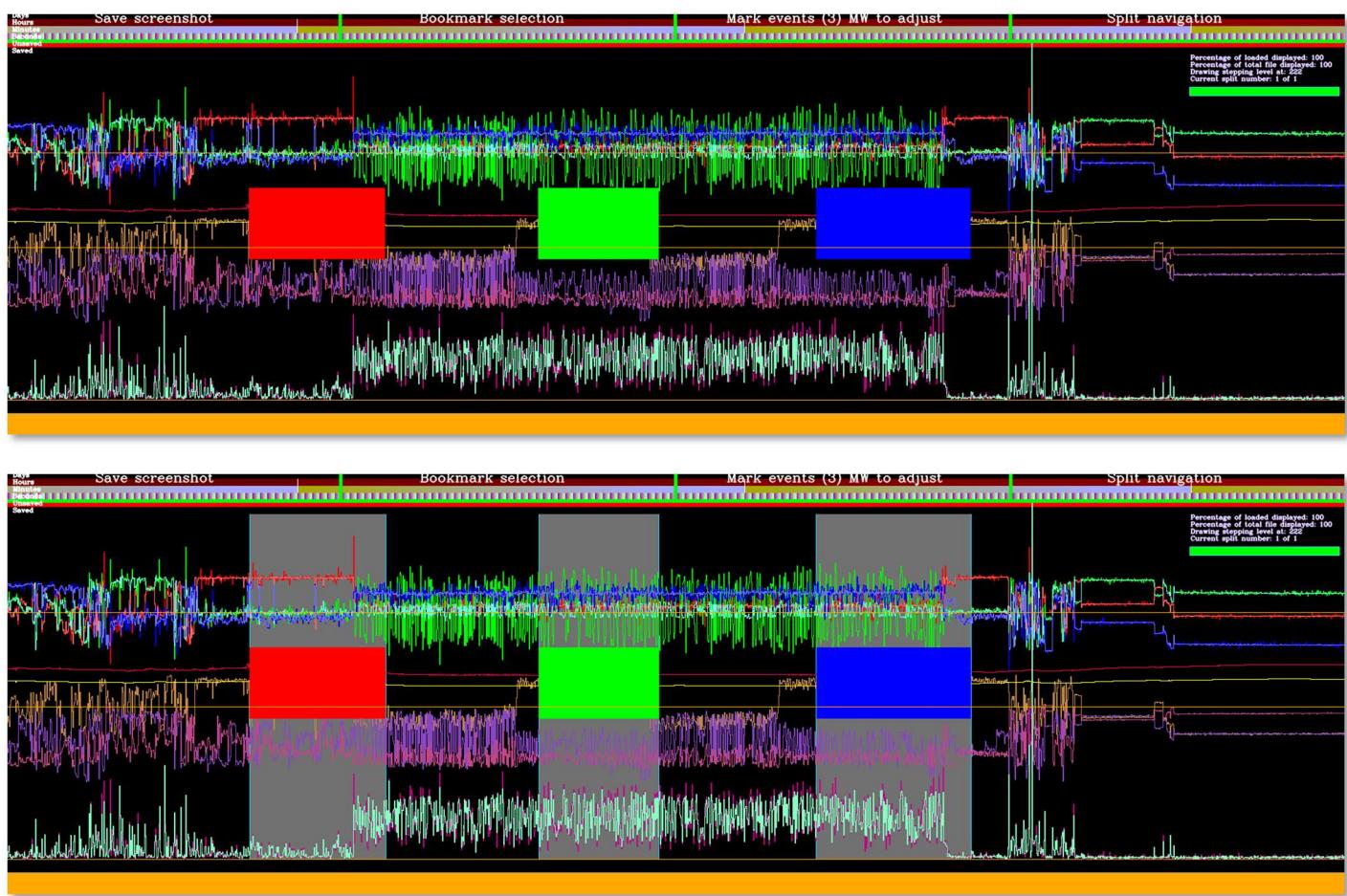


Figure 6.17 Autocreate bookmarks by converting Marked Events to Bookmarks



7. Tilt-Compensated Compass

The magnetometer and accelerometer on the Daily Diary can be combined to create a tilt-compensated compass. This means the accelerometer corrects the magnetic compass heading when the Daily Diary is tilted from the horizontal such that, when turning in a circle, from a plan-view, the magnetometer data would generate a circle as opposed to an ellipse. Circular magnetic data allows the generation of heading, while an ellipse would be a distortion of this.

Knowing the actual orientation of the Daily Diary relative to the ground can be helpful, but not essential, in informing DDMT which axis should be used to calculate pitch and roll respectively, and which axes are the most important to determine compass heading.

To correctly determine magnetic heading, your data must first be offset-corrected, achieved by following the steps outlined in chapter 3. This brings the bulk of the magnetic data to sit around the 0,0,0 origin.

Once the magnetic offsets have been corrected, click the **Show 3D** button on the right side of the control panel and right-click in the view and select **Device Orientation** near the bottom of the list (inset). This will then show a simplified image of the ground, and a forward pointing light-green triangle. The data logger hovering above here can be switched between the Square (as below), the Elongated, the Thumb, or the Alice.

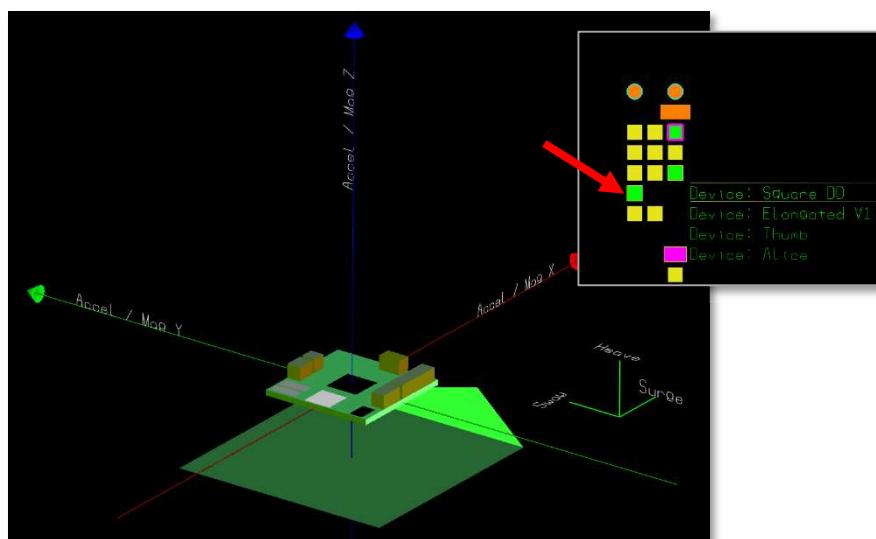


Figure 7.1 Setting Daily Diary orientation, device selection

The 5th button down in the third menu column allows the user to flick with the mouse wheel through all 24 possible orientations (perpendicular to the ground). The visual of the selected logger will orient itself accordingly. The user must left click this button to confirm the selection.



Figure 7.2 Logger orientation selection



Once confirmed, the magnetometer axes adjustments and pitch/roll accelerometer axes selections will be input into the DD Orientation tab interface, and heading, pitch, and roll etc. will all be calculated based on this preference.

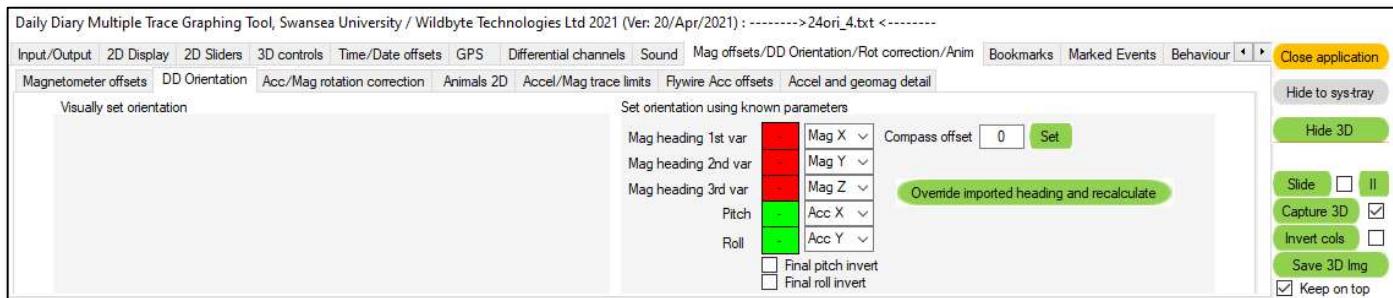


Figure 7.3 Manual orientation controls. The settings are automatically set through the 3D Device Orientation visual, but can also be used to set the parameters required for other manufacturer's device's data streams.

Any visualisations using heading, pitch, or roll, including the 2D graphs will update after they've been recalculated.

To save your magnetic heading data, click on the **Input/Output** and go to **File output/Save**. Ensure **Magnetic heading** is selected on the **What to save/Additional / metrics channels** sub-tab.

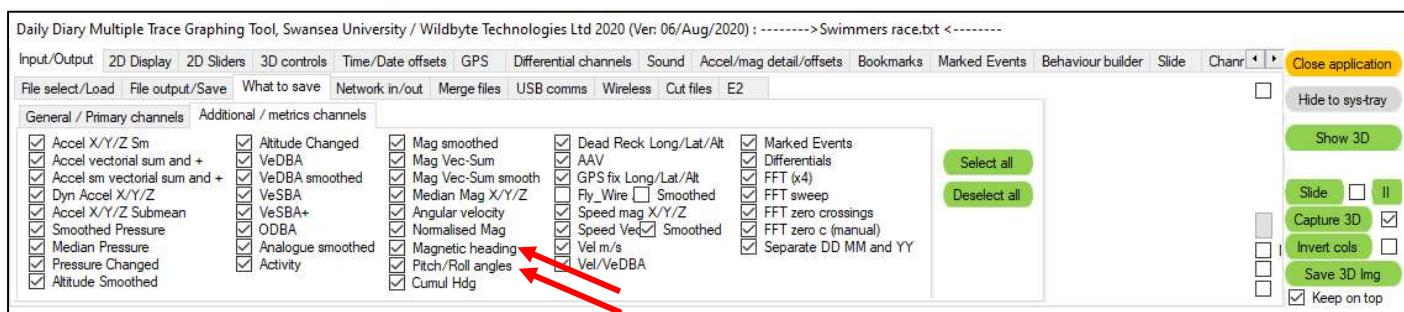


Figure 7.4 Ensure **Magnetic heading** is selected under **What to save**

Setting the device orientation manually by defining **Heave, Surge, and Sway**

Alternatively, to visually setting the logger orientation, one can observe the heave, surge, and sway in the 2D graph acceleration traces.

In the 2D graphs below, a motion test is performed with a logger where it is gently pitched up, and down, and then rolled left, and right, at all 4 compass points, before finally being held level, and spun around 360 degrees (see YouTube video: <https://www.youtube.com/watch?v=uDzJ6SQFeoQ>).

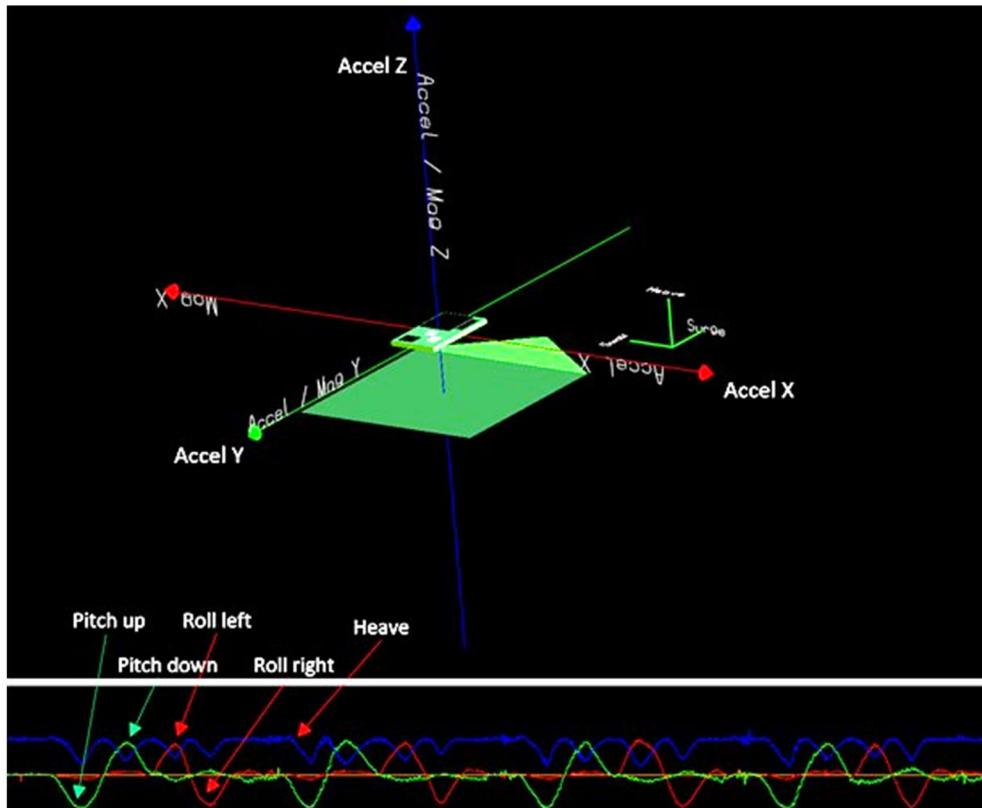


Figure 7.5 Observing acceleration traces to identify Heave, Surge, and Sway

It is best to start by setting Heave as, in the above, this is obviously, Z acceleration (blue), at 1g.

From here, you have 4 possible combinations; surge is +ve or -ve on one of the two remaining axes, and the same for sway.

The next step would be to define Surge as, perhaps with a dive behaviour acceleration (not pitch) data increases or decreases on the trace, when the animal pitches down. If the animal were to pitch up, but the accel trace goes negative, then a negative acceleration axis should be selected. In the above, surge is green (Accel Y), and so, for surge, select Accel -Y as it goes negative for positive pitch and vice versa. Sway can now only be one possible value, in this case Accel -X, as +ve Accel X is not possible due to the previous 2 selections.

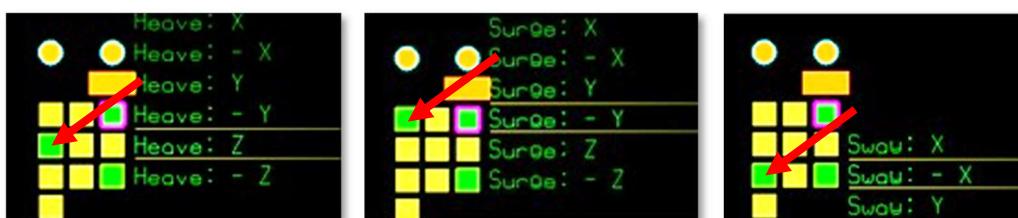


Figure 7.6 Setting the Heave, Surge, and Sway axes via the menu structure

Once Heave, Surge, and Sway have been selected and define a possible orientation, the ortho' lines all green (at the right side of Figure 7.5 will turn green, and the device orientation selection, as shown in Figure 7.2 and the visual orientation of the logger will switch to that defined by these parameters. Left clicking on the device orientation button (5th down) then causes the magnetometer and pitch/roll parameters to be correctly set on the DD Orientation tab (Figure 7.3).



Pitch and Roll

Pitch and Roll are determined from the acceleration data. Ideally, pitch and roll should be determined from static Daily Diary data but dynamic acceleration can be dealt with using the **Acc smoothing for compass** on the **2D Sliders/Channel smoothing** tab. The level of smoothing here is independent of any other acceleration smoothing in the software. The smoothed acceleration channels for compass heading can be viewed on the **drawing** window by selecting **Graph acceleration, for compass**.

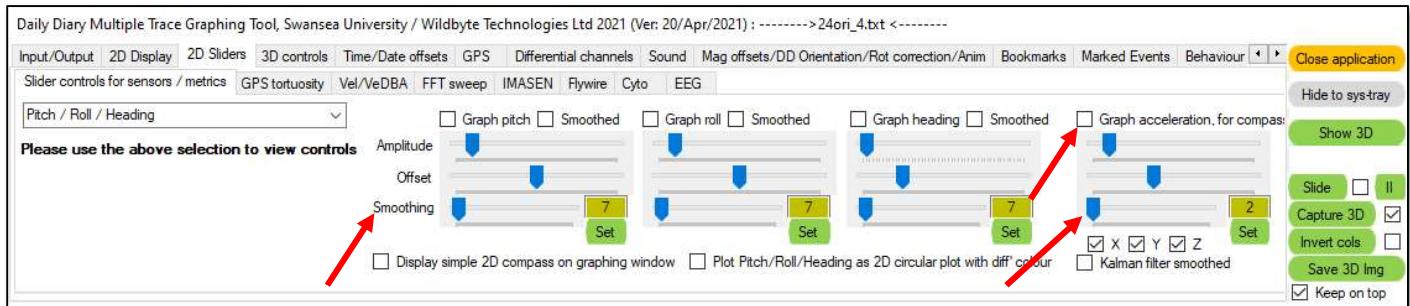


Figure 7.7 Set the acceleration smoothing level to remove dynamic acceleration from the signal

8. Fast Fourier Transform (FFT)

Fast Fourier Transform is a method which can be used to determine frequencies of regular oscillatory signals in accelerometry and geomagnetism data. For oscillating data on either acceleration or geomagnetism, the FFT function can be used to determine both the frequency and period of oscillation. In Figure 8.1, the strongest oscillating signal is in the **AccY** (green) channel:

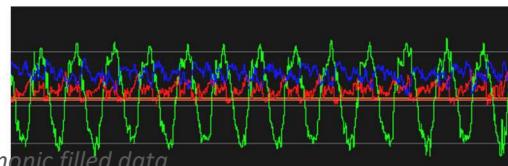


Figure 8.1 Example of oscillating, harmonic filled data

There are several FFT functions built into the software to analyse the data. The first, and easiest to use, is **FFT on the fly**, enabled by a checkbox on the **FFT** tab (Figure 8.2). The small scroll buttons (indicated with 3rd arrow below) may have to be clicked to scroll along to the **FFT** tab.

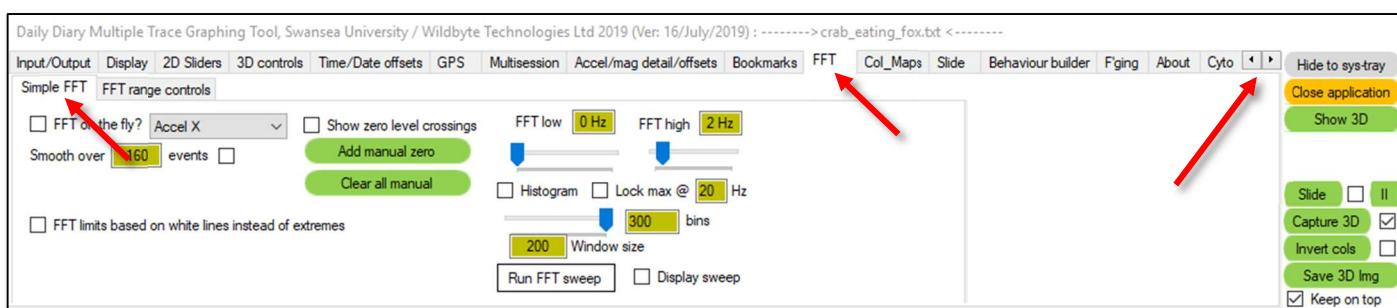


Figure 8.2 Select the **FFT/Simple FFT** tab for access to the **FFT on the fly** function

Note: The **FFT on the fly** function will be disabled if the extremes width (the number of events on display in the **drawing** window) exceeds 50,000 events.

After enabling the **FFT on the fly** function, the next option is to select the appropriate channel. Use the drop-down menu (default **Accel X**) to choose from any of the accelerometer or magnetometer, raw or smooth channels, or the magnetometer's **Vec-sum** or **Vec-sum smoothed** channels. The separate FFT window displays the frequencies and their amplitude components (Figure 8.3). The **FFT** routine automatically picks out the frequency with the greatest amplitude and displays the frequency and period based on this peak for the data currently displayed on-screen.

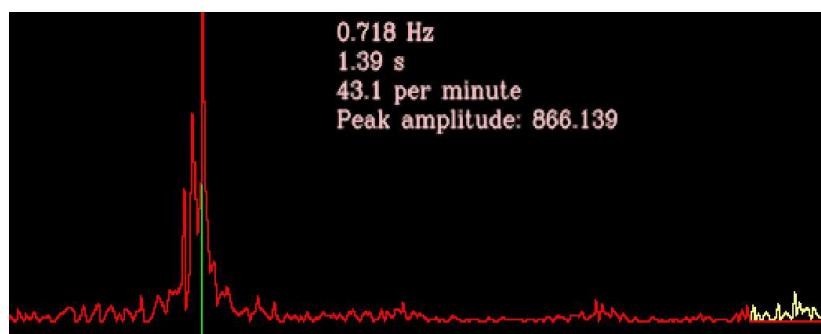


Figure 8.3 The frequency spectrum window showing a real-time adjusted spectrum based on the channel/data in-view on the drawing window



The controls on the next tab, **FFT range controls**, allow the user to zoom in on a particular range of frequencies within this spectrum (Figure 8.4). If the zoom excludes the current main peak, then a new main peak will be selected. This feature might be useful if there is more than one dominant frequency within the data.

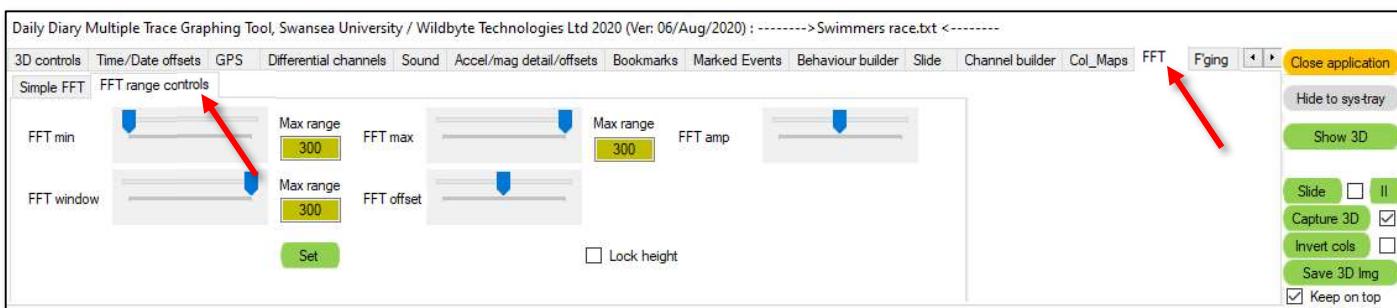


Figure 8.4 The **FFT/FFT range controls** tab, controls to zoom into the frequency spectrum

The **drawing** window will now display a yellow trace. This trace is a slightly cleaned (i.e. less higher harmonics) inverse FFT of the source trace (Figure 8.5).

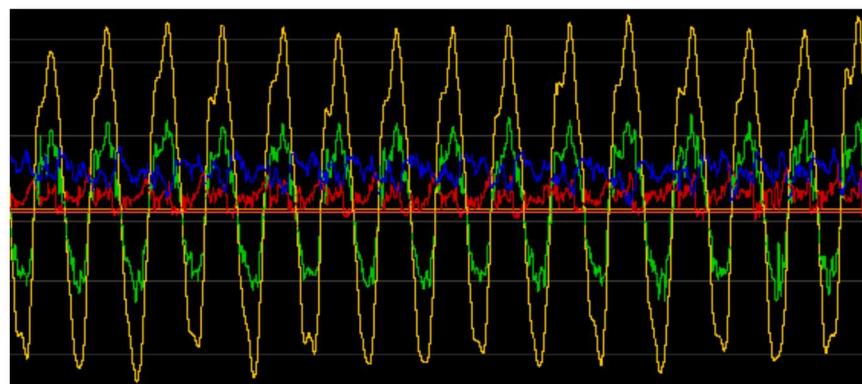


Figure 8.5 Yellow trace, the inverse FFT of the source waveform; in this instance, Accel Y (green)

The controls are simply **FFT min** and **FFT max** frequencies (the limits visible in the separate FFT window shown in Figure 8.3), **FFT window** which defines the spectral width centred around the currently auto-selected main frequency peak, and finally **FFT amp** and **FFT offset** which set the amplitude and offset of the yellow trace on the **drawing** window.

Reducing the **FFT window** slider eliminates harmonics from the inverse FFT waveform, thus eventually reducing it to that of the fundamental (main peak) frequency (Figure 8.6). The FFT window shows the same trace in yellow with the 'windowed' trace overlaid in red. Winding this 'FFT win' to the left reduces the width of the window and thus, at the low end, will result in only the **fundamental** being passed back i.e. the single highest peak. The frequency and period measures will remain the same, as these are based on the dominant peak value.

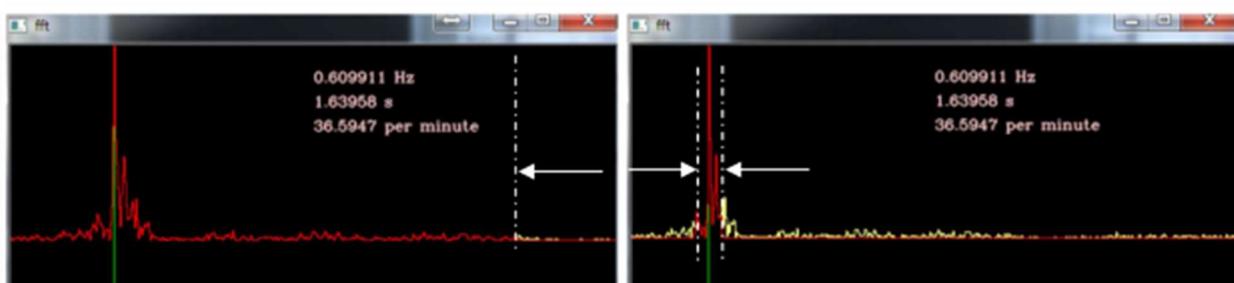


Figure 8.6 A frequency spectrum on the left, with a window narrowing down what will be presented to the **inverse FFT** function to then generate the resulting yellow waveform on the **drawing** window

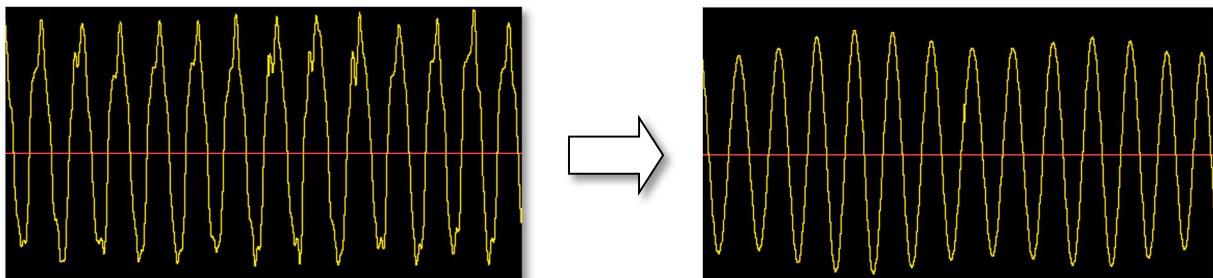


Figure 8.7 Cleaning up the inverse FFT waveform

This second waveform (Figure 8.8) is almost purely sinusoidal as there are very few harmonics present, and will remain in-phase with the source signal:

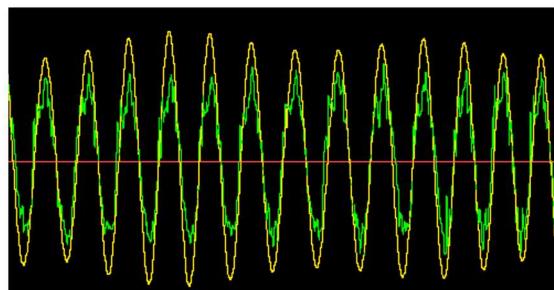


Figure 8.8 In-phase source and inverse FFT waveforms

FFT on the fly is calculated in real-time. The user therefore can slide through the data either by adjusting the extreme limits of the orange bar at the bottom of the **drawing** window, zooming with the mouse-wheel, or holding down Left-alt and left clicking in the **drawing** window to relocate focus on some other section of data. As long as the extremes width is less than 50,000 events, the function should continue to work without issue.

This inverse FFT waveform is also stored in memory and is saved along with all other data. In the example below, a **Selection** of data has been made using the two white extreme lines:

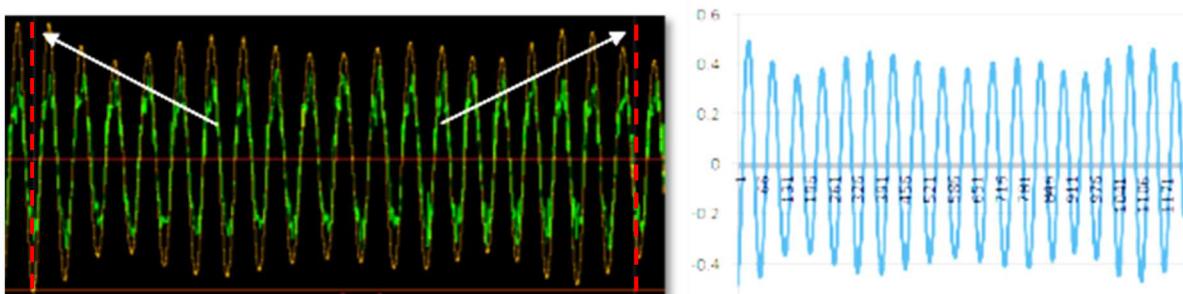


Figure 8.9 Inverse FFT data is exported, in alignment with the source channel data

The dominant frequency and period of the data viewed on-screen can be displayed on the overlay in the main graph window by adjusting the **Overlay selection** within the **Display** tab. These values are the same as those displayed in the separate **FFT** window.

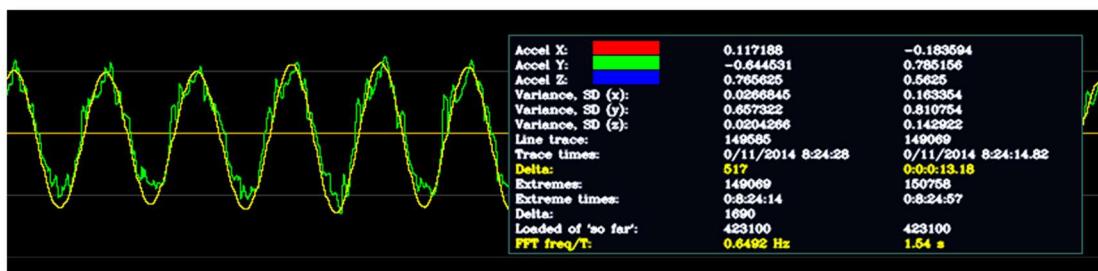


Figure 8.10 FFT peak frequency and period displayed in the overlay table



9. Importing GPS Data (Pre-Loading)

GPS data files containing parameters **Date**, **Time**, **Latitude**, **Longitude**, and **Altitude** can be imported into DDMT by saving any GPS data as a .csv file, i.e. comma-delimited, with the column headers as shown in Figure 9.1. Note all other columns will simply be ignored during the import process. If altitude data is not included, ensure there is still an altitude column with zero values or the file will not load correctly.

The **Date** column must be in DD/MM/YYYY format, whilst the **Time** column must be in a 24 hour clock with formatting as HH:MM:SS. **Latitude** and **Longitude** will be assumed to be in decimal degrees, and **Altitude** in metres. The ordering of the columns is equally important as the header row is not read.

	A	B	C	D	E	F	G
1	Date	Time	Latitude	Longitude	Altitude	Speed	Co
2	31/10/2014	16:05:42	54.58593	-5.94168	36.52	720	
3	31/10/2014	16:07:41	54.58593	-5.94168	36.52	2196	
4	31/10/2014	16:11:36	54.58591	-5.9419	13.55	1296	
5	31/10/2014	16:13:53	54.58599	-5.94227	20.57	3312	
6	31/10/2014	16:16:01	54.58576	-5.94168	84.01	1800	
7	31/10/2014	16:17:44	54.58604	-5.94221	63.82	468	
8	31/10/2014	16:19:19	54.58596	-5.9421	43.73	1044	
9	31/10/2014	16:20:51	54.58575	-5.94092	18.54	828	
10	31/10/2014	16:22:20	54.58578	-5.94121	50.45	1620	

Figure 9.1 Format required of GPS data to be imported

To import a GPS .csv file (this must be imported prior to selecting / loading the main data file) click **Preload GPS data file** at the top left of the **Input / Output** tab. Once the GPS data has been imported, the user is then free to select a Daily Diary data file and load the first split.

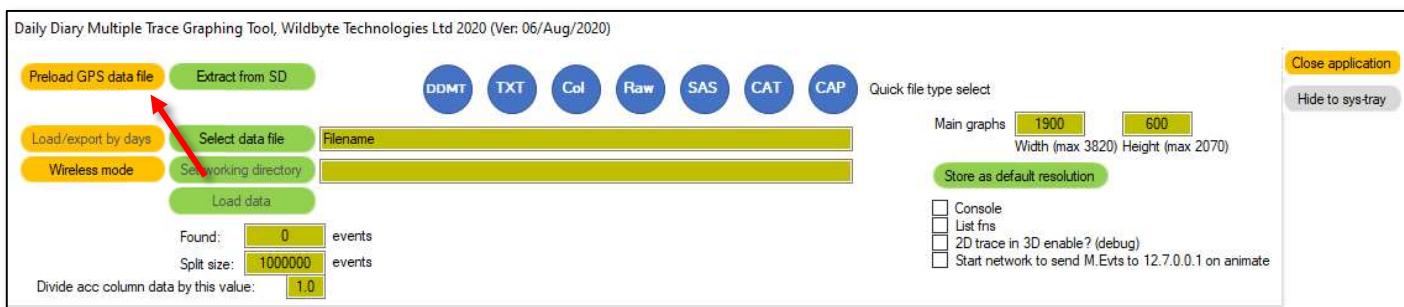


Figure 9.2 Preload GPS data file button must be selected before clicking Select data file



Syncing GPS data to Daily Diary data

By synchronising the GPS with DD data sets, DD data becomes a great deal more powerful due to the ability to couple environmental attributes with behaviours. GPS data will be assumed to have the correct time and date information. To set or correct the Daily Diary data time and date, refer to Chapter 3.

Before synchronising GPS data to Daily Diary (DD) logger data, it is important to ensure that any date/time corrections are performed on the DD data, as the sync function looks for matches in both date and time information. The maximum time resolution DDMT performs is 1 Hz so GPS data must be prefiltered to 1 Hz prior to presenting it to DDMT.

Synchronisation will occur from either the start of the Daily Diary dataset, or, from the event number where time/date corrections have been applied, if at all.

To sync the two data sets, click the **Sync GPS file to DD** button on the **GPS** tab. DDMT may be non-responsive for a few seconds or more, depending on how much data it has to assess; please be patient with large datasets. A red banner on the **GPS** tab will state **Syncing GPS** and switch to **Sync complete** when finished.

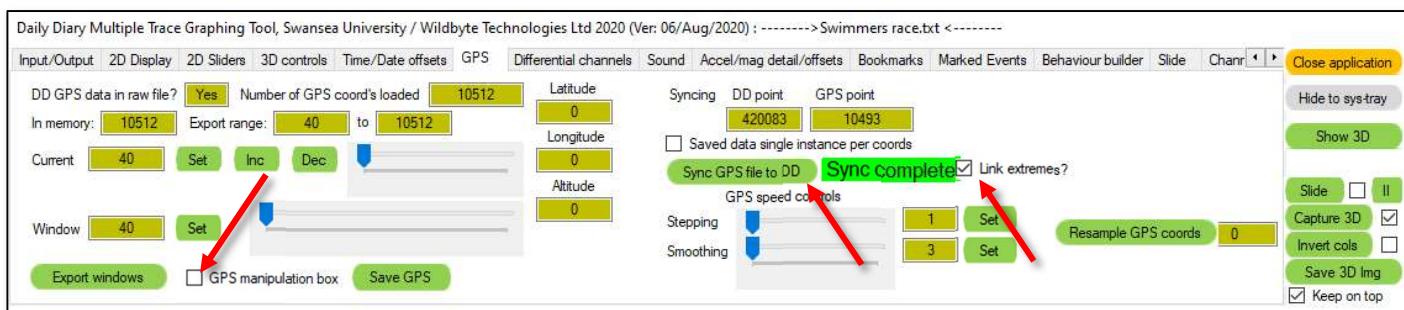


Figure 9.3 Sync GPS file to DD on the GPS tab

Link extremes tickbox must then be clicked to ensure that GPS sync keeps up when the user zooms in/out of the data for other parts of the software.

To show the synchronicity of the DD and GPS data sets, close the **GPS manipulation box** (untick if selected) and enable the 2D graphic view of the GPS fixes by ticking the **Show GPS** checkbox on the **2D Sliders/GPS** sub-tab (Figure 9.4). This is also where the amplitude and offset of the GPS trace can be adjusted on the **drawing** window. It is also possible to add a velocity graph (measured in metres per second, m/s) calculated from GPS fixes and adjust its amplitude and offset.

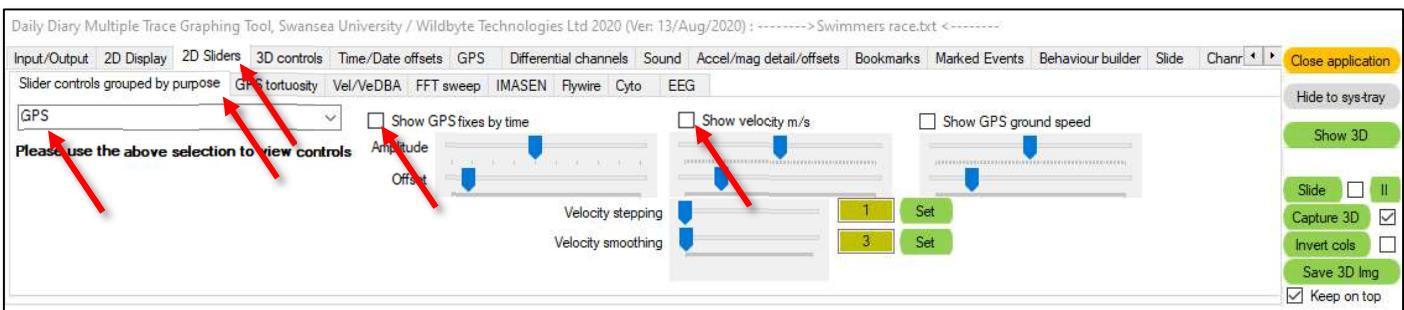


Figure 9.4 Show GPS and Show velocity m/s settings



The **drawing** window will now show a vertical blue line for every DD data point and GPS fix with the same time and data attributes. 1 Hz is the maximum frequency for syncing data, therefore with Daily Diary data at 40 Hz, there will be 40 blue lines for every GPS fix match.

The first blue Daily Diary event fix synchronised to a GPS fix will also have a small red dot (the other blue lines are considered in-sync due to sub-second sampling frequency).

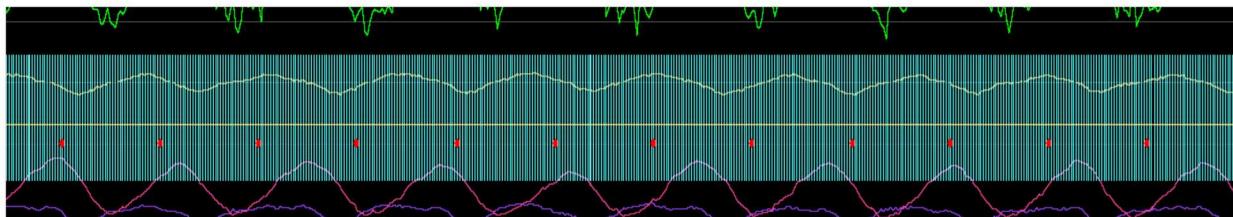


Figure 9.5 Blue lines indicating a DD to GPS fix match

Figure 9.6 shows us DD to GPS fix matches in blue, with a blue base line where a DD to GPS fix match did not occur.



Figure 9.6 Blue lines indicating a DD to GPS fix match while blue dots indicate a lack of a fix match

Once the data has been successfully synchronised, **GPS Long/Lat/Alt** and **GPS mean velocity m/s** for events that have a GPS fix match can be seen in the **Overlay** on the **Display/Overlay selection** sub-tab.

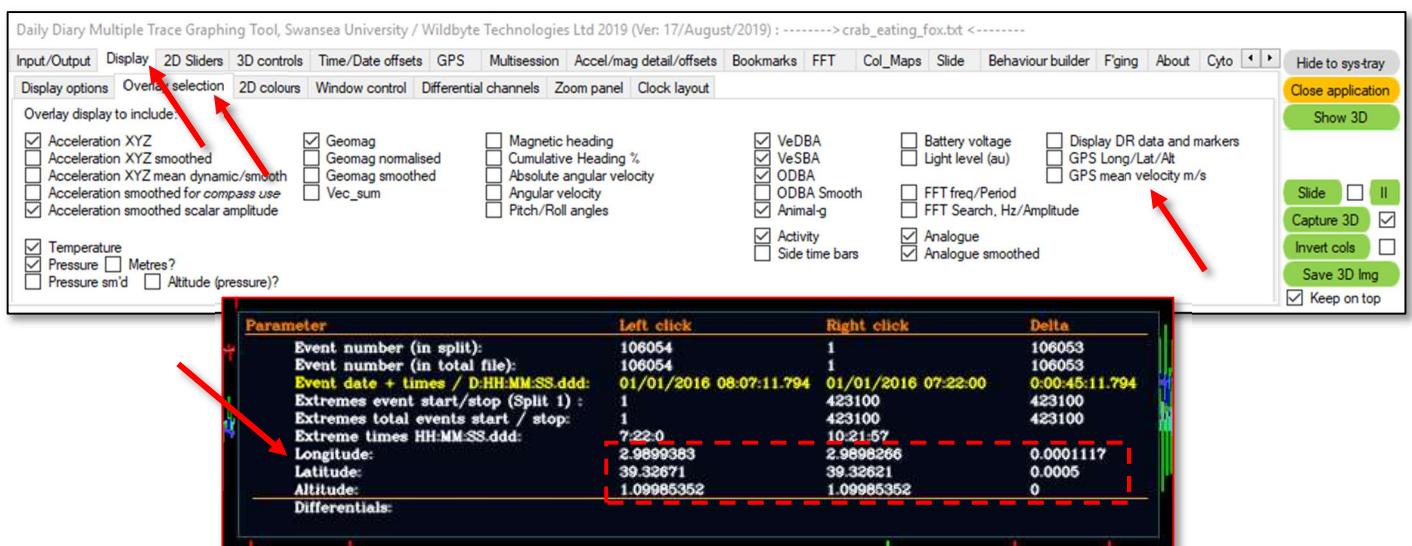


Figure 9.7 Enabling GPS information on the Overlay

The values will then appear at the bottom of the **Overlay**, for data points that have a DD to GPS fix match. GPS parameters are selected for export by default on the **Input/Output/What to save/Additional / metrics channels** sub-tab.


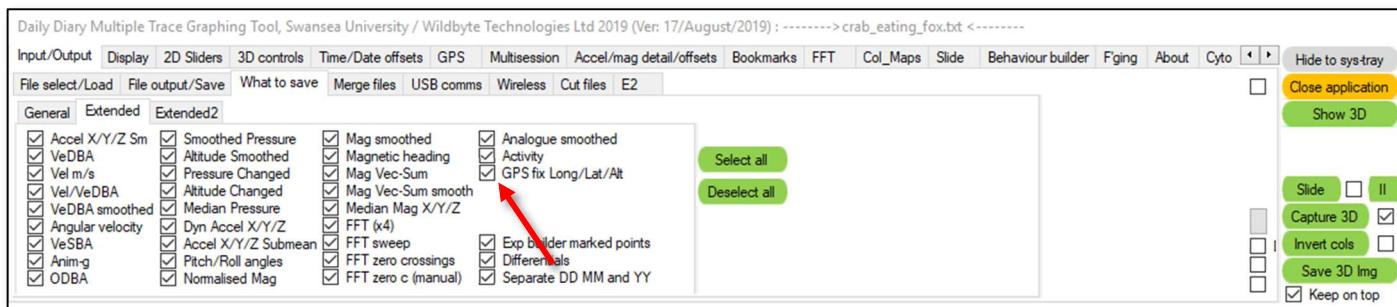


Figure 9.8 Exporting the DD to GPS fix match GPS coordinates

GPS straight-line interpolation

It is possible to fill in gaps between GPS points to provide a second-by-second fix, simply derived by straight-line interpolation between two consecutive GPS fixes. This should be used with care as it does not consider heading or potential changes in velocity but fills data gaps. Click **Resample GPS coords** on the **GPS** tab and click the **Sync GPS file to DD** again to account for interpolation.

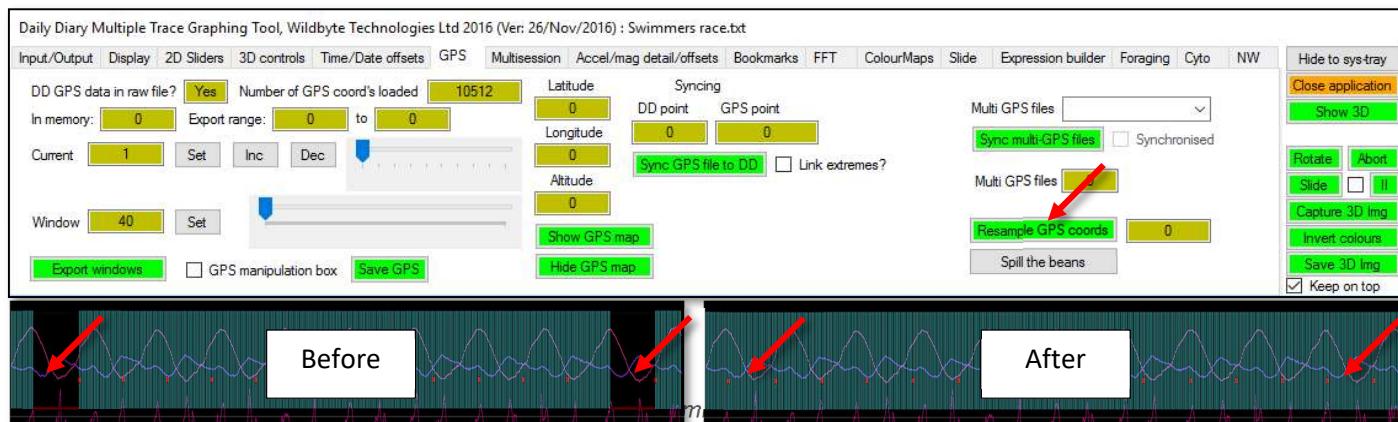

Bracketing

With the DD and GPS datasets synchronised, the software can slide through the data, centring on GPS points one after another with a user-defined window width of events around that point. To do this, open the **GPS** tab.

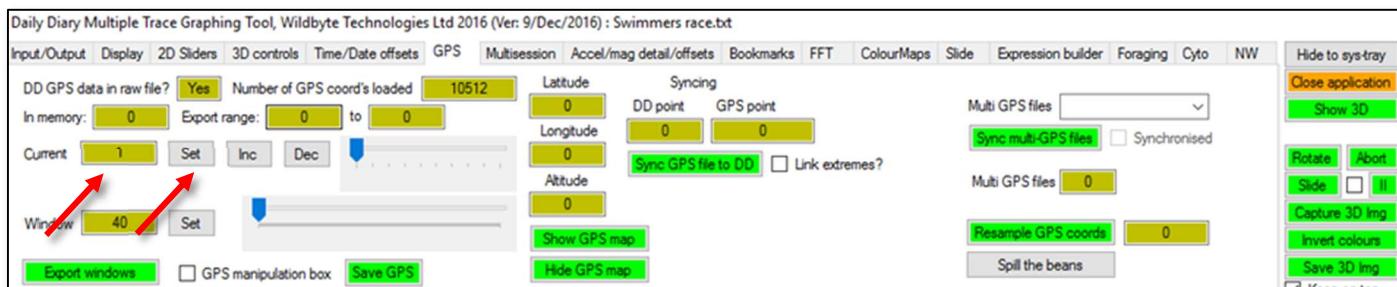


Figure 9.10 GPS data bracketing controls

Current is the currently displayed GPS point with a window of 1000 events around the point (+/- 500 points, limited at either extreme of the data set). To enable this function, simply click **Set** next to **Current**'s value-box. The **drawing** window will immediately lock down and display only 1000 events with this GPS data point at the centre. Clicking the green **Inc** or **Dec** buttons will move the window forward or backward to the next GPS fix match. Slider bars are also provided to the right of the GPS point window width buttons. As each GPS point is shifted into position, the **Longitude**, **Latitude** and **Altitude** value boxes will



automatically update. The centred points are the red dots described earlier in this GPS section, which denote the first event to have a time/date match, between the Daily Diary and GPS data sets.

In memory tells the user how many GPS points are in memory (within the currently loaded split of the full dataset) and the **Export range** allows the user to specify which points will be used in the data export process described below. If all the Daily Diary data has been loaded as a single split, then this will be from 1 to n . The **Current** slider and **Inc/Dec** buttons will be confined to the allowed range.

The **Export windows** function makes the software cycle through all the in-memory GPS points (within the user specified **Export range**, and export all the windows of Daily Diary data surrounding them as a single file; the data channels exported is as per the **What to save** tab tick-boxes. Please note that this can create enormous files! If you consider having 10,000 GPS data points each having a Daily Diary data window of 1,000 events – it'll create a single file (with a single line break between each window) of 1,000x10,000 events, of n channels of data.



10. Behaviour Builder

Behaviour builder is a Boolean conditional search function, used to find data that fits user-specified criteria. This can be in the form of a single statement that will generate **Marked Events** if true, or a double statement that, when found true, will create a new **Bookmark**. The strength of this utility lies both in the quality of the data, and in the user's ability to identify essential criteria.

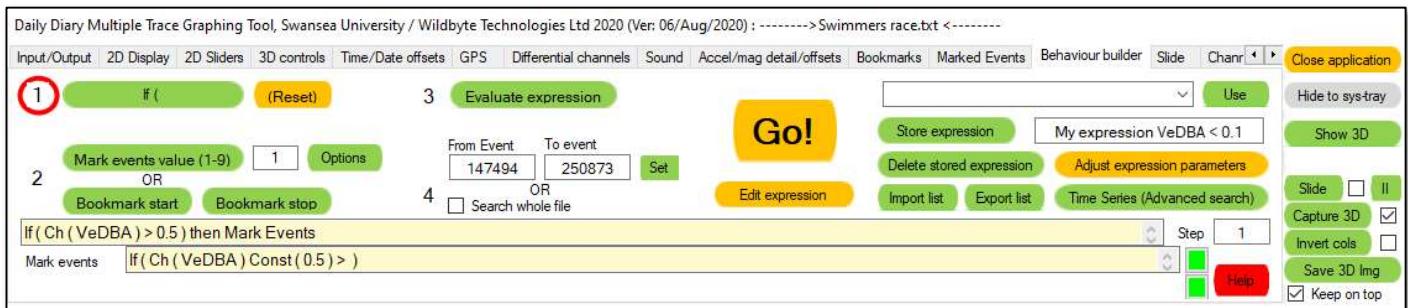


Figure 10.1 The initial behaviour builder tab

Some examples below of Behaviour Builder expressions (notice the brackets around the inequalities before and after the **AND** function – required (also for OR and NOT) to ensure correct parsing of these types of double-barrelled expressions):

- `If((SM(Acc X, 50) > 0.5) AND (Temperature < 28)) then mark events`
- `If((VeDBA Smoothed > 0.5) AND (ABS(Acc X Smoothed) < 0.1)) then mark events`
- `If((Temperature < 25) AND () AND (GPS Lat < 4.231) AND (GPS Long > 52.514) AND (GPS Long < 53.154)) then mark events`



Using Behaviour Builder to create marked events

To build our first expression, perform the following:

1. Click **If(** button

The view will now change to allow us to enter our expression

Daily Diary Multiple Trace Graphing Tool, Swansea University / Wildbyte Technologies Ltd 2020 (Ver: 06/Aug/2020) : ----->Swimmers race.txt <-----

Input/Output 2D Display 2D Sliders 3D controls Time/Date offsets GPS Differential channels Sound Accel/mag detail/offsets Bookmarks Marked Events Behaviour builder Slide Chann [] Close application

Standard Channel VeDBA ▾ Pre-smooth? (by:) 2 Enter * + < > = 0.5 Enter

Differential Channel ▾ Pre-smooth? (by:) 2 Enter / - () AND OR NOT

ABS (Vector FFT Square root

GRADIENT(# Variance Blindsight

If(Next step Undo last entry Step 1

Mark events If(Ch(VeDBA) Const(0.5) >) Help

Show 3D

Slide Capture 3D Invert cols Save 3D Img Keep on top

2. Now click the drop-down box indicated in Figure 10.3 and change it to VeDBA and click **Enter** to its right

Daily Diary Multiple Trace Graphing Tool, Swansea University / Wildbyte Technologies Ltd 2020 (Ver: 06/Aug/2020) : ----->Swimmers race.txt <-----

Input/Output 2D Display 2D Sliders 3D controls Time/Date offsets GPS Differential channels Sound Accel/mag detail/offsets Bookmarks Marked Events Behaviour builder Slide Chann [] Close application

Standard Channel VeDBA ▾ Pre-smooth? (by:) 2 Enter * + < > = 0.5 Enter

Differential Channel ▾ Pre-smooth? (by:) 2 Enter / - () AND OR NOT

ABS (Vector FFT Square root

GRADIENT(# Variance Blindsight

If(Ch(VeDBA) Next step Undo last entry Step 1

Mark events If(Ch(VeDBA) Const(0.5) >) Help

Show 3D

Slide Capture 3D Invert cols Save 3D Img Keep on top

3. Now click **>**, enter 0.1 into the values box and click its **Enter** button, followed by **Next step**

Daily Diary Multiple Trace Graphing Tool, Swansea University / Wildbyte Technologies Ltd 2020 (Ver: 06/Aug/2020) : ----->Swimmers race.txt <-----

Input/Output 2D Display 2D Sliders 3D controls Time/Date offsets GPS Differential channels Sound Accel/mag detail/offsets Bookmarks Marked Events Behaviour builder Slide Chann [] Close application

Standard Channel VeDBA ▾ Pre-smooth? (by:) 2 Enter * + < > = 0.1 Enter

Differential Channel ▾ Pre-smooth? (by:) 2 Enter / - () AND OR NOT

ABS (Vector FFT Square root

GRADIENT(# Variance Blindsight

If(Ch(VeDBA) > 0.1) Next step Undo last entry Step 1

Mark events If(Ch(VeDBA) Const(0.5) >) Help

Show 3D

Slide Capture 3D Invert cols Save 3D Img Keep on top

4. Now click **Next step** button, and the previous view will be returned

Daily Diary Multiple Trace Graphing Tool, Swansea University / Wildbyte Technologies Ltd 2020 (Ver: 06/Aug/2020) : ----->Swimmers race.txt <-----

Input/Output 2D Display 2D Sliders 3D controls Time/Date offsets GPS Differential channels Sound Accel/mag detail/offsets Bookmarks Marked Events Behaviour builder Slide Chann [] Close application

1 If((Reset) 3 Evaluate expression Go! Use

2 Mark events value (1-9) 1 Options Store expression My expression VeDBA < 0.1

OR Delete stored expression Adjust expression parameters

Bookmark start Bookmark stop 4 From Event To event Set

From Event To event 147494 250873 Set

OR Edit expression Import list Export list Time Series (Advanced search)

Search whole file

If(Ch(VeDBA) > 0.1) Step 1


Mark events If(Ch(VeDBA) Const(0.5) >) Help

Show 3D


Slide Capture 3D Invert cols Save 3D Img Keep on top




5. Now there is the choice of defining either a single expression that will evaluate single events as true or false, or double expressions that will define the start and stop of a **Bookmark**. In this first example, click **Mark events value (1-9)** to generate **Marked Events**

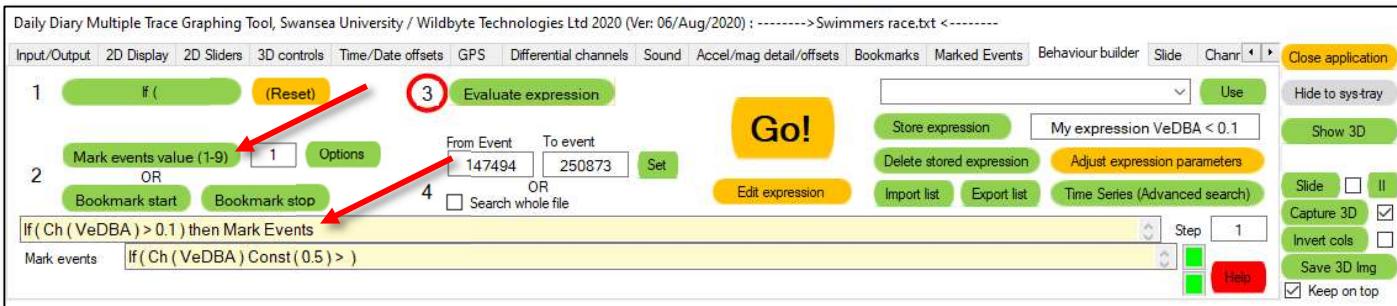

6. Now we need to tell the algorithm that we've completed the expression and request it convert it to Reverse Polish Notation (RPN), a format more easily parsed rapidly by a search algorithm. Click **Evaluate expression** and the translated expression appears in the lower box as **If (Ch (VeDBA) Const (0.1) >)** (RPN)

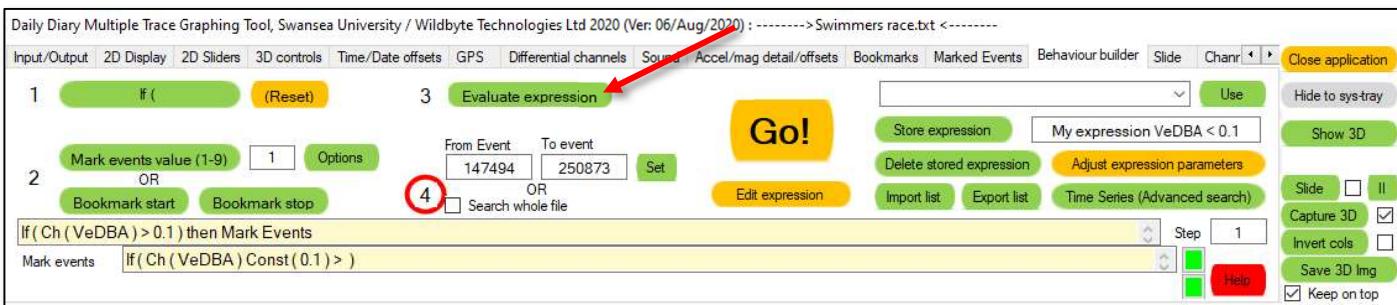

7. Next, define what data we want to evaluate, so place white lines using the left/right mouse buttons on the **2D graphing window**, and click **Set**, or click **Search whole file**. Note - if **Search whole file** is ticked, this will override the start/stop event markers in the two boxes.

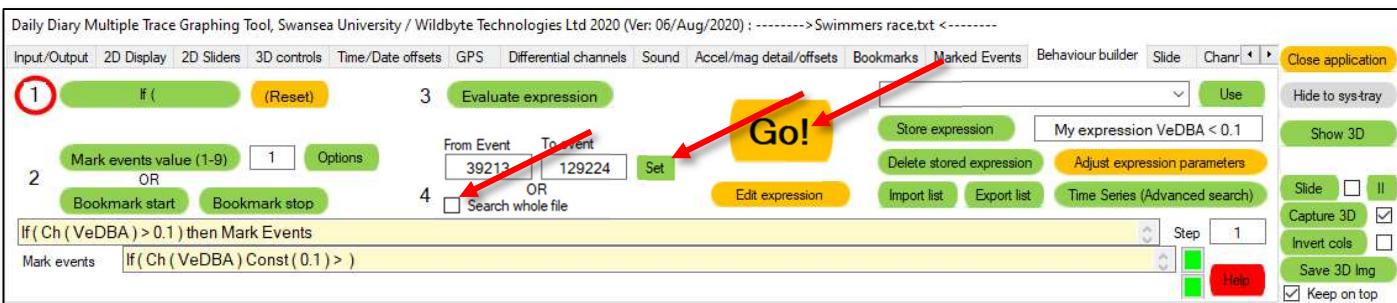

8. Finally, click **Go!** Because a value of 1 was in the box adjacent to the button labelled **Mark Events value (1-9)**, this will place a **Marked Event** value of 1 on any event where the VeDBA is less than 0.1 g. By default, **Marked Events** value 1 is red. The colour scheme can be adjusted on the **Marked Events** tab.

Here, some data with low (< 0.1 g) has been “Marked” as value 1 (red):

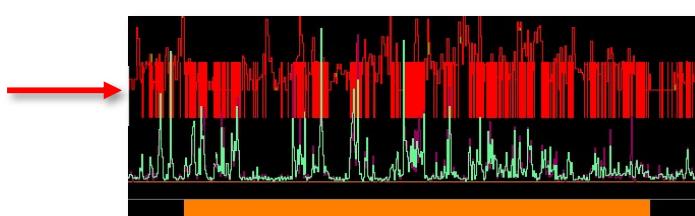


Figure 10.2 Marked events on events where VeDBA < 0.1 g



This expression can be stored in a list by clicking **Store expression**. Before storing, a name can be entered into the box to the right of this button “Name your expression”

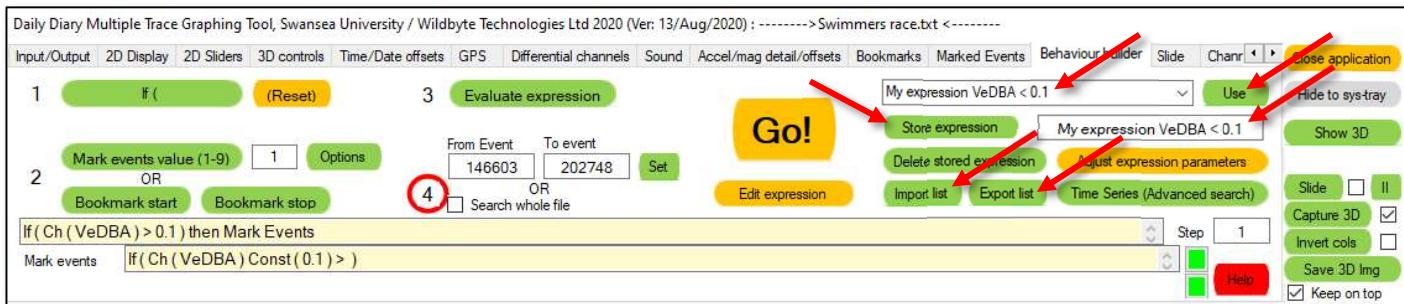


Figure 10.3 Storing this first expression

Multiple expressions of this **Marked Event** type can be stored in the list. There are also buttons for exporting and importing these expressions lists **Export list** and **Import list**. Expressions can also be deleted from the list by clicking **Delete stored expression**.

To recall an expression from the list, select it from the drop-down box, and click **Use**, and its expression and RPN conversion will be loaded into the two boxes at the bottom.

The alternative to the single **Marked Event** expression, as discussed briefly above, is to create two expressions that define the start and stop points for a bookmark. To use this feature, instead of clicking **Mark events value (1-9)** button (step 5 above), click either **Bookmark start** or **Bookmark stop**. Below, Figure 10.5, are two expressions using the start/stop pair, in RPN.



Figure 10.4 Creating an expression pair that will define the start/stop for a Bookmark

Clicking the **Go!** button results in multiple bookmarks where the two expressions are found true. Figure 10.6 illustrates this. Note, the start of the bookmark is where $\text{VeDBA} < 0.1 \text{ g}$, and the stop if where $\text{VeDBA} > 0.1 \text{ g}$.

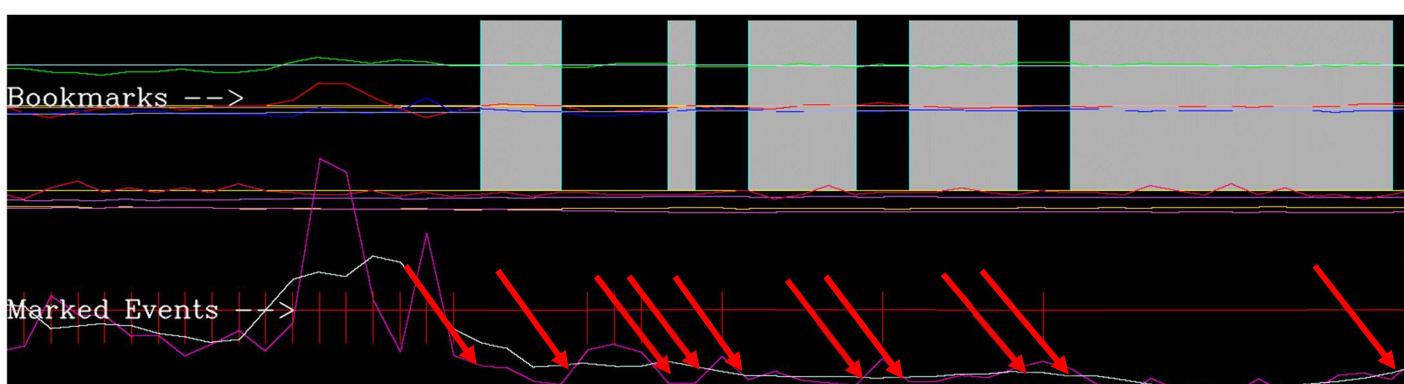


Figure 10.5 Creating an expression pair that will define the start/stop for a Bookmark



Behaviour builder pre-smoothing

Channels may be pre-smoothed on-the-fly if required by ticking the **Pre-smooth? (by:)** checkbox and entering a value in the adjacent box (the averaging window) before clicking the **Enter** button to add that data channel to the expression. See Figure 10.6 for an example where **Accel X** is used in an expression with no pre-smoothing and then the behaviour is instead added with a pre-smoothing function **If (SM(Acc X, 25)...** telling the software to smooth with a centred averaging window of 25 events.

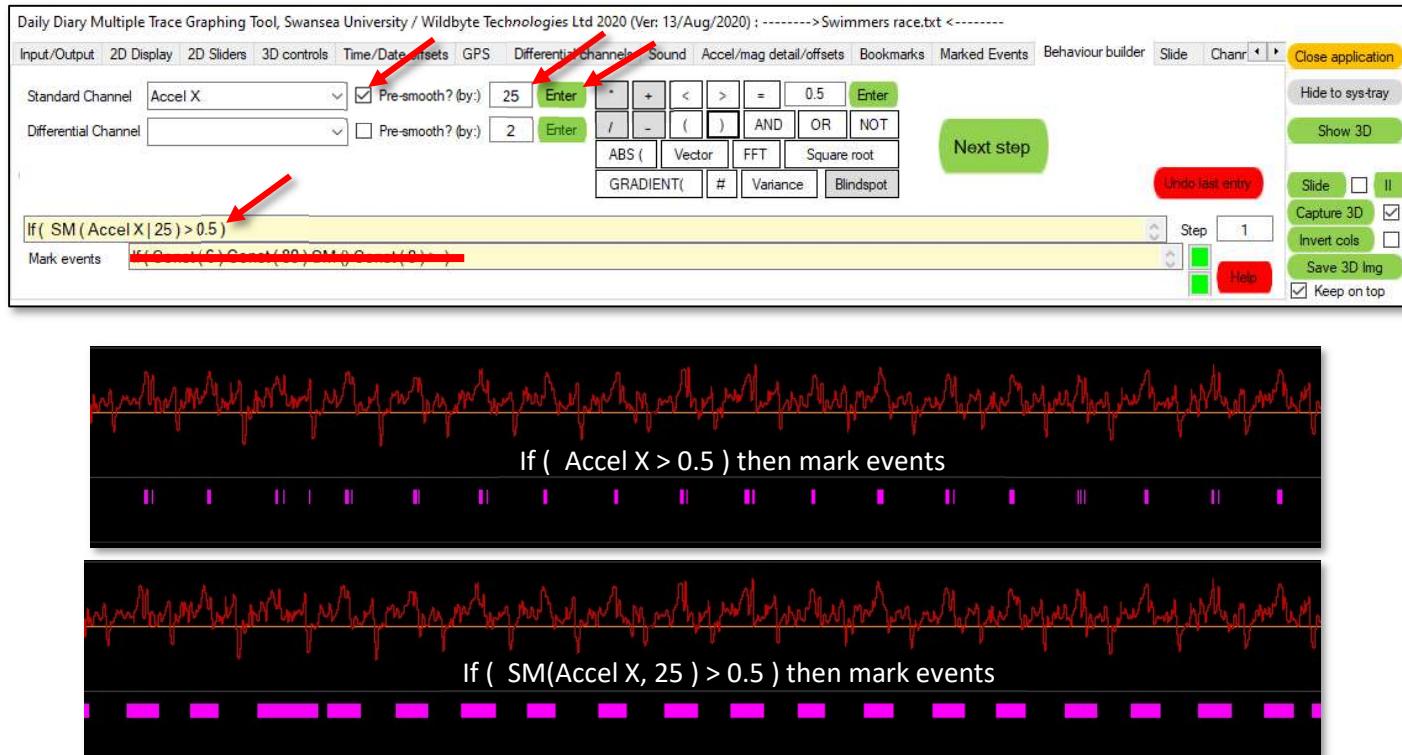


Figure 10.6 Marked events with and without pre-smoothing

Behaviour Builder function - Variance

Variance of any channel can be used as part of an expression. An example might be:

If (Variance(Accel X smoothed | 1000) > 2.5) then Mark Events

By starting a **Behaviour Builder** expression with **If (**, clicking the **Variance** button will show its controls; Figure 10.7.

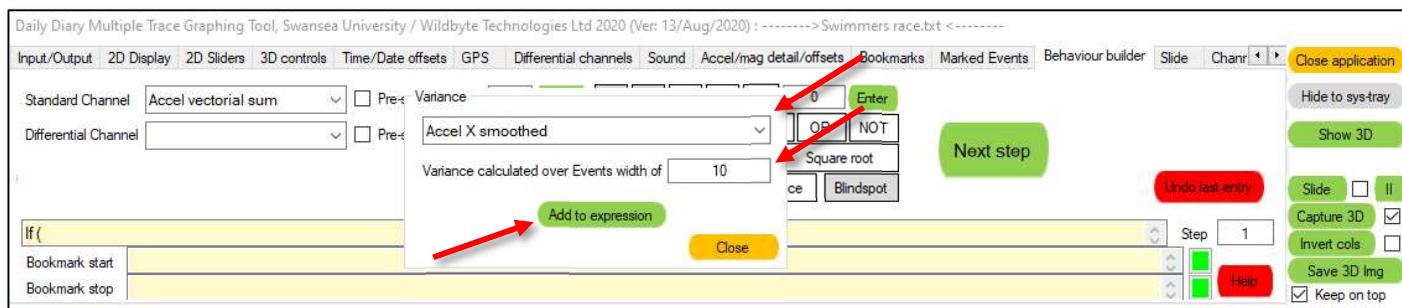


Figure 10.7 Marked events with and without pre-smoothing

Variance parameters, such as over how wide a window to use (per event), and the resulting variance value, can both be predetermined by generating a Variance channel as discussed in Chapter 2 (Metric – Variance controls / Figure 2.28), and looking at the Variance values on the **Values Overlay**.



Gradient function

The **Gradient** function allows one to determine the gradient across X number of any points, for any channel listed (**Standard Channels** only, not **Differential Channels**). To use this, select the **Gradient** button, then select a channel and click the green **Enter** button to push the channel into the behaviour box. Click the # button to add in a vertical line separator. Next, enter a value for the number of events over which to determine the gradient into the numeric box and click the green **Enter** button, then close the bracket with the) button. The user can then add an inequality (< or >) and another value. The **Gradient** function can be used for both **Marked Events** or **Bookmarks**. Note that the gradient is calculated at point n such that it is “*Gradient value(n + step) – Gradient value(n)*”.

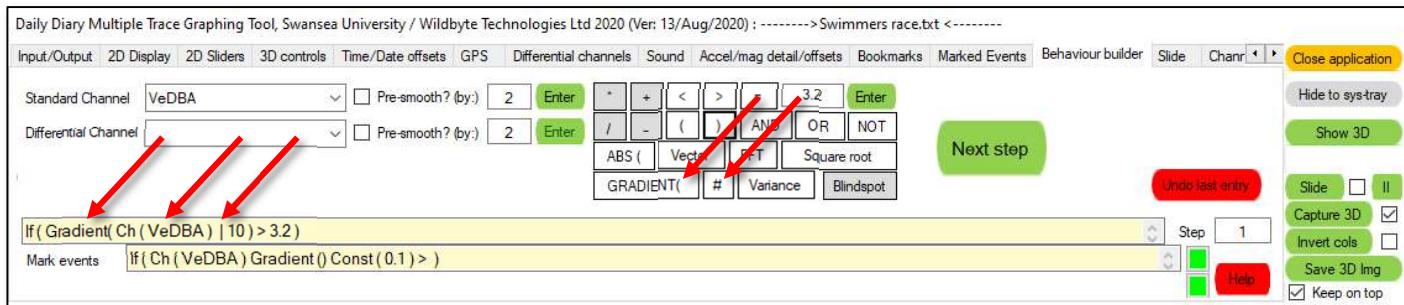


Figure 10.8 Entering a function using **Gradient**

An example expression using the **Gradient** would be:

If (Gradient(Ch (Mag Z) | 10) > 0.22) then mark events

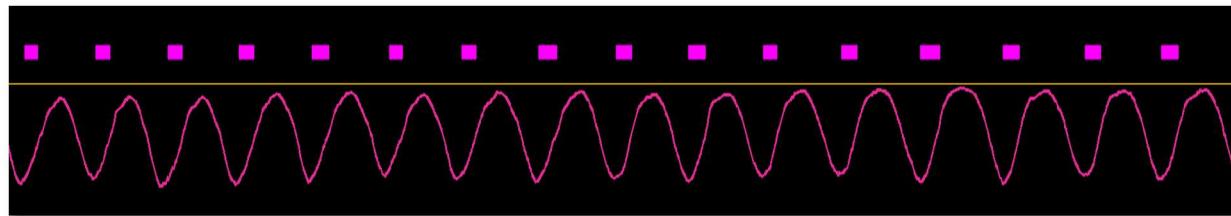


Figure 10.9 Gradient function with marked events used to pick out the positive gradients

Figure 10.9 shows how the **Gradient** function can be used to pick out the positive-going slopes in an oscillatory signal. Here, we get more than one marked event per slope.

Blindspot

The **Blindspot** function can be added to the end of an expression so that DDMT stops looking for any match of the behaviour for a defined number of events after it finds one; essentially stepping over n events. Knowing the approximate period of the signal above, it is possible to set the **Blindspot** range so that there is only one **Marked Event** per region of **Events**. This is useful for functions such as **Gradient** that might generate multiple **Marked Events** where only one is required.

If (Gradient(Ch(Mag Z) | 10) > 0.22) Blindspot(30) then mark events

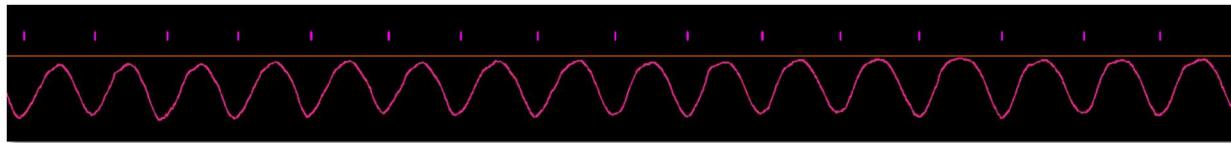


Figure 10.10 Gradient function with marked events used to pick out the positive gradients including **Blindspot**

The positive slope of the waveform in Figure 10.10 is approximately 30 events, so setting the **Blindspot** to 30 events ensures that the next marked event will occur at the next trough. Using the **Blindspot** function, this allows a single marker per oscillation of the waveform. Data exported with the checkbox **Export behaviour builder marked events only** on the **Input/Output / File output/Save** tab will also, by default, export the timing of single events.



Differential channels in the Behaviour Builder

Below the **Standard Channel** selection on the Behaviour Builder control panel, other user generated channels can be found such as the **Differential channels**. These are used in the same way as **Standard Channels**; any differential channels that are created by the user will appear in this list. Pre-smoothing is also available for these channels. Note that differentials can also be smoothed as part of their construction; see Chapter 4.

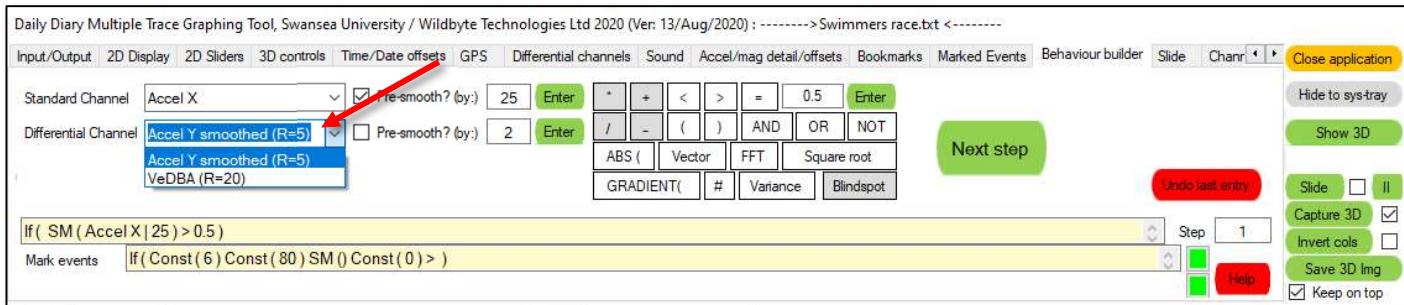


Figure 10.11 Marked events with and without pre-smoothing

Categories

When bookmarking data via the **Behaviour Builder**, **Bookmarks** will be assigned to a category. The currently selected category will be assigned to each bookmark that is created.

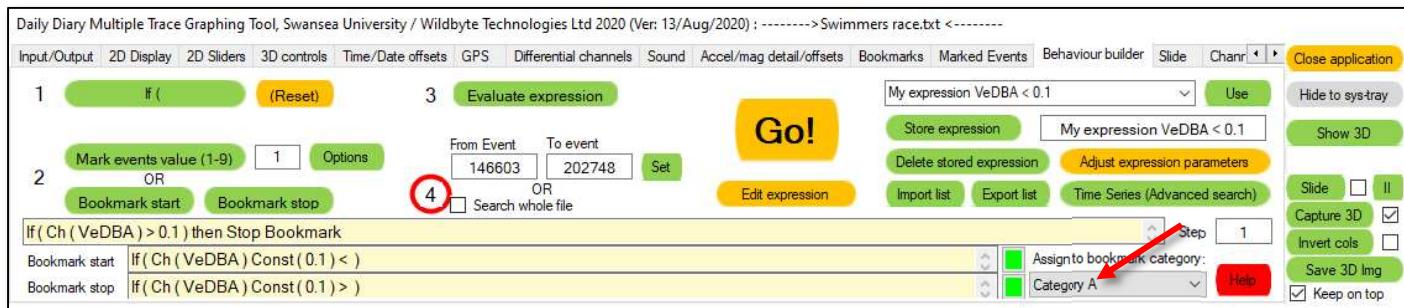


Figure 10.12 Assigning categories to Bookmarks generated by the Behaviour Builder

When category names are updated on the **Bookmarks** tab, this down-down list of categories will also be updated. Default category names are **Category A**, **Category B...**

Summary

In summary, the basic **Behaviour Builder** evaluates all **Events** within the selected range of events, or the whole file, for the currently selected expression and, either marks it as a **Marked Event** or as a start/stop marker to generate a **Bookmark**.

Other functions available are *, /, +, -, Absolute (**ABS**), =, inequalities < and >, logical expressions **AND**, **OR**, and **NOT**, **Square root**, and some special functions including **FFT**, **Gradient** and **Vector**

Time series – Advanced Behaviour Builder

A series of Behaviour Builder expressions can be linked in time to define a multi-step behaviour. This works by defining multiple individual Behavioural Builder expressions, that stipulate levels of VeDBA, or the rate of change of Accel Y for instance. These are then placed in a time series, as a list of elements in a sequence, each having a defined number of events over which its condition must be true, a rule that defines when the next successive element can begin with its own search, and a specified range for



flexibility. This flexibility value means that each element is not fixed rigidly in time, rather they can adjust automatically to the natural variability of animal movement data.

The expressions passed to the Time series part of the **Behaviour Builder** must be of the **Marked Events** type, not the **Bookmark Start/Stop** type.

To use this, first create some equations and store them so they appear within the drop-down boxes shown in Figure 10.13.

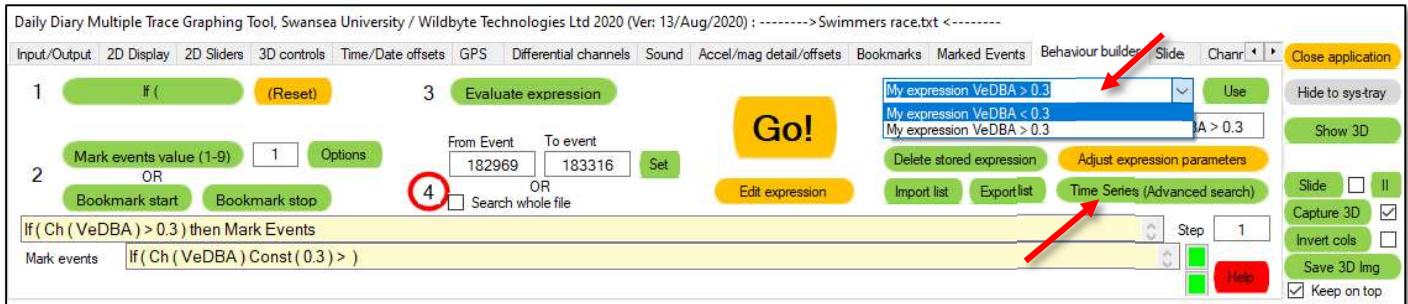


Figure 10.13 Creating and storing behaviours for use within the time series behavioural identification algorithm

Now click the button labelled **Time Series (Advanced search)** to open the time series control panel which will have these same expressions listed (Figure 10.14).

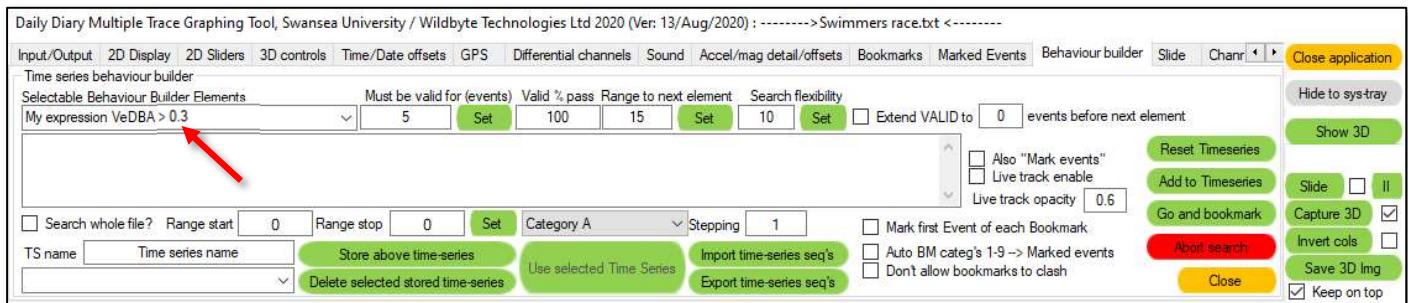


Figure 10.14 Time series control panel; using expressions previously defined and saved within the Behavioural Builder

The aim of this part of the algorithm is to allow previously defined expressions to be placed along a time line at points of “interest”, enveloping a behaviour i.e. a series of changes in the data that express a known pattern(s) or movement(s).



First example of Time series Behaviour builder

Our first simple case below is an oscillation in the smoothed Acceleration Y channel. We wish to identify all the sections where the data is > 0 . The horizontal orange line defines where $Y = 0$.

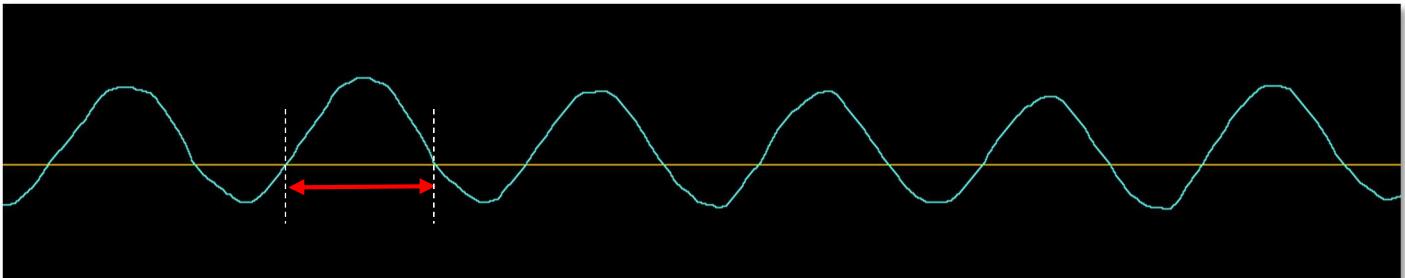


Figure 10.15 Acceleration Y smoothed; horizontal orange line is $Y = 0$

In the **Behaviour Builder**, we create two expressions; the first defines $\text{AccYsm} > 0$ and the second $\text{AccYsm} < 0$, Figure 10.16.



Figure 10.16 Two Behaviour Builder expressions

To identify a full-width peak as in Figure 10.15, we need to identify the crossings. There are several ways, perhaps better than what we'll do here, but the following method will illustrate the simplicity of this advanced search.

What we're going to do is add an expression, which has several parameters associated.

Must be valid (events)

Number of consecutive events that expression n is valid (true) for

Valid % pass (0-100%)

% of the above consecutive events that must be valid for expression n to be valid (true)

Range to next element (events)

Number of events from expression n to expression $n+1$

Search flexibility (events)

Number of events within which expression $n+1$ can be searched for / found to be valid i.e. expression $n+1$ must exist within the flexibility region of expression n else expression $n+1$ will be found invalid

To begin, we select our first expression $\text{AccYsm} > 0$, change **Must be valid** to 1, and **Search flexibility** to 100 (for this particular dataset) and click **Add to Timeseries**. This first **Element** will be added to the central behaviour definition box.

Next, we select our second expression $\text{AccYsm} < 0$, and then just click **Add to Timeseries**; the **Range to next element** and **Search flexibility** parameters will be ignored for the 2nd **Element** as there is no 3rd **Element** to which they will be applied. See Figure 10.17.



Expressions may be used more than once in a behaviour definition.

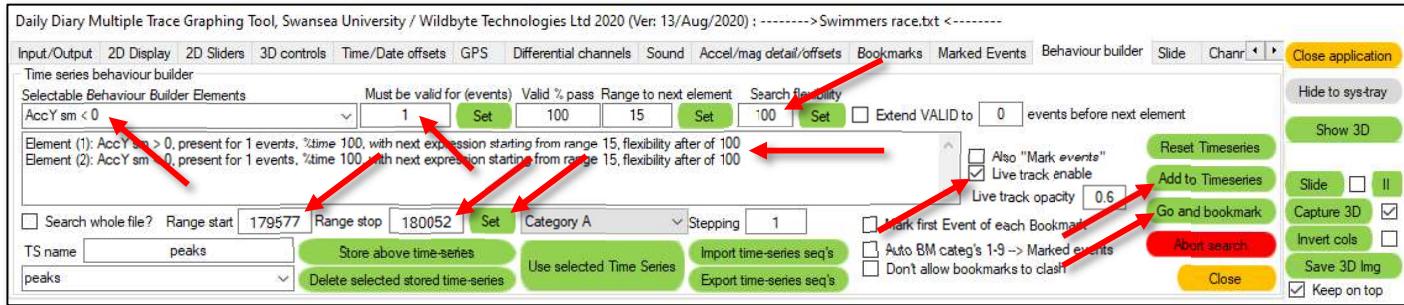


Figure 10.17 Element parameters when added to the **Behaviour definition**

Highlight a region covering the data to search with the mouse's left/right white lines, and click **Set** and those start/stop **Event** values will be entered into the **Range start** and **Range stop** boxes.

Now, to use this **Behaviour**, we must first store it and tell the software we wish to use it. First assign a name to the **Behaviour**; in Figure 10.17, we called this **peaks**. Clicking **Store above time-series** adds this **Behaviour** to a buffer. Finally, click **Use selected Time Series** and the data within the search range just defined will be validated for the defined **Elements**. Clicking **Live track enable** tells the software to validate the **Behaviour** for wherever the user moves the **Elements**.

Live tracking is a feature that allows the user to "test" a section of data and perfect the **Behaviour** definition.

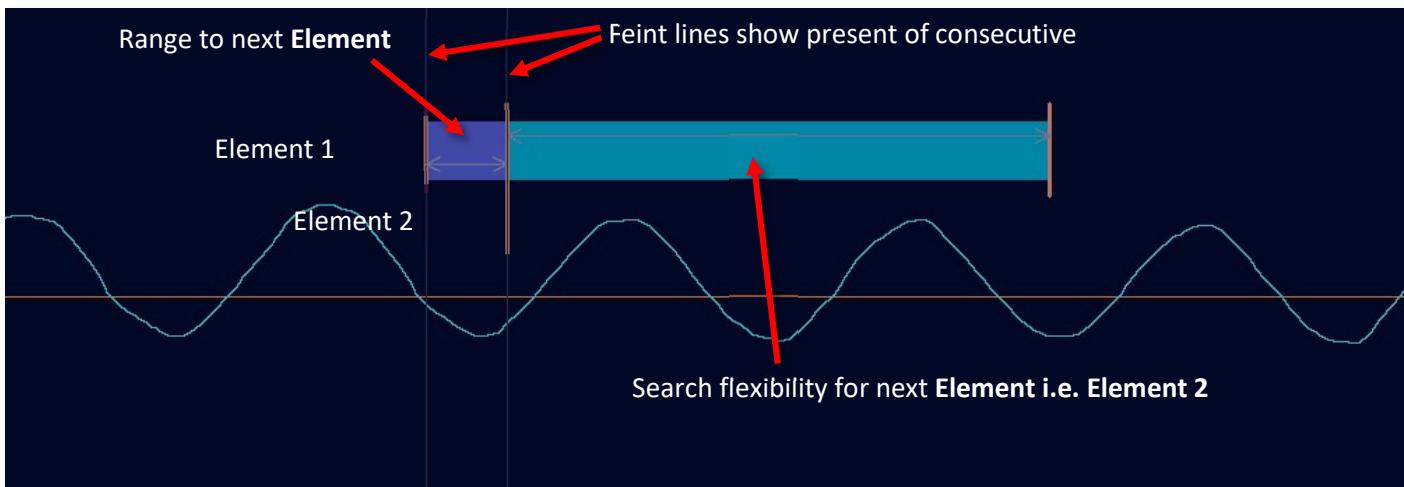


Figure 10.18 Layout of **Elements** on the **Live tracking** graphic

Once **Live track enable** is ticked, the user can hold down **left ctrl** and **left shift** and left-click the mouse and drag the graphic shown in Figure 10.18 through the data.

By default, when **Live tracking** is enabled, all **Elements** will move to sit under last invalid **Element**. Once the first **Element** finds a point where it is valid, its **Element** line will turn green and the subsequent **Element** (in our example here, **Element 2**) will move to the start of the **Search flexibility** region (i.e. the end of the previous **Element's Range to next element** point). In Figure 10.19, the second **Element** moved to the start of the **Search flexibility point** of the first **Element**, and then searched for the first point where it too is valid i.e. $\text{AccYsm} < 0$. Once it was valid, the system determined that all (2) **Elements** were valid and so the **Behaviour** was valid, and placed a horizontal green line spanning the first to last **Elements**. Now, if a search (as opposed to **Live tracking** – a test feature) were performed, a **Bookmark** would have been created here.

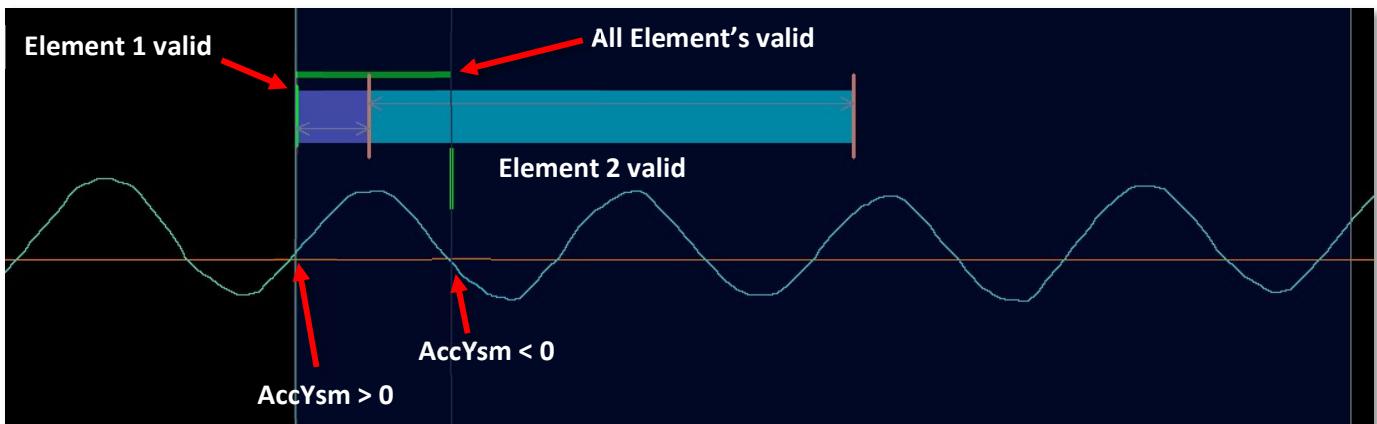


Figure 10.19 Elements 1 and 2 valid. Horizontal green band along the top spans the full distance between the first and last Elements that will define the width of the **Bookmark** that will be created in a search

To test the actual search function of the Time series, click **Go and bookmark**. There will be a pause while DDMT validates the **Elements** for every event currently in memory. The search is performed immediately after.

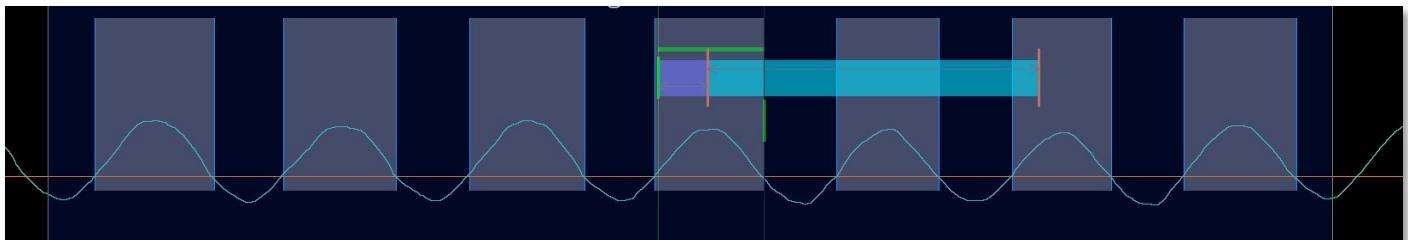


Figure 10.20 Bookmarks created in the defined search range

Second example of Time series Behaviour builder

As another example, here is some pressure data from a diving animal. The pressure has been inverted such that greater pressures are at the bottom of the graph.

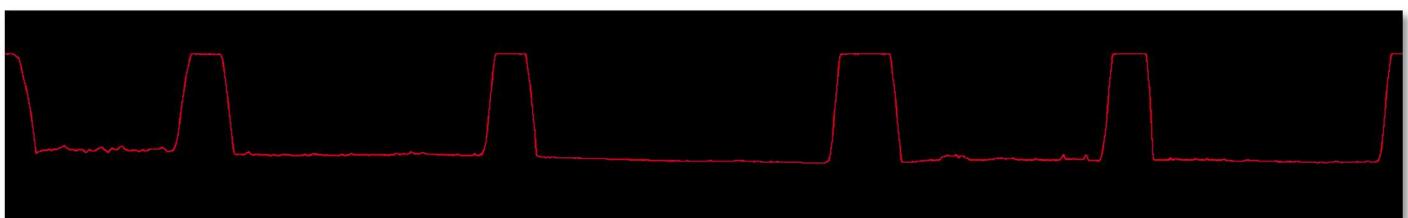


Figure 10.21 Second data example of diving behaviour

Here, we would like to **Bookmark** each of the dives, and subsequently use the power of the **Bookmarking** functionality of DDMT to export all the stats of each individual dive.

First, we define expressions for the low pressure and high pressure.

If (Ch (Pressure) > 0.15) then mark events

If (Ch (Pressure) < 0.15) then mark events

We measure the number of **Events** to span the dive here as ~64,000 events. So we'll set the **Search flexibility** of the first **Element** to 100,000 events.

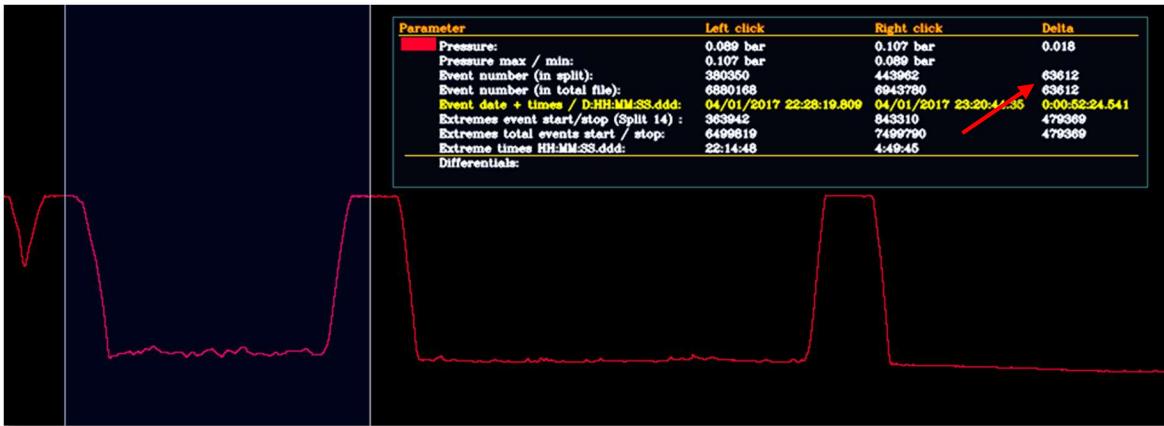


Figure 10.23 The span of one of the dives is approximately 64,000 events, and so the first **Element's Search flexibility** needs to be at least this for the 2nd **Element** to fit

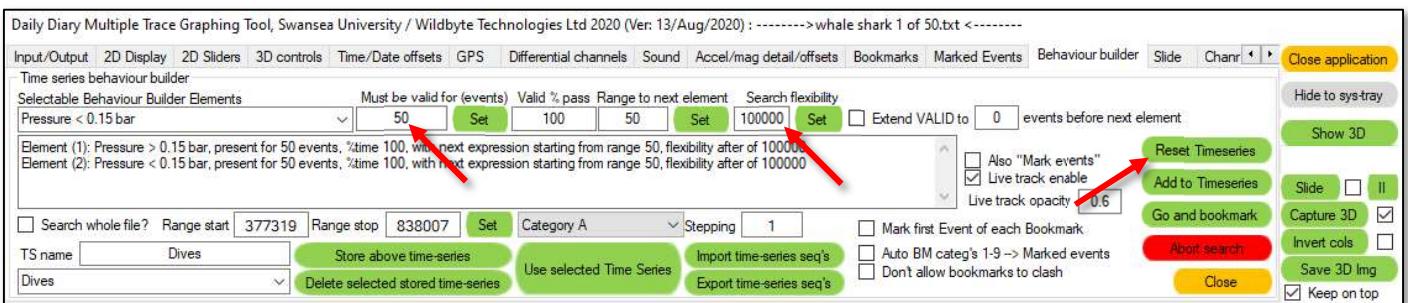


Figure 10.23 The two **Elements** with **Must be valid** values of 50 to negate a hit on any small fluctuations i.e. ensure it is a real dive

Saving these two **Elements** as "Dives", setting the data range to cover a few of these dives, clicking **Use selected Time Series** ("Dives"), and clicking **Live track enable** we get Figure 10.24.

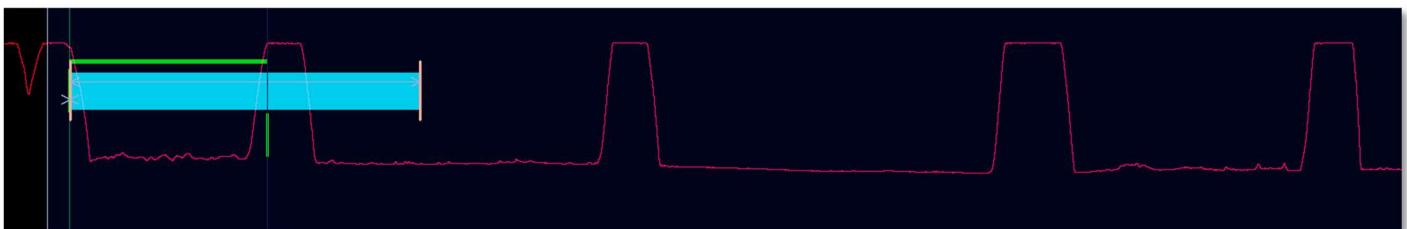


Figure 10.24 The two **Elements** both pass and so the overall **Behaviour** is valid – horizontal green line

Moving the live tracker further along, the longer dives are too wide for the flexibility. To save resetting the **Behaviour** by clicking **Reset Timeseries**, we can dynamically adjust the **Search flexibility** by moving mouse over the first **Element's** vertical bar – a label will appear showing the flex, range, and valid values; Figure 10.25.

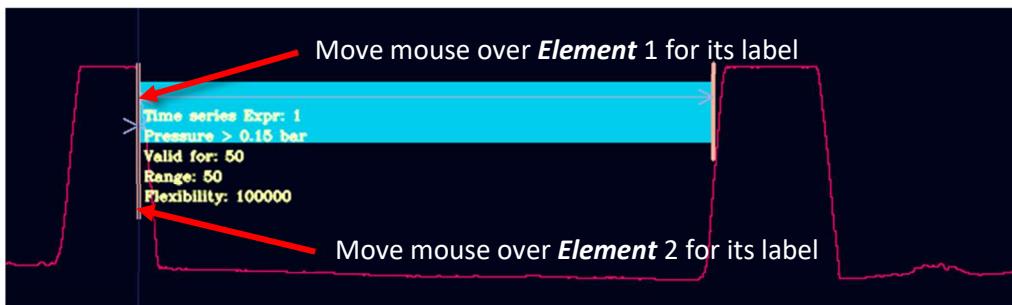


Figure 10.25 Element values labels. In the above, the **flex** value is not large enough for it to span the gap between the dive down and surface aspects of the “behaviour”

When the **Element**’s label is visible, changes made as:

Left shift + Left ctrl + Left alt	Mousewheel to adjust Range to next element
Left shift + Left alt	Mousewheel to adjust Search flexibility
Left shift + Left ctrl	Mousewheel to adjust Must be valid for

With each change, the **Live tracker** will dynamically check the validity of all **Elements** both visually on the **2D graphing window**, and also numerically on the **Behaviour definition** window.

Note that the mouse-wheel only adjusts these values by +/- 1, so if these values are drastically different from that required, one would need to reset the time series expression window by clicking the **Reset** button and redefining from **Element 1**.

Adjusting these values only affects the **Live track** behaviour i.e. the **Valid**, **Flex**, and **Range** values in the **Behaviour** definition window (where the n elements are listed along with their individual parameter values). The updated time-series definition must first be saved with the **Store above time-series** button, before clicking the **Go and bookmark** button before it can be used in a search. The **Live track** function allows the user to make adjustments and determine optimal limits for behaviour definition(s).

Time series definitions can be exported by clicking **Export time-series seq’s** which will create a .csv file. Opening this will reveal a simple format that allows the user to combine multiple exports for the purpose of creating a library of time-series definitions. Within the .csv file, the first line contains a value defining the number of time-series definitions within the file. If multiple files are combined, simply count the number of lines (time series definitions) and enter this number at the top line to tell the software how many it should expect within the file.

Extend VALID to n events before next element

If **Extend VALID to n events before next element** is ticked prior to clicking the **Add to Timeseries** button, this will require that **Element n** is valid from after its initial **Valid** window to the start of where the next **Element** is valid minus n **events**. Note that the **Must be valid for** events has its own **Valid % pass** validation i.e. 80% pass rate, while the **Extend VALID to n events before next element** is validated separately from this i.e. 80% pass rate. Do not tick this and add it with the final **Element** in a **Behaviour** definition.

Additional Time series functionality

- Ticking **Also “Mark events”** will generate **Marked Events** throughout any **Bookmarks** created during a search.
- To search an entire data file (all splits), tick the box **Search whole file ?** and the search algorithm will search from **Split 1** to **Split 2**.
- If a search is taking a long time to complete for a **Split**, click **Abort search** to step out of the search algorithm



11. 3D Data Visualisation

Enable 3D, and syncing with changes in data view on 2D graph window

To view data in 3D, simply click the **Show 3D** button as indicated. The control panel will be extended with the new window appearing beneath. The single grey, unlabelled, button toggles the 3D drawing window between the old and the new system. The new 3D system is enabled by default.

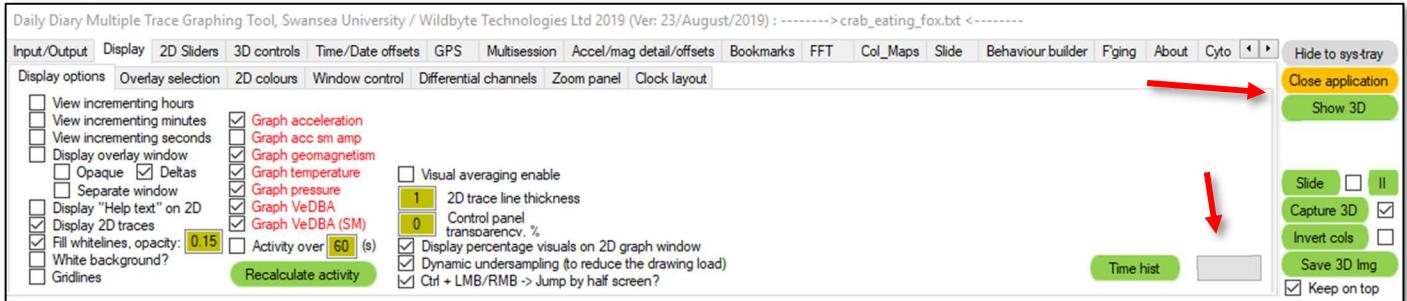


Figure 11.1 3D display options

Moving the mouse over the new area will show a series of buttons in the top left corner.

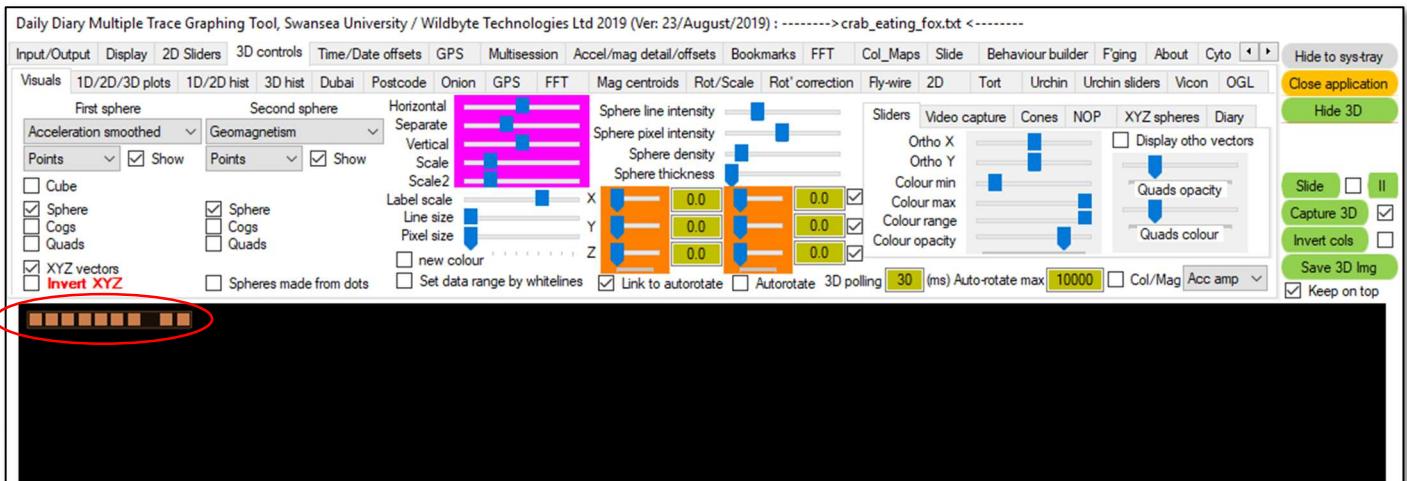


Figure 11.2 Initial 3D display

These buttons, from left to right are (more will be added over time):

1. **Toggle 3D full window** – full view will extend the 3D window vertically up to cover the control panel. Moving the mouse out of the window, or to a border will force the control panel to reappear.
2. **Toggle colour** – provides a choice between five colour schemes, with orange as the default.
3. **Colour map Generator** – discussed later.



4. **Show 3D help** brings up help text outlining the mouse controls.

Right-click on background to obtain vis-menu
Right-click on Data (reset orientation) button for 3 sub-buttons
Right-click on top sub-button for more menu items
To autorotate a vis, left Alt and slow-drag the vis
To cancel an autorotate, right-click on its centre
To move a visual colourbar around, left Ctrl and drag it
To move a manual-entry values box, left Shift and drag it
Hover mouse over a button to show a graphic of its mouse interactions
(some buttons are a combination of left, right, and/or mousewheel)
To unhook from a vis, left-click then right-click on backdrop
Unhook from current vis to interact with the colourbar generator
Visuals can be rotation-linked by setting them to the same vis channel
(bottom menu button to enable, mousewheel to select channel)
Hold down left Alt to rotate Plan-view with dead-reckoning visualisations
Hold down CTRL and SHIFT and rotate mouse wheel to adjust vis axes font size

Figure 11.3 Mouse controls in the 3D window

5. **Toggle background** alters between honeycomb/patchwork/no background. When a background is present, the mouse-wheel, while situated over this button, can be used to modulate the contrast of the background colour.
6. **Capture screen** – once clicked, will cause the 3D window to be captured to memory once the mouse is moved off the window area; it waits for the mouse to leave so that the visualisation's menu is removed before capture. Note that this image capture will follow the configuration of the tick-boxes on the RHS of the control panel. The checkbox adjacent to the **Capture 3D** button tells DDMT to immediately save the 3D screen capture to the user's working directory. The checkbox adjacent to the **Invert cols** button will request the image's RGB colours be inverted automatically upon capture. When used with the direct save above, this will result in an inverted image being saved to the user's working directory. The **Save 3D Img** button is used to save the captured image if auto-save is not selected. This allows the user to continuously attempt to get the correct image, and only then save the required image.

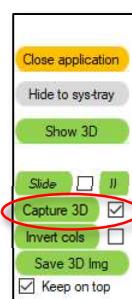


Figure 11.4 Create screen captures of the 3D window and save them to the working directory

7. **Initiate magnetometer correction algorithm** – this will take the user through correcting the magnetometer data. This is outlined in Chapter 4 and written instructions are available on screen.
8. **Create Smoothed Accel sphere** – brings a spherical acceleration plot onto the screen. This sphere corresponds to the data currently on view in the **2D graphing window**.
9. **Create Smoothed Mag sphere** – brings a spherical magnetometry plot onto the screen. This sphere also corresponds to the data currently on view in the **2D graphing window**.

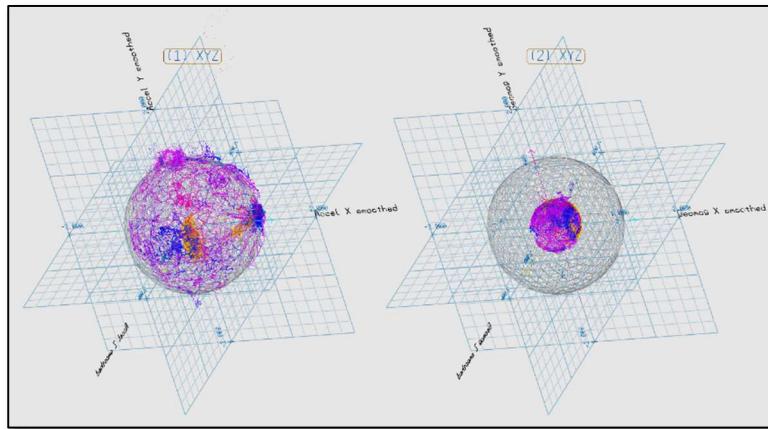


Figure 11.5 Smoothed Acceleration (left) and Magnetometry (right) spheres

Creating a visualisation

Right-clicking on the 3D window will show a menu from which one of many visualisation types can be selected:

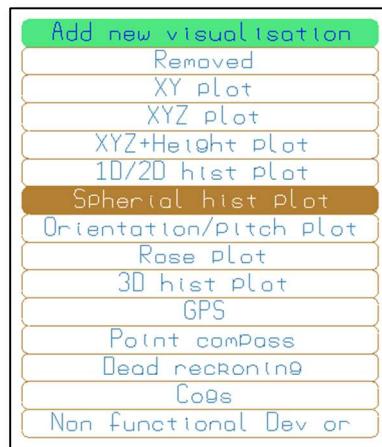


Figure 11.6 Visualisation selection menu

- XY plot** Plot of two variables, X vs Y.
- XYZ plot** Plot of three variables, X vs Y vs Z.
- XYZ+Height plot** Plot of acceleration or geomagnetism, normalised, and then modulated by a 4th channel such as pressure or temperature.
- 1D/2D hist plot** Generates either a 1D or 2D (if a second axis is selected) histogram of channel data.
- Spherical hist plot** A spherical histogram plot of data, using either acceleration or geomagnetism to determine from where the point data is collated from.
- Orientation/pitch plot** Plot (O-sphere) of data points generated by heading being the circumference, with pitch (up to +/- 90 degrees) varying the point between the two poles.
- Rose plot** Display a rose plot of the data, i.e. a histogram of heading, segmented to 5 degrees resolution.
- 3D hist plot** Generates a 3D histogram from 3 selected channels.
- GPS** Plots the GPS track currently loaded into memory.

**Point compass**

Creates a red arrow with a compass that aligns itself to the currently calculated heading value according to the position of the left white line in the **2D graphing window**. Tilting according to pitch is also included.

Dead reckoning #1

First attempt at a dead-reckoning algorithm whereby a track is generated according to levels of VeDBA above a user-determined threshold.

Visualisation controls

Most of the visualisations can be rotated using the mouse via what is known as the “arc-ball” rotation method, i.e. as though there were a ball on the screen and you grab the surface with the left mouse-button and drag it around. Some visualisations, such as the 1D/2D histogram, GPS, Dead reckoning, etc., contain data that do not make sense when fully inverted by rotation and therefore use a different rotation algorithm. In this case, the mouse can be used to rotate around the vertical (Y) axis, and also pitch the visual up/down. Holding down the left Alt key lets the user rotate the visual on the Z axis (the axis coming out of the screen). The user will find this makes a lot more sense for these visuals. Note that a visualisation must be selected in order to rotate it at all. Select or deselect a visualisation by clicking the centre square with a white cross (green when selected, red when deselected).

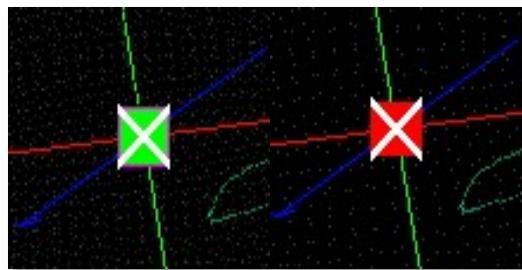


Figure 11.7 Centre square when deselected (left) and selected (right)

The white cross indicates that the data currently being displayed on the 3D visualisation will automatically update when the user moves to a different section of data on the **2D graphing window**. A notice will appear at the top left of the viewing area to inform the user that data tracking is on (as opposed to the data being locked). Moving the mouse over this red square will expand the notice information to provide a more complete explanation.



Figure 11.8 Data tracking notice

Position a visual by either left-clicking and dragging the centre square, or holding down shift and left clicking to force the currently selected visual to move to the mouse position. Right-click on the viewing area to deselect the visual. Once deselected, right-clicking again on the area will show the **Add new visualisation** menu.

When a visualisation is created, it becomes the focus. This means that any left-click and drag of the mouse will cause the visualisation to rotate, while the mouse-wheel will change the scale/zoom.

All visualisations, when selected, will show a menu of buttons when the mouse is within the viewing area. Moving the mouse over the buttons will provide a written description of the button’s action.

The top two buttons, and the buttons highlighted in red below are identical for all visualisations, while for the remaining buttons, if present, will vary in function with visualisation type.

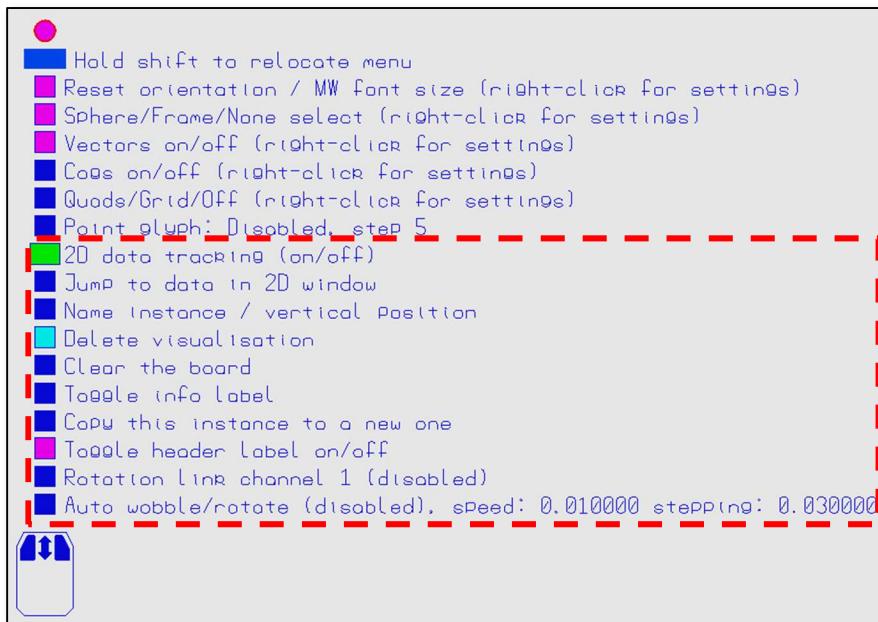


Figure 11.9 Visualisation menu controls

The top button in the menu (coloured orange to make it stand out – blue in the inverted image above!) can be used to adjust the location of the menu within the 3D window (sometimes it gets in the way or is not entirely visible due to screen dimensions); hold down the left shift, and “then” left click the mouse and drag the menu where required. The second button down has 3 mouse-interactions, as indicated by the mouse image at the bottom (left, right, and mouse-wheel). The left mouse button will reset the rotational orientation of the visualisation, the mouse-wheel (MW) will adjust the font size of the text (sometimes useful on larger screens), while the right mouse button will produce another column of buttons. More on this later. The circular button at the top of the menu structure toggles all button labels below it, else the button name will appear when the mouse passes over it.

Some visualisations will have an associated “ortho-viewer” appear than can be dragged into any location on the 3D area, which best shows the 3D orientation of the current visualisation. Right clicking on this also cycles the rotation of the visualisation through a number of set orientations, including X-Y, and X-Z etc.

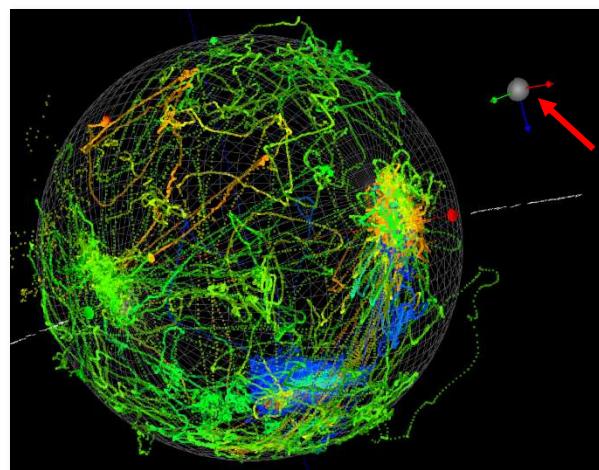


Figure 11.10 Ortho-viewer for visualisation orientation, and also for auto-switching the rotation through one of a number of preset orientations

The other buttons, highlighted with the dashed red line, have a variety of useful function as described below, and are common to all visualisations.

**Data tracking**

When a visualisation is created, this is enabled by default (pink to emphasize it is enabled). When data tracking is enabled, any change to the viewed data in the 2D viewing area (or perhaps changes in smoothing level of various data channels) will update this data-tracking enabled visual. The visual doesn't need to be selected at the time for it to be updated by a change in the viewing on the 2D window. When data tracking is disabled, the visual is data-frozen; moving to another split will not update the visual with the newly loaded data channels.

Jump to data in 2D window

If the currently viewed split contains the data of the selected visual, then the 2D window will jump to / zoom in on the data used to create it. Only to be used when **Data tracking** is enabled for the visualisation in use

Name visualisation instance

This lets the user set the name of the instance (visual).

Delete visualisation

Does what it says on the tin. Note there is currently no confirmation so do not accidentally click it. Coloured red to emphasize its use/importance.

Clear the board

This will scale and position all the currently generated visuals to the top of the screen.

Toggle info label

Moving the mouse over the centre square of a visual will show a small info label detailing the split origin of the visual, its start/stop data points, visual type etc. Enabling the info label will make this show permanently, for this selected visual.

Copy this instance to a new one

Will simply create a copy. Useful if you wish to perhaps show the same data sets, coloured by two different properties such as temperature, and pressure.

Toggle visual label

This will toggle the presence of the name of the visual above it. The mouse-wheel can be used to adjust the vertical position of the name label.

Rotation link / channel

All visualisations can be linked via a rotation "channel" (1-9). When enabled, any other visualisations that are enabled (and on the same rotation "channel") will rotate at the same time. Note that clicking the orientation reset button will also orientation-reset any other rotational linked visuals. – useful if the user would like to match rotation orientations. Two or more visualisations do not need to begin with the same rotation offset, one might be rotated by 40 deg in the X axis, and 27 degrees in the Z axis, but will still move the same relative amount indicated by the user's mouse

Switch visualisation type

Will reassign the visual to a different type, while retaining the same data start/stop limits



With the **Data reset orientation** button at the top of the menu (2nd button down), right-clicking on this with the mouse will show a further 3 (or 4) buttons to the left of it/ The reset button will be outlined in purple to emphasize that these 3 new buttons relate to the reset button as can be seen below (just below the orange rectangular “menu move” button):

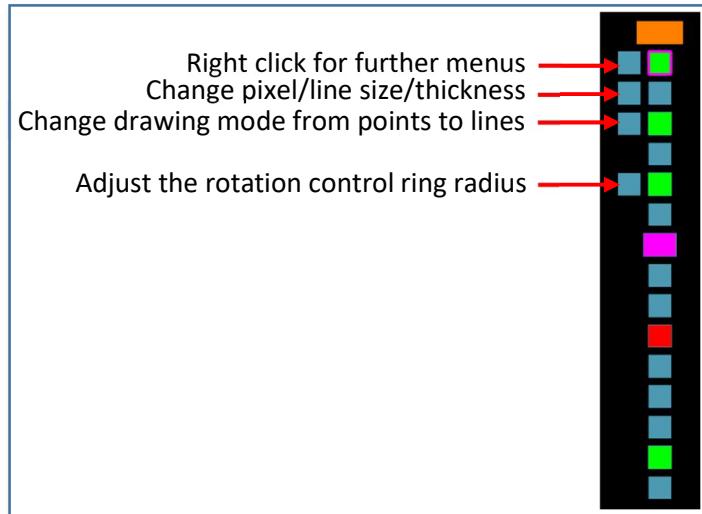


Figure 11.11 First extended menu

The second button from the top (2nd column from right) will adjust the point or line thickness (depending on drawing type) of the data in the visualisation; adjust with the mouse-wheel.

The third button toggles the drawing mode between points and lines.

The top is right-click only, providing access to the remainder of the controls for each visual, the button-content of which will be visualisation-type dependent. Additional round buttons appear above each column of controls to allow all button labels to be simultaneously displayed. Some visualisations will have 3, 4, or maybe even 5 columns in total, depending on the complexity of the algorithm. Enabling one round button for button labels will toggle off all other round buttons.

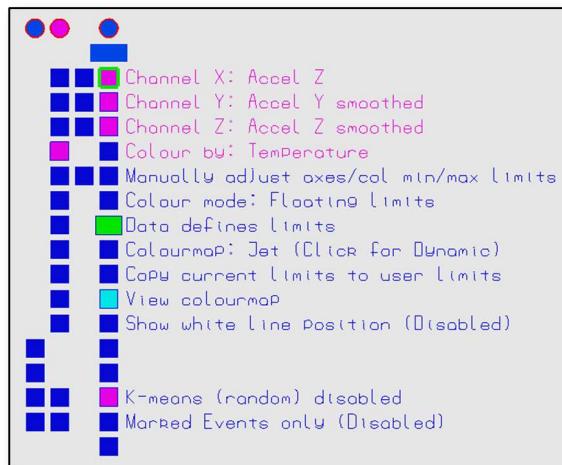


Figure 11.12 Second extended menu

The 3rd column of buttons (from the right) above allows selection of channels for the various axes and colour schemes of the visualisations, and then more complex adjustments depending on the visualisation type.



Several of the visualisations have some other common controls, with buttons placed in the first (RHS) column:

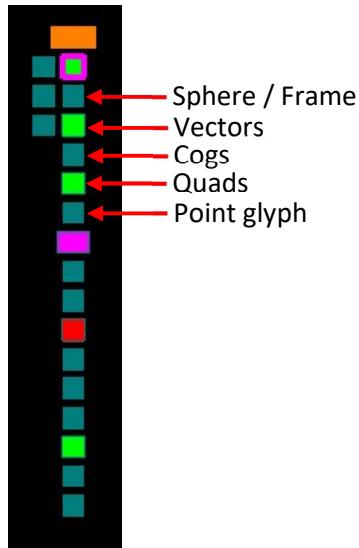


Figure 11.13 Customisable visualisation options

Sphere/Frame

This will place a sphere of points, or a sphere of lines, or a box frame around the visual; left mouse click cycles from one mode to the next. Right mouse-click will shift the above purple box to over the *Sphere* button and the 4 buttons in the 2nd column will then be attributes for sphere (Intensity, Density, dot/line thickness, and colour) -use the mouse-wheel to adjust these attributes.

Vectors

This button places the X, Y, and Z vectors onto the graphic along with labels indicating which channel data is responsible for each axis. Right mouse-click will shift the above purple box to over the *Vectors* button and the three buttons in the 2nd column will then be attributes for vectors (line thickness, invert colours, and vectors magnitude) -use the mouse-wheel to adjust these attributes.

Cogs

This button places 3 orthogonal intersecting cogs over the visual. Right mouse-click will shift the above purple box to over the *Cogs* button and the three buttons in the 2nd column will then be attributes for cogs (opacity, divisions, and cog thickness) -use the mouse-wheel to adjust these attributes.

Quads

This button will place 3 orthogonal quad planes on the visual; left mouse button cycles between simple planes, and gridlines with numbered scales. Right mouse-click will shift the above purple box to over the *Quads* button and the three buttons in the 2nd column will then be attributes for quads (opacity, colour, and size) -use the mouse-wheel to adjust these attributes.

Point glyph

This button will place small glyphs placed at every X points of data. Their size is fixed, but the value of X can be adjusted with the mouse-wheel. The point glyph will assume the colour of the data point it sits over.



Visualisation controls – by visualisation “type”

XY plot

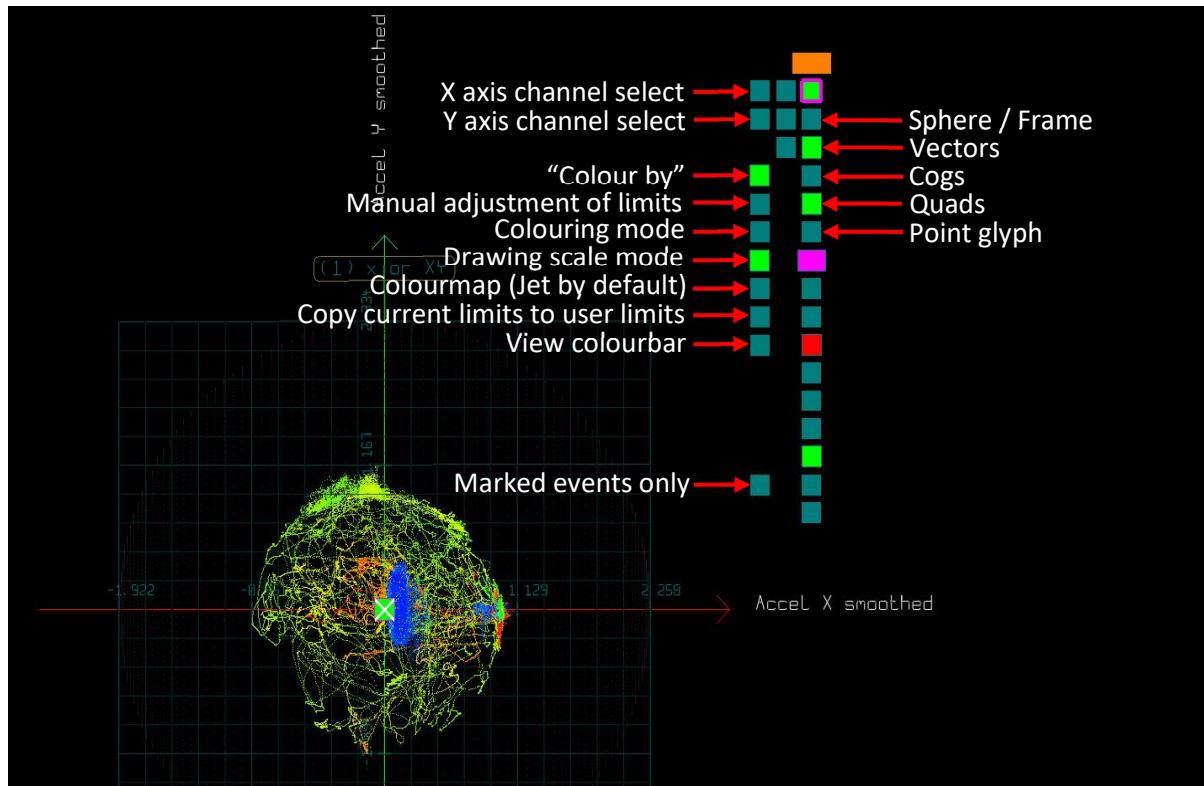


Figure 11.14 XY plot data selection and visualisation options

X Axis channel select

Use the mouse-wheel to select the appropriate data channel for the x axis.

Y Axis channel select

Use the mouse-wheel to select the appropriate data channel for the y axis.

Colour by select

Use the mouse-wheel to select the appropriate data channel to colour the data points. Left-clicking this button switches from colouring by channel data, to a single colour; use the mouse-wheel to adjust the actual colour.

Manual adjustment of limits

This manifests a control panel whereby the user can adjust the X/Y limits of the above grid scale. Note that it is the halfway points that will be defined:

20/30 min/max limits			
Extremes range X	-0. 961	1. 5256	
Extremes range Y	-1. 101	3. 9663	
Not used	-1. 000	1. 0000	
Extremes range Col	23. 600	38. 200	
Set range X	-1. 000	[1. 0000]	
Set range Y	-1. 000	[1. 0000]	
Not used	-1. 000	[1. 0000]	
Set Col range	-1. 000	[1. 0000]	

Figure 12.14 Control panel allowing X and Y channel limits to be defined

Hold down the shift key and mouse click and drag the above window around to better position it for viewing the values. The top few rows of data indicate the currently selected channels’ min/max values. Click on any of the ringed values to set the limits for the halfway +/- values for each axis. A numeric touchpad will appear. Enter the



required value and then click anywhere else on the viewing area to accept this value and the graphic will update (if **Drawing scale mode** is set to **Axes: User limits**). For example, to double the size of the drawing set the min/max values for the two axes to +/- 0.5, rather than the default +/- 1.0. As the data, in the above graphic, is approximately +/- 1.0, it will double the size of the sphere. Note that this values panel is used in other visuals so some of the rows are “Not used”. In the above example, the graphic is coloured by **Temperature**, and the values in the panel are min 23.6C and max 38.2C, which, by the “Jet” colour scheme, is blue through to red.

Click the **Manual adjustments of limits** button again to hide this values panel.

Colouring mode

When **Colour by select** is using channel data to colour the data points, it can either use the colour limits as defined by the user in the above **Manual adjustment of limits** for the min/max colour values, or use floating limits i.e. scale the colour to the min/max values in the currently selected data.

Drawing scale mode

There are 3 drawing scale modes to display the data:

1. Auto-normalise/centralise on the origin.
2. By the user-defined limits through the **Manual adjustment of limits** control above.
3. Data defines the limits; centre 0.0,0.0, with grid max at +/- 2.0.

Colourmap

Select the colourmap to use via the mouse-wheel (only Jet available by default). Left mouse-click to switch to dynamic colouring. The graphic will update in real-time while adjusting the colourmap in the colourmap generator app (3rd button from the left on the top left of the viewing area).

Copy current limits to user limits

This copies the values for the X and Y min/max extremes to the **Set range X/Y** user limits boxes:

2D/3D min/max Limits		
Extremes range X	-0. 961	1. 1291
Extremes range Y	-1. 094	1. 1674
Not used	-1. 000	1. 0000
Extremes range Col	23. 600	35. 700
Set range X	-0. 961	1. 1291
Set range Y	-1. 094	1. 1674
Not used	-1. 000	1. 0000
Set Col range	-1. 000	1. 0000

Figure 12.15 Min/max limits control panel

View colourbar

This toggles the view of the colourbar, showing the **colour by** channel name, and min/max limits currently being applied.

Marked events only

By default, this is disabled (value 0). Use the mouse-wheel to cycle this through values 1-9 in order to filter by data that has been marked with events values 1-9 either manually, or via the Behaviour Builder facility.



XYZ plot

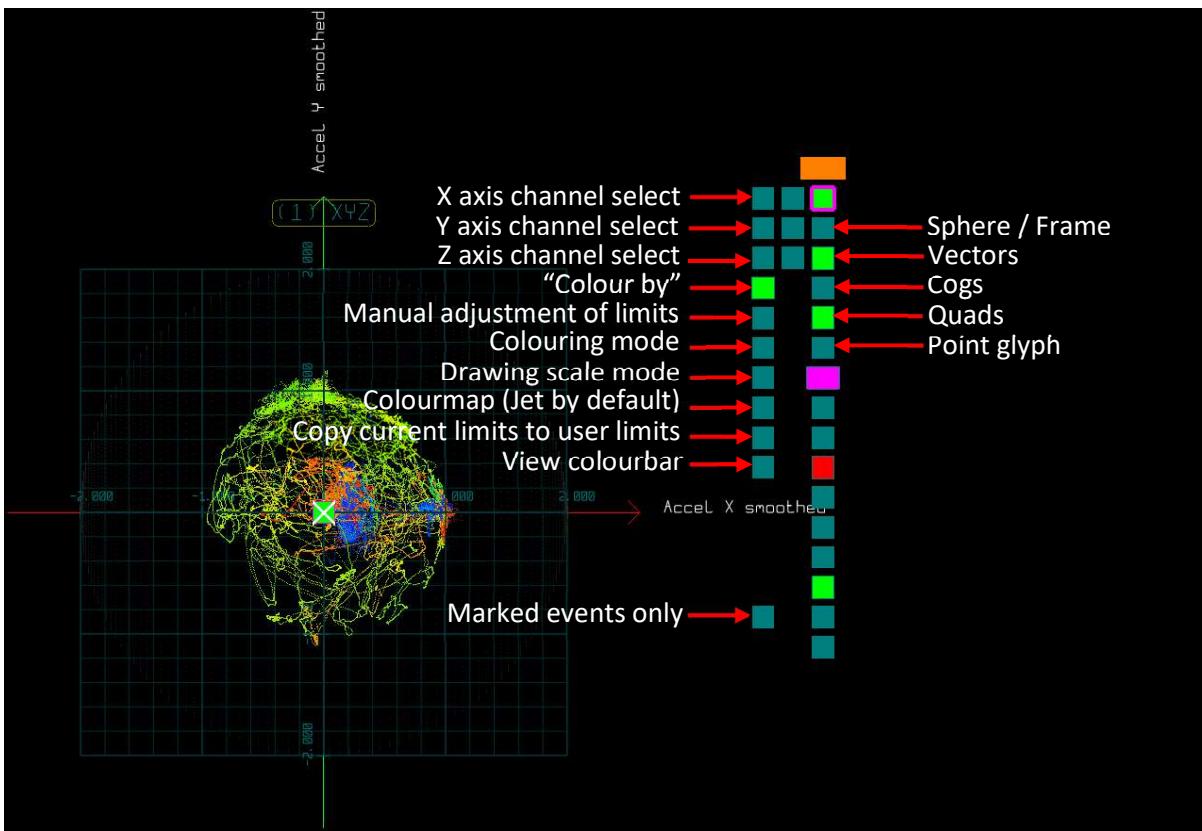


Figure 12.16 XYZ plot data selection and visualisation options

X Axis channel select

Use the mouse-wheel to select the appropriate data channel for the x axis.

Y Axis channel select

Use the mouse-wheel to select the appropriate data channel for the y axis.

Z Axis channel select

Use the mouse-wheel to select the appropriate data channel for the z axis.

Colour by select

Use the mouse-wheel to select the appropriate data channel to colour the data points. Left-clicking this button switches from colouring by channel data, to a single colour; use the mouse-wheel to adjust the actual colour.

Manual adjustment of limits

This manifests a control panel whereby the user can adjust the X/Y limits of the above grid scale. Note that it is the halfway points that will be defined:

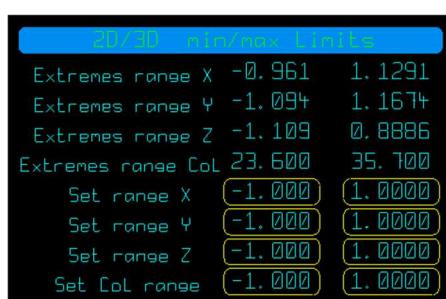


Figure 12.17 Control panel allowing channel limits to be defined



Hold down the shift key and mouse click and drag the above window around to better position it for viewing the values. The top few rows of data indicate the currently selected channels' min/max values. Click on any of the ringed values to set the limits for the halfway +/- values for each axis. A numeric touch-pad will appear. Enter the required value and then click anyway else on the viewing area to accept this value, and the graphic will update (if **Drawing scale mode** is set to **Axes: User limits**). For example, to double the size of the drawing in the image below set the min/max values for the three axes to +/- 0.5, rather than the default +/- 1.0. As the data, in the above graphic, is approximately +/- 1.0, it will double the size of the sphere. In the above example, the graphic is coloured by **Temperature**, and the values in the panel are min 23.6C and max 35.7C, which, by the "Jet" colour scheme, is blue through to red. Click the **Manual adjustments of limits** button again to hide this values panel.

Colouring mode

When **Colour by select** is using channel data to colour the data points, it can either use the colour limits as defined by the user in the above **Manual adjustment of limits** for the min/max colour values, or use floating limits i.e. scale the colour to the min/max values in the currently selected data.

Drawing scale mode

There are 3 drawing scale modes to display the data:

1. Auto-normalise/centralise on the origin.
2. By the user-defined limits through the **Manual adjustment of limits** control above.
3. Data defines the limits; centre 0.0,0.0, with grid max at +/- 2.0.

Colourmap

Select the colourmap to use via the mouse-wheel (only Jet available by default). Left mouse-click to switch to dynamic colouring. The graphic will update in real-time while adjusting the colourmap in the colourmap generator app (3rd button from the left on the top left of the viewing area).

Copy current limits to user limits

This copies the values for the X, Y, and Z min/max extremes to the **Set range X/Y/Z** user limits boxes:

2D/3D min/max limits		
Extremes range X	-0. 961	1. 1291
Extremes range Y	-1. 094	1. 1674
Extremes range Z	-1. 109	0. 8886
Extremes range Col	23. 600	35. 700
Set range X	-0. 961	1. 1291
Set range Y	-1. 094	1. 1674
Set range Z	-1. 109	0. 8886
Set Col range	-1. 000	1. 0000

Figure 12.18 Control panel allowing channel limits to be defined or copied from current extremes

View colourbar

This toggles the view of the colourbar, showing the **colour by** channel name, and min/max limits currently being applied.

Marked events only

By default, this is disabled (value 0). Use the mouse-wheel to cycle this through values 1-9 in order to filter by data that has been marked with events values 1-9 either manually, or via the Behaviour Builder facility.



XYZ + Height plot

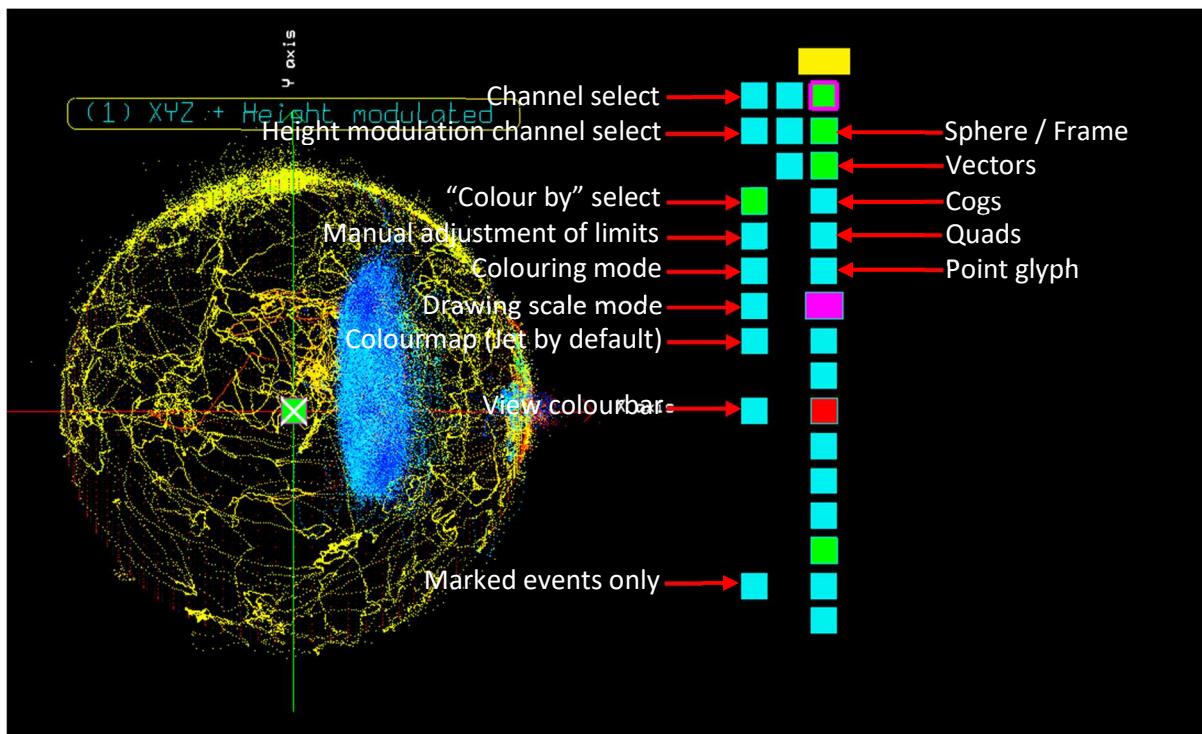


Figure 12.19 XYZ plot data selection and visualisation options

Channel select

Use the mouse-wheel to select the appropriate data channel for the x axis.

Height modulation channel select

Use the mouse-wheel to select the appropriate data channel for the y axis.

Colour by select

Use the mouse-wheel to select the appropriate data channel to colour the data points. Left-clicking this button switches from colouring by channel data, to a single colour; use the mouse-wheel to adjust the actual colour.

Manual adjustment of limits

This manifests a control panel whereby the user can adjust the height channel limit.

Hold down the shift key and mouse click and drag the above window around to better position it for viewing the values. Clicking on the height min/max ringed values and a numeric touch-pad will appear. Enter the required value and then click anyway else on the viewing area to accept this value, and the graphic will update (if **Drawing scale mode** is set to **Axes: User limits**). If the **Drawing scale mode** is set to floating limits, the height modulation will auto-scale to the min/max data values in the height modulation channel selected. Click the **Manual adjustments of limits** button again to hide this values panel.

Colouring mode

When **Colour by select** is using channel data to colour the data points, it can either use the colour limits as defined by the user in the above **Manual adjustment of limits** for the min/max colour values, or use floating limits i.e. scale the colour to the min/max values in the currently selected data.

Drawing scale mode

There are 2 drawing scale modes to display the data:

1. By the user-defined limits through the **Manual adjustment of limits** control above.
2. Floating limits whereby the min/max values within the height modulation channel selected will be used to auto-scale this effect.

**Colourmap**

Select the colourmap to use via the mouse-wheel (only Jet available by default). Left mouse-click to switch to dynamic colouring. The graphic will update in real-time while adjusting the colourmap in the colourmap generator app (3rd button from the left on the top left of the viewing area).

Copy current limits to user limits

This copies the values for the X, Y, and Z min/max extremes to the **Set range X/Y/Z** user limits boxes.

View colourbar

This toggles the view of the colourbar, showing the **colour by** channel name, and min/max limits currently being applied.

Marked events only

By default, this is disabled (value 0). Use the mouse-wheel to cycle this through values 1-9 in order to filter by data that has been marked with events values 1-9 either manually, or via the Behaviour Builder facility.



1D/2D Histogram

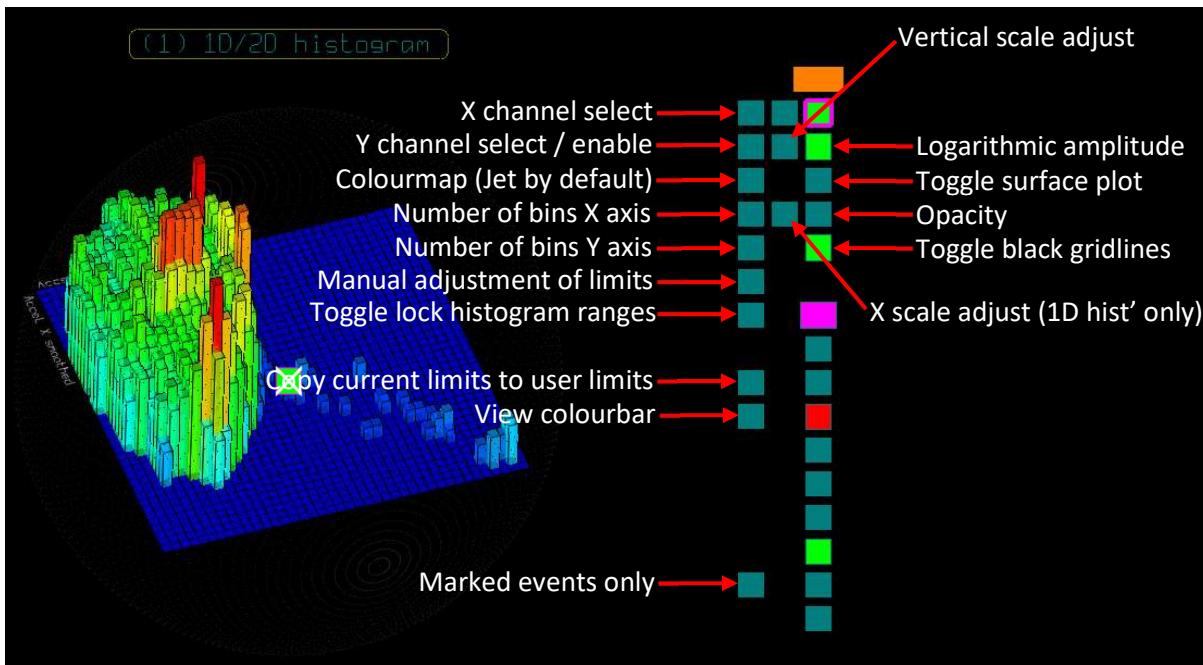


Figure 12.20 Histogram data selection and visualisation options

Vertical scale adjust

Adjust the vertical height of the histogram.

Logarithmic amplitude

Rescale the vertical height of the histogram logarithmically.

Toggle surface plot

Toggle between blocks to a surface plot.

Opacity

Alter the opacity of the graphic.

Toggle black gridlines

Outline the blocks with black lines.

X scale adjust (1D histogram only)

Stretch out the 1D histogram.

X channel select

Use the mouse-wheel to select the appropriate data channel for the x axis.

Y channel select / enable

Use the mouse-wheel to select the appropriate data channel for the y axis; left mouse-click to toggle between 1D and 2D histogram i.e. disable the Y channel.

Colourmap

Select the colourmap to use via the mouse-wheel (only Jet available by default). Left mouse-click to switch to dynamic colouring. The graphic will update in real-time while adjusting the colourmap in the colourmap generator app (3rd button from the left on the top left of the viewing area).

Number of bins X axis

Use the mouse-wheel to adjust the number of bins for the X axis.

Number of bins Y axis

Use the mouse-wheel to adjust the number of bins for the Y axis.

Manual adjustment of limits

This manifests a control panel whereby the user can adjust the X/Y limits of the histogram. Hold down the shift key and mouse click/drag the above window around to better position it for viewing the values. Click on any of the ringed values to set the limits for the halfway +/- values for each axis. A numeric touch-pad will appear. Enter the required value and then click anyway else on the viewing area to accept this value, and the graphic will update (if "Lock histogram ranges" is enabled). Click the **Manual adjustments of limits** button again to hide this values panel.

***Toggle lock histogram ranges***

Linked with the above button; when enabled, the X (and Y if 2D) axes are limited to those min/max values set by the user.

Copy current limits to user limits

This copies the values for the X, and Y, min/max extremes to the ***Set range X/Y*** user limits boxes.


View colourbar

This toggles the view of the colourbar, showing the ***colour by*** channel name, and min/max limits currently being applied.

Marked events only

By default, this is disabled (value 0). Use the mouse-wheel to cycle this through values 1-9 in order to filter by data that has been marked with events values 1-9 either manually, or via the Behaviour Builder.

Spherical Histogram

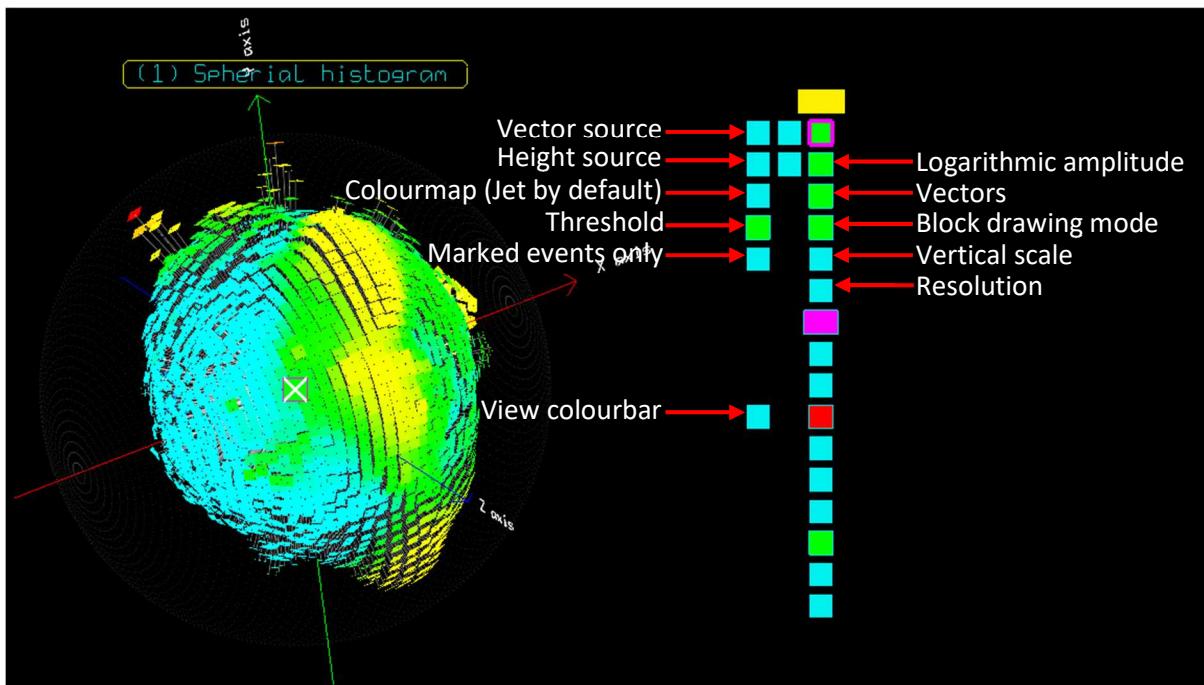



button to disable the threshold, then it switches to a mode whereby any spikes/blocks that fall below the threshold will not be drawn at all.

Marked events only

By default, this is disabled (value 0). Use the mouse-wheel to cycle this through values 1-9 in order to filter by data that has been marked with events values 1-9 either manually, or via the Behaviour Builder.

View colourbar

This toggles the view of the colourbar, showing the **colour by** channel name, and min/max limits currently being applied.



Orientation sphere (o-sphere)

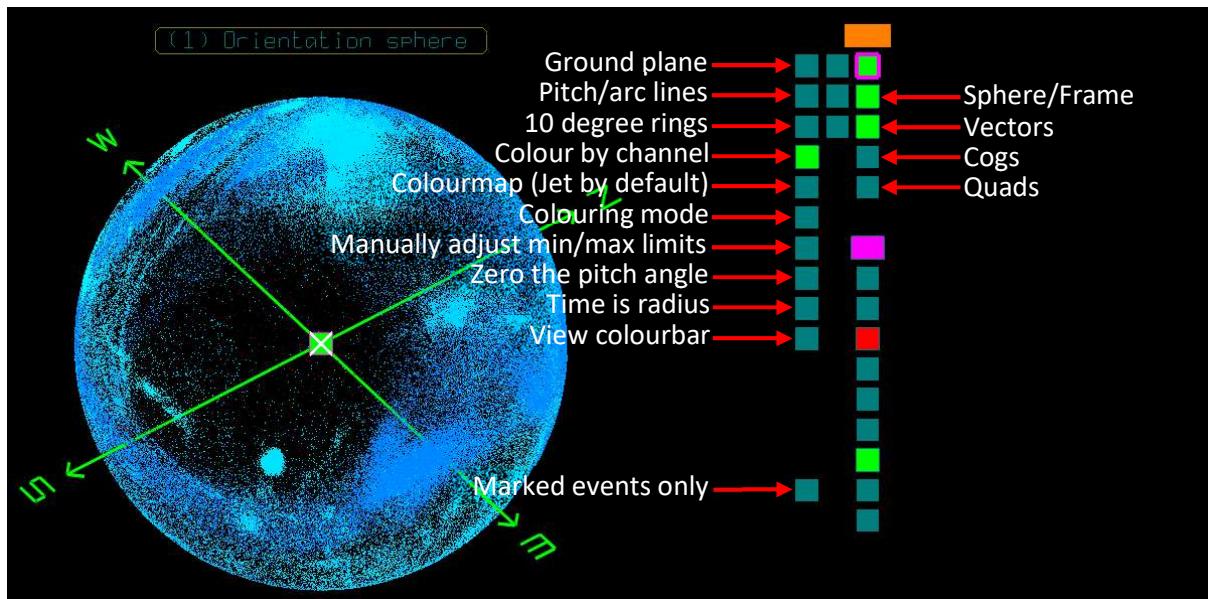


Figure 12.22 Orientation sphere data selection and visualisation options

Ground plane

This simply toggles the presence of a green square beneath the graphic.

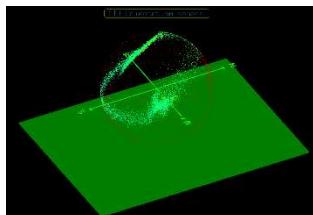


Figure 12.23 Spherical histogram with **ground plane** enabled

Pitch/arc lines

Generates lines from the origin out to the sphere in 10-degree increments from +90 to -90 degrees. The heading vector that these all point at can be adjusted by using the mouse-wheel on this button to rotate through 0-360 degrees.

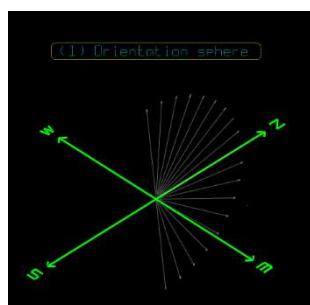


Figure 12.24 Spherical histogram with **pitch/arc lines** enabled



10-degree rings

Generates circles at 10-degree increments from +90 to -90 degrees that draw around the full 360-degree heading.

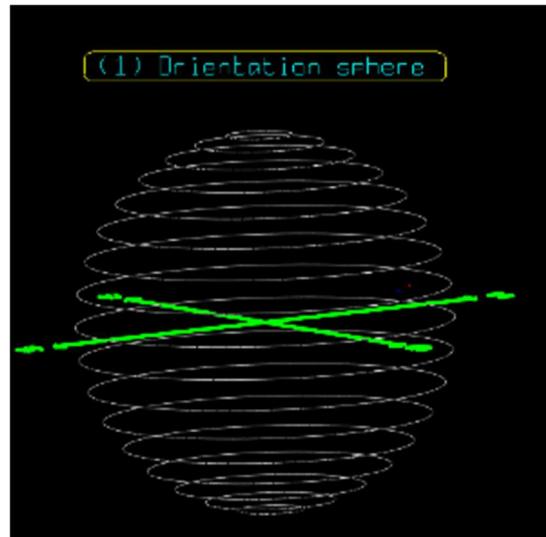


Figure 12.25 Spherical histogram with 10-degree rings enabled

Colour by channel

Select the data channel to be used to colour each data point.

Colourmap

Select the colourmap to use via the mouse-wheel (only Jet available by default). Left mouse-click to switch to dynamic colouring. The graphic will update in real-time while adjusting the colourmap in the colourmap generator app (3rd button from the left on the top left of the viewing area).

Colouring mode

When **Colour by select** is using channel data to colour the data points, it can either use the colour limits as defined by the user in the above **Manual adjustment of limits** for the min/max colour values, or use floating limits i.e. scale the colour to the min/max values in the currently selected data.

Manually adjust min/max limits

This manifests a control panel whereby the user can adjust the min/max limits of the colour channel. Hold down the shift key and mouse click/drag the above window around to better position it for viewing the values. Click on any of the ringed values to set the limits for the colour min/max values. A numeric touch-pad will appear. Enter the required value and then click anyway else on the viewing area to accept this value, and the graphic will update (if the colour mode is set to **user limits**).

Click the **Manual adjustments of limits** button again to hide this values panel.

Zero the pitch angle

This reduces all pitch data to zero so that the visual is simply one of heading only. This button can be used together with the **Time is radius** button below.

Time is radius

If selected, instead of the data always sitting on the surface of a sphere, it instead originates at the origin of the sphere, and radiating out to the sphere surface. This can be used with or without zeroing the pitch angle.

View colourbar

This toggles the view of the colourbar, showing the **colour by** channel name, and min/max limits currently being applied.

Marked events only

By default, this is disabled (value 0). Use the mouse-wheel to cycle this through values 1-9 in order to filter by data that has been marked with events values 1-9 either manually, or via the Behaviour Builder.



Rose plot

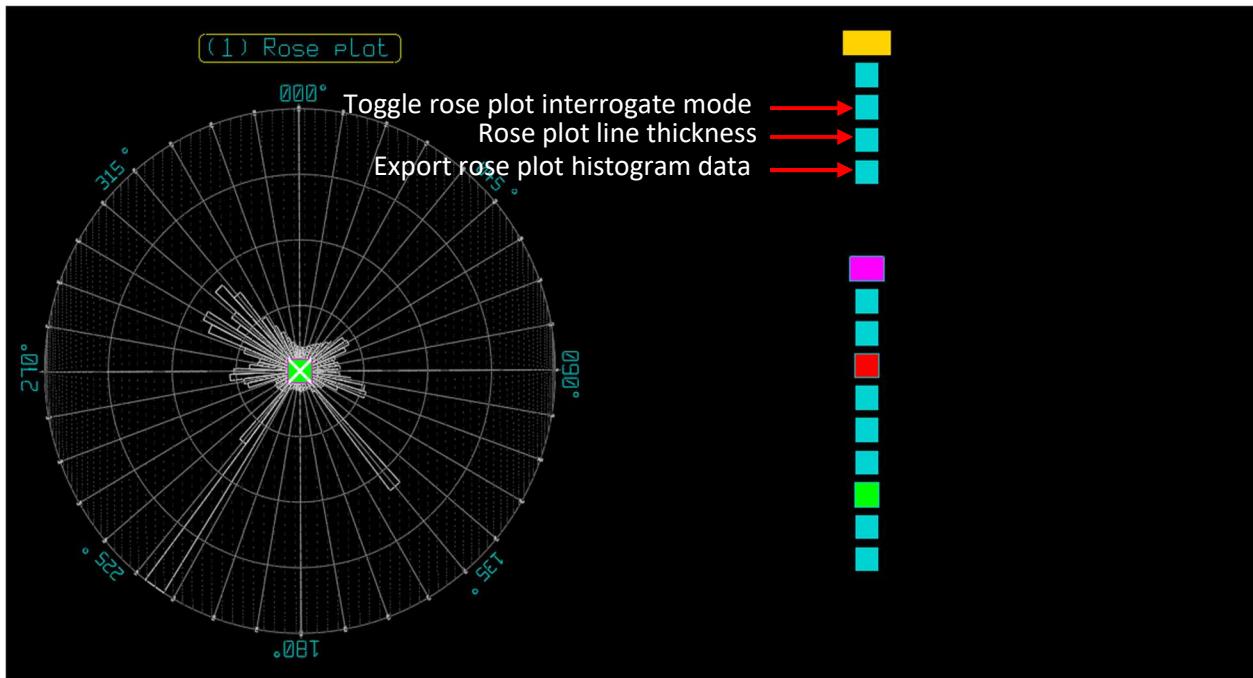


Figure 12.26 Rose plot and associated options

Toggle rose plot interrogate mode

Resets the orientation of the visual. Any mouse movements over the visual then shows the histogram value for that sector.

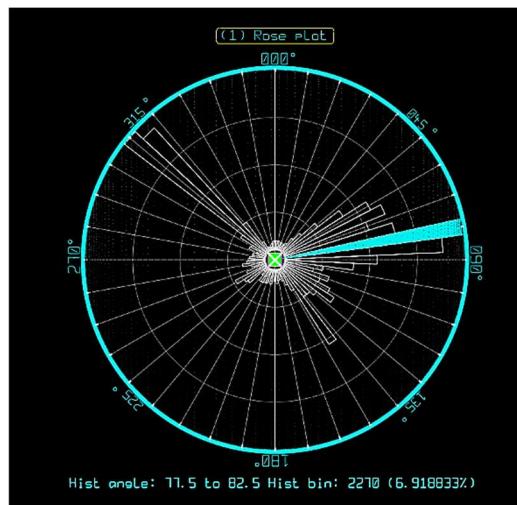


Figure 12.27 Rose plot with interrogate mode enabled

Rose plot line thickness

Simply adjusts the line thickness of the graphic; sometimes necessary for screenshots and image-size reduction.

Export rose plot histogram data

Creates a csv file containing the histogram data by sector.

3D histogram

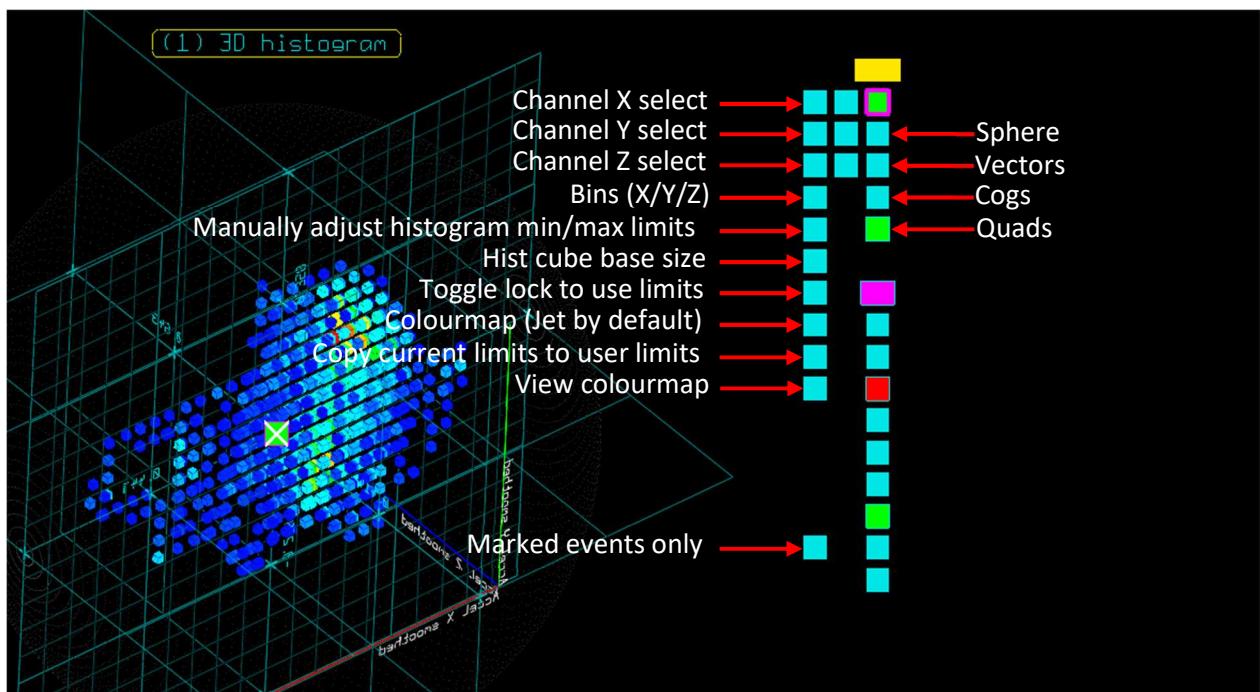


Figure 12.28 3D histogram with data selection and visualisation options

Channel X select

Use the mouse-wheel to select the appropriate data channel for the x axis.

Channel Y select

Use the mouse-wheel to select the appropriate data channel for the y axis.

Channel Z select

Use the mouse-wheel to select the appropriate data channel for the z axis.

Bins (X/Y/Z)

Define the bin resolution of the histogram. All axes have the same number of bins. There is a 50-bin limit for each axis i.e. 50x50x50 (= 125,000) total bins.

Manual adjustment histogram min/max limits

This manifests a control panel whereby the user can adjust the X/Y/Z limits of the histogram data:

2D/3D min/max LIMITS		
Extremes range x	0. 1362	0. 4470
Extremes range y	-0. 229	0. 3588
Extremes range z	0. 4145	0. 6430
Not used	-1. 000	1. 0000
Set range X	-1. 000	1. 0000
Set range Y	-1. 000	1. 0000
Set range Z	-1. 000	1. 0000
Not used	-1. 000	1. 0000

Figure 12.29 Setting min/max limits for 3D histograms

Hold down the shift key and use the mouse to click and drag the above window around to better position it for viewing the values. The top few rows of data indicate the currently selected channels' min/max values. Click on any of the ringed values at the end of this list to set the limits for each axis. A numeric touchpad will appear. Enter the required value and then click anywhere else on the viewing area to accept this value, and the graphic will update (if **Toggle lock to use user limits** is set).

Click **Manual adjustments of limits** again to hide this values panel.

**Hist cube base size**

This defines the base size of the cubes, to which is added the value of the histogram bin. This ensures that small bin data is always visible.

Toggle lock to use limits

This tells the visual to use the user-defined limits as set in **the Manual adjustment histogram min/max limits**.

Colourmap

Select the colourmap to use via the mouse-wheel (only Jet available by default). Left mouse-click to switch to dynamic colouring. The graphic will update in real-time while adjusting the colourmap in the colourmap generator app (3rd button from the left on the top left of the viewing area).

Copy current limits to user limits

This copies the values for the X, Y, and Z min/max extremes to the **Set range X/Y/Z** user limits boxes (as shown in Figure 10.29).

View colourbar

This toggles the view of the colourbar, showing the min/max limits currently being applied.

Marked events only

By default, this is disabled (value 0). Use the mouse-wheel to cycle this through values 1-9 in order to filter by data that has been marked with events values 1-9 either manually, or via the Behaviour Builder facility.



Point compass

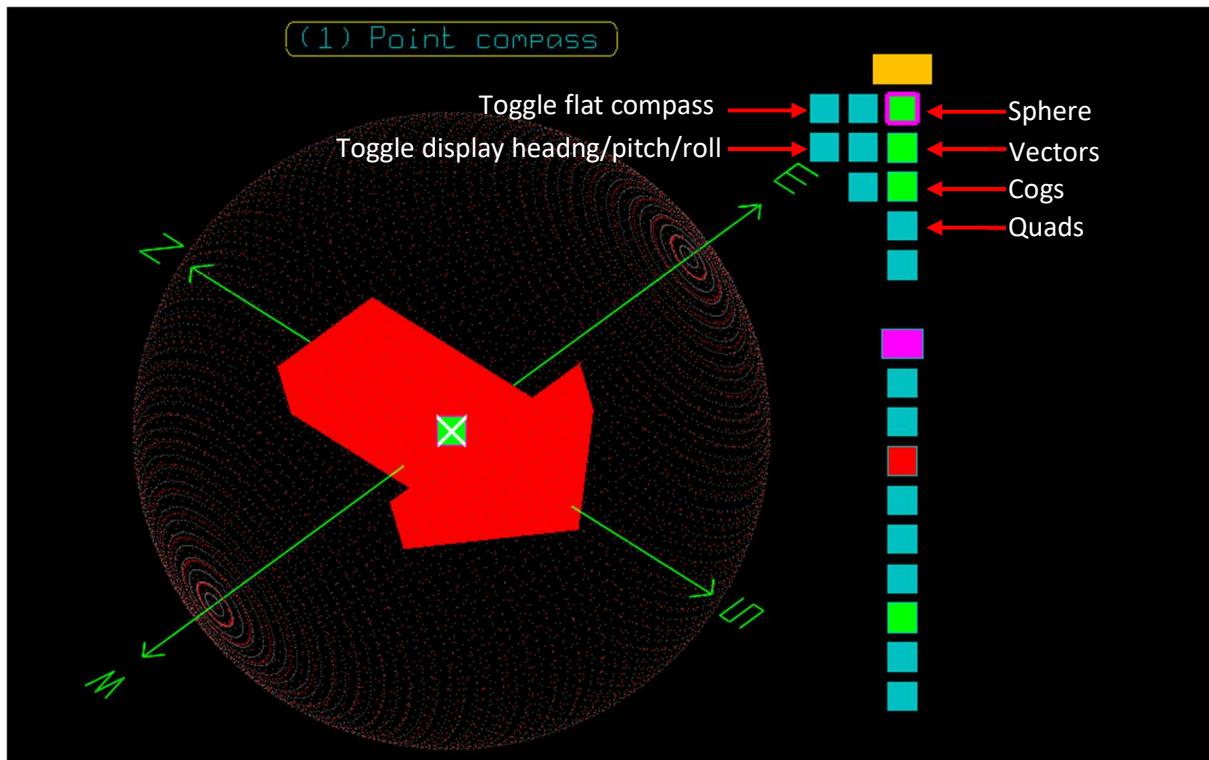


Figure 12.32 Point compass with visualisation options

Toggle flat compass

Normally the compass will tilt up/down by pitch value and rotate by heading. If flat compass is enabled, then pitch is always zero.

Toggle display heading/pitch/roll

Display heading, pitch, and roll values for current left white line position.

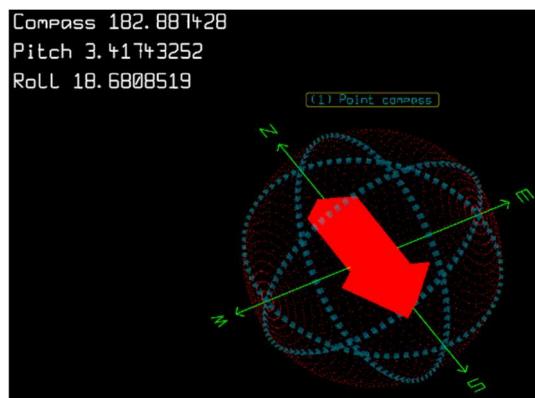


Figure 12.33 Point compass with display values enabled



Dead reckoning plot

The dead reckoning algorithm creates a 2- or 3-dimensional path, recreated using all or some of the sensor data presented to it. Initially, it takes a combination of accelerometry and magnetometry data to derive heading, a technique known as a tilt-compensated compass. Together with heading information, pressure, or some other altimeter measure, can be used to determine the 3rd dimension. The quality of the path generated will be a function of the data and the way in which it is selectively applied by the user. A dead-reckoned path has a tendency to become more erroneous over time due to accumulation of errors in metrics such as speed, or incorrect association of behaviours to particular parts of data such as indicating that an animal is moving in an easterly direction whereas it is in fact stationary, and scratching. GPS data can be synchronised to the dataset to provide intermittent ground truthing, back correcting both range and bearing of data between fixes.

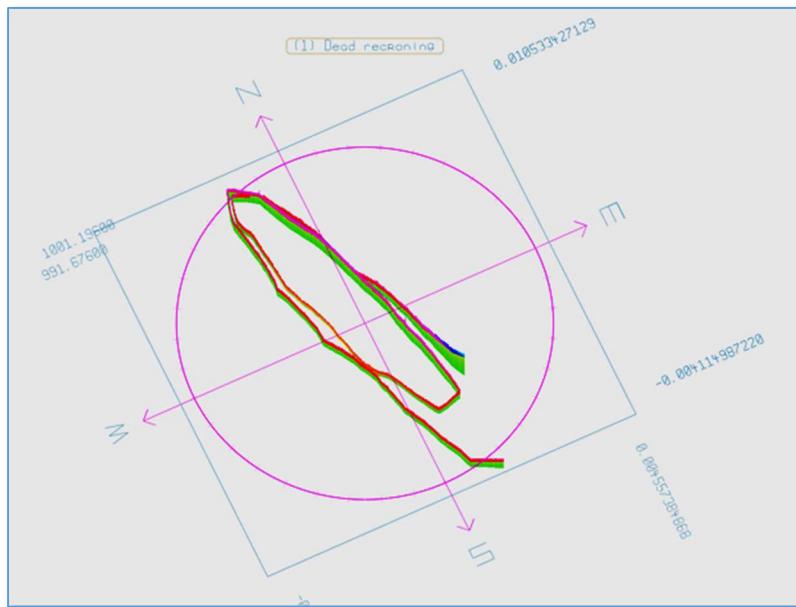


Figure 11.30 Dead reckoning track

In the above image we can see an oval track, circled twice. This was derived from heading data, with pressure (barometric) providing some height information.



The controls for dead reckoning are (3rd column):

		Toggle 3 dimensional (Enabled)
		Colour by altitude/depth/time/marred events (altitude/depth)
		Toggle ground grid (Disabled)
		Toggle solid ground grid, intensity 4 (Disabled)
		Toggle centralise on left white line (Disabled)
		DR using Pressure (MW to toggle invert) not pitch for altitude
		Pitch offset: 0 Deg. (Disabled)
		Vertical scale: 16%
		Toggle show left white line position on trace
		View colourmap
		Toggle show vertical stack: By altitude
		VeDBA limit set: (Disabled) 1.500000 e
		Set VeDBA movement threshold / coefficient
		Switch graph modes (off)
		Dead-reckon by Marked Events only (Disabled)
		Export dead-reckoned trace
		Set initial Long/Lat/Alt coords
		Toggle Long/Lat/Alt coords
		Toggle vectors green ring/cross
		DR by smoothed or raw VeDBA (raw selected)

...and (4th column):

		Auto push DR data to main array
		Clear DR main array
		Overlay GPS onto DR trace (disabled)
		Continuous calc of DR from this point and export to split 1
		Heading correction value: 269
		Altitude in metres (disabled)
		Inherit starting GPS coords (disabled)
		GPS overlay colour
		GPS undersample: 1
		GPS line thickness: 2
		GPS correction (disabled)
		Show VeDBA vs GPS speed graph (disabled)
		Scale bar (non-func)
		Graph Bearing / Range error (disabled)
		Apply wind vector (disabled)
		Dead-reckon using first/last GPS points only (Disabled)
		Ground truth using GPS with Marked Events only (Disabled)

**Definitions for these controls are as follows:****Toggle 3D**

Switches between dead reckoning on a flat plane, or including height/depth

Colour by altitude/depth/time

Toggle between colouring the track based on altitude or depth, or, on time

Toggle ground grid

Toggle view of a ground grid

Toggle solid ground grid

Toggle view of a solid ground grid

Toggle centralise on left white line

Centralise on the position of the left white line in the 2D graphing window

Dead-reckon using pressure, not pitch, for altitude change

Use the smoothed pressure channel to determine depth/altitude, instead of using pitch to determine rate of change of depth/altitude. Mouse wheel inverts the vertical scale

Pitch offset

Apply an offset to the pitch angle to compensate for placement on the animal

Vertical scale

Adjust the vertical scaling of the visual with the mouse wheel

Toggle show left white line as white block

Show the position of the left white line in the track with a white block. Useful when determining where in the 2D data window something occurs in the dead reckoned track

View colourmap

Show the currently selected colourmap (currently redundant; cannot yet apply to the visual)

Toggle show vertical stack

Show the vertical stack – vertical lines from each point down to the lowest in the current track. Adjusting this button with the mouse-wheel switches vertical stack height to be a function of either altitude or VeDBA (smoothed)

VeDBA limit set

When dead reckoning using VeDBA as a proxy for step size, it is useful to provide a limit / cap to this variable. Let click to enable/disable. Mouse wheel to adjust limit

Set VeDBA movement threshold / coefficient

Here the user can check a box to use VeDBA as the step size for movement, and also provide a threshold value (VeDBA) which must be exceeded for movement to occur by an event, and a coefficient that is multiplied by the value of VeDBA that exceeds this threshold
i.e. $(VeDBA - Threshold) * Coefficient$ is the distance moved in metres for that event (data point)

DR VeDBA Threshold/Coefficient
Threshold/Coeff [0.10000000] [0.05000000] <input checked="" type="checkbox"/> By VeDBA coeff?

Alternatively, with “By VeDBA coeff?” deselected, a unit movement of $coefficient * 1.0$ metres occurs for every event (data point) i.e. always moving, unless a data point is not selected due to a **Marked Events** filter (see below)

The user can either use the mouse wheel to adjust the threshold and coefficient, or click on the value and enter a new value using the keypad that appears



Switch graph modes

Clicking this button cycles the visual through having either one of the two graphs above, or none. The two graphs below represent either a linear scale, or proportional to the amount of movement, per unit time on the x axis.

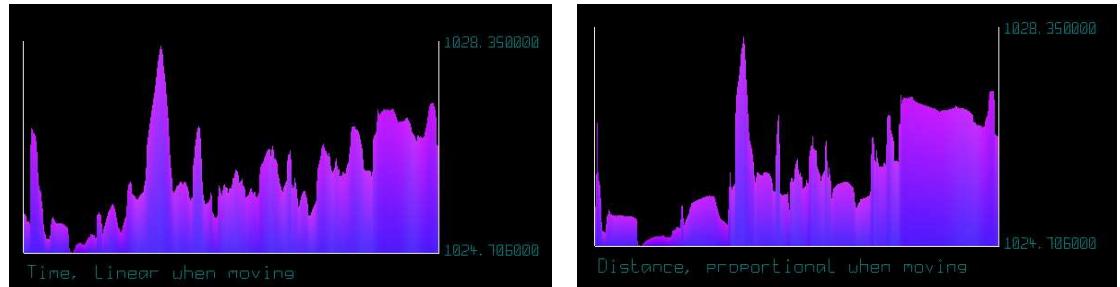


Figure 12.31 Graph mode options, with time on the left and distance on the right

Dead reckon by **Marked events** only

The dead reckoned track is generated from events selected by **Marked Events** only. Enabling this will modify the layout of the DR VeDBA Threshold / Coefficient settings panel:

DR VeDBA Threshold/Coefficient			
Walking	(0.1000000)	(0.1000000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Sleeping	(0.1000000)	(0.0500000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Running	(0.1000000)	(0.2000000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Running fast	(0.1000000)	(0.0500000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Foaging	(0.1000000)	(0.0500000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Category F	(0.1000000)	(0.0500000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Category G	(0.1000000)	(0.0500000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Category H	(0.1000000)	(0.0500000)	<input checked="" type="checkbox"/> By VeDBA coeff?
Category I	(0.1000000)	(0.0500000)	<input checked="" type="checkbox"/> By VeDBA coeff?
(Individual Marked Events settings)			

Here, the **Marked Events** value (1-9) inherits the **Bookmarks** category names (A-I) i.e. Walking, Running, etc. Now that “by **Marked Events**” is selected for the movement step scale, the user can select different threshold and coefficient values for different categories of behaviours. As an example, the user might highlight a selection of data they believe is running and manually set this to **Marked Events** value 3 for *Running*, and another section of data to 1 for *Walking*. Which this **Marked Events** filter enabled, only those events that have **Marked Events** of a non-zero value will contribute to the dead reckoned track. The user would then “adjust up” the coefficient for the *Running* to perhaps 0.2, and “adjust down” the third coefficient to perhaps 0.1. Note that if there is data marked with different values (than 1 or 3 in this example), and the dead reckoned track should not move i.e. the animal is scratching, or sleeping (non-translational events), then set the coefficient to 0 (zero). The dead reckoned path is dynamically updated when a change is made to the track either by manually entering a value, or by adjustment with the mouse wheel

Export dead-reckoned track

Export the date/time/longitude/latitude/altitude data of the dead-reckoned track to a csv file



Set initial long/lat/alt coords

Set the starting longitude/latitude/altitude coordinates of the dead reckoned track. Click the yellow ringed values and use the numerical pad that appears to set each value. Click anywhere else in the 3D window to accept the value

The dialog box is titled "DR Starting Coordinates". It contains three rows of input fields: "Longitude" (two fields), "Latitude" (two fields), and "Altitude" (two fields). Each field has a numerical input area with a small yellow ring and a "Next" button.

Toggle long/lat/alt coords

Toggle the view of the longitude/latitude/altitude min/max limits on the visual

Toggle vectors green ring/cross

Toggle the view of the green ring/central cross on the graphic as these might obscure the track

DR by smoothed or raw VeDBA

Left click to toggle between using raw or smoothed VeDBA as the step size for dead reckoned movement

Auto push DR data to main array

Enable this control to have the algorithm push the dead reckoned track to the system main array so that the DR longitude/latitude etc. is synchronised with the main sensor/timing data and can therefore be exported using the normal save/export functions. Under the **What to save** tab, there is an option to export these pseudo longitude/latitude/altitude data

The screenshot shows the "What to save" tab of the software interface. It lists numerous checkboxes under two categories: "General / Primary channels" and "Additional / metrics channels". A red arrow points to the "Dead Reck Long/Lat/Alt" checkbox in the "Additional / metrics channels" section.

Clear DR main array

Clears the DR data in the main sensor/timing data array. If selected on the **What to save** tab, only zeros will be output

Overlay GPS onto DR trace

If GPS data has been preloaded into DDMT and synchronised timewise to the main array, then selecting this option will show the GPS coordinates within the time window selected in the **2D graphing window**

Continuous calc of DR from this point and export to split 1

This enables a mechanism whereby the DR track on display export up to around 75% of the current split, and then load the next alternate split i.e. either even or odd, auto-generating the DR track, and will continue from 25% into the data, and export to the 75% point, alternating split to split. The last split it will export is that shown on the button label. The user can use the mouse wheel, while hovering on the button, to adjust this up or down. The minimum split value the user can select will be the current split. The data exported is determined by the user selection on the **What to save** tab. Separate files are output per split. These can either be merged later in R, or using the **Merging csv / tab delimited files** function, discussed in the first chapter, **Importing and Exporting Data with DDMT**, to create a single large file for the total DR track

Heading correction value: xxx

This value synchronises with the **Compass offset** value on the **Accel/mag detail/offsets / DD Orientation** tab. It allows the user to manually account for any heading offset in the logger's magnetometer, possibly due to collar fitting issues. Any rotation correction will be redundant when using GPS for ground truthing, as heading errors are automatically corrected when the algorithm rotates DR data between GPS ground truthing fixes. The user



can adjust this value using the mouse wheel while hovering over the button. The dead reckoned track (when not using GPS) will dynamically update

Altitude in metres

When enabled, the altitude values are in metres; disabled, the values are in pressure. The vertical scale, in either case, can be inverted (toggled) using the mouse wheel on the **Dead-reckon using pressure, not pitch, for altitude change control** (adjacent button in 3rd column)

Inherit starting GPS coords

When the user has selected the **Overlay GPS onto DR trace** control, selecting this control will force the dead reckoning to start from the first GPS (left most) fix in time, within the current time-window on the **2D graphing window**. The first GPS coords (in time) will be loaded into the **Set initial long/lat/alt coords** boxes. Dead reckoning continues on from these starting coordinates. This allows the GPS and dead reckoning data to be simultaneously displayed, else their coordinates may be too far apart to be viewed. Note that only the first GPS coordinates are used. Once this has been enabled, it is recommended that the user adjust any heading offsets, VeDBA thresholds/coefficients to achieve a best match of the dead reckoned path with the currently viewed GPS coordinates. Once this is done, the user can then enable the synchronisation of the remainder of the GPS coordinates with the **GPS correction** control below

GPS overlay colour


The user can adjust the colour of the GPS fixes using the mouse wheel

GPS undersample


When using high frequency GPS data, it might be desirable to under-sample the GPS data points i.e. 1 Hz GPS could be under-sampled to a point every 10 seconds, or every 10 minutes, by using the mouse wheel to set the interval in points here

**GPS line thickness** (incl. point size)

The GPS points' size, blue in the image below, are set to the "line thickness +2", so when the visual's points/lines mode is set to lines (2nd column, 3rd button down), the GPS points will be joined up with this line thickness

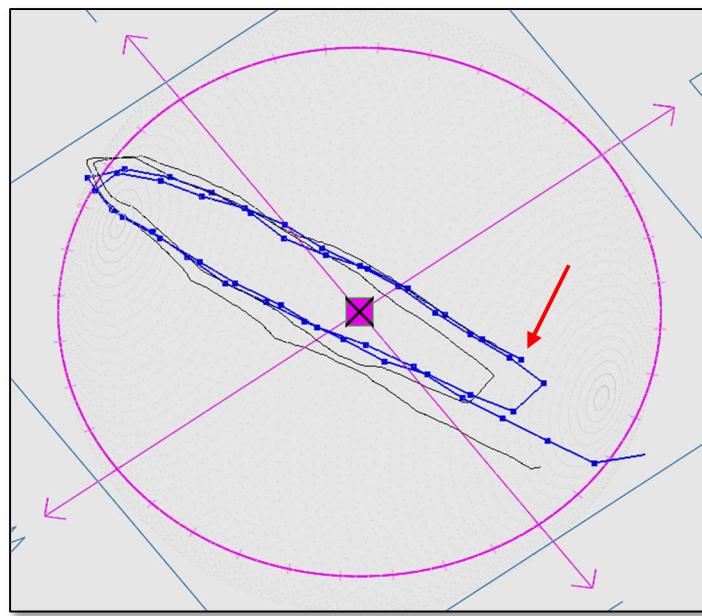**GPS correction**

It may not be obvious in the image above, but the point marked with the red arrow is both the first GPS point, and the start of the DR track, due to the "**Inherit starting GPS coords**" option above. Selecting **GPS correction** reruns the dead reckoning algorithm and uses all available GPS fixes for ground truthing. See image below

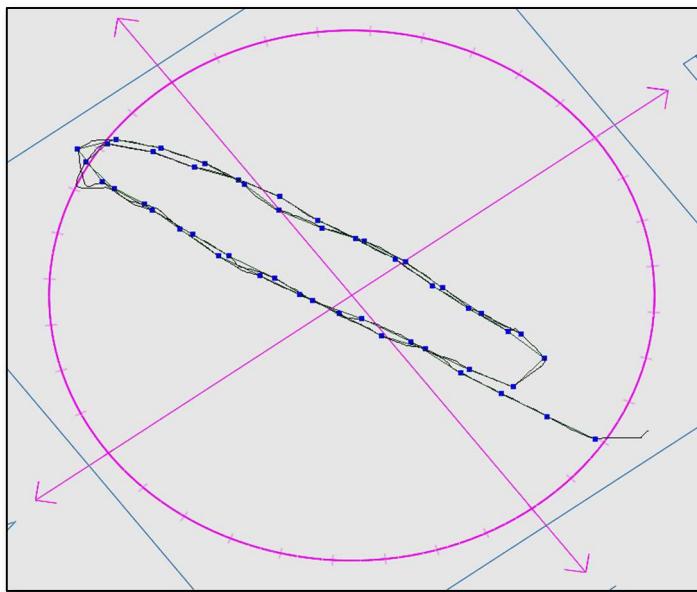**Show VeDBA vs GPS speed graph**

Still in development

Scale bar

Still in development


Graph Bearing / Range error

Still in development

Apply wind vector

This option allows a constant vector to be applied to all data points used in the creation of the dead reckoned track.

The original track, without a wind vector applied, inset, and with the wind vector applied, main:

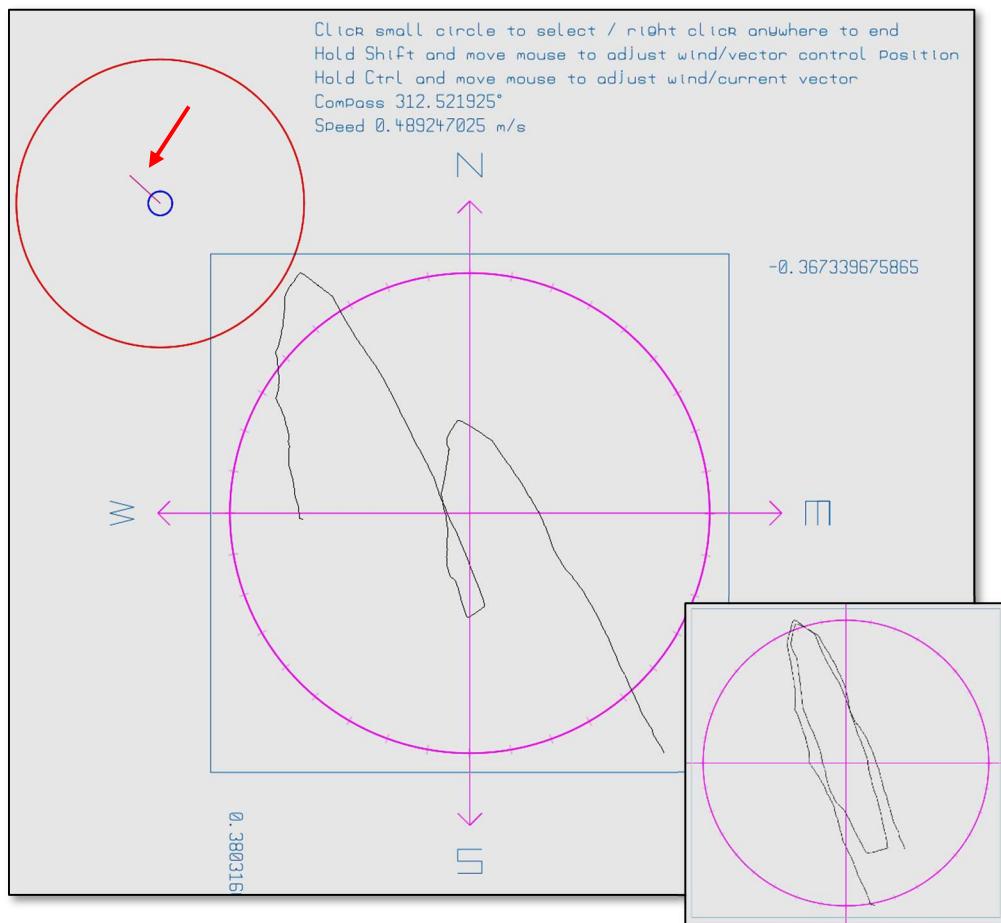

When the wind vector control is enabled, the button's label switches to **Apply wind vector (enabled) Multiplier at :xx** where xx is a coefficient. When enabled, there is no wind vector. The user must first click on the small circle, within the larger circle, top left of the above image. Then, by holding down the CTRL button and moving the mouse around this larger circle, the bearing of the wind is determined, and also the magnitude (multiplied by the button's coefficient) being the length of the vector from the centre of the circle, is the speed in m/s. In the above image (main), a vector has been applied in a northwest direction with a strength of approximately 5 m/s. This has caused the dead reckoned track to stretch in that direction. The dead reckoned track is updated dynamically with any change in the wind vector direction or magnitude. Note that with GPS ground truthing enabled, applying a wind vector may have limited results due to the autocorrection with the GPS fixes.

Dead reckon using first/last GPS points only

Still in development

Ground truth using GPS with Marked Events only

Still in development

Using the dead reckoning algorithm

By default, when generating a dead reckoning visual in DDMT, the user is presented with a first approximation of a track using heading data, and step length per data point will be proportional to the value of VeDBA at each.

Assuming the DD orientation settings are correct to generate heading data, along with pitch and roll from the accelerometer, then, if GPS data has been preloaded and synchronised with the DD data, and the GPS data is believed to be low in positional error, it is recommended to enable the following, to bring GPS ground truthing into play:


GPS overlay



Inherit starting GPS coords (enabled automatically when **GPS overlay** enabled)

Do not enable GPS correction yet (this should be enabled later, once the VeDBA threshold(s) / coefficient(s) have been determined and applied)

The next step would be to either manually apply **Marked Events** values to sections of data such that, perhaps, sections of data where an animal is running would have a **Marked Events** value of 1, walking 2, sleeping or scratching (non-translational) 3 etc. In the image below, the first section has been “marked” with value 1 (red), the second with value 2 (green), and the third value 1 (red). The other data before and after this data is value 0 i.e. “not marked”.

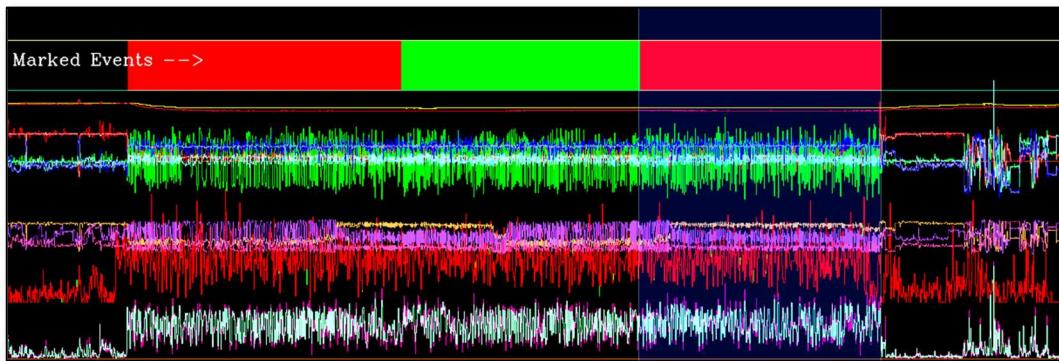

Enabling **Dead reckon by Marked Events only** (3rd column, near the bottom) and selecting **Set VeDBA movement threshold / coefficient**, the user should proceed to adjust the threshold and coefficient values for the various **Marked Events** values. In the settings below, the middle section of data, value 2, I have set the coefficient to 0.005, while the other sections value 1, I have left the coefficient at the default value of 0.05, but depending on the animal under mapping, these values may well be significantly different.


DR VeDBA Threshold/Coefficient		
Category A	0.10000000	0.05000000
Category B	0.10000000	0.05000000
Category C	0.10000000	0.05000000
Category D	0.10000000	0.05000000
Category E	0.10000000	0.05000000
Category F	0.10000000	0.05000000
Category G	0.10000000	0.05000000
Category H	0.10000000	0.05000000
Category I	0.10000000	0.05000000
(Individual Marked Events settings)		

Another useful option to try is to switch the colouring of the dead reckoned track to using the **Marked Events** colouring scheme i.e. 1 = red, 2 = green, 3 = blue etc. by rolling the mouse wheel over the **Colour by...** button (2nd column, 2nd button down) until “marked events” appears in the brackets:






This results in the following track (with the data used here):

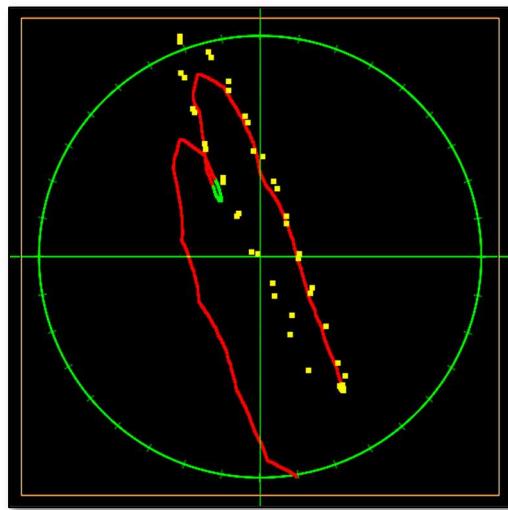

One can see that the dead reckoned track is locked to the first GPS point (bottom right), but from then on, there is no synchronisation, as this is not yet selected. The dead reckoned track is generated as a function of VeDBA, with the first third and last third of the track's VeDBA coefficient set to 0.05, while the middle section has a 1/10th coefficient of 0.005, resulting in a much smaller green track-length.

Enabling **GPS correction** results in the following track, where the middle 3rd of data, while having a low VeDBA coefficient, is compensated for due to the GPS correction algorithm and is stretched to approximately recreate the original track (slightly distorted due to the highly under-sampled GPS fix array).

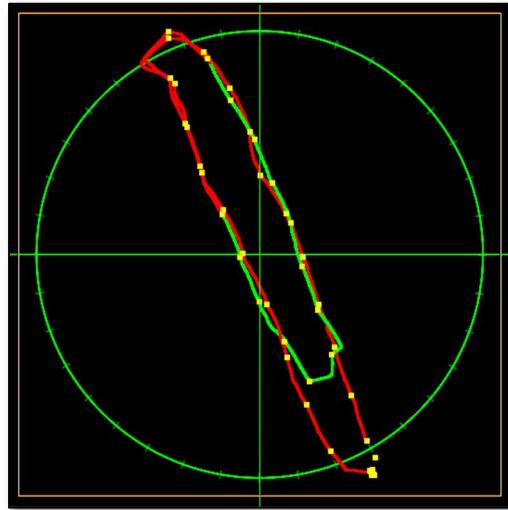



Colourmap Generator



Store current map

Stores the current map using the set name (default is **User defined name**) and makes it available to all visualisation types.

Import colour map

Imports a .csv file containing colour map(s) and adds them to the list for all visualisation types.

Export colour map

Exports all currently in-memory colour maps to a .csv file.

Set map name

Reveals a virtual keyboard that the user can use to enter the name for the current dynamic colourmap.

Show custom maps

Under development



Figure 12.34 Creating a custom colourmap

Notice the 7x triple buttons. These are the red/green/blue components (top/middle/bottom, respectively). Positioning the mouse over one of these many smaller buttons will reveal a small tag showing the red/green/blue values for that “slider”. The mouse-wheel can be used to adjust the value for these 3 component colours (0-255). Left click on any of the three small buttons for each of the middle 5 sliders, and its position will then track with the mouse. Left click again, and the slider’s position will be released. In the above, several sliders have been adjusted horizontally.

Important note – In order to interact with the colourmap generator, disconnect from all visualisations by right-clicking on nothing in the 3D area i.e. turn a currently selected visual (green square at the centre of any currently generated visualisation to red)